This product (collection kit in combination with the authorized test) has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by the authorized laboratories.

This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

For use by individuals age 18 years or older (self-collected), 14 years and older (self-collected under adult supervision), or 2 years and older (collected with adult assistance). For Prescription Use Only.

IVD	2°C		8°C	EUA ONLY
		-		

Distributed by DolNeedaCOVID19Test.com If you have questions during your specimen collection, please contact COVID19@eTrueNorth.com or Call (800) 635-8611.

DolNeedaCOVID19Test.com

Self-Collection Kit

INSTRUCTIONS FOR USE

PRECAUTION - TIME SENSITIVE TEST - READ BELOW BEFORE STARTING

- Follow the sample and packaging instructions closely to prevent inaccurate results or rejected samples. *******PLEASE ACTIVATE YOUR KIT******* To activate, collect and package your sample, please visit: www.DolNeedaCOVID19Test.com
- To get a valid diagnostic test result, it is important to return your sample to the store for pickup on the SAME DAY OF COLLECTION.


DolNeedaCOVID19Test.com Self-Collection Kit

* **PRECAUTION** - TIME SENSITIVE TEST - READ ALL INSTRUCTIONS BEFORE STARTING For use by individuals age 18 years or older (self-collected), 14 years and older (self-collected under adult supervision, i.e., individuals 14-17 years old and adults that require supervision), or 2 years and older (collected with adult assistance, i.e., individuals 2-13 years old and adults that require assistance).

TEST KIT CONTENTS

1

2

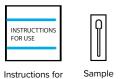

swab

Use Manual



Tube

label






Saline filled Biohazard bag with absorbent pad



Activate your kit before starting the sample collection. Scan the QR code above or go to: www.DolNeedaCOVID19Test.com. You MUST log on and enter the tube number located below the barcode on the sticker to receive a test result. If you are under the age of 18 please ask your parent or guardian for help with this kit. Children 14-17 can self-collect under adult supervision, while children ages 2-13 require adult assistance.




Place kit contents on a clean, dry surface with enough room to collect your sample.

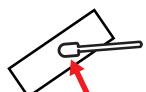



Wash your hands for at least 20 seconds and drv.

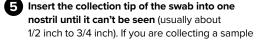


from a younger child, insert swab less than 1/2 inch). Rotate swab firmly against the inside wall of the nostril 4 times taking about 15 seconds to collect the sample.

Place same tip into liquid

using the same tip.

Use same tip

Repeat step 5 in your other nostril,

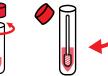


6



Open pre-filled saline collection tube provided and place sample swab in tube, tip first into the liquid. Screw cap back on.



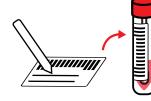

Wash your hands for at least 20 seconds and dry.



Place tube into biohazard bag with absorbent pad. Seal the bag.

Using a pen, neatly write your last, first name and date of birth on the tube label. Peel off the sticker back and place the label length-wise on the tube.



On the day you take the test, return it to the same location where you picked it up.

Follow location instructions for dropoff (i.e., hand to location staff, place in collection bin, etc.). Be sure to check the retail location's current hours of operation.

Remove sample swab from its packaging. Do not touch the tip and do not set the swab down.

*FAILURE TO DROP OFF THE SAMPLE SAME DAY OF COLLECTION WILL RESULT IN SAMPLE REJECTION

This product (collection kit in combination with the authorized test) has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA for use by the authorized laboratories.

This product has been authorized only for the detection of nucleic acid from SARS-CoV-2, not for any other viruses or pathogens; and

The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19 under Section 564(b)(1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb-3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

For use by individuals age 18 years or older (self-collected), 14 years and older (self-collected under adult supervision), or 2 years and older (collected with adult assistance). For Prescription Use Only.

IVD	2°C	P	8°C	EUA ONLY
		\sim		

Distributed by DolNeedaCOVID19Test.com If you have questions during your specimen collection, please contact COVID19@eTrueNorth.com or Call (800) 635-8611.

DolNeedaCOVID19Test.com

Self-Collection Kit

INSTRUCTIONS FOR USE

PRECAUTION - TIME SENSITIVE TEST - READ BELOW BEFORE STARTING

- Follow the sample and packaging instructions closely to prevent inaccurate results or rejected samples.
 *****PLEASE ACTIVATE YOUT KIT*****
 To activate, collect and package your sample, please visit:
 www.DolNeedaCOVID19Test.com
- To get a valid diagnostic test result, shipment of your sample must take place on the SAME DAY OF COLLECTION.
- Return your sample to a USPS Blue Collection Box, USPS Post Office, or you may schedule a pickup online at www.USPS.com and use your home mailbox. Be sure to check the USPS location's current hours of operation.

DolNeedaCOVID19Test.com Self-Collection Kit

* **PRECAUTION** - TIME SENSITIVE TEST - READ ALL INSTRUCTIONS BEFORE STARTING For use by individuals age 18 years or older (self-collected), 14 years and older (self-collected under adult supervision, i.e., individuals 14-17 years old and adults that require supervision), or 2 years and




Activate your kit before starting the sample collection. Scan the QR code above or go to: www.DolNeedaCOVID19Test.com. You MUST log on and enter the tube number located below the barcode on the sticker to receive a test result. If you are under the age of 18 please ask your parent or guardian for help with this kit. Children 14-17 can self-collect under adult supervision, while children ages 2-13 require adult assistance.

5

SCAN ME



2 Place kit contents on a clean, dry surface with enough room to collect your sample.

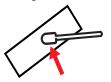

4 Remove sample swab from its packaging. Do not touch the tip and do not set the swab down.



3 Wash your hands for at least 20 seconds and dry.



Insert the collection tip of the swab into one nostril until it can't be seen (usually about 1/2 inch to 3/4 inch). If you are collecting a sample from a younger child, insert swab less than 1/2 inch). Rotate swab firmly against the inside wall of the nostril 4 times taking about 15 seconds to collect the sample.



6 Repeat step 5 in your other nostril, using the same tip.




Wash your hands for at least 20 seconds and dry.



Place same tip

into liquid


Open pre-filled saline collection tube provided and place sample swab in tube, tip first into the liquid. Screw cap back on.

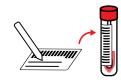

Using a pen, neatly write your last, first name and date of birth on the tube label. Peel off the sticker back and place the label length-wise on the tube.



Place tube into biohazard bag with absorbent pad. Seal the bag.

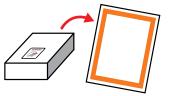

Place lab sample box into shipping envelope and seal.



Place biohazard bag into lab sample box and insert tabs to seal box.





On the day you take the test, place in an outbound USPS mailbox. Acceptable outbound mailboxes include: USPS Blue Collection Box, USPS Post Office or you may schedule a pickup at www.USPS.com and use your home mailbox. Be sure to check the USPS location's current hours of operation.

*FAILURE TO DROP OFF THE SAMPLE SAME DAY OF COLLECTION WILL RESULT IN SAMPLE REJECTION