




Quality in Mammography

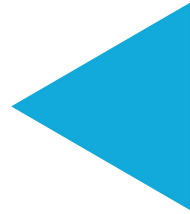


Julian Marshall

Chief Marketing Officer
Volpara Solutions, Inc.

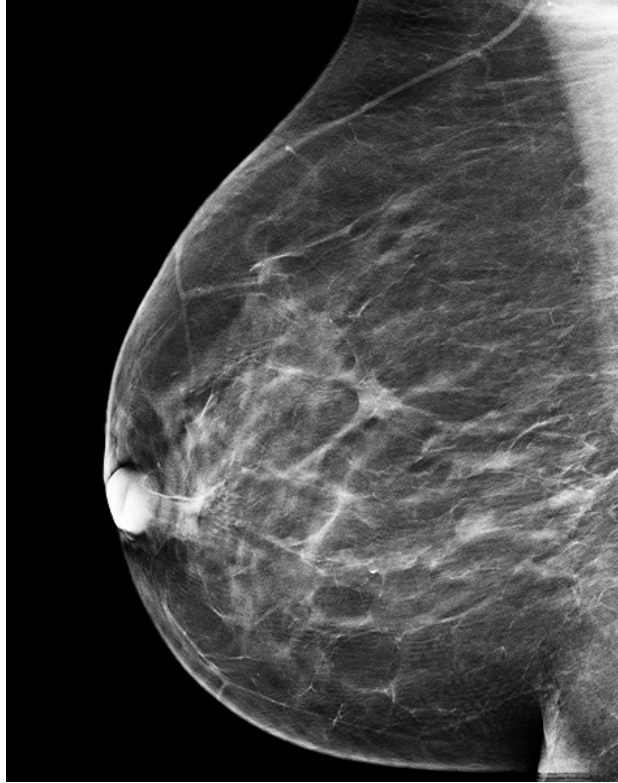

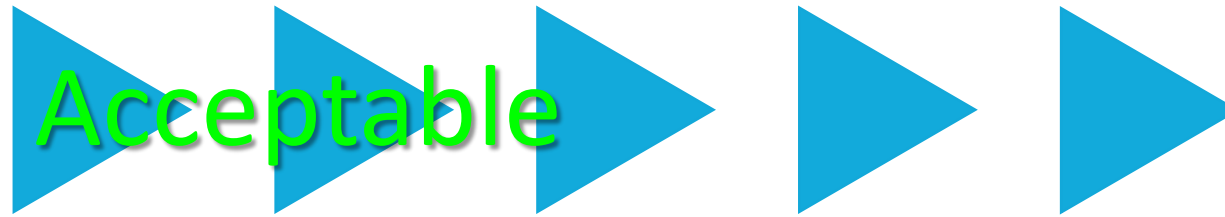


- PACS
- Workstations
- Density
- CAD

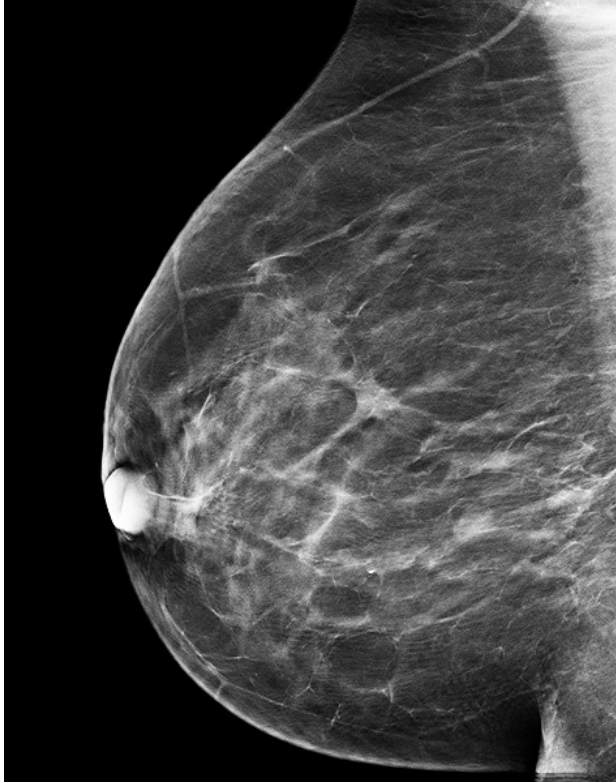


Not Acceptable

Reject Reason



- Typical gantry cannot send rejected images to a different destination
- Digital mammography has removed our ability to retrospectively analyze rejected images:
 - Were the rejections legitimate?
 - Is re-training called for?
- We must change that

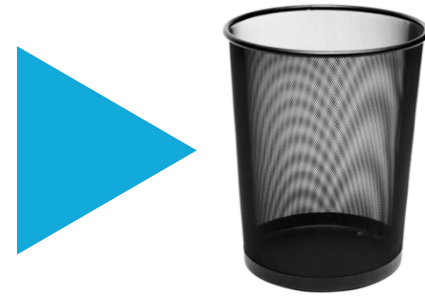



Training

40 hr didactic + 25 supervised studies

On-going CE

15 units every 36 mo

Triennial Accreditation

Their skill is judged

Technologist



Weekly QC

Bi-annual Physicist Visit

Triennial Accreditation

Mammo System





Technologist

Training

40 hr didactic + 25 supervised studies

On-going CE

15 units every 36 mo

Triennial Accreditation

Their skill is judged



Mammo System

Weekly QC

Bi-annual Physicist Visit

Triennial Accreditation

The current process
relies on “spot checks”





Technologist

Training

40 hr didactic + 25 supervised studies

On-going CE

15 units every 36 mo

Triennial Accreditation

Their skill is judged



Mammo System

Weekly QC

Bi-annual Physicist Visit

Triennial Accreditation

If **spot checks** are sufficient,
why do sites
fail accreditation?



Poor Positioning Responsible For Most Clinical Image Deficiencies, Failures

Mammography combines "the science of imaging and the art of positioning" [1]. Although there have been many significant and exciting changes to the technology of mammography since the passage of MQSA in 1999, including the introduction of full-field digital mammography (FFDM) and digital breast tomosynthesis (DBT), one aspect of mammography that remains unchanged and critically important is proper patient positioning.

Positioning is so important because only those portions of the breast which are included on the mammographic image can be evaluated for signs of cancer. Any portion of the breast which is not imaged cannot be evaluated, and cancers in those portions of the breast can be missed. In a 2002 study, the "[s]ensitivity [of mammography] dropped from 84.4% among cases with passing positioning to 66.3% among cases with failed positioning" [2].

Poor positioning has been found to be the cause of most clinical image deficiencies and most failures of accreditation. In 2015, the American College of Radiology (ACR), the largest FDA-approved accreditation body (AB), found that of all clinical images which were deficient on the first attempt at accreditation, 92% were deficient in positioning. Also, in ACR-accredited facilities, 79% of all unit accreditation failures in 2015 were due to positioning. Similar results were noted by the Iowa and Texas state ABs: in 2015, positioning was a cause of 91% of clinical image failures in Iowa and 100% of clinical



92% of 1st attempt
ACR (2015): clinical image deficiencies
were due to positioning

*[images submitted were the result of
the facility selecting high quality images
specifically for accreditation]*






Constant Quality Metrics

- Assess positioning of every mammography image
- Provide daily statistics to chief technologist to drive additional training, as needed
- Form basis for new, stronger quality programs




Image Positioning Score per Tech

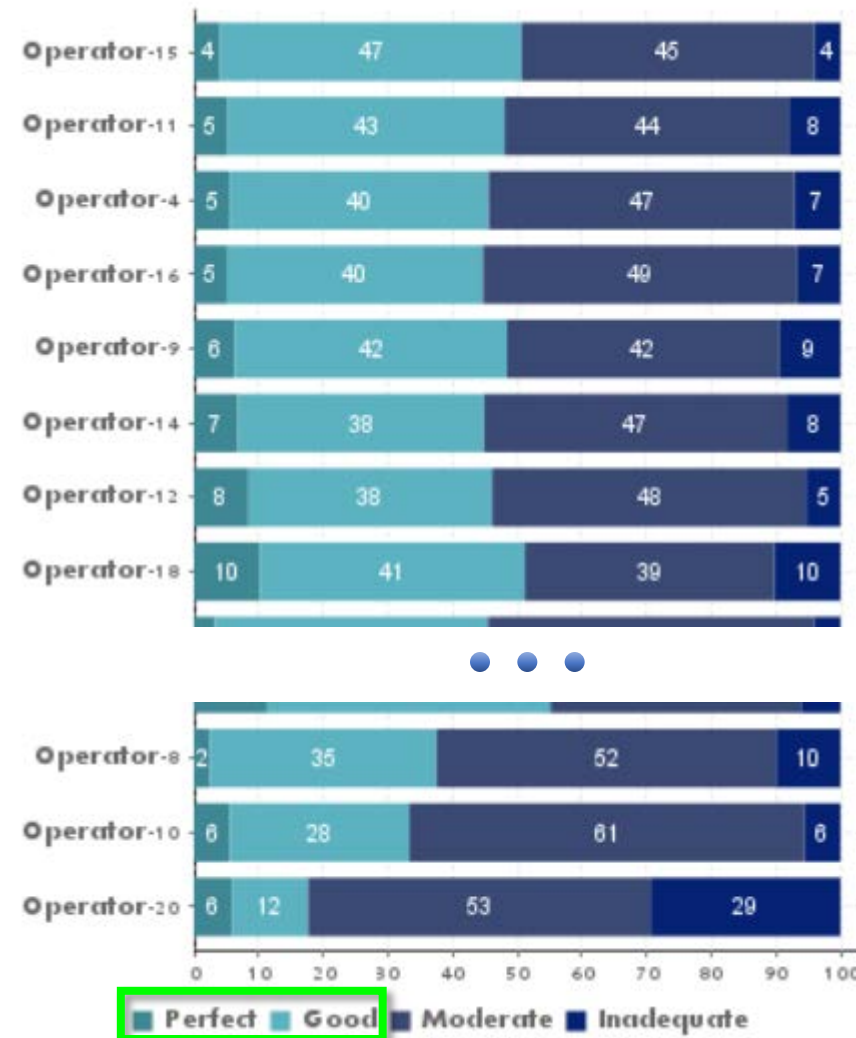

Limited Diagnostic Value

Image Positioning Score per Tech

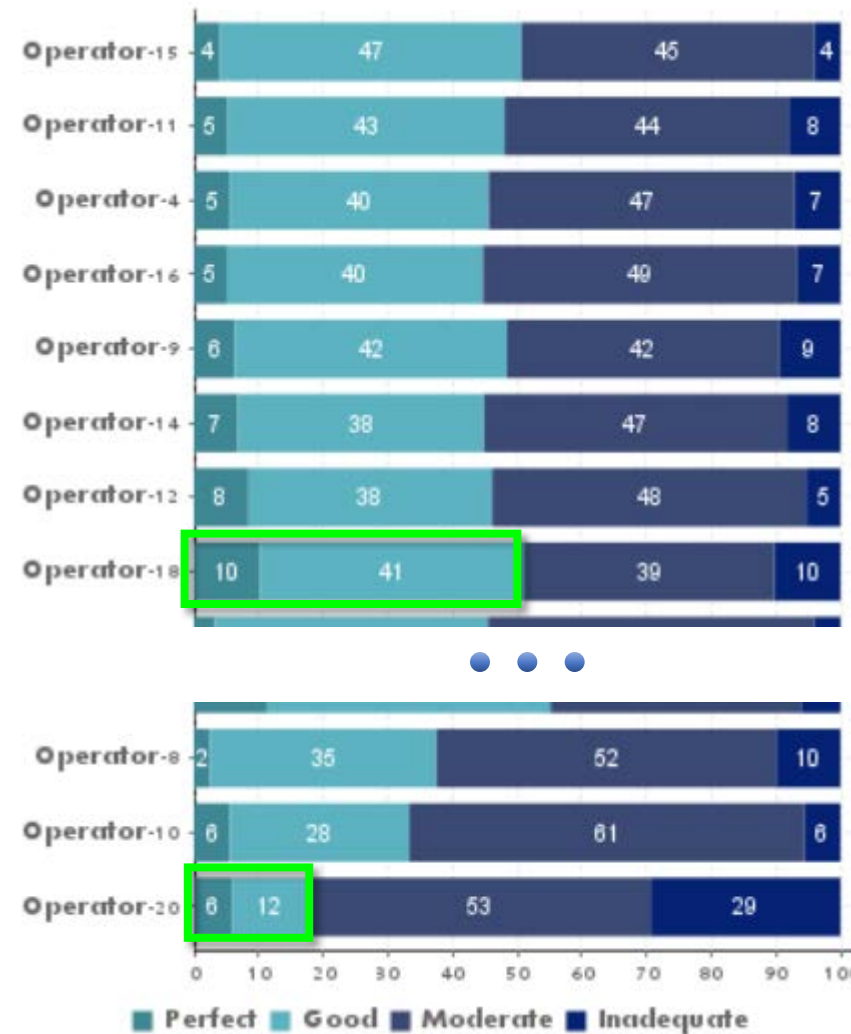

51%

Good or Perfect

18%

Good or Perfect

Significant
operator
variability

Learning from Issues


Its critical to understand **why** an image or study is **Inadequate**

Looking at statistics and trends lets us re-train techs on the **exact issues** with which they struggle

Performed on 7/23/2016 10:23:38 AM
Performed by Operator-1
Quality Score: Inadequate
Nipple is not in profile.

volpara[®]enterprise

Image Result
Anonymized Patient ID: 2027

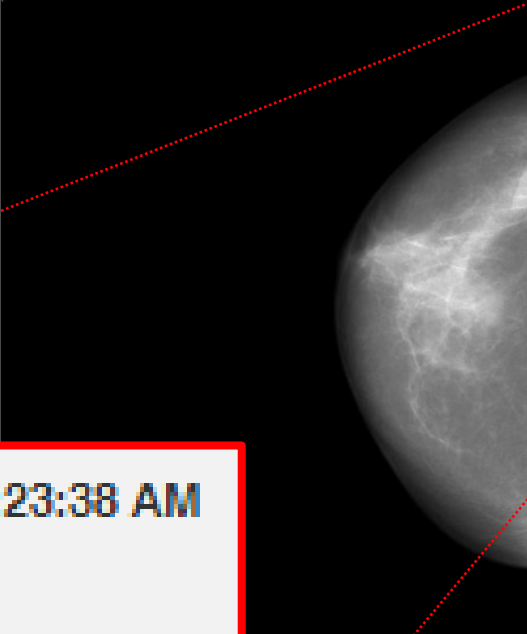

Performed on 7/23/2016 10:23:38 AM
Performed by Operator-1
Quality Score: Inadequate
Nipple is not in profile.
Nipple to chestwall: 65.5 mm
Breast Density: 8.8 %
Breast Volume: 318.3 cm³
Breast Thickness: 43.7 mm
Compression Pressure: 34.2 kPa
Compression Force: 160 N
Volpara Dose: 1.2 mGy
Tube Voltage: 28 kVp
HVL: 0.4 mm
Exposure: 54 mAs
Filter Material: Rhodium
Anode Target Material: Rhodium


For diagnostic use or density assessment.

© 2016. Volpara Health Technologies Limited. All rights reserved.


Remote System Monitoring

[Dose](#)
[Dose Search](#)
[Radiation Parameters](#)
[Compression](#)
[Machine Monitoring](#)
[EUREF Dose Audit](#)
[Export Table](#)
[Image Quality](#)


Tube Voltage, Total Exposure, HVL and Entrance Dose

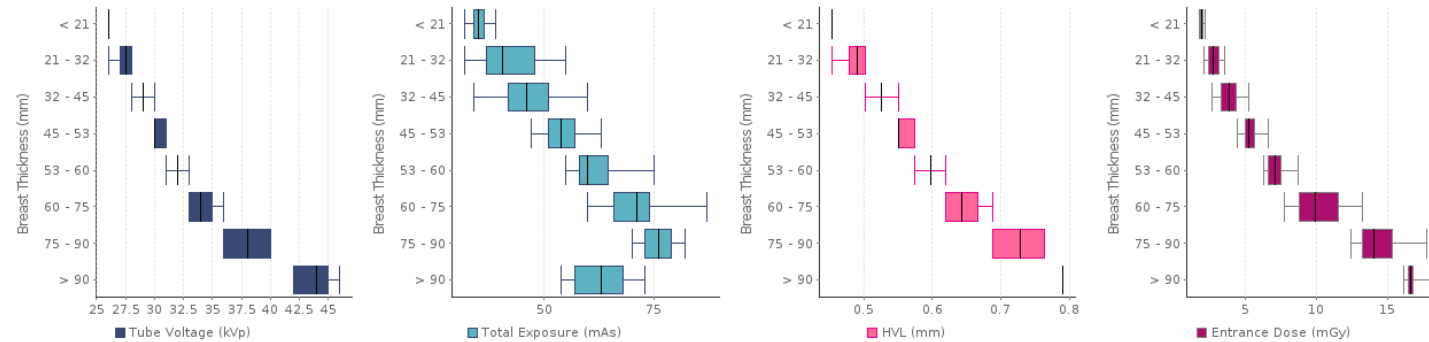

Anode/Filter Combination

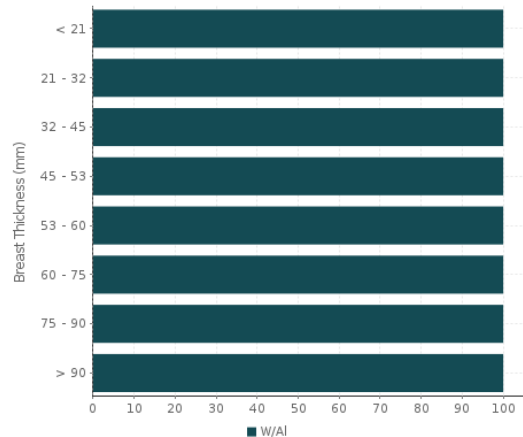

VolparaDose vs. Breast Thickness

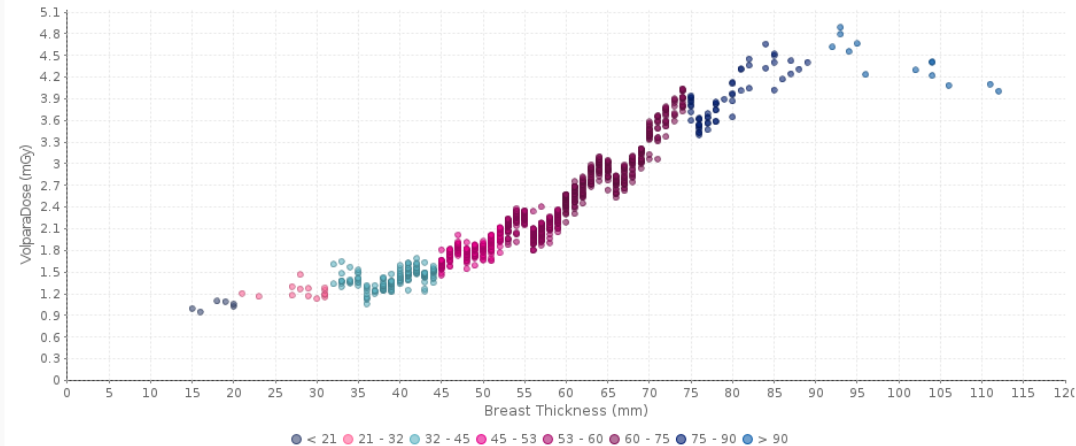


Dynamic monitoring of equipment performance:

- Tube voltage
- Total exposure
- HVL
- Entrance dose
- Anode/filter combination
- Dose vs. breast thickness


Discovery of an Issue

VolparaDose by Mammography System

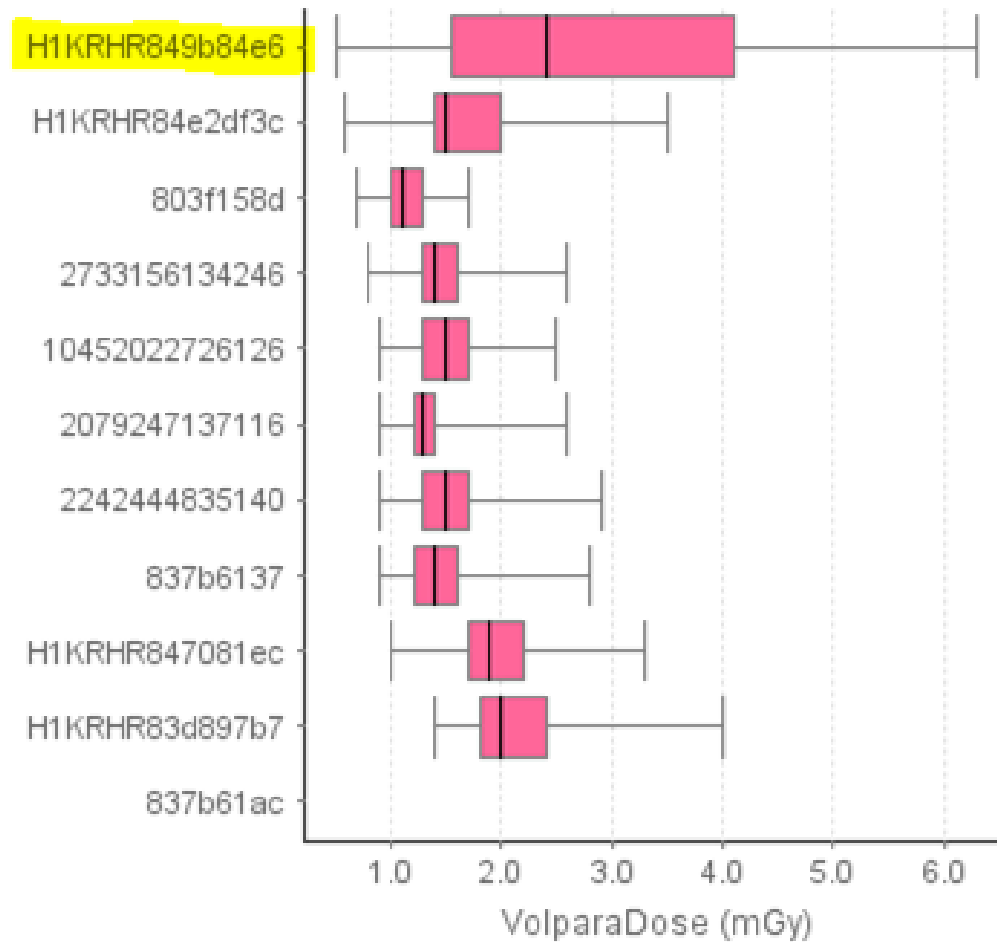

Discovery of an Issue

Average VolparaDose vs. Breast Thickness

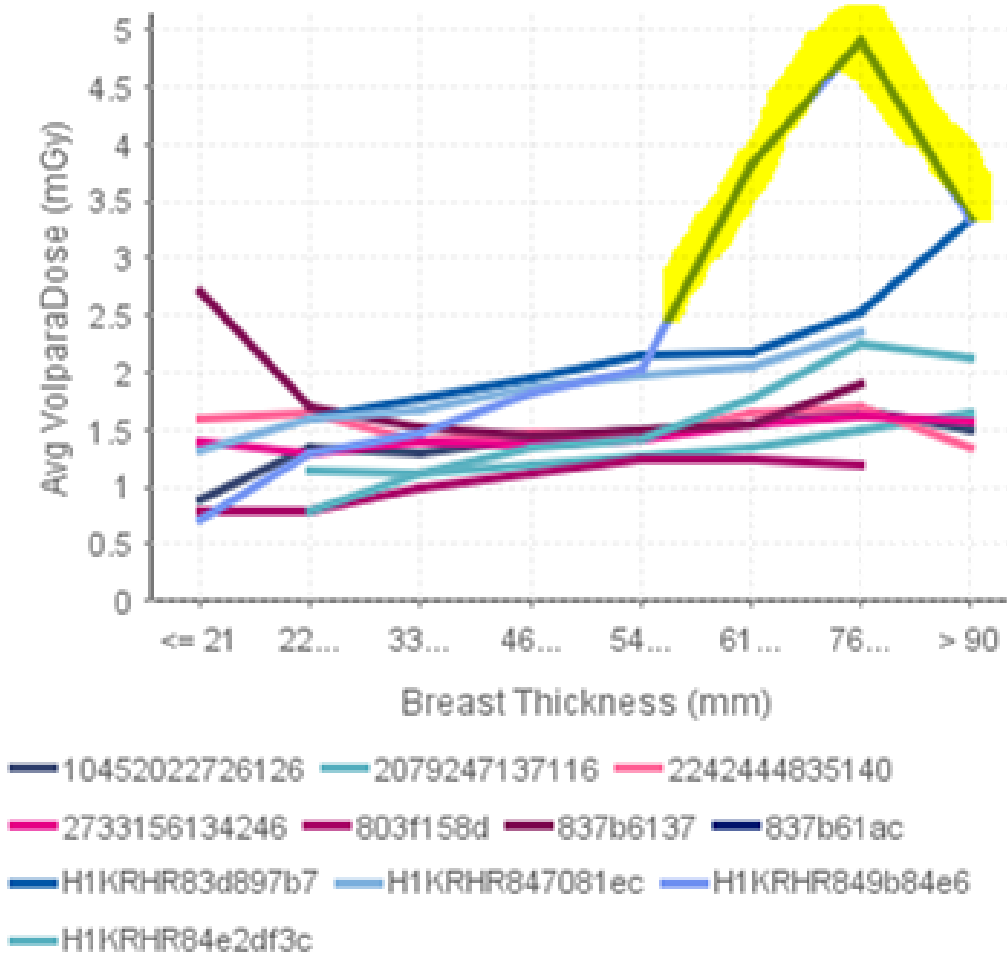

Solution?

- Recalibrate system
- HVL incorrectly entered

In the end:

- Patients doses were normal
- Reported dose was excessive

Demand ConstantQuality!

- Every woman deserves a diagnostic quality mammogram
 - Why wait for the next “spot check”?
 - Check every mammogram and every system – every time!



Thank You

