

MAY 23 2000

K001024

510(k) Summary of Safety and Effectiveness

- (1) **Submitter's name:** Encore Orthopedics, Inc.
Submitter's address: 9800 Metric Blvd, Austin, TX 78758
Submitter's telephone number: (512) 834-6237
Contact person: Debbie De Los Santos
Date summary prepared: March 29, 2000
- (2) **Trade or proprietary device name:** PASS Spinal System
Common or usual name: Pedicle screw spinal system
Classification name: Class II
- (3) **Legally marketed predicate device:** ISOBAR (Scient'x) and Mirage (Alphatec)
- (4) **Subject device description:**

The PASS Spinal System consists of pedicle screws, rods, nuts and crosslink members. It can be used for single or multiple level fixation. All components are manufactured from titanium alloy (Ti-6Al-4V) that conforms to ASTM F136.

This system consists of pedicle screws and clamps that allow the rod to be offset from the screw.

- (5) **Subject device intended use:**

The PASS Spinal System consists of pedicle screws, rods, nuts and crosslink members utilized to provide immobilization and stabilization of spinal segments in skeletally mature patients as an adjunct to fusion in the treatment of the following acute and chronic instabilities or deformities of thoracic, lumbar, and sacral spine: degenerative spondylolisthesis with objective evidence of neurological impairment, fracture, dislocation, scoliosis, kyphosis, spinal tumor, and failed previous fusion (pseudarthrosis).

As a pedicle screw system, the PASS Spinal System is intended for the treatment of severe spondylolisthesis (Grade 3 and 4) of the L5-S1 vertebrae in skeletally mature patients receiving fusion by autogenous bone graft having implants attached to the lumbar and sacral spine (L3 to sacrum) with removal of the implants after the attainment of a solid fusion.

- (6) **Performance data:**
The Food and Drug Administration have established no performance standards applicable to pedicle screw spinal systems. However, static and fatigue compression testing of the PASS Screws with Crosslinking was performed according to ASTM F1717-96.
- (7) **Basis for substantial equivalence:**

The PASS Spinal System has similar design characteristics, i.e., material, screw size, lateral screw/rod connection, crosslinking connectors and indications, as the ISOBAR (K991889/K000020) and Mirage (K951846) systems distributed by Scient'x and Alphatec Manufacturing, respectively.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

MAY 23 2000

Ms. Joanna Monahon
Regulatory Specialist
Encore Orthopedics, Inc.
9800 Metric Boulevard
Austin, Texas 78758

Re: K001024/S1
Trade Name: PASS Spinal System
Regulatory Class: II
Product Code: MNH and MNI
Dated: May 15, 2000
Received: May 17, 2000

Dear Ms. Monahon:

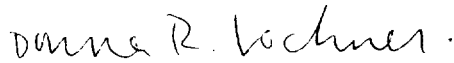
We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general control provisions of the Act. The general control provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the current Good Manufacturing Practice requirement, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic (QS) inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4659. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or at (301) 443-6597, or at its Internet address "<http://www.fda.gov/cdrh/dsmamain.html>".

Sincerely yours,

A handwritten signature in dark ink, appearing to read "Celia M. Witten".

Celia M. Witten, Ph.D., M.D.
Director
Division of General, Restorative and
Neurological Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

510(k) Number (if known): K001024

Device Name: PASS Spinal System

Indications For Use:

PASS Spinal System

Indications For Use

The PASS Spinal System consists of pedicle screws, rods, nuts and crosslink members utilized to provide immobilization and stabilization of spinal segments in skeletally mature patients as an adjunct to fusion in the treatment of the following acute and chronic instabilities or deformities of thoracic, lumbar, and sacral spine: degenerative spondylolisthesis with objective evidence of neurological impairment, fracture, dislocation, scoliosis, kyphosis, spinal tumor, and failed previous fusion (pseudarthrosis).

As a pedicle screw system, the PASS Spinal System is intended for the treatment of severe spondylolisthesis (Grade 3 and 4) of the L5-S1 vertebrae in skeletally mature patients receiving fusion by autogenous bone graft having implants attached to the lumbar and sacral spine (L3 to sacrum) with removal of the implants after the attainment of a solid fusion.

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Donner R. Lockner
(Division Sign-Off)
Division of General Restorative Devices
510(k) Number K001024

Prescription Use ✓
(per 21 CFR 801.109)

OR

Over-The-Counter Use _____

(Optional Format 1-2-96)_