

SEP 27 2000

510(k) Summary

K002146
Page 1 of 2

July 10, 2000

Eastman Kodak Company
343 State Street
Rochester NY 14650

Contact: Stephen G. Slavens
1 Imation Way, 304-3B-61
Oakdale, MN 55128
Phone: 651-393-1395
FAX: 651-393-1025

Device: Trade name: KODAK DRYVIEW 8610 Laser Imager / for
mammography
Common name: Laser Printer
Classification name: Medical Image Hardcopy Device
21 CFR 892.2040

Predicate devices: 3M (Imation™) 8600 Laser Imager (K972822)

Description And Intended Use of Device:

The KODAK DRYVIEW 8610 Laser Imager / for mammography is intended use as a high resolution hard copy device for output from digital imaging source modalities for use in medical imaging diagnosis and referral. Electronic image information signals are managed in the 8610 and transformed optically to expose KODAK DRYVIEW mammography imaging media. The system is intended for use with a variety of digital modalities, including, but not limited to digital radiography and full field digital mammography for diagnostic use by medical radiologists and communications to referring physicians and their patients.

Technological Characteristics:

The subject device and predicate devices use the same technical design base. The printers receive image data from the modality. Modality data and printing functions are performed by the IMS (Image management System). User control is performed by a keypad or directly by the modality through the host control. KODAK DRYVIEW mammography imaging media is removed from a daylight cartridge and transported to the laser imaging station. Image data and media merge at the laser station and the film is scanned. The exposed media is transported through the integrated processor and exits the printer.

Software is used to control the image management and machine functions. AIQC (Automated Image Quality Control) matches printing power with film characteristics to provide consistently high image quality.

Performance Data:

K002146
Page 2 of 2

Safety and effectiveness are assured via meeting voluntary standards, including IEC601-1, IEC601-1-1, IEC601-1-2 and UL2601.

Conclusion:

The subject device, like the predicate, has no patient contact. The devices also do not control, monitor or otherwise affect any devices directly connected to or affecting the patient. Medical personnel review images displayed by the subject device and its predicates. This offers ample opportunity for competent human intervention in case of a malfunction or other failure.

The subject KODAK DRYVIEW 8610 laser imager and predicate device Imation DryView™ Model 8600 have both been designed to the same safety standard. As with this predicate device, a test pattern generator and automatic image quality control (AIQC) system are incorporated to assure consistency between input signals and output density. The subject device has been designed to have the same resolution as the predicate Imation DryView™ Model 8600, which has the indications for use of full view digital mammography and high resolution computed radiography.

Eastman Kodak therefore concludes that the KODAK DRYVIEW 8610 Laser Imager / for mammography is as safe and effective as the predicate device.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

SEP 27 2000

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

Eastman Kodak Company
C/o Stephen G. Slavens
Regulatory Affairs Manager
Kodak Health Imaging
1 Imation Way, 304-4B-75
Oakdale, MN 55128-3414

Re: K002146
KODAK DRYVIEW 8610 Laser Imager/
for mammography
Dated: July 10, 2000
Received: July 17, 2000
Regulatory class: II
21 CFR 892.2040/Procode: 90 LMC

Dear Mr. Slavens:

We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the Current Good Manufacturing Practice requirements, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic QS inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4639. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its internet address "<http://www.fda.gov/cdrh/dsma/dsmamain.html>".

Sincerely yours,

Daniel G. Schultz, M.D.
Captain, USPHS
Acting Director, Division of Reproductive,
Abdominal, and Radiological Devices
Office of Device Evaluation
Center for Devices and Radiological Health

Enclosure (s)

ATTACHMENT 2

Statement of Indications for Use:

510(K) Number (if known): K002146

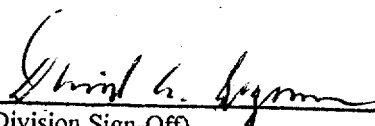
Device Name: KODAK DRYVIEW 8610 Laser Imager / for Mammography

Indications for Use:

The KODAK DRYVIEW 8610 laser Imager / for mammography is intended to provide high-resolution hard copy images from digital imaging source output signals. The device is intended for use with a variety of digital modalities, including, but not limited to digital radiography (DR) and full field digital mammography. The images are to be used for medical diagnosis and referral to physicians and patients.

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use X OR Over the Counter Use _____
(Per 21 CFR 801.109)


(Division Sign-Off)
Division of Reproductive, Abdominal, ENT,
and Radiological Devices
510(k) Number K002146