



JAN 16 2001

K003991

**Special 510(k) Summary**  
**CD-Cal Plus Hematology Calibrator**

Date of Summary:	December 21, 2000
Company Name:	R&D Systems, Inc. 614 McKinley Place N.E. Minneapolis, MN 55413
Contact name:	Kenneth T. Edds, Ph.D. 612-379-2956, FAX 612-379-6580
Classification name:	Hematology calibrator
Classification code:	81KRX Calibrator for cell indices 81KRY Calibrator for platelet counting 81KRZ Calibrator for hemoglobin and hematocrit measurement 81KSA Calibrator for red and white cell counting
Product name:	CD-Cal Plus™ Hematology Calibrator
CFR section:	864.8150, 864.8165, 864.8175, 864.8185
Device Class:	Class II

Device to which substantial equivalence is claimed:  
CD-Cal Plus Hematology Calibrator, manufactured by R&D Systems, Inc. 510(k)  
number: K955925

The product is an *in vitro* diagnostic reagent composed of human erythrocytes, mammalian leukocytes and platelets in a plasma-like fluid with preservatives. CD-Cal Plus is composed of stable materials that provide a means of calibrating Abbott CELL DYN hematology systems. CD-Cal Plus allows the calibration of multiple parameters including, with this device modification, the White cell Impedance Count (WIC) on the CELL DYN 3500 and 3700 instruments.

Intended use: CD-Cal Plus™ is a whole blood calibrator for use in calibration of CELL-DYN® hematology instruments. Refer to the assay table for specific instrument models. Values are provided for WBC, RBC, HGB, PLT, MCV, and MPV.

CD-Cal Plus Hematology Calibrator has an intended use that is identical to the predicate device. The technologies of the two devices are identical.

Nonclinical testing of 3 validation lots centered on the performance attributes of stability and precision. CD-Cal Plus Hematology Calibrator passed the acceptance criteria of remaining within the acceptable range over the life of the product. CD-Cal Plus calibrator also demonstrated precision as indicated by the small standard deviations and %CVs obtained during testing. Expiration dating has been established at 30 days in the customers hands (closed vial) and 7 days open vial when stored at 2-8°C and handled according to instructions for use.

R & D Systems, Inc.  
614 McKinley Place N.E.  
Minneapolis, MN 55413

PHONE: (612) 379-2956  
FAX: (612) 379-6580

WATS: (800) 343-7475  
E-MAIL: info@rndsystems.com



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
2098 Gaither Road  
Rockville MD 20850

JAN 16 2001

Kenneth T. Edds, Ph.D.  
Director, RA/QA  
R & D Systems, Inc.  
614 McKinley Place N.E.  
Minneapolis, Minnesota 55413

Re: K003991

Trade Name: CD-Cal Plus™ Hematology Calibrator

Regulatory Class: II

Product Code: KRX

II

KRY

II

KRZ

II

KSA

Dated: December 21, 2000

Received: December 26, 2000

Dear Dr. Edds:

We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895.

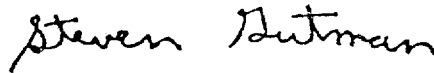
A substantially equivalent determination assumes compliance with the Current Good Manufacturing Practice requirements, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic QS inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

Page 2

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4588. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its internet address "<http://www.fda.gov/cdrh/dsma/dsmamain.html>".

Sincerely yours,

A handwritten signature in black ink that reads "Steven Gutman". The signature is written in a cursive, slightly slanted style.

Steven I. Gutman, M.D., M.B.A.  
Director  
Division of Clinical Laboratory Devices  
Office of Device Evaluation  
Center for Devices and Radiological Health

Enclosure

510(k) Number: K003991

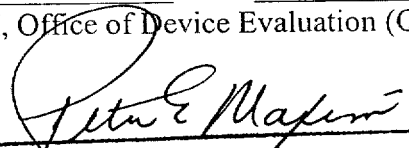
Device Name: CD-Cal Plus Hematology Calibrator

Indications for Use:

CD-Cal Plus is a whole blood calibrator for use in calibration of CELL-DYN® hematology instruments. Refer to the assay tables for specific instrument models. Values are provided for WBC, RBC, HGB, PLT, and MPV.

(PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

  
(Division Sign-Off)  
Division of Clinical Laboratory Devices  
510(k) Number K003991

Prescription Use ✓  
(Per 21 CFR 801.109)

OR

Over-The-Counter Use \_\_\_\_\_

(Optional Format 1-2-96)