

510(k) Summary

As Required by 21 section 807.92 (c)

- 1-Submitter Name:** A & A Medical, Inc.
2-Address: 9370 Industrial Trace
Alpharetta, GA 30004
3-Phone: (770) 343- 8400
4-Fax: (770) 343- 8985
5-Contact Person: Jihad Mansour
6-Date summary prepared: January 4th, 2001
7-Device Trade or Proprietary Name: Uterine Manipulator Injector Cannula
8-Device Common or usual name: Uterine Manipulator Injector Cannula
9-Device Classification Name: Cannula, Manipulator/Injector, Uterine
10-Substantial Equivalency is claimed against the following device:
- Zumi-4.5TM- Zinnanti Uterine Manipulator Injector

11-Description of the Device:

The device is to be used by physicians in hospitals

The Uterine Manipulator Injector Cannula is both a uterine manipulator and a uterine injector for single use. This is a sterile (by ethylene oxide) disposable product made out of clear plastic, which meets USP recommendations for implant testing. This product is designed with a double lumen, one for inflation of a 10cc intrauterine cuff and the other for injection of fluid through a distal endport. The product is curved to facilitate forward uterine manipulation. The product features an inflation valve and pilot balloon assembly, an endport, an inflatable cuff, centimeters depth markings, cervical stop, a removable rigid plastic handle and a luer fitting to accommodate a syringe. The instrument has a length of 13 inches (33cm) and an outer diameter of 4.5mm (.18 inches)

12-Intended use of the device:

This device is indicated for use in Diagnostic Laparotomy, Minilaparotomy, Fertility, Examinations, and Salpingoplastic procedures where manipulation of the uterus is required. This product also facilitates the sealing of the cervical os while providing a fluid or air injection port.

13-Safety and Effectiveness of the device:

This device (**Uterine Manipulator Injector Cannula**) is safe and effective as the other predicate device cited above.

This is better expressed in the tabulated comparison (Paragraph 14 below)

14-Summary comparing technological characteristics with other predicate device:

Please find below a tabulated comparison supporting that **Uterine Manipulator Injector Cannula** is substantially equivalent to other medical devices in commercial distribution. Also, Equivalency overview chart path is attached.

FDA file reference number	510k 941458
Attachments inside notification submission file	REFER TO TABLE ON PAGE 11 OF 12 FOR DETAILS
TECHNOLOGICAL CHARACTERISTICS	Comparison result
Indications for use	Equivalent
Target population	Equivalent
Design	Equivalent
Materials	Equivalent
Performance	Equivalent
Sterility	Similar (Ethylene Oxide but different parameters)
Biocompatibility	Equivalent
Mechanical safety	Equivalent
Chemical safety	Equivalent
Anatomical sites	Equivalent
Human factors	Equivalent
Energy used and/or delivered	Equivalent
Compatibility with environment and other devices	Equivalent
Where used	Equivalent
Standards met	Equivalent
Electrical safety	Equivalent (not applicable)
Thermal safety	Equivalent (not applicable)
Radiation safety	Equivalent (not applicable)



Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

FEB - 1 2001

Mr. Jihad Mansour
Quality and Regulatory Manager
A&A Medical, Inc.
9370 Industrial Trace
ALPHARETTA GA 30004

Re: K010056
Uterine Manipulator Injector Cannula
Dated: January 6, 2001
Received: January 8, 2001
Unclassified
Procode: 85 LKF

Dear Mr. Mansour:


We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the Current Good Manufacturing Practice requirements, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic QS inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note: this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4639. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its internet address "<http://www.fda.gov/cdrh/dsma/dsmamain.html>".

Sincerely yours,

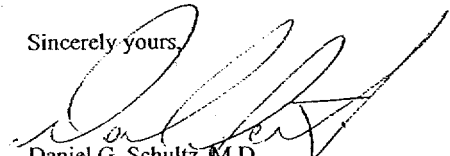

Daniel G. Schultz, M.D.
Captain, USPHS
Acting Director, Division of Reproductive,
Abdominal, and Radiological Devices
Office of Device Evaluation
Center for Devices and Radiological Health

Enclosure (s)

510(k) Number (if known): K010056

Device Name: UTERINE MANIPULATOR INJECTOR CANNULA

Indications For Use:

THIS DEVICE IS INDICATED FOR USE IN
DIAGNOSTIC LAPAROTOMY, MINILAPAROTOMY,
FERTILITY, EXAMINATIONS, AND SALPINGOPLASTIC
PROCEDURES WHERE MANIPULATION OF THE UTERUS
IS REQUIRED. THIS DEVICE ALSO FACILITATES
THE SEALING OF THE CERVICAL OS WHILE
PROVIDING A FLUID OR AIR INJECTION PORT

(PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ☒
(Per 21 CFR 801.109)

OR

Over-The-Counter Use _____

(Optional Format 1-2-96)

David L. Segura
(Division Sign-Off)
Division of Reproductive, Abdominal, ENT,
and Radiological Devices

Page 4 of 12

510(k) Number K010056