



K010499

MAY - 3 2001

### 3.0 Summary of Safety and Effectiveness Information

**SPONSOR:** Synthes (USA)  
1690 Russell Road  
Paoli, PA 19301  
(610) 647-9700  
Contact: Thomas M. Maguire

**DEVICE NAME:** Midface Distractor

**CLASSIFICATION:** Class II, 21 CFR 872.4760: Bone plate  
**PREDICATE DEVICE:** Cohen Distractor (Howmedica Leibinger Inc.)

**DEVICE DESCRIPTION:** The Synthes Midface Distractor is a craniofacial distraction device consisting of two telescoping components with attached footplates. The device is intended to be placed subcutaneously, with an anterior footplate fastened to the lateral orbital rim, extending down to the maxilla and spanning the zygomaticomaxillary suture; and a posterior footplate fastened to the temporal region of the cranium. The plates are fixed to the bone through unthreaded screw holes using 1.5 mm or 2.0 mm Cortex screws.

**INTENDED USE:** For use in adult and pediatric populations for the treatment of cranial or midface conditions for which reconstructive osteotomy and segment advancement are indicated. This includes conditions such as syndromic craniosynostosis and midfacial retrusion. The device is intended to provide temporary stabilization and gradual lengthening of the cranial or midfacial bones.

**MATERIAL:** Titanium Alloy, Titanium, Chromium Cobalt



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

MAY - 3 2001

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

Mr. Thomas M. Maguire  
Project Leader of Regulatory Affairs  
Synthes (USA)  
1101 Synthes Avenue  
Monument, Colorado 80132

Re: K010499  
Trade/Device Name: Synthes (USA) Midface Distractor  
Regulation Number: 872.4760  
Regulatory Class: II  
Product Code: JEY  
Dated: February 20, 2001  
Received: February 21, 2001

Dear Mr. Maguire:

We have reviewed your Section 510(k) notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.


If your device is classified (see above) into either class II (Special Controls) or class III (Premarket Approval), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 895. A substantially equivalent determination assumes compliance with the Current Good Manufacturing Practice requirements, as set forth in the Quality System Regulation (QS) for Medical Devices: General regulation (21 CFR Part 820) and that, through periodic QS inspections, the Food and Drug Administration (FDA) will verify such assumptions. Failure to comply with the GMP regulation may result in regulatory action. In addition, FDA may publish further announcements concerning your device in the Federal Register. Please note:

this response to your premarket notification submission does not affect any obligation you might have under sections 531 through 542 of the Act for devices under the Electronic Product Radiation Control provisions, or other Federal laws or regulations.

This letter will allow you to begin marketing your device as described in your 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4692. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its internet address "<http://www.fda.gov/cdrh/dsma/dsmamain.html>".

Sincerely yours

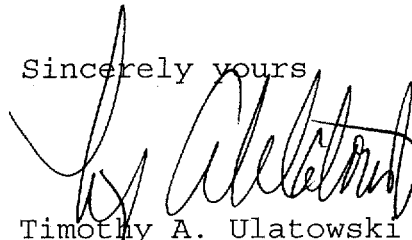

Timothy A. Ulatowski  
Director  
Division of Dental, Infection Control  
and General Hospital Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure



K010499

2.0 Indications for Use Statement

Page 1 of 1

510(k) Number (if known):

K010499

Device Name: Synthes (USA) Midface Distractor

Indications/Contraindications:

For use in adult and pediatric populations for the treatment of cranial or midface conditions for which reconstructive osteotomy and segment advancement are indicated. This includes conditions such as syndromic craniosynostosis and midfacial retrusion. The device is intended to provide temporary stabilization and gradual lengthening of the cranial or midfacial bones.

(PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ☒  
(Per 21 CFR 801.109)

OR

Over-The-Counter Use ☐

*Susan Purges*

(Official Sign-Off)  
Division of Dental, Infection Control,  
General Hospital Devices  
Number K010499