

OCT 10 2002

**510(k) Summary of Safety & Effectiveness**

---

<b>Submitter</b>	Vanguard Medical Concepts, Inc. 5307 Great Oak Drive Lakeland, FL 33815
<b>Contact</b>	Mike Sammon, Ph.D. Director, Research and Development (863) 904-1628 (801) 327-3339 (facsimile) msammon@safe-reuse.com
<b>Date</b>	July 15, 2002
<b>Device</b>	<ul style="list-style-type: none"><li>• Trade Names: Vanguard Reprocessed Diagnostic Electrophysiology Catheters ⇒ EP Technologies™ Explorer 360™ Diagnostic Electrophysiology Catheters</li><li>• Common Name: Electrode Recording Catheter, Diagnostic Electrophysiology (EP) Catheter</li><li>• Classification: 21 CFR 870.1220 – Class II – Catheter, Electrode Recording, or Probe, Electrode Recording</li><li>• Product Code DRF</li></ul>
<b>Predicate Devices</b>	EP Technologies™ Explorer 360™ diagnostic EP catheters legally marketed under various 510(k) premarket notifications.
<b>Indications for Use</b>	This catheter is intended for temporary intracardiac pacing and recording during electrophysiology studies, e.g. evaluation of arrhythmias or cardiac mapping.

---

*Continued on next page*

## 510(k) Summary of Safety & Effectiveness, Continued

### Contra- indications

- Patients with active systemic infection.
- Patients with prosthetic valves.
- Retrograde approach in patients with aortic valve replacement.
- Transseptal approach in patients with left atrial thrombus or myxoma, or interatrial baffle patch.
- Diagnostic EP catheters are not intended for electrical ablation.

### Device Description

Explorer 360™ diagnostic electrophysiology catheters are constructed of a 6 French hollow polymer shaft 100 cm in length. The quadripolar catheters are available with various fixed distal curves and electrode spacing configurations. The polymer shaft is manufactured with additives (typically barium sulfate) that enhance its radiopacity during fluoroscopic positioning of the catheter. The platinum alloy electrodes are sealed to the distal catheter and internally wired to a proximal shielded Nexus connector for bi-directional transmission of electrical signals (pacing and recording). The connector is designed to be attached to a compatible instrument cable that interfaces with various standard types of sensing, recording, stimulation, and pacing equipment.

Vanguard receives previously used diagnostic EP catheters from healthcare facilities; cleans, inspects, tests, refurbishes, applies a unique serial number, repackages, and sterilizes the devices; and returns them to the healthcare facility.

*Continued on next page*

## 510(k) Summary of Safety & Effectiveness, Continued

---

<b>Technological Characteristics</b>	The Vanguard reprocessed Explorer 360™ diagnostic catheters are essentially identical to the currently marketed OEM devices. No changes are made to the currently marketed device's specifications and they possess the same technological characteristics. Biocompatibility and performance/functional testing demonstrate that the devices are equivalent and continue to be safe and effective for their intended use.
<b>Test Data</b>	Cleaning, sterilization, and packaging validations; and functional/performance and biocompatibility testing demonstrates that the reprocessed devices perform as intended and are safe and effective.
<b>Conclusion</b>	Based on the information provided herein and the 510(k) "Substantial Equivalence" Decision Making Process Chart, we conclude that the Vanguard Reprocessed Diagnostic EP Catheters (EP Technologies™ Explorer 360™ Catheters) are substantially equivalent to the predicate devices under the Federal Food, Drug and Cosmetic Act.

---



Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

OCT 1 0 2002

Vanguard Medical Concepts, Inc.  
c/o Mike Sammon, Ph.D.  
Director, Research and Development  
5307 Great Oak Drive  
Lakeland, FL 33815

Re: K022316

Trade Name: Vanguard Reprocessed Diagnostic Electrophysiology Catheters  
Regulation Number: 21 CFR 870.1220  
Regulation Name: Electrode Recording Catheter  
Regulatory Class: Class II (two)  
Product Code: DRF  
Dated: July 15, 2002  
Received: July 17, 2002

Dear Dr. Sammon:

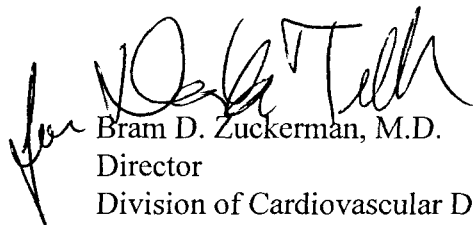
We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050. This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 21 CFR Part 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4646. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read "Bram D. Zuckerman". The signature is stylized with a large "Z" and "U".

Bram D. Zuckerman, M.D.  
Director  
Division of Cardiovascular Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

## Indications for Use

510(k) Number: K022316

Device Name: Vanguard Reprocessed Diagnostic Electrophysiology Catheters

### Indications for Use:

This catheter is intended for temporary intracardiac pacing and recording during electrophysiology studies, e.g. evaluation of arrhythmias or cardiac mapping.

(PLEASE DO NOT WRITE BELOW THIS LINE – CONTINUE ON ANOTHER PAGE IF NEEDED.)


Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use ✓

OR

Over-the-Counter Use \_\_\_\_\_

(Per 21 CFR 801.109)

  
Division of Cardiovascular & Respiratory Devices  
510(k) Number K022316

(Optional Format 1-2-96)