

**510(k) Summary of Safety & Effectiveness**

---

<b>Submitter</b>	Vanguard Medical Concepts, Inc. 5307 Great Oak Drive Lakeland, FL 33815
------------------	---

---

<b>Contact</b>	Mike Sammon, Ph.D. Director, Research and Development (863) 904-1628 (801) 327-3339 (facsimile) msammon@safe-reuse.com
----------------	--

---

<b>Date</b>	January 10, 2002
-------------	------------------

---

<b>Device</b>	<ul style="list-style-type: none"><li>• Trade Names: Vanguard Reprocessed Diagnostic Electrophysiology (EP) Catheters</li><li>• Common Name: Electrode Recording Catheter, Diagnostic Electrophysiology (EP) Catheter</li><li>• Classification: 21 CFR 870.1220 – Class II – Catheter, Electrode Recording, or Probe, Electrode Recording</li><li>• Product Code DRF</li></ul>
---------------	--

---

<b>Predicate Devices</b>	<ul style="list-style-type: none"><li>• Daig™ Response™ Fixed Curve Diagnostic Catheters (K002976)</li><li>• Daig™ Supreme™ Fixed Curve Diagnostic Catheters (K002976)</li><li>• Vanguard Reprocessed Diagnostic EP Catheters (K012687, K012688, K022316, K023180)</li></ul>
--------------------------	--

---

<b>Indications for Use</b>	This catheter is intended for temporary intracardiac pacing and recording during electrophysiology studies, e.g. evaluation of arrhythmias or cardiac mapping.
----------------------------	--

<b>Contra-indications</b>	<ul style="list-style-type: none"><li>• Patients with active systemic infection.</li><li>• Patients with prosthetic valves.</li><li>• Retrograde approach in patients with aortic valve replacement.</li><li>• Transseptal approach in patients with left atrial thrombus or myxoma, or interatrial baffle patch.</li><li>• Diagnostic EP catheters are not intended for electrical ablation.</li><li>• Diagnostic EP catheters are not intended for coronary artery mapping.</li></ul>
---------------------------	---

---

*Continued on next page*

## 510(k) Summary of Safety & Effectiveness, Continued

<b>Device Description</b>	<p>Reprocessed diagnostic electrophysiology catheters are constructed of a hollow polymer shaft approximately 65 to 125 cm in length that terminates with a hand piece or connector. A range of diameters is available; the most clinically utilized sizes are 4 – 7 French. Various configurations of distal platinum alloy electrodes are wired to a proximal connector for bi-directional transmission of electrical signals (pacing and recording). The connector is attached to an interconnecting cable that interfaces with various standard types of sensing, recording, stimulation and pacing equipment. The catheters are available with various distal curves, either fixed or steerable for remote manipulation of the distal tip segment that facilitates precise positioning of the electrode array.</p> <p>The catheters are also available in a variety of electrode configurations, connector compatibility and torque-transmitting properties that are selected by the clinician based on preference and/or indication. The shaft polymer is manufactured with additives (typically barium sulfate) that enhance the catheter's radiopacity to enable positioning under fluoroscopic guidance. No lumens of catheters reprocessed by Vanguard are open to the patient bloodstream.</p> <p>Vanguard collects previously used diagnostic EP catheters from healthcare facilities; cleans, inspects, tests, applies a unique serial number, packages, labels, and sterilizes each device for return to the healthcare facility for an additional clinical use.</p>
<b>Technological Characteristics</b>	<p>Vanguard Reprocessed Diagnostic EP Catheters are essentially identical to the Original Equipment Manufacturer (OEM) devices. No changes are made to the device materials or specifications and the reprocessed catheters possess identical technological characteristics.</p>
<b>Test Data</b>	<p>Cleaning, sterilization, and packaging validations; and performance and biocompatibility testing demonstrate that the reprocessed devices perform as intended and are safe and effective.</p>
<b>Conclusion</b>	<p>Based on the information provided herein and the 510(k) "Substantial Equivalence" Decision Making Process Chart, we conclude that the Vanguard Reprocessed Diagnostic EP Catheters are substantially equivalent to their predicate devices under the Federal Food, Drug and Cosmetic Act.</p>



Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

MAR 21 2003

Vanguard Medical Concepts Inc.  
c/o Mike Sammon, Ph.D.  
Director, Research and Development  
5307 Great Oak Drive  
Lakeland, FL 33815

Re: K030114

Trade Name: Vanguard Reprocessed Diagnostic Electrophysiology (EP) Catheters  
Regulation Number: 21 CFR 870.1220  
Regulation Name: Electrode Recording Catheter or Electrode Recording Probe  
Regulatory Class: Class II (two)  
Product Code: DRF  
Dated: January 10, 2002  
Received: January 13, 2003

Dear Dr. Sammon:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

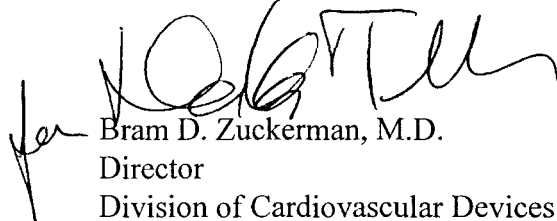
Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act

or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (301) 594-4646. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read "Bram D. Zuckerman", is written over the typed name. To the left of the signature is a small, stylized mark that looks like a checkmark or the letter "f".

Bram D. Zuckerman, M.D.

Director

Division of Cardiovascular Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

## Indications for Use

**510(k) Number:**

**Device Name:** Vanguard Reprocessed Diagnostic Electrophysiology Catheters

### Indications for Use:

This catheter is intended for temporary intracardiac pacing and recording during electrophysiology studies, e.g. evaluation of arrhythmias or cardiac mapping.

(PLEASE DO NOT WRITE BELOW THIS LINE – CONTINUE ON ANOTHER PAGE IF NEEDED.)


Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use X

OR

Over-the-Counter Use \_\_\_\_\_

(Per 21 CFR 801.109)

  
(Division Sign-Off)  
Division of Cardiovascular Devices

(Optional Format 1-2-96)

510(k) Number K030114

iv