# 

FEB 2 5 2004

510(k) SUMMARY

(Per 21 CFR 807.92)

page 1 of 2

## General Company Information

Name:

Alveolus Inc.

Contact:

Howard Schrayer

Regulatory Affairs Consultant

Address:

401 North Tryon Street, 10th Floor

Charlotte, NC 28202

Telephone:

(704) 998 - 6009

čax:

(704)998 - 5699

Date Prepared

February 17, 2004

#### General Device Information

Trade Name:

TB-STS™ Tracheobronchial Stent System

Classification:

"Tracheal Prosthesis", Product code: JCT

21 CFR 878.3720 - Class II

#### **Predicate Devices**

Alveolus, Inc. TB-STS<sup>TM</sup> Tracheobronchial Stent System

[510(k) Number K030947]

Rusch International Polyflex Stent Kit

[510(k) Number K013266]

Novatech S.A. Endoxane® Stent

[510(k) Number K971509]

### Description

The Alveolus Tracheobronchial Stent Technology System is comprised of two components: the radiopaque stent and the delivery system. The nitinol stent is completely covered with a biocompatible polyurethane (ChronoFlex™) membrane and is self-expanding. The stent expansion results from the mechanical properties of the metal and the proprietary geometry. The stent is designed with a slightly larger diameter near the distal and proximal ends to minimize the possibility of migration. The stent ends are slightly vauited inwardly in order to minimize possible airway injury from the stent edges. The overall stent geometry is designed to maintain a constant length over the entire range of possible diameters. As a result of this

K033023

page 2012.

unique design the stent has virtually no foreshortening, thus facilitating the selection of the appropriate stent length.

# Intended Use (Indications)

The Alveolus TB-STS<sup>TM</sup> Tracheobronchial Stent System is indicated for use in the treatment of tracheobronchial strictures and airway compression (stenosis) produced by malignant neoplasms. Because the device is removable it may also be used to treat benign conditions such as tracheo-esophageal fistulae and strictures.

# Substantial Equivalence


This submission supports the position that the Alveolus Tracheobronchial Stent is substantially equivalent to a number of previously cleared devices, including the Boston Scientific Corp. Inc. Ultraflex<sup>TM</sup> Tracheobronchial Stent System [501(k) Number K963241], the Rusch International Polyflex Stent Kit [510(k) Number K013266] and the Novatech S.A. Endoxane® Stent [510(k) Number K971509].

The 510(k) Notice contains a report of an *in vivo* study that demonstrates that the completely covered stent is removable like the silicone stent predicates and is thus suitable for use in the treatment of benign conditions in addition to the original indications for use in treating malignant disease.

The single-patient-use components of the TB-STS™ Tracheobronchial Stent System are provided sterile.

#### Conclusions

Alveolus Inc. believes that the information provided establishes that similar legally marketed devices have been used for the same clinical applications as the Alveolus Tracheobronchial Stent. The materials from which the Alveolus device is fabricated have an established history of use in clinical applications, and the devices produced by Alveolus have been tested in accordance with applicable FDA guidelines.

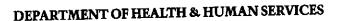



Food and Drug Administration 9200 Corporate Boulevard Rockville MD 20850

FEB 2 5 2004

Mr. Howard L. Schrayer Alveolus, Inc. 401 North Tryon Street, 10<sup>th</sup> Floor Charlotte, North Carolina 28202

Re: K033053

Trade/Device Name: Alveolus TB-STS™ Tracheobronchial Stent System

Regulation Number: 21 CFR 878.3720 Regulation Name: Tracheal Prosthesis

Regulatory Class: II Product Code: JCT Dated: January 7, 2004 Received: January 8, 2004

Dear Mr. Schrayer:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the <u>Federal Register</u>.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (301) 594-4659. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address http://www.fda.gov/cdrh/dsma/dsmamain.html

Sincerely yours,

Manh M. Mulkers.

Celia M. Witten, Ph.D., M.D.

Director

Division of General, Restorative and Neurological Devices Office of Device Evaluation

Center for Devices and Radiological Health

Enclosure

# Indications for Use

510(k) Number (if known): K033053
Device Name: <u>Alveolus, TB-STS™ Tracheobronchial Stent System</u>
Indications For Use:
The Alveolus TB-STS™ Tracheobronchial Stent System is indicated for use in the treatment of tracheobronchial strictures and airway compression (stenosis) produced by malignant neoplasms. Because the device is removable it may also be used to treat benign conditions such as tracheo-esophageal fistulae and strictures resulting from surgical anastomosis of the airway.
Prescription Use AND/OR Over-The-Counter Use (Part 21 CFR 801 Subpart D) (21 CFR 807 Subpart C)  (PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE IF NEEDED)
Concurrence of CDRH, Office of Device Evaluation (ODE)
Division Sign-Off)  Division of General, Restorative,  and Metarological Devices

510(\* Tumber K033053