

FDA
Electronic Submission Gateway (ESG) NextGen
AS2 Guide for Industry Users

Version 2.2
July 2026

Version History

Version Number	Implemented By	Revision Date	Approved By	Approval Date	Description of Change
2.1	Nicholas Oja	01/30/2026	Jessica Bernhardt	01/30/2026	Updated section 2.3 'Additional Steps for Password Protected Endpoints' to indicate user should email credentials.
2.1.1	Amy McKeehan	06/15/2026	Jessica Bernhardt	06/15/2026	Section 2 added link to AS2 Set-up Instructions. Section 2.4 links for AS2 routing IDs updated
2.2	Felicity Thomas	7/1/2026			Removed references to ZZFDATST certificate. Removed "metadata"

					from X- Cyclone-True- Receiver
--	--	--	--	--	--------------------------------------

Contents

Version History.....	1
1. Introduction	3
1.1 Purpose of the Guide	3
1.2 Target Audience	3
2. Getting Started with AS2	3
2.1 AS2 Submission Process.....	4
2.2 Steps to Update AS2 Configuration	4
2.3 Additional Steps for Password Protected AS2 Endpoints	6
2.4 Send a Submission	6
2.5 Confirm Receipt of MDN and ACKs	8
3. Unchanged from Legacy ESG.....	9
4. Changes from Legacy ESG	9
5. Common Issues and Trouble Shooting	10
6. References	11

1. Introduction

1.1 Purpose of the Guide

This document provides guidelines for submitting files and folders through Applicability States 2 (AS2) with the FDA's ESG NextGen system.

AS2 is a point-to-point protocol for securely exchanging data. AS2 allows businesses to transmit data, such as electronic data interchange (EDI) messages, securely and reliably. AS2 is considered a legacy technology. The modernization of AS2 submissions is to submit via API. If you are looking to create a new point-to-point submission to connect to ESG NextGen, API is the recommended method. For more information on API submissions, refer to the ESG NextGen API Guide which can be found [here](#) under User Guides.

1.2 Target Audience


This guide is intended for industry participants who want to utilize AS2 to submit regulatory information to the FDA. The audience includes current AS2 submitters and user technical support resources.

2. Getting Started with AS2

All legacy pre-production and production ESG AS2 accounts have been migrated to ESG NextGen. Please review the entire AS2 Guide to understand how to integrate with ESG NextGen. If you have never used AS2 with legacy ESG, please contact support at: ESGNGSupport@fda.hhs.gov and review the [ESG NextGen AS2 Setup Instructions](#).

2.1 AS2 Submission Process

Below is the high level AS2 submission process in ESG NextGen

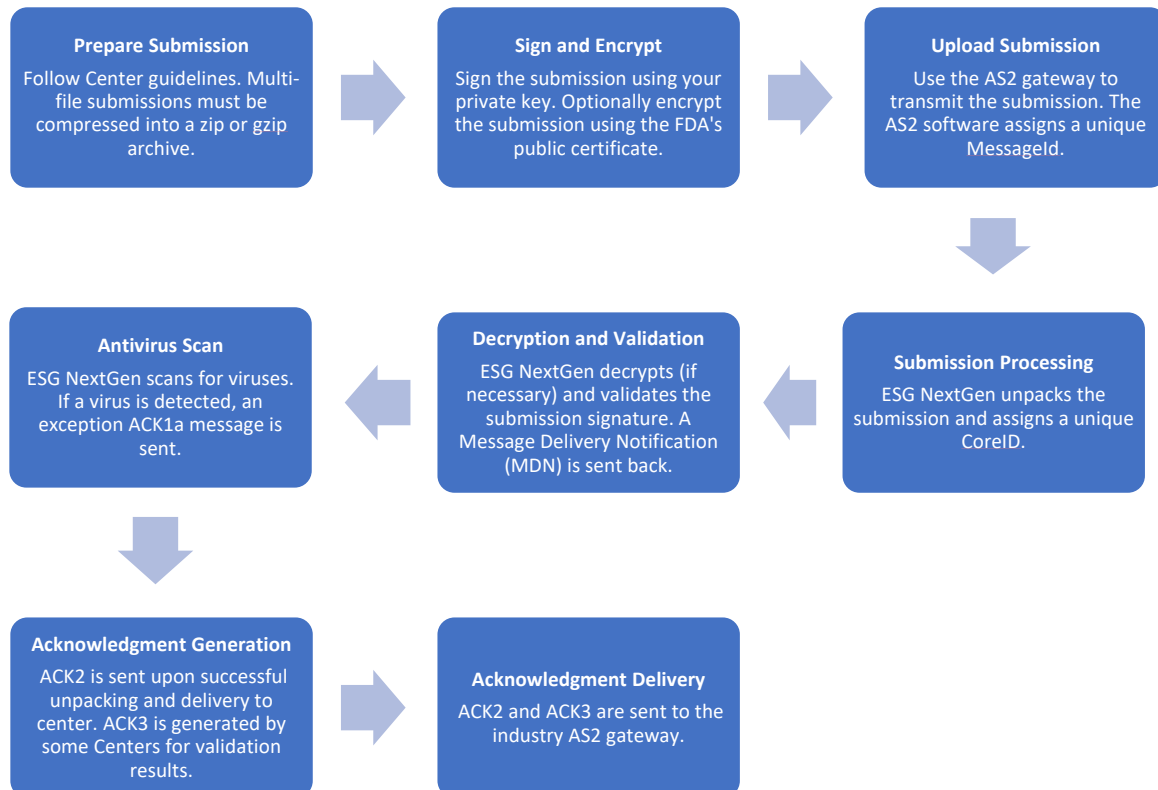

2.2 Steps to Update AS2 Configuration

1. Access AS2 Configuration:

- Navigate to your AS2 gateway configuration settings.

2. Update AS2 Endpoint:

- If your company limits its AS2 communications via a firewall IP whitelist they will need to add new ranges for ESG Next Gen. These are different than the addresses and/or ranges used for legacy ESG. Retain all legacy ESG firewall rules.
 - Inbound ranges (Industry to FDA) - 15.205.247.22 and 3.31.183.245
 - Outbound ranges (FDA to Industry) - 150.148.0.0/16
- If your company did not need to make changes to its whitelist in legacy ESG, then it should not be required for ESG NextGen.
- Modify the AS2 endpoint to deliver submissions based on production or test submission context.
 - For **production submissions**, update the AS2 endpoint to: <https://upload-api-esgng.fda.gov:4080/as2/receive>
 - For **test submissions**, update the AS2 endpoint to: <https://upload-api-esgng.fda.gov:4080/as2/receive/test>
- Advanced Encryption Standard (AES) encryption is recommended for submissions in ESG NextGen
 - The following encryption methods are supported but AES is preferred:
 - AES-128
 - AES-192
 - **AES-256**
 - AES-256-GCM is currently unsupported.
 - RC2
 - Triple DES (DESEde)
 - Cast5
 - Idea

3. Save and Confirm Changes

- Save the updated configuration and ensure all changes are correctly applied.

2.3 Additional Steps for Password Protected AS2 Endpoints

***This step is only required if the industry ACK AS2 endpoint is secured with a username and password. If it is not secured in this fashion, please skip to section 2.4.**

1. Prepare JSON Payload:

- Construct a JSON document with the following structure:

```
{  
  "username": "Replace this with the username which ESG NextGen  
should use to authenticate to your gateway. An empty string can be sent  
to remove this if desired.",  
  "password": "Replace this with the password which ESG NextGen  
should use to authenticate to your gateway. An empty string can be sent  
to remove this if desired.",  
  "url": "Replace this with the URL of your gateway."  
}
```

2. Email the JSON payload to ESG NextGen Support:

- Email the JSON payload to ESGNGSupport@fda.hhs.gov. This will create a ticket automatically.

3. Receive Confirmation:

- The ESG NextGen support team will update your endpoint credentials and respond to your ticket with confirmation of the update.

2.4 Send a Submission

1. Prepare either a Production or Test Submission:

- Follow the standard procedure to prepare a submission. This will be specific to your organization.

2. Send Submission:

- Transmit the submission to ESG NextGen using the updated AS2 configuration. Use either a routing ID or headers to indicate the receiving center and submission type. Please reference [Center Submission Types | FDA](#) for list of AS2 routing IDs.

An industry partner submitting via AS2 must indicate to ESG NextGen on how to route the submission via one of 3 methods:

Option 1: Center / submission type specific AS2 routing ID in the AS2-To header.

- Valid AS2 center / submission type routing ID's for industry can be found at [Center Submission Types | FDA](#) (scroll down to second table).
- Messages sent using this option should be **encrypted using the ZZFDA certificate** or left unencrypted.

Option	-	Production Submission - https://upload-api-esgng.fda.gov/as2/receive	Test Submission - https://upload-api-esgng.fda.gov/as2/receive/test
Option 1	Submission Header: AS2-To	FDA_AERS	FDA_AERS
	Submission Header: X-Cyclone-True-Receiver	DO NOT USE	DO NOT USE
	Submission Header: X-Cyclone-Metadata-FDASubmissionType	DO NOT USE	DO NOT USE
	Submission Header: X-Cyclone-Metadata-FDACenter	DO NOT USE	DO NOT USE
	Submission encryption certificate:	ZZFDA / esg.fda.gov	ZZFDA / esg.fda.gov
	ACK signing certificate:	ZZFDA / esg.fda.gov	ZZFDA / esg.fda.gov
	ACK Header: AS2-From	FDA_AERS	FDA_AERS

Option 2: Custom X-Cyclone-Metadata-FdaSubmissionType and X-Cyclone-Metadata-FdaCenter headers paired with **ZZFDA for production submissions** or **ZZFDATST for test submissions** in the AS2-To header.

- Valid custom header options can be found at [Center Submission Types | FDA](#) (scroll down to second table).
- Production submissions addressed to ZZFDA should be encrypted using the ZZFDA certificate or left unencrypted.
- Test submissions addressed to ZZFDATST should be encrypted using the ZZFDA certificate or left unencrypted.

Option	-	https://upload-api-esgng.fda.gov/as2/receive	https://upload-api-esgng.fda.gov/as2/receive/test
Option 2	Submission Header: AS2-To	ZZFDA	ZZFDATST
	Submission Header: X-Cyclone-True-Receiver	DO NOT USE	DO NOT USE
	Submission Header: X-Cyclone-Metadata-FDACenter	CDER	CDER
	Submission Header: X-Cyclone-Metadata-FDASubmissionType	AERS	AERS
	Submission encryption certificate:	ZZFDA / esg.fda.gov	ZZFDA / esg.fda.gov
	ACK signing certificate:	ZZFDA / esg.fda.gov	ZZFDA / esg.fda.gov
	ACK Header: AS2-From	ZZFDA	ZZFDATST

Option 3: Center / submission type specific AS2 routing ID in the X-Cyclone-True-Receiver header paired with **ZZFDA for production submissions** or **ZZFDATST for test submissions** in the AS2-To header.

- Valid AS2 center / submission type routing ID's for industry can be found at [Center Submission Types | FDA](#) (scroll down to second table).
- Production submissions addressed to ZZFDA should be encrypted using the ZZFDA certificate or left unencrypted.
- Test submissions addressed to ZZFDATST should be encrypted using the ZZFDA certificate or left unencrypted.

Option	-	https://upload-api-esgng.fda.gov/as2/receive	https://upload-api-esgng.fda.gov/as2/receive/test
Option 3	Submission Header: AS2-To	ZZFDA	ZZFDATST
	Submission Header: X-Cyclone-True-Receiver	FDA_AERS	FDA_AERS
	Submission Header: X-Cyclone-Metadata-FDACenter	DO NOT USE	DO NOT USE
	Submission Header: X-Cyclone-Metadata-FDASubmissionType	DO NOT USE	DO NOT USE
	Submission encryption certificate:	ZZFDA / esg.fda.gov	ZZFDA / esg.fda.gov
	ACK signing certificate:	ZZFDA / esg.fda.gov	ZZFDA / esg.fda.gov
	ACK Header: AS2-From	FDA_AERS	FDA_AERS

ESG NextGen utilizes the same encryption certificate as legacy ESG.

2.5 Confirm Receipt of MDN and ACKs

1. Monitor Acknowledgements:

- Verify the receipt of the following acknowledgements via AS2:
 - MDN (Message Disposition Notification) (optional)
 - ACK2 (if applicable to your submission type)
 - ACK3/4 (if applicable to your submission type)

3. Unchanged from Legacy ESG

The following aspects of the ESG NextGen system remain unchanged from the previous ESG system:

- FDA and Industry Certificates: All certificates loaded in legacy ESG Preprod and Prod will be migrated to ESG NextGen Production.
- FDA Routing IDs
- FDA Custom Routing Headers
- Industry AS2 Contact Information
- Industry Routing IDs
- ACK2 Format
- ACK3 Format

4. Changes from Legacy ESG

- Core ID Format: ci<timestamp>.<GUID> (47 characters).
- ESG NextGen AS2 Message ID follows RFC 4130 guidelines.
- IP Ranges:
 - Inbound range (industry to FDA) used for synchronous MDN's and submissions:
Prod: 15.205.247.22 and 3.31.183.245
 - Outbound range (FDA to industry) used for asynchronous MDN's, ACK2's, and ACK3's: 150.148.0.0/16
- AS2 URL
 - For **production submissions**: <https://upload-api-esgng.fda.gov:4080/as2/receive>

- For **test submissions**: <https://upload-api-esgng.fda.gov:4080/as2/receive/test>
- Industry Credentials must be resubmitted via JSON payload or encrypted email. For more information on the JSON file, refer to [Section 2.3 Additional Steps for Password Protected AS2 Endpoints](#).
- ESG NextGen only accepts one signing/encryption certificate per industry routing ID regardless of which environment the routing ID is utilized. If multiple signing certificates are required by your organization, then multiple industry routing IDs are required.
 - Ex: Same routing ID is used in the industry pre-prod/test environment as their production environment and a different certificate is used in each environment, this is not allowed. A new routing ID would need to be generated for one of the environments.
- Please contact ESGNGSupport@fda.hhs.gov for further assistance, if needed.

5. Common Issues and Trouble Shooting

The following table describes behaviors and possible causes:

Symptom	Potential Cause	Potential Action Items
Not receiving MDN	Security blocks	Verify industry firewalls are allowing FDA IP's to connect
		Verify certificates in place on both sides
	Bad request	Ensure header being used is a supported header

Symptom	Potential Cause	Potential Action Items
	Sender has requested not to receive MDN in header	Reconfigure sending software to request an MDN
	AS2 gateway issues	MDN's are being delivered to industry AS2 gateway but rejected with a 400 or 500 error. Contact ESG NextGen Support for help troubleshooting.
Not receiving ACKs	Credentials	Verify if credentials are in use and if ESG NextGen has them
	Firewall	Verify industry firewalls are allowing FDA IP's to connect
	AS2 gateway issues	ACK's are being delivered to industry AS2 gateway but rejected with a 400 or 500 error. Contact ESG NextGen Support for help troubleshooting.

6. References

- FDA ESG NextGen Support: ESGNGSupport@fda.hhs.gov

- For ESG NextGen encryption and signing certificates, email ESGNGSupport@fda.hhs.gov or download them from the ESG NextGen Gateway Configuration Information document linked in Step 3 of [ESG NextGen AS2 Account Set-Up Steps | FDA](#).
- For additional information pertaining to AS2, please reference [User Guides | FDA](#) and the AS2 User Guide section.