

BIONETICS

Summary of mutagenicity screening studies, host-mediated assay cytogenetics dominant
lethal assay-Contract FDA 71-268 & Compound FDA 71-8 (Potassium Nitrate)
11/24/72

P23

SUMMARY OF MUTAGENICITY
SCREENING STUDIES
CONTRACT FDA 71-268
COMPOUND FDA 71-8
POTASSIUM NITRATE
HOST-MEDIATED ASSAY
CYTOGENETICS
DOMINANT LETHAL ASSAY

1401 Wisconsin Avenue
Baltimore, Maryland
21204

LBI PROJECT #2311

SUMMARY OF MUTAGENICITY
SCREENING STUDIES
CONTRACT FDA 71-268
COMPOUND FDA 71-8
POTASSIUM NITRATE
HOST-MEDIATED ASSAY
CYTOGENETICS
DOMINANT LETHAL ASSAY

SUBMITTED TO

FOOD & DRUG ADMINISTRATION
DEPARTMENT OF HEALTH, EDUCATION AND WELFARE
ROCKVILLE, MARYLAND

SUBMITTED BY

LITTON BIONETICS, INC.
7315 WISCONSIN AVENUE
BETHESDA, MARYLAND

NOVEMBER 24, 1972



BIONETICS



BIONETICS

7315 Wisconsin Avenue, Bethesda, Maryland 20014 301 881-5600

November 24, 1972

Mr. Leonard Appleby, Contracting Officer
Department of Health, Education and Welfare
Public Health Service
Food and Drug Administration, CA-212
5600 Fishers Lane, Room 5C-13
Rockville, Maryland 20852

Reference: Contract FDA 71-268; LBI Project #2311

Dear Mr. Appleby:


Litton Bionetics, Inc. is pleased to submit a report for the referenced contract entitled "Mutagenicity Screening Studies" for compound FDA 71-8, Potassium Nitrate.

Included in this report are the results and raw data of the three tests conducted: Host-Mediated Assay; Cytogenetic Studies; and Dominant Lethal Assay. Eight (8) copies are being submitted for your review.

If there are any questions concerning this report, or, if additional information is required, please do not hesitate to contact us.

Sincerely yours,

LITTON BIONETICS, INC.

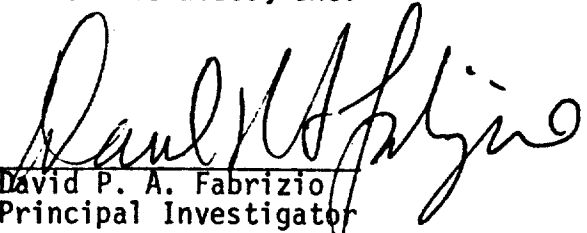
David P. A. Fabrizio
Principal Investigator

DPAF:11s
Enclosures (8)

TABLE OF CONTENTS

	Page No.
I. REPORT	1
A. Introduction	1
B. Objective	2
C. Compound	3
1. Test Material	3
2. Dosages	3
D. Methods	4
E. Summary	4
1. Host-Mediated Assay	4
2. Cytogenetics	4
a. <u>In vivo</u>	4
b. <u>In vitro</u>	4
3. Dominant Lethal Assay	5
F. Results and Discussion	5
1. Toxicity	5
a. <u>In vivo</u>	5
b. <u>In vitro</u>	6
c. Toxicity data sheets	7
2. Host-Mediated Assay	11
a. Host-mediated assay summary sheets	12
b. Host-mediated assay data sheets	17
3. Cytogenetics	42
a. <u>In vivo</u>	42
b. <u>In vitro</u>	42
c. Cytogenetics summary sheets	43
4. Dominant Lethal Study	47
a. Acute study	47
b. Subacute study	47
c. Dominant lethal assay summary tables	48
II. MATERIALS AND METHODS	65
A. Animal Husbandry	65
1. Animals (Rats and Mice)	65
2. Preparation of Diet	65
3. Husbandry	65
B. Dosage Determination	65
1. Acute LD ₅₀ and LD ₅ Determination	65
2. Subacute Studies	67
C. Mutagenicity Testing Protocols	67
1. Host-Mediated Assay	67
a. Acute study	68
b. Subacute study	70
c. <u>In vitro</u> study	70
2. Cytogenetic Study	71
a. <u>In vivo</u> study	71
b. <u>In vitro</u> study	73



BIONETICS

TABLE OF CONTENTS (continued)

Page No.

II. MATERIALS AND METHODS (continued)

C.	Mutagenicity Testing Protocols (continued)	
3.	Dominant Lethal Assay	75
D.	Supplementary Materials and Methods	76
1.	Host-Mediated Assay <u>In Vitro</u> and Formulae	76
a.	Bacterial <u>in vitro</u> plate tests ...	76
b.	<u>In vitro</u> for mitotic recombination	76
c.	Minimal medium (bacteria)	77
d.	Complete medium (bacteria)	78
e.	Complete medium (yeast)	78
2.	Cytogenetics <u>In Vitro</u> Preparation of Anaphase Chromosomes	79
3.	Statistical Analyses of Dominant Lethal Studies	80
a.	The fertility index	80
b.	Total number of implantations	80
c.	Total number of <u>corpora lutea</u>	80
d.	Preimplantation losses	80
e.	Dead implants	81
f.	One or more dead implants	81
g.	Two or more dead implants	81
h.	Dead implants per total implants .	81
E.	References	84
1.	Host-Mediated Assay	84
2.	Cytogenetics	84
3.	Dominant Lethal	85
F.	Abbreviations	86



BIONETICS

I. REPORT

A. Introduction

Litton Bionetics, Inc. (LBI) has investigated the possible mutagenicity of compounds selected and provided by the Food and Drug Administration under Contract 71-268. LBI's investigation utilized the three mammalian test systems herein described -- Host-Mediated Assay, Cytogenetic Studies and Dominant Lethal Assay. These tests provide information as to the types of genetic damage caused by environmental compounds -- pesticides, chemicals, food additives, drugs and cosmetics.

The Host-Mediated Assay is based upon the assumption that the action of a mutagen on the genetics of bacteria is similar to that in man. This is further strengthened by the use of an eukaryotic organism (Saccharomyces cerevisiae). Since the mutation frequencies are well established for the indicator organism, any deviation due to the action of the test compound is readily detectable. As some compounds are mutagenic in bacteria and not in the host animal, and vice versa, this test is able to differentiate an action which may have been due to hosts' ability to detoxify or potentiate a suspected mutagen. This action is dependent upon the ability of the compound to gain access to the peritoneal cavity. Coupled with the direct action of the compound on the indicator organism in vitro, the assay provides a clear insight into host-mediation of mutagenicity.

Cytogenetics provides a valuable tool for the direct observation of chromosomal damage in somatic cells. Alteration of the chromosome number and/or form in somatic cells may be an index of mutation. These studies utilized examination of bone marrow cells arrested in C-metaphase from rats exposed to the test compound as compared to positive and negative control animals. If mutational

changes occur, the types of damage expected due to the action of chemicals are structural rearrangements, breaks and other forms of damage to the chromosomal complement of the cells exposed.

For the in vitro cytogenetic studies, we have a more rapid and inexpensive means of determining chromosomal damage. This is accomplished by observing cells in anaphase. As the chromatids separate and move along the spindle, aberrations may occur. Chromatids which do not migrate to the daughter cells may lead to uneven distribution of parts or of entire chromatids (mitotic nondysjunction). These give rise to "side arm" bridges which have been interpreted as point stickiness or localized failures of chromosome duplication point errors. These aberrations (bridges, pseudochiasmata, multipolar cells, acentric fragments, etc.) are extremely sensitive indicators of genetic damage.

The Dominant Lethal Test is an accurate and sensitive measure of the amount and type of fetal wastage which may occur following administration of a potential mutagen. Dominant lethal mutations are indicators of lethal genetic lesions. The effects of mutagens on the chromosomal complement of the spermatozoa of treated males results in alterations of form and number of chromosomes. Structural rearrangements and aneuploidy may lead to the production of non-viable zygotes, early and late fetal deaths, abortions and congenital malformations. In addition, aberrations could lead to sterility or reduced reproductive capacity of the F_1 generation. The action of a mutagen on specific portions of spermatogenesis is also apparent in this test.

B. Objective

The purpose of these studies is to determine any mutagenic effect of the test compound by employing the Host-Mediated Assay, Cytogenetic Studies



and the Dominant Lethal Assay, both in vivo and in vitro tests are employed with the cytogenetic and microbial test systems. These tests and their descriptions are referenced in the Appendices A through F.

C. Compound

1. Test Material

Compound FDA 71-8, Potassium Nitrate, as supplied by the Food and Drug Administration.

2. Dosages

The animals employed, the determination of the dosage levels and the route of administration are contained in the technical discussion.

The dosage levels employed for compound FDA 71-8 are as follows for Cytogenetics Studies in vivo in rats.

Low Level	3.0 mg/kg
Intermediate Level	30 mg/kg
LD ₅	300 mg/kg
Negative Control	Saline
Positive Control (TEM*)	0.3 mg/kg

The dosage levels employed for compound FDA 71-8 are as follows for the Host-Mediated Assay in vivo in mice.

Low Level	3.0 mg/kg
Intermediate Level	30 mg/kg
LD ₅	300 mg/kg
Negative Control	Saline
Positive Control (EMS**)	350 mg/kg
(DMN***)	100 mg/kg

- * Triethylene Melamine
- ** Ethyl Methane Sulfonate
- *** Dimethyl Nitrosamine



The dosage levels employed for compound FDA 71-8 are as follows for the Dominant Lethal Assay in vivo in rats.

Low Level	3.0 mg/kg
Intermediate Level	30 mg/kg
LD ₅	300 mg/kg
Negative Control	Saline
Positive Control (TEM*)	0.5 mg/kg

The in vitro cytogenetics studies were performed employing three logarithmic dose levels.

Low Level	1 mcg/ml
Medium Level	10 mcg/ml
High Level	100 mcg/ml
Negative Control	Saline
Positive Control (TEM*)	0.1 mcg/ml

*Triethylene Melamine

The discussion of this test is contained in the technical discussion.

D. Methods

The protocols employed are explained in Appendices C and D.

E. Summary

1. Host-Mediated Assay

This compound was non-mutagenic at the dose levels tested in this study.

2. Cytogenetics

a. In vivo

The compound produced no detectable significant aberration of the bone marrow metaphase chromosomes of rats when administered orally at the dosage levels employed in this study.

b. In vitro

The compound produced no significant aberration



in the anaphase chromosomes of human tissue culture cells when tested at the dosage levels employed in this study.

3. Dominant Lethal Assay

Compound FDA 71-8 is considered to be non-mutagenic in the Dominant Lethal Study in rats employing the dosage levels used in this study.

F. Results and Discussion

1. Toxicity

a. In vivo

Compound FDA 71-8 was suspended in 0.85% saline and administered to 10 male rats by gastric intubation. The average weight of 26 animals was 340 grams and each received a dose of 5,000 mg/kg. All animals were found dead after 24 hours. Findings at necropsy indicated patchy stomach mucosa.

Dose levels of 100, 250, 500, 1,000, 2,000, and 5,000 mg/kg were selected to determine an acute LD₅₀.

The toxicity data is presented on the LD₅₀ reporting form using the Litchfield-Wilcoxon method with the enclosed graph. These are indicated in the report under Toxicity Data Sheets. The LD₅₀ was determined to be 650 mg/kg. The LD₅ dose level was derived from the raw data LD₅₀ probit line (uncorrected). The LD₅₀ derived from both the connected probit line and the uncorrected probit line were within confidence limits. The acute and sub-acute doses used were LD₅ - 300 mg/kg, intermediate level - 30 mg/kg, and the low level - 3 mg/kg. The data on the dose levels, number of animals and the necropsy findings are presented in the Toxicity Data Sheets.



b. In vitro

The compound was suspended in 0.85% sterile saline at the concentrations listed below. It was introduced into the tubes containing the WI-38 cells in a logarithmic phase of growth. The cells were observed for the presence of CPE and mitoses with the following results:

<u>Tube Number</u>	<u>No. of Cells</u>	<u>Conc. mcg/ml</u>	<u>CPE</u>	<u>Mitoses</u>
1	5×10^5	1000	+	-
2	5×10^5	1000	+	-
3	5×10^5	500	+	+
4	5×10^5	500	+	+
5	5×10^5	100	-	+
6	5×10^5	100	-	+
7	5×10^5	10	-	+
8	5×10^5	10	-	+
9	5×10^5	0.1	-	+
10	5×10^5	0.1	-	+

Since a CPE and also inhibition of mitoses were observed, a closer range of concentrations was employed.



c. TOXICITY DATA SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-8

POTASSIUM NITRATE



BIONETICS

TOXICITY DATA
CONTRACT FDA 71-268
COMPOUND FDA 71-8
POTASSIUM NITRATE

Solvent: 0.85% saline suspension

Animals: Male rats with an average body weight of 340 grams. All animals were observed for 10 days.

Range Finding:

	Dose mg/kg	No. Dead No. Animals	Necropsy and Day of Death
	5000	10/10	Found dead within 24 hours. Patchy stomach mucosa.
LD ₅₀	100	0/5	Day 5, reddened intestinal mucosa.
	250	1/5	Day 4, reddened intestinal mucosa.
	500	3/5	Day 3 (1), Day 4 (2), reddened intestinal mucosa.
	1000	3/5	Day 2 (3), reddened intestinal mucosa.
	2000	4/5	Day 1 (4), reddened intestinal mucosa.
	5000	4/5	Day 1 (4), reddened intestinal mucosa and stomach mucosa vascular.



LD₅₀ REPORTING FORM USING LITCHFIELD-WILCOXON METHOD

DOSE EFFECT CURVE FOR FDA 71-8

Potassium Nitrate

DOSE	# dead / # tested PROPORTION	OBSERVED PERCENT	EXPECTED PERCENT	OBS-EXPT PERCENT	CONTRIB. TO (CHI) ²
100	.5 / 5	.100	.138	-.038	.061
250	1 / 5	.200	.288	-.088	.189
500	3 / 5	.600	.436	.164	.547
1000	3 / 5	.600	.595	.005	.001
2000	4 / 5	.800	.742	.058	.088
5000	4 / 5	.800	.881	-.081	.313

Total animals = 30

Number Doses, K = 6

Animals/Dose = 5

(CHI)² for n of k-2 = 9.49

1 OBS-EXPT

Total = .434

(CHI)² = 1.199

Degrees of Freedom, n=k-2= 4

since 1.199 is less than 9.49,
therefore data not significantly
heterogeneous

LD₈₄ = 3,600

LD₅₀ = 650

LD₁₆ = 115

fLD₅₀ = $S \frac{2.77}{\sqrt{N!}} = \frac{5.595}{\sqrt{N!}} = \frac{5.595}{\sqrt{20}} = \frac{2.77}{\sqrt{20}} = 2.905$


LD₅₀ × fLD₅₀ = 1888.3

LD₅₀ = 223.7

fLD₅₀

LD₅₀ and 19/20 Confidence Limits = P(224 ≤ LD₅₀ ≤ 1888) = .95

Attached should be a plot of the dose-effect curve on log-probit paper.

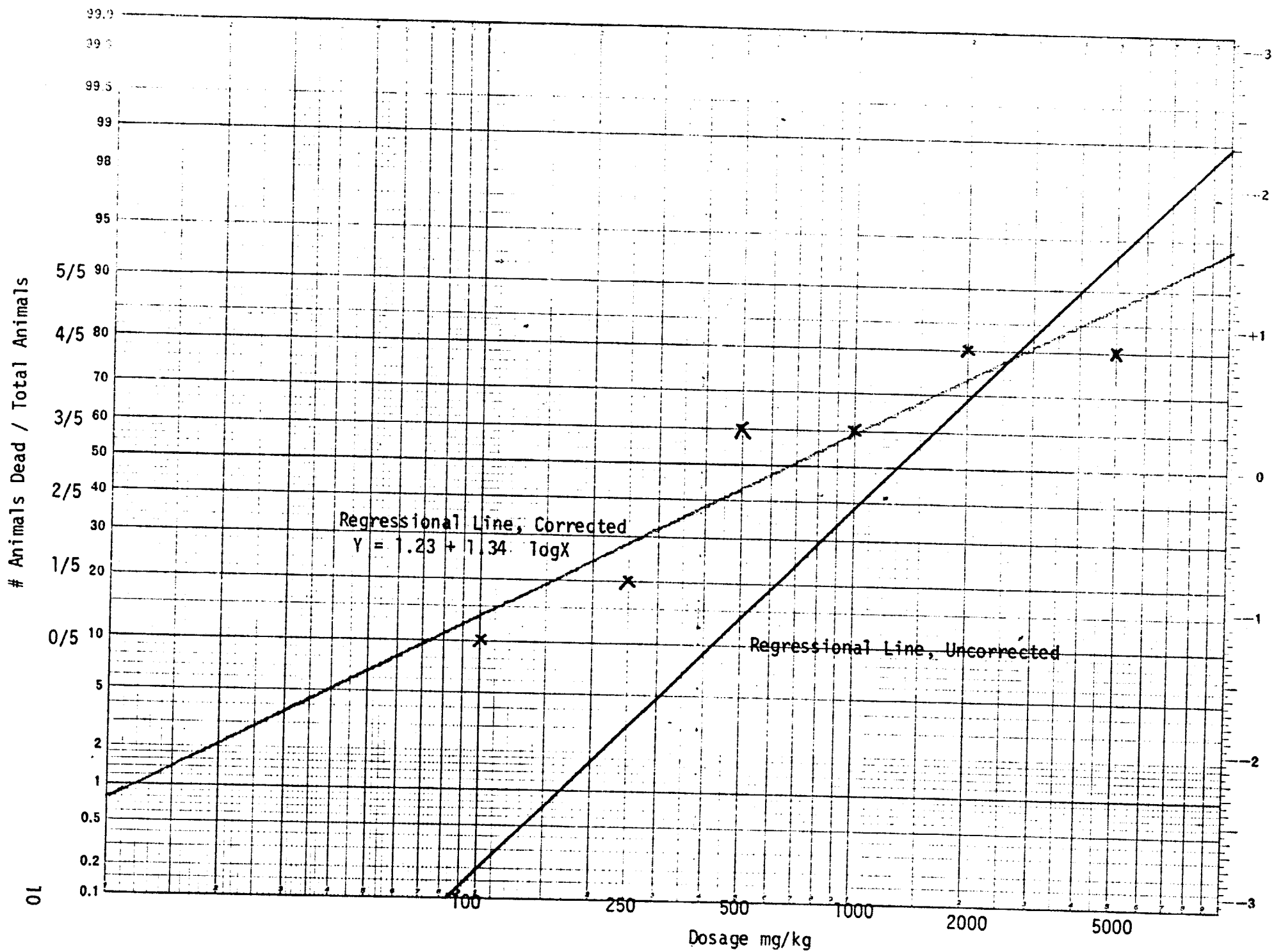

2. Host-Mediated Assay

Compound FDA 71-8 showed no significant increases in mutation frequencies when tested in vivo against Salmonella G-46 and TA-1530 and Saccharomyces D-3. The in vitro studies were also negative when tested against these organisms.



a. HOST-MEDIATED ASSAY SUMMARY SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-8.

POTASSIUM NITRATE



BIONETICS

HOST MEDIATED ASSAY

SUMMARY SHEET

OUTLIERS REMOVED

COMPOUND: FDA 71-8

	SALMONELLA		G-46		SACCHAROMYCES D-3	
	TA1530					
	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MRF (X 10E-5)	MRT/MRC
ACUTE						
NC	3.50		1.45		3.92	
PC	21.48	6.14	124.19	85.65	59.79	15.25
AL	2.96	.85	1.34	.92	4.03	1.03
AI	4.79	1.37	1.77	1.22	7.54	1.92
ALD5	2.17	.62	.95	.66	8.25	2.10
SUBACUTE						
NC	3.50		1.45		3.92	
SL	3.18	.91	1.99	1.37	5.65	1.44
SI	2.93	.84	3.19	2.20	5.45	1.39
SLD5	3.12	.89	3.44	2.37	4.56	1.16

IN VITRO

TA1530

G-46

% CONC

D-3
% SURVIVAL

R X 10E5

NC
PC

CSCX CSC85F 04 DEC 72 13:44:45 USER CFU007 190

CARDS IN 411 OUT 0 LINES 473 PROCESSING TIME 16.38 SECONDS

HOST MEDIATED ASSAY

SUMMARY SHEET

OUTLIERS INCLUDED

COMPOUND: FDA 71-8

	SALMONELLA		SACCHAROMYCES D-3			
	TA1530	G-46				
	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MRF (X 10E-5)	MRT/MRC
ACUTE						
NC	3.50		1.72		3.92	
PC	21.48	6.14	145.81	84.77	59.79	15.25
AU	4.00	1.14	1.85	.96	4.95	1.26
AI	4.79	1.37	2.13	1.24	7.54	1.92
ALDS	2.17	.62	1.05	.61	10.14	2.59
SUBACUTE						
NC	3.50		1.72		3.92	
SU	3.18	.91	1.99	1.16	6.60	1.68
SI	2.93	.84	3.19	1.85	5.45	1.39
SLDS	3.89	1.11	3.44	2.00	5.06	1.29

IN VITRO	TA1530	G-46	D-3
			% SURVIVAL
			R X 10E5
NC			
PC			

HOST MEDIATED ASSAY

SUMMARY SHEET

COMPOUND: FDA 71-8

SALMONELLA

SACCHAROMYCES D-3

TA1530

G-46

	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-5)	MRT/MRC
ACUTE						
NC	3.50		1.72		5.57	
PC	21.48	6.14	145.81	84.77	80.80	14.51
AU	4.00	1.14	1.65	.96	6.86	1.23
AI	4.79	1.37	2.13	1.24	10.14	1.82
ALD5	2.17	.62	1.05	.61	12.32	2.21
SUBACUTE						
NC	3.50		1.72		5.57	
SU	3.18	.91	1.99	1.16	3.51	1.53
SI	2.93	.84	3.19	1.85	7.00	1.26
SLD5	3.89	1.11	3.44	2.00	6.49	1.17

IN VITRO

TA1530

G-46

D-3

	% CONC	% SURVIVAL	R X 10E5
TCPD	-	-	-
NC	25	67	17
PC	10	68	341

CSCX CSC85F 22 NOV 72 21:40:59 USER CFU007 200

CARDS IN 75 OUT 0 LINES 50 PROCESSING TIME 2.96 SECONDS

HOST MEDIATED ASSAY (OUTLIERS REMOVED)

SUMMARY SHEET

COMPOUND: FDA 71-8

	TA1530		SALMONELLA G-46		SACCHAROMYCES D-3	
	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MRF (X 10E-5)	MRT/MRC
ACUTE						
NC	3.50		1.45		5.57	
PC	21.48	6.14	124.19	85.65	80.80	14.51
AU	2.96	.85	1.34	.92	5.09	.91
AI	4.79	1.37	1.77	1.22	10.14	1.82
ALD5	2.17	.62	.95	.66	12.32	2.21
SUBACUTE						
NC	3.50		1.45		5.57	
SU	3.18	.91	1.99	1.37	6.70	1.20
SI	2.93	.84	3.19	2.20	7.00	1.26
SLD5	3.12	.89	3.44	2.37	5.45	.98

IN VITRO

TA1530

G-46

% CONC

D-3
% SURVIVAL

R X 10E5

NC
PC

SAME AS ON PRECEDING SUMMARY SHEET

DATA CARDS ENCOUNTERED BY SYSTEM - IGNORED

CSCX CSC85F 24 NOV 72 15:31:46 USER CFU007 190
CARDS IN 74 OUT 0 LINES 48 PROCESSING TIME 2.89 SECONDS

READY

b. HOST-MEDIATED ASSAY DATA SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-8

POTASSIUM NITRATE



BIONETICS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: NEGATIVE CONTROL - WATER

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	36.10	6.02	14.00	2.33
2	24.20	4.03	13.00	3.22
3	35.10	5.52	5.00	.91
4	37.20	6.20	6.00	.97
5	12.00	2.00	15.00	7.50
6	30.10	5.02	8.00	1.59
7	9.00	1.50	12.00	8.00

NO. OF ANIMALS EQUALS 7
TOTAL CFU OUT OF RANGE EQUALS 3

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	4.33	10.43	3.50
RANGE	4.70	10.00	7.09
MAX	6.20	15.00	8.00
MIN	1.50	5.00	.91

NO OUTLIERS

CSCX CSC85F 21 NOV 72 17: 1:27 USER CFU007 200

CARDS IN 236 OUT 0 LINES 63 PROCESSING TIME 6.6 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: POSITIVE CONTROL + DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	30.00	5.00	72.00	14.40
2	12.90	2.15	24.00	11.16
3	15.20	2.53	109.00	43.03
4	33.10	5.52	112.00	20.30
5	57.60	9.60	22.00	2.29
6	12.00	2.00	72.00	36.00
7	15.90	2.65	32.00	12.03
8	19.50	3.25	106.00	32.61

NO. OF ANIMALS EQUALS 8
TOTAL CFU OUT OF RANGE EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	4.09	60.62	21.48
RANGE	7.60	90.00	40.73
MAX	9.60	112.00	43.03
MIN	2.00	22.00	2.29

NO OUTLIERS

CSCX CSC85F 21 NOV 72 17: 1:39 USER CFU007 200

CARDS IN 236 OUT 0 LINES 64 PROCESSING TIME 5.80 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: LOW - 3.0 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	7.30	1.22	2.00	1.64
2	49.80	8.30	10.00	1.20
3	30.00	5.00	14.00	2.80
4	12.10	2.02	13.00	6.45
5	16.10	3.02	12.00	3.98
6	33.10	5.52	12.00	2.18
7	14.20	2.37	5.00	2.11
8	18.00	3.00	10.00	3.33
9	7.30	1.22	15.00	12.33

NO. OF ANIMALS EQUALS 9
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.52	10.33	4.00
RANGE	7.08	13.00	11.12
MAX	8.30	15.00	12.33
MIN	1.22	2.00	1.20

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.80	9.75	2.96
RANGE	7.08	12.00	5.24
MAX	8.30	14.00	6.45
MIN	1.22	2.00	1.20

CSCX CSCBSF 21 NOV 72 17: 1:49 USER CFU007 200

CARDS IN 236 OUT 0 LINES 76 PROCESSING TIME 5.92 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: INTERMEDIATE - 30 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-6
1	20.50	3.42	15.00	4.39
2	11.00	1.83	12.00	6.55
3	9.70	1.62	11.00	6.80
4	11.00	1.83	10.00	5.45
5	9.90	1.65	8.00	4.85
6	10.30	1.72	12.00	6.99
7	33.00	5.50	5.00	.91
8	17.40	2.90	7.00	2.41

NO. OF ANIMALS EQUALS 8
NO. OF DEAD ANIMALS EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-6)
MEAN	2.56	10.00	4.79
RANGE	3.88	10.00	6.08
MAX	5.50	15.00	6.99
MIN	1.62	5.00	.91

NO OUTLIERS

CSCX CSC85F 21 NOV 72 17: 1159 USER CFU007 200

CARDS IN 232 OUT 0 LINES 64 PROCESSING TIME 6. 2 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: LD5 - 300 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10L-8
1	24.00	4.00	10.00	2.50
2	35.20	6.03	9.00	1.49
3	22.90	3.82	8.00	2.10
4	16.10	2.68	8.00	2.98
5	25.90	4.32	10.00	2.32
6	16.40	2.73	6.00	2.20
7	20.10	3.35	8.00	2.39
8	30.60	5.10	7.00	1.37

NO. OF ANIMALS EQUALS 8

NO. OF DEAD ANIMALS EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10L-8)
MEAN	4.00	8.25	2.17
RANGE	3.35	8.00	1.01
MAX	6.03	10.00	2.98
MIN	2.68	6.00	1.37

NO OUTLIERS

CSCX CSC85F 21 NOV 72 171 21 9 USER CFU007 200

CARDS IN 232 OUT 0 LINES 64 PROCESSING TIME 5.92 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-B

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: LOW - 3.0 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/D) X 10E-8
1	57.90	9.65	4.00	.41
2	9.50	1.58	11.00	6.95
3	43.30	7.22	8.00	1.11
4	35.10	5.05	9.00	1.54
5	49.80	8.30	10.00	1.20
6	11.10	1.85	12.00	6.49
7	9.70	1.62	8.00	4.95
8	13.30	2.22	10.00	4.51
9	55.20	9.70	14.00	1.44

NO. OF ANIMALS EQUALS 9

TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	5.33	9.56	3.18
RANGE	8.12	10.00	6.53
MAX	9.70	14.00	6.95
MIN	1.58	4.00	.41

NO OUTLIERS

CSCX CSCB5F 21 NOV 72 17: 4:35 USER CFU007 200

CARDS IN 236 OUT 0 LINES 65 PROCESSING TIME 6.11 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: INTERMEDIATE - 30 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	24.10	4.02	10.00	2.49
2	30.20	5.03	9.00	1.79
3	24.10	4.02	7.00	1.74
4	14.30	2.38	14.00	5.87
5	33.00	6.33	11.00	1.74
6	12.20	2.03	3.00	1.48
7	0.60	1.43	8.00	5.88
8	8.70	1.45	4.00	2.70

NO. OF ANIMALS EQUALS 8

NO. OF DEAD ANIMALS EQUALS 1

TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.34	6.25	2.93
RANGE	4.90	11.00	4.40
MAX	6.33	14.00	5.87
MIN	1.43	3.00	1.48

NO OUTLIERS

SCX CSC85F 21 NOV 72 17: 6:12 USER CFU007 200

CARDS IN 234 OUT 0 LINES 65 PROCESSING TIME 5.95 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: F0A 71-6

ORGANISM: SALMONELLA 141530

DOSE LEVEL: LD5 - 300 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 10, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E8/1.0ML	D MUTATION FRE (C/D) X 10E-6
1	16.20	2.70	2.00	.74
2	16.00	3.00	6.00	2.00
3	25.20	4.20	12.00	2.86
4	24.10	4.02	15.00	3.73
5	6.00	1.00	10.00	10.00 *
6	54.10	9.02	16.00	1.77
7	12.10	2.02	11.00	5.45
8	9.80	1.65	3.00	1.84
9	9.10	1.52	10.00	6.59

NO. OF ANIMALS EQUALS 9
NO. OF DEAD ANIMALS EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-6)
MEAN	3.23	9.44	3.89
RANGE	8.02	10.00	9.26
MAX	9.02	16.00	10.00
MIN	1.00	2.00	.74

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-6)
MEAN	3.51	9.38	3.12
RANGE	7.50	14.00	5.85
MAX	9.02	16.00	6.59
MIN	1.52	2.00	.74

CSCX C5C85F 22 NOV 72 18:23:19 USER CFU007 200

CARDS IN 234 OUT 0 LINES 76 PROCESSING TIME 5.99 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-B

ORGANISM: SALMONELLA G-46

DOSE LEVEL: NEGATIVE CONTROL - WATER

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/H) X 10E-8
1	19.50	3.25	5.00	1.54
2	20.50	3.42	8.00	2.34
3	12.50	2.08	8.00	3.84
4	13.60	2.27	6.00	2.65
5	47.50	7.92	2.00	.25
6	12.00	2.00	2.00	1.00
7	24.00	4.00	4.00	1.00
8	13.70	2.28	5.00	2.19
9	37.60	6.27	4.00	.64

NO. OF ANIMALS EQUALS 9
NO. OF DEAD ANIMALS EQUALS 0

1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.72	4.89	1.72
RANGE	5.92	6.00	3.59
MAX	7.92	8.00	3.84
MIN	2.00	2.00	.25

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.93	4.50	1.45
RANGE	5.92	6.00	2.39
MAX	7.92	8.00	2.65
MIN	2.00	2.00	.25

CSCX C5CB5F 21 NOV 72 16:53: 0 USER CFU007 200

CARDS IN 234 OUT 0 LINES 76 PROCESSING TIME 5.05 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA G-46

DOSE LEVEL: POSITIVE CONTROL -- DMN = 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	24.70	4.12	670.00	162.75
2	13.20	2.20	290.00	131.82
3	12.00	2.00	320.00	160.00
4	14.50	2.42	666.00	275.58
5	13.60	3.10	300.00	96.77
6	21.40	3.57	256.00	71.77
7	11.90	1.98	242.00	122.01

NO. OF ANIMALS EQUALS 7

NO. OF DEAD ANIMALS EQUALS 3

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	2.77	392.00	145.81
RANGE	2.13	428.00	203.81
MAX	4.12	670.00	275.58
MIN	1.98	242.00	71.77

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	2.63	346.33	124.19
RANGE	2.13	428.00	90.98
MAX	4.12	670.00	162.75
MIN	1.98	242.00	71.77

CSCX CSC85F 21 NOV 72 16147147 USER CFU007 200

CARDS IN 230 OUT 0 LINES 74 PROCESSING TIME 6.4 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LOW - 3.0 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	43.00	7.17	5.00	.70
2	25.30	4.22	5.00	1.19
3	6.30	1.05	4.00	3.81
4	13.50	2.25	5.00	2.22
5	53.60	8.93	10.00	1.12
6	13.60	2.27	4.01	1.77
7	49.70	8.28	10.00	1.21
8	10.00	1.67	2.00	1.20

NO. OF ANIMALS EQUALS 8
TOTAL CFU OUT OF RANGE EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	4.48	5.63	1.65
RANGE	7.88	8.00	3.11
MAX	8.93	10.00	3.81
MIN	1.05	2.00	.70

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	4.97	5.86	1.34
RANGE	7.27	8.00	1.52
MAX	8.93	10.00	2.22
MIN	1.67	2.00	.70

CSCX C5C85F 21 NOV 72 16:52:50 USER CFU007 200

CARDS IN 236 OUT 0 LINES 75 PROCESSING TIME 5.75 SECONDS

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA G-46

DOSE LEVEL: INTERMEDIATE - 30.0 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.0ML	B TOTAL CFU X 10E6/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	11.90	1.98	3.00	1.51
2	11.80	1.97	3.00	1.53
3	25.00	4.17	11.00	2.64
4	15.70	2.52	6.00	2.29
5	23.20	3.87	5.00	1.29
6	23.70	3.95	6.00	1.52
7	15.70	2.28	3.00	1.31
8	25.80	4.30	10.00	2.33
9	13.40	2.23	12.00	5.37
10	11.70	1.95	3.00	1.54

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	2.93	6.20	2.13
RANGE	2.35	9.00	4.08
MAX	4.30	12.00	5.37
MIN	1.95	3.00	1.29

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.01	5.56	1.77
RANGE	2.35	8.00	1.35
MAX	4.30	11.00	2.64
MIN	1.95	3.00	1.29

CSCX CSC85F 21 NOV 72 16:53:10 USER CFU007 200

CARDS IN 236 OUT 0 LINES 76 PROCESSING TIME 5.87 SECONDS

DOSE LEVEL: HIGH + 300 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E6/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	31.90	5.32	6.00	1.13
2	47.00	7.83	15.00	1.91
3	36.30	6.05	8.00	1.32
4	27.30	4.55	2.00	.44
5	17.70	2.95	2.00	.68
6	23.30	3.88	5.00	1.29
7	44.70	7.45	6.00	.61
8	39.00	6.50	8.00	1.23
9	25.00	4.17	5.00	1.20
10	49.30	8.22	4.00	.49

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E6)	COL. D (X 10E-8)
MEAN	5.69	6.10	1.05
RANGE	5.27	13.00	1.48
MAX	8.22	15.00	1.91
MIN	2.95	2.00	.44

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E6)	COL. D (X 10E-8)
MEAN	5.45	5.11	.95
RANGE	5.27	6.00	.88
MAX	8.22	8.00	1.32
MIN	2.95	2.00	.44

CSCX C5C85F, 21 NOV 72 16:53:56 USER CFU007 200

CARDS IN 236 OUT 0 LINES 76 PROCESSING TIME 6.21 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LOW - 3.0 mg/kg

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	15.40	2.57	9.00	3.51
2	15.20	2.53	6.00	2.37
3	27.90	4.65	2.00	.43
4	21.90	3.65	10.00	2.74
5	30.00	5.00	4.00	.80
6	171.20	28.53	12.00	.42
7	122.00	20.33	10.00	.49
8	8.40	1.40	5.00	3.57
9	13.00	2.17	9.00	4.15
10	21.50	3.58	5.00	1.40

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.44	7.20	1.99
RANGE	27.13	10.00	3.73
MAX	28.53	12.00	4.15
MIN	1.40	2.00	.42

NO OUTLIERS

ICX CSC85F 21 NOV 72 16:54:27 USER CFU007 200

CARDS IN 236 OUT 0 LINES 65 PROCESSING TIME 5.87 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-B

ORGANISM: SALMONELLA G-46

DOSE LEVEL: INTERMEDIATE - 30.0 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E6/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	49.00	8.17	19.00	2.33
2	23.70	3.95	20.00	5.06
3	65.10	10.85	40.00	3.69
4	14.70	2.45	6.00	2.45
5	20.40	3.40	10.00	4.71
6	17.00	2.83	7.00	2.47
7	11.80	1.97	4.00	2.03
8	14.40	2.40	9.00	3.75
9	13.70	2.28	5.00	2.19

NO. OF ANIMALS EQUALS 9

NO. OF DEAD ANIMALS EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	4.26	14.00	3.19
RANGE	8.88	30.00	3.03
MAX	10.85	40.00	5.06
MIN	1.97	4.00	2.03

NO OUTLIERS

CSCX CSC85F 21 NOV 72 16:54:37 USER CFU007 200

ARDS IN 234 OUT 0 LINES 65 PROCESSING TIME 5.86 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-3

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LD5 - 300 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 31, 1972

ANIMAL NUMBER	A RAE CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	15.00	2.50	10.00	4.00
2	11.40	1.90	7.00	3.68
3	22.00	3.67	13.00	3.55
4	19.10	3.18	15.00	4.71
5	78.00	13.00	18.00	1.33
6	15.00	2.50	6.00	2.40
7	12.50	2.08	9.00	4.52

NO. OF ANIMALS EQUALS 7

NO. OF DEAD ANIMALS EQUALS 3

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	4.12	11.14	3.44
RANGE	11.10	12.00	3.33
MAX	13.00	18.00	4.71
MIN	1.90	6.00	1.33

NO OUTLIERS

CSCX CSC85F 21 NOV 72 16:54:34 USER CFU007 200

CARDS IN 230 OUT 0 LINES 63 PROCESSING TIME 5.93 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: NEGATIVE CONTROL - WATER

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	130.00	.13	1.00	7.69
2	190.00	.19	1.00	5.26
3	591.00	.59	2.00	3.33
4	500.00	.50	1.00	2.00
5	200.00	.20	2.00	10.00
6	351.00	.35	1.00	2.85
7	142.00	.14	2.00	14.08
8	750.00	.75	2.00	2.67
9	461.00	.46	1.00	2.17

TOTAL 3.31 13.00

NO. OF ANIMALS EQUALS 9
TOTAL SCREENED OUT OF RANGE EQUALS 1

MEAN C/MEAN B = 3.92

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.37	1.44	5.57
RANGE	.62	1.00	12.08
MAX	.75	2.00	14.08
MIN	.13	1.00	2.00

NO OUTLIERS

CSCX CSC85F 21 NOV 72 17: 5: 5 USER CFU007 200

CARDS IN 236 OUT 0 LINES 70 PROCESSING TIME 6.6 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: POSITIVE CONTROL - EMS - 350 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	500.00	.50	14.00	28.00
2	211.00	.21	8.00	37.91
3	160.00	.16	16.00	100.00
4	143.00	.14	12.00	63.92
5	252.00	.25	11.00	43.85
6	141.00	.14	19.00	134.75
7	133.00	.13	18.00	135.34
8	521.00	.52	10.00	19.19
9	180.00	.18	26.00	144.44
TOTAL		2.24	134.00	

NO. OF ANIMALS EQUALS 9
TOTAL SCREENED OUT OF RANGE EQUALS 1

MEAN C/MEAN B = 59.79

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.25	14.89	80.80
RANGE	.39	18.00	125.25
MAX	.52	26.00	144.44
MIN	.13	8.00	19.19

NO OUTLIERS

CSCX CSC85F 21 NOV 72 17: 5:15 USER CFU007 200

CARDS IN 236 OUT 0 LINES 70 PROCESSING TIME 6.8 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LOW ~ 3.0 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	261.00	.26	2.00	7.66
2	143.00	.14	3.00	20.98
3	284.00	.28	.00	.00
4	354.00	.35	2.00	5.65
5	152.00	.15	1.00	6.58
6	613.00	.61	1.00	1.64
7	140.00	.14	1.00	7.14
8	200.00	.20	2.00	10.00
9	481.00	.48	1.00	2.05
TOTAL		2.62	13.00	

NO. OF ANIMALS EQUALS 9
NO. OF DEAD ANIMALS EQUALS 1

MEAN C/MEAN B = 4.95

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.29	1.44	6.86
RANGE	.47	3.00	20.98
MAX	.61	3.00	20.98
MIN	.14	.00	.00

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 4.03

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.31	1.25	5.09
RANGE	.47	2.00	10.00
MAX	.61	2.00	10.00
MIN	.14	.00	.00

CSCX CSC85F 21 NOV 72 17: 5:26 USER CFU007 200

CARDS IN 234 OUT 0 LINES 04 PROCESSING TIME 6.28 SECONDS 36

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: INTERMEDIATE - 30 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	460.00	.46	2.00	4.35
2	200.00	.20	3.00	15.00
3	100.00	.10	1.00	10.00
4	190.00	.19	4.00	21.05
5	300.00	.30	.00	.00
6	490.00	.49	2.00	4.08
7	100.00	.10	2.00	20.00
8	150.00	.15	1.00	6.67
TOTAL		1.99	15.00	

NO. OF ANIMALS EQUALS 8

TOTAL SCREENED OUT OF RANGE EQUALS 2

MEAN C/MEAN B = 7.54

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.25	1.88	10.14
RANGE	.39	4.00	21.05
MAX	.49	4.00	21.05
MIN	.10	.00	.00

NO OUTLIERS

SCX CSC85F 21 NOV 72 17: 6:22 USER CFU007 200

CARDS IN 236 OUT 0 LINES 69 PROCESSING TIME 6.40 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LD5 - 300 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	191.00	.19	1.00	5.24
2	313.00	.31	.00	.00
3	402.00	.40	5.00	12.44
4	224.00	.22	2.00	8.93
5	500.00	.50	3.00	6.00
6	130.00	.13	2.00	15.38
7	200.00	.20	2.00	10.00
8	252.00	.25	1.00	3.97
9	211.00	.21	4.00	18.96
10	143.00	.14	6.00	42.25
TOTAL		2.56	26.00	

NO. OF ANIMALS EQUALS 10

MEAN C/MEAN B = 10.14

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.26	2.60	12.32
RANGE	.37	6.00	42.25
MAX	.50	6.00	42.25
MIN	.13	.00	.00

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 8.25

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.27	2.22	8.99
RANGE	.37	5.00	18.96
MAX	.50	5.00	18.96
MIN	.13	.00	.00

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LOW - 3.0 MG/KG

TREATMENT: IN VIVO, ORAL, SUBCUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	280.00	.28	2.00	7.14
2	650.00	.65	1.00	1.54
3	121.00	.12	3.00	24.79 *
4	143.00	.14	1.00	6.99
5	330.00	.33	2.00	6.06
6	260.00	.26	2.00	7.69
7	150.00	.15	1.00	6.67
8	134.00	.13	1.00	7.46
9	182.00	.18	2.00	10.99
10	173.00	.17	1.00	5.76

TOTAL 2.42 16.00

NO. OF ANIMALS EQUALS 10

MEAN C/MEAN B = 6.60

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.24	1.60	6.51
RANGE	.53	2.00	23.25
MAX	.65	3.00	24.79
MIN	.12	1.00	1.54

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 5.65

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.26	1.44	6.70
RANGE	.52	1.00	9.45
MAX	.65	2.00	10.99
MIN	.13	1.00	1.54

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-8

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: INTERMEDIATE - 30 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMB/ANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	370.00	.37	2.00	5.41
2	360.00	.36	1.00	2.78
3	560.00	.56	1.00	1.79
4	162.00	.16	1.00	6.17
5	363.00	.36	2.00	5.51
6	304.00	.30	1.00	3.29
7	164.00	.16	2.00	12.20
8	222.00	.22	2.00	9.01
9	113.00	.11	1.00	8.85
10	133.00	.13	2.00	15.04
TOTAL		2.75	15.00	

NO. OF ANIMALS EQUALS 10

MEAN C/MEAN B = 5.45

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.28	1.50	7.00
RANGE	.45	1.00	13.25
MAX	.56	2.00	15.04
MIN	.11	1.00	1.79
NO OUTLIERS			

CSCX CSCBSF 21 NOV 72 17: 6143 USER CFU007 200

CARDS IN 236 OUT 0 LINES 70 PROCESSING TIME 6. 3 SECONDS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-3

ORGANISM: SACCHAROMYCES D-5

DOSE LEVEL: LDS - 300 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 3, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	360.00	.36	3.00	8.33
2	391.00	.39	1.00	2.56
3	650.00	.65	5.00	7.69
4	700.00	.70	2.00	2.56
5	220.00	.22	2.00	9.09
6	202.00	.20	3.00	14.85
7	741.00	.74	2.00	2.70
8	753.00	.75	2.00	2.66
9	130.00	.13	1.00	7.69

TOTAL 4.15 21.00

NO. OF ANIMALS EQUALS 9
NO. OF DEAD ANIMALS EQUALS 1

MEAN C/MEAN B = 5.06

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.46	2.33	5.49
RANGE	.62	5.00	12.29
MAX	.75	5.00	14.85
MIN	.13	1.00	2.56

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 4.56

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.49	2.25	5.45
RANGE	.62	5.00	6.53
MAX	.75	5.00	9.09
MIN	.13	1.00	2.56

CSCX CSC85F 21 NOV 72 17: 6:53 USER CF0007 200

CARDS IN 234 OUT 0 LINES 64 PROCESSING TIME 6.10 SECONDS

3. Cytogenetics

a. In vivo

(1) Acute study

The negative control groups were within normal control values. Two groups of the experimental compound dose levels were somewhat higher than the negative control group. These were the 48-hour low level with 10% breaks and the 48-hour high level with 8% breaks. These are probably not significant and in the absence of an appropriate statistical analysis can be regarded only as "high". The lack of any observed reunions in these groups indicates that the effect, if any, of the compound is not severe. The positive control group produced the severe chromosomal damage expected. Mitotic indices were within normal values.

(2) Subacute study

The negative control group was within normal limits. The low dosage group contained 4% cells with breaks as did the high dosage level group. The medium level dosage groups contained a slightly elevated percentage of cells with breaks - 8%, but is not considered to be significant, especially with the absence of reunions.

b. In vitro

The negative control groups contained 1% cells with a bridge. All other groups were negative except the high level group, which contained 4% of cells with acentric fragments. While this is higher than the negative control group, it is within normal control values as observed in this laboratory. The positive control should expect severe damage.

c. CYTOGENETICS SUMMARY SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-8

POTASSIUM NITRATE



BIONETICS

FDA 71-8
ACUTE STUDY
METAPHASE SUMMARY SHEET

<u>Compound</u>	<u>Dosage (mg/kg)</u>	<u>Time*</u>	<u>No. of Animals</u>	<u>No. of Cells</u>	<u>Mitotic Index %</u>	<u>% Cells with Breaks</u>	<u>% Cells with Reunions</u>	<u>% Cells other Aber.**</u>	<u>% Cells with Aber.</u>
Negative Control	Saline	6	3	150	9	3	0	0	3
	Saline	24	3	150	6	2	0	0	2
	Saline	48	3	150	10	6	0	0	6
Low Level	3	6	5	250	9	6	0	0	6
	3	24	5	250	12	8	0	0	8
	3	48	5	250	10	10	0	0	10
Intermediate	30	6	5	250	10	4	0	0	4
	30	24	5	250	11	6	0	0	6
	30	48	5	250	12	3	0	0	3
LD ₅ Level	300	6	5	250	6	2	0	0	2
	300	24	5	250	8	6	0	0	6
	300	48	5	250	8	8	0	0	8
Positive Control (TEM)***	0.30	48	5	250	4	24	12	6 (a)	31

*Time of sacrifice after injection (hours).

**Cells that have polyploidy (P), pulverization (pp), or greater than 10 aberrations (a).

***Acute dose only one time. Sample taken at 48 hours.

FDA 71-8
SUBACUTE STUDY
METAPHASE SUMMARY SHEET

<u>Compound</u>	<u>Dosage*</u> <u>(mg/kg)</u>	<u>No. of</u> <u>Animals</u>	<u>No. of</u> <u>Cells</u>	<u>Mitotic</u> <u>Index %</u>	<u>%</u> <u>Cells</u> <u>with</u> <u>Breaks</u>	<u>%</u> <u>Cells</u> <u>with</u> <u>Reunions</u>	<u>% Cells</u> <u>other</u> <u>Aber.**</u>	<u>%</u> <u>Cells</u> <u>with</u> <u>Aber.</u>
Negative Control	Saline	3	150	14	6	0	0	6
Low	3	5	250	6	4	0	0	4
Medium	30	5	250	6	8	0	0	8
LD ₅	300	5	250	5	4	0	0	4

*Dosage 1x/day x 5 days

**Cells that have polyploidy (P), pulverization (pp), or greater than 10 aberrations (a).

FDA 71-8
ANAPHASE SUMMARY SHEET

<u>Compound</u>	<u>Dosage** (mcg/ml)</u>	<u>Mitotic Index</u>	<u>No. of Cells</u>	<u>% Cells with Acentric Frag.</u>	<u>% Cells with Bridges</u>	<u>% Multipolar Cells</u>	<u>% Cells Other Aber.*</u>	<u>% Cells with Aber.</u>
Low Level	1	1	100	0	0	0	0	0
Medium Level	10	2	100	0	0	0	0	0
High Level	100	1	100	4	0	0	0	4
Negative Control	Saline	3	100	0	1	0	0	1
Positive Control (TEM)	0.1	1	100	10	11	0	2 (pp)	22

*Cells that have polyploidy (P), pulverization (pp), or greater than 10 aberrations (a).

**Cells harvested 24 hours after addition of the compound.

4. Dominant Lethal Study

a. Acute study

In general, significant differences between the negative control and experimental groups were shown in a few instances, but no strong indications of change were seen.

b. Subacute study

The results were similar to those found in the acute study.

c. DOMINANT LETHAL ASSAY
SUMMARY TABLES
CONTRACT FDA 71-268
COMPOUND FDA 71-8
POTASSIUM NITRATE



BIONETICS

TABLE 1
COMPOUND 8 STUDY ACUTE
FERTILITY INDEX

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
		1	43/ 60=0.72	13/20=0.65	11/20=0.55	14/20=0.70	13/20=0.65	15/20=0.75
		2	47/ 60=0.79	12/20=0.60	10/20=0.50 *	19/20=0.95**	16/20=0.80	17/20=0.85
		3	53/ 60=0.89	14/20=0.70	13/20=0.65 *	17/20=0.85	18/20=0.90	16/20=0.80
		4	55/ 60=0.92	12/20=0.60 **	11/20=0.55 **	14/20=0.70 *	12/20=0.60 **	17/20=0.85
		5	52/ 60=0.87	15/20=0.75	15/20=0.75	17/20=0.85	18/20=0.90	17/20=0.85
		6	51/ 60=0.85	15/20=0.75	15/20=0.75	17/20=0.85	17/20=0.85	19/20=0.95
		7	52/ 60=0.87	15/20=0.75	15/20=0.75	15/20=0.75	14/20=0.70	19/20=0.95
		8	52/ 60=0.87	16/20=0.80	15/19=0.79	15/19=0.79	17/20=0.85	17/20=0.85

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE II
COMPOUND 8 STUDY ACUTE

AVERAGE NUMBER OF IMPLANTATIONS PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
		1	517/ 43=12.0	158/13=12.2	126/11=11.5	160/14=11.4	159/13=12.2	194/15=12.9
		2	547/ 47=11.6	155/12=12.9	121/10=12.1	241/19=12.7	194/16=12.1	197/17=11.6
		3	624/ 53=11.8	181/14=12.9	149/13=11.5	231/17=13.6 **@I	213/18=11.8	215/16=13.4 **@
		4	642/ 55=11.7	135/12=11.3	133/11=12.1	163/14=11.6	156/12=13.0 @I	197/17=11.6
		5	619/ 52=11.9	182/15=12.1	192/15=12.8	206/17=12.1	215/18=11.9	194/17=11.4
		6	608/ 51=11.9	179/15=11.9	179/15=11.9	200/17=11.8	187/17=11.0	228/19=12.0
		7	634/ 52=12.2	182/15=12.1	174/15=11.6	175/15=11.7	171/14=12.2	223/19=11.7
6 II 6611 6611		8	605/ 52=11.6	197/16=12.3	178/15=11.9	170/15=11.3	232/17=13.7 **@I	196/17=11.5

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

@ AND * = TWO-TAILED TEST
I AND @ = ONE-TAILED TEST

ONE I, @, *, @ = SIGNIFICANT AT P LESS THAN 0.05
TWO I, @, *, @ = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
@, I SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE III
COMPOUND 8 STUDY ACUTE

AVERAGE CORPORA LUTEA PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
6	1	1	546/ 43=12.7	188/13=14.5 *@I	161/11=14.6 *@@I	189/14=13.5	188/13=14.5 **@@I	204/15=13.6
		2	593/ 47=12.6	175/12=14.6 *@I	149/10=14.9 *@@I	265/19=14.0 *@I	222/16=13.9	237/17=13.9
6	1	3	673/ 53=12.7	188/14=13.4	164/13=12.6	244/17=14.4 **@@I	251/18=13.9 @I	224/16=14.0 *@I
		4	689/ 55=12.5	160/12=13.3	155/11=14.1	174/14=12.4	165/12=13.8	203/17=11.9
		5	666/ 52=12.8	190/15=12.7	198/15=13.2	207/17=12.2	224/18=12.4	203/17=11.9
		6	647/ 51=12.7	179/15=11.9	202/15=13.5@I	219/17=12.9	201/17=11.8	229/19=12.1
		7	664/ 52=12.8	188/15=12.5	181/15=12.1	185/15=12.3	178/14=12.7	225/19=11.8 @D
6	11	8	660/ 52=12.7	202/16=12.6	184/15=12.3	173/15=11.5 *@D	232/17=13.7	202/17=11.9

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

@ AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, @, *, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, @, *, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
@, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (READING OF COLUMN)

TABLE IV
COMPOUND 8 STUDY ACUTE

AVERAGE PREIMPLANTATION LOSSES PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
8811	1	1	29/ 43= 0.7	30/13= 2.3 *01	35/11= 3.2	29/14= 2.1 **001	29/13= 2.2 *001	10/15= 0.70D
		2	46/ 47= 1.0	20/12= 1.7	28/10= 2.8 *01	24/19= 1.3	28/16= 1.8 *01	40/17= 2.4
8811 8811 8 1 8811		3	49/ 53= 0.9	7/14= 0.5	15/13= 1.2	13/17= 0.8	38/18= 2.1**001 *001	9/16= 0.6
		4	47/ 55= 0.9	25/12= 2.1	22/11= 2.0	11/14= 0.8	9/12= 0.8	6/17= 0.4*00
!		5	47/ 52= 0.9	8/15= 0.5	6/15= 0.4 0D	1/17= 0.1 **00D	9/18= 0.5	9/17= 0.5
		6	39/ 51= 0.8	0/15= 0.0 **00D	23/15= 1.5**001	19/17= 1.1**001	14/17= 0.801	1/19= 0.1 **01
		7	30/ 52= 0.6	6/15= 0.4	7/15= 0.5	10/15= 0.7	7/14= 0.5	2/19= 0.10D **01
8811 8811		8	55/ 52= 1.1	5/16= 0.3 *00D	6/15= 0.4 *0D	3/15= 0.2 **00D	0/17= 0.00D **00D	6/17= 0.4 *001

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

8 AND * = TWO-TAILED TEST
! AND 0 = ONE-TAILED TEST

ONE !, 0, 0, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, 0, 0, * = SIGNIFICANT AT P LESS THAN 0.01

*, 0 SIGNIFICANTLY DIFFERENT FROM CONTROL
8, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (READING OF COLUMN)

TABLE V
COMPOUND 8 STUDY ACUTE

AVERAGE RESORPTIONS (DEAD IMPLANTS) PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
		1	9/ 43=0.21	4/13=0.31	5/11=0.46	9/14=0.65	7/13=0.54	3/15=0.20
E I		2	20/ 47=0.43	7/12=0.59	6/10=0.60	16/19=0.85 *@I	14/16=0.88	26/17=1.53**@I **@I
I E I		3	25/ 53=0.48	17/14=1.22 **@I	10/13=0.77	7/17=0.42*@@D	20/18=1.12 @I	23/16=1.44 *@I
I		4	27/ 55=0.50	4/12=0.34	7/11=0.64	16/14=1.15	10/12=0.84	50/17=2.95**@I **@I
		5	28/ 52=0.54	10/15=0.67	4/15=0.27	9/17=0.53	9/18=0.50	30/17=1.77*@I **@I
E I E I		6	27/ 51=0.53	7/15=0.47	9/15=0.60	14/17=0.83	2/17=0.12 **@@D	4/19=0.22 *@D
E E I I E I I		7	32/ 52=0.62	2/15=0.14 **@@D	3/15=0.20 *@D	6/15=0.40	0/14=0.0 **@@D	2/19=0.11 **@@D
I I E I I E I		8	30/ 52=0.58	8/16=0.50	12/15=0.80	13/15=0.87	25/17=1.48@I @I	13/17=0.77

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

E AND * = TWO-TAILED TEST
I AND @ = ONE-TAILED TEST

ONE !, E, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, E, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
E, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VI
COMPOUND 8 STUDY ACUTE

PROPORTION OF FEMALES WITH ONE OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
		1	9/ 43=0.21	4/13=0.31	5/11=0.46	4/14=0.29	3/13=0.24	3/15=0.20
1		2	14/ 47=0.30	5/12=0.42	4/10=0.40	12/19=0.64	8/16=0.50	14/17=0.83*
1						*		**
		3	16/ 53=0.31	10/14=0.72	7/13=0.54	6/17=0.36*	10/18=0.56	9/16=0.57
				**				
		4	21/ 55=0.39	4/12=0.34	6/11=0.55	8/14=0.58	7/12=0.59	16/17=0.95**
								**
		5	18/ 52=0.35	6/15=0.40	4/15=0.27	9/17=0.53	6/18=0.34	13/17=0.77*
								**
1		6	21/ 51=0.42	5/15=0.34	5/15=0.34	9/17=0.53	2/17=0.12	4/19=0.22
1							*	
11	11	7	22/ 52=0.43	2/15=0.14	2/15=0.14	4/15=0.27	0/14=0.0	1/19=0.06
11	11			*	*		**	**
		8	20/ 52=0.39	6/16=0.38	9/15=0.60	9/15=0.60	10/17=0.59	8/17=0.48

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VII
COMPOUND 8 STUDY ACUTE

PROPORTION OF FEMALES WITH TWO OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG	POSITIVE CONTROL
11	1	1	0/ 43=0.0	0/13=0.0	0/11=0.0	1/14=0.08	2/13=0.16	0/15=0.0
11	1						**	
		2	6/ 47=0.13	2/12=0.17	1/10=0.10	4/19=0.22	3/16=0.19	9/17=0.53*
								**
		3	7/ 53=0.14	5/14=0.36	3/13=0.24	1/17=0.06*	5/18=0.28	6/16=0.38
								*
		4	6/ 55=0.11	0/12=0.0	1/11=0.10	1/14=0.08	2/12=0.17	11/17=0.65**
								**
		5	8/ 52=0.16	3/15=0.20	0/15=0.0	0/17=0.0	3/18=0.17	7/17=0.42
								*
		6	6/ 51=0.12	2/15=0.14	3/15=0.20	3/17=0.18	0/17=0.0	0/19=0.0
		7	6/ 52=0.12	0/15=0.0	1/15=0.07	2/15=0.14	0/14=0.0	1/19=0.06
		8	8/ 52=0.16	1/16=0.07	3/15=0.20	4/15=0.27	6/17=0.36*	2/17=0.12

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE 1,* = SIGNIFICANT AT P LESS THAN 0.05
TWO 1,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

1 SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 8

TABLE VIII
STUDY ACUTE

WEEK	DEAD IMPLANTS / TOTAL IMPLANTS						POSITIVE CONTROL
	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG		
1	9/ 517=0.02	4/158=0.03	5/126=0.04	9/160=0.06	7/159=0.05		3/194=0.02
2	20/ 547=0.04	7/155=0.05	6/121=0.05	16/241=0.07	14/194=0.08		26/197=0.14
3	25/ 624=0.05	17/181=0.10	10/149=0.07	7/231=0.04	20/213=0.10		23/215=0.11
4	27/ 642=0.05	4/135=0.03	7/133=0.06	16/163=0.10	10/156=0.07		50/197=0.26
5	28/ 619=0.05	10/182=0.06	4/192=0.03	9/206=0.05	9/215=0.05		30/194=0.16
6	27/ 608=0.05	7/179=0.04	9/179=0.06	14/200=0.07	2/187=0.02		4/228=0.02
7	32/ 634=0.06	2/182=0.02	3/174=0.02	6/175=0.04	0/171=0.0		2/223=0.01
8	30/ 605=0.05	8/197=0.05	12/178=0.07	13/170=0.08	25/232=0.11		13/196=0.07

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE HISTORICAL CONTROL GROUP

* = TWO-TAILED TEST

@ = ONE-TAILED TEST

ONE *.@ = SIGNIFICANT AT P LESS THAN 0.05

TWO *.@ = SIGNIFICANT AT P LESS THAN 0.01

*,@ SIGNIFICANTLY DIFFERENT FROM CONTROL

TABLE I
COMPOUND 8 STUDY SUBACUTE

FERTILITY INDEX

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
		1	44/ 60=0.74	14/20=0.70	14/20=0.70	9/20=0.45 *	14/20=0.70
		2	44/ 60=0.74	16/20=0.80	16/20=0.80	16/20=0.80	12/20=0.60
		3	48/ 60=0.80	15/20=0.75	19/20=0.95	16/20=0.80	16/20=0.80
		4	48/ 60=0.80	15/20=0.75	16/20=0.80	18/20=0.90	17/20=0.85
		5	48/ 60=0.80	15/20=0.75	15/20=0.75	16/20=0.80	17/20=0.85
		6	50/ 60=0.84	17/20=0.85	18/20=0.90	18/19=0.95	15/20=0.75
!!	!!	7	49/ 58=0.85	19/20=0.95	16/20=0.80	17/20=0.85	10/20=0.50**
!!	!!						**

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE II
COMPOUND 8 STUDY SUBACUTE

AVERAGE NUMBER OF IMPLANTATIONS PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
		1	492/ 44=11.2	173/14=12.4	169/14=12.1	118/ 9=13.1 @I	176/14=12.6
		2	540/ 44=12.3	205/16=12.8	182/16=11.4	208/16=13.0	145/12=12.1
		3	580/ 48=12.1	175/15=11.7	246/19=13.0	195/16=12.2	176/16=11.0
		4	561/ 48=11.7	198/15=13.2 *@I	199/16=12.4	229/18=12.7	195/17=11.5@D
		5	579/ 48=12.1	185/15=12.3	181/15=12.1	203/16=12.7	207/17=12.2
6 1 6 11 6 11 6@11		6	610/ 50=12.2	215/17=12.7	214/18=11.9	228/18=12.7	211/15=14.1@I **@@I
6 1 6 11		7	545/ 49=11.1	249/19=13.1 **@@I	215/16=13.4 **@@I	214/17=12.6 *@I	111/10=11.1*@D

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

@ AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, @, *, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, @, *, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

!, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE III
COMPOUND 8 STUDY SUBACUTE

AVERAGE CORPORA LUTEA PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
EEII	EEII	1	523/ 44=11.9	201/14=14.4 * $\omega\omega$ I	184/14=13.1 * ω I	133/ 9=14.8 ** $\omega\omega$ I	204/14=14.6 ** $\omega\omega$ I
		2	566/ 44=12.9	221/16=13.8	215/16=13.4	226/16=14.1 ω I	159/12=13.3
I	I	3	612/ 48=12.8	205/15=13.7	261/19=13.7	212/16=13.3	194/16=12.1 ω D
I	I	4	594/ 48=12.4	198/15=13.2	206/16=12.9	232/18=12.9	199/17=11.7
		5	605/ 48=12.6	196/15=13.1	185/15=12.3	208/16=13.0	211/17=12.4
I	E II	6	641/ 50=12.8	220/17=12.9	217/18=12.1	233/18=12.9	211/15=14.1 * ω I
EEII	EEII	7	583/ 49=11.9	258/19=13.6 * $\omega\omega$ I	218/16=13.6 ** $\omega\omega$ I	225/17=13.2 ** $\omega\omega$ I	112/10=11.2** $\omega\omega$ D

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST

I AND ω = ONE-TAILED TEST

ONE I, E, ω , * = SIGNIFICANT AT P LESS THAN 0.05

TWO I, E, ω , * = SIGNIFICANT AT P LESS THAN 0.01

*. ω SIGNIFICANTLY DIFFERENT FROM CONTROL

E, I SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE IV
COMPOUND 8 STUDY SUBACUTE

AVERAGE PREIMPLANTATION LOSSES PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
6 11 8 1		1	31/ 44= 0.7	28/14= 2.0 * ϕ I	15/14= 1.1	15/ 9= 1.7	28/14= 2.0 ϕ I
		2	26/ 44= 0.6	16/16= 1.0	33/16= 2.1 ϕ I ** ϕ ϕ I	18/16= 1.1	14/12= 1.2
		3	32/ 48= 0.7	30/15= 2.0 * ϕ ϕ I	15/19= 0.8 ϕ D	17/16= 1.1	18/16= 1.1
6 1		4	33/ 48= 0.7	0/15= 0.0 ** ϕ ϕ D	7/16= 0.4	3/18= 0.2 * ϕ ϕ D	4/17= 0.2 * ϕ D
		5	26/ 48= 0.5	11/15= 0.7	4/15= 0.3	5/16= 0.3	4/17= 0.2
6 11 8 1		6	31/ 50= 0.6	5/17= 0.3	3/18= 0.2 * ϕ D	5/18= 0.3 ϕ D	0/15= 0.0 ** ϕ ϕ D
6 1		7	38/ 49= 0.8	9/19= 0.5	3/16= 0.2 ** ϕ ϕ D	11/17= 0.7	1/10= 0.1 ** ϕ ϕ D

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

6 AND * = TWO-TAILED TEST
1 AND ϕ = ONE-TAILED TEST

ONE 1,6, ϕ ,* = SIGNIFICANT AT P LESS THAN 0.05
TWO 1,6, ϕ ,* = SIGNIFICANT AT P LESS THAN 0.01

*. ϕ SIGNIFICANTLY DIFFERENT FROM CONTROL
6,1 SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE V
COMPOUND 8 STUDY SUBACUTE

AVERAGE RESORPTIONS (DEAD IMPLANTS) PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
8 1	8 1	1	12/ 44=0.28	11/14=0.79 @I	11/14=0.79 @I	3/ 9=0.34	15/14=1.08 @I
		2	21/ 44=0.48	4/16=0.25	15/16=0.94* @I	7/16=0.44	11/12=0.92
		3	31/ 48=0.65	7/15=0.47	14/19=0.74	17/16=1.07** **@I	12/16=0.75
		4	20/ 48=0.42	4/15=0.27	9/16=0.57	9/18=0.50	7/17=0.42
6 11	8 1	5	34/ 48=0.71	9/15=0.60	14/15=0.94	5/16=0.32 @D	4/17=0.24 * @D
8 11	8 1	6	25/ 50=0.50	11/17=0.65	11/18=0.62	10/18=0.56	9/15=0.60
		7	36/ 49=0.74	5/19=0.27	16/16=1.00* @I	5/17=0.30	2/10=0.20 @D

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

8 AND * = TWO-TAILED TEST

1 AND @ = ONE-TAILED TEST

ONE !, @, *, = SIGNIFICANT AT P LESS THAN 0.05

TWO !, @, *, = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

6, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VI
COMPOUND 8 STUDY SUBACUTE

PROPORTION OF FEMALES WITH ONE OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
		1	12/ 44=0.28	7/14=0.50	7/14=0.50	3/ 9=0.34	7/14=0.50
		2	16/ 44=0.37	4/16=0.25	10/16=0.63*	6/16=0.38	5/12=0.42
		3	20/ 48=0.42	6/15=0.40	9/19=0.48	14/16=0.88**	8/16=0.50
		4	13/ 48=0.28	3/15=0.20	6/16=0.38	8/18=0.45	5/17=0.30
!	!	5	23/ 48=0.48	8/15=0.54	10/15=0.67	4/16=0.25	4/17=0.24
		6	19/ 50=0.38	5/17=0.30	8/18=0.45	6/18=0.34	7/15=0.47
		7	15/ 49=0.31	5/19=0.27	10/16=0.63*	4/17=0.24	2/10=0.20

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VII
COMPOUND 8 STUDY SUBACUTE

PROPORTION OF FEMALES WITH TWO OR MORE DEAD IMPLANTATIONS

LOG ARITH DOSE DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
	1	0/ 44=0.0	3/14=0.22 **	3/14=0.22 **	0/ 9=0.0	3/14=0.22 **
	2	3/ 44=0.07	0/16=0.0	3/16=0.19	1/16=0.07	2/12=0.17
	3	7/ 48=0.15	1/15=0.07	3/19=0.16	3/16=0.19	2/16=0.13
	4	6/ 48=0.13	1/15=0.07	3/16=0.19	1/18=0.06	2/17=0.12
!	5	9/ 48=0.19	1/15=0.07	3/15=0.20	1/16=0.07	0/17=0.0
!	6	4/ 50=0.08	3/17=0.18	2/18=0.12	2/18=0.12	2/15=0.14
	7	10/ 49=0.21	0/19=0.0 *	4/16=0.25*	1/17=0.06	0/10=0.0

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VIII
COMPOUND 8 STUDY SUBACUTE

WEEK	DEAD IMPLANTS / TOTAL IMPLANTS				
	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 3.000 MG/KG	DOSE LEVEL 30.000 MG/KG	DOSE LEVEL 300.000 MG/KG
1	12/ 492=0.03	11/173=0.07	11/169=0.07	3/118=0.03	15/176=0.09
2	21/ 540=0.04	4/205=0.02	15/182=0.09	7/208=0.04	11/145=0.08
3	31/ 580=0.06	7/175=0.04	14/246=0.06	17/195=0.09	12/176=0.07
4	20/ 561=0.04	4/198=0.03	9/199=0.05	9/229=0.04	7/195=0.04
5	34/ 579=0.06	9/185=0.05	14/181=0.08	5/203=0.03	4/207=0.02
6	25/ 610=0.05	11/215=0.06	11/214=0.06	10/228=0.05	9/211=0.05
7	36/ 545=0.07	5/249=0.03	16/215=0.08	5/214=0.03	2/111=0.02

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE HISTORICAL CONTROL GROUP

* = TWO-TAILED TEST

@ = ONE-TAILED TEST

ONE *,@ = SIGNIFICANT AT P LESS THAN 0.05

TWO *,@ = SIGNIFICANT AT P LESS THAN 0.01

*,@ SIGNIFICANTLY DIFFERENT FROM CONTROL

APPENDICES

II. MATERIALS AND METHODS

A. Animal Husbandry

1. Animals (Rats and Mice)

Ten to twelve week old rats (280 to 350 g) and male mice (25 to 30 g) were fed a commercial 4% fat diet and water ad libitum until they were put on experiment. Flow Laboratories random-bred, closed colony, Sprague-Dawley CD strain rats were used in the cytogenetic studies. Flow Laboratories ICR male mice were employed in the Host-Mediated Assay.

2. Preparation of Diet

A commercial 4% fat diet was fed to all animals. Periodic tests to verify the absence of coliforms, Salmonella and Pseudomonas sp. were performed.

3. Husbandry

Animals were held in quarantine for 4-11 days. Mice were housed five to a cage and rats one to five to a cage. Animals were identified by ear punch. Sanitary cages and bedding were used, and changed two times per week, at which time water containers were cleaned, sanitized and filled. Once a week, cages were repositioned on racks; racks were repositioned within rooms monthly. Personnel handling animals or working within animal facilities wore head coverings and face masks, as well as suitable garments. Individuals with respiratory or other overt infections were excluded from the animal facilities.

B. Dosage Determination

1. Acute LD₅₀ and LD₅ Determination

Since the compounds proposed for testing are included in



the food additive regulations as "generally recognized as safe" (GRAS), it was expected that a large number of them would be sufficiently non-toxic so that determination of a LD_{50} or a LD_5 would be of no practical value. In fact, this has been our experience with previously tested compounds from this list. In the case of these relatively non-toxic compounds, attempts were made to assure that the amounts to be administered would not affect the animals by means (mechanical, physical, etc.) related to their bulk rather than to their toxicity. In the cases of certain compounds where a LD_{50} or a LD_5 could not be determined, an exceedingly high concentration, 5 g/kg, was employed and accepted as the LD_5 level. In cases where the toxicity was high enough to allow determination of a LD_5 , the following protocol was used.

Thirty rats of the strain chosen for studies described below and of approximately the age and weight specified were assigned at random to six groups. Each group was then given, using the chosen route of administration, one of a series of dosages of the test compound following a logarithmic dosage scheme. The series of dosages were derived from a consideration of whatever toxicity information was available for the particular test compound. The objective in selecting dosages was to choose values which would cause mortalities between 10% and 90%.

When information was inadequate to derive a suitable series of dosages, five rats were used to identify the proper range. Each of these was given one of a widely spaced (differing by 10X) series of doses. This was confidently expected to suffice for derivation of the series of dosages to be used in the LD_{50} determination.

The mortalities observed when the series of dosages were given to the 30 rats were then subjected to a probit analysis and calculation of LD_{50} , LD_5 , slope and confidence limits by the method of Litchfield and Wilcoxon. The highest dose level used was either a finite LD_5 or 5000 mg/kg. The intermediate level used was either 1/10 of the finite LD_5 or 2500 mg/kg. The low level used was either 1/100 of the finite LD_5 or 30 mg/kg.

2. Subacute Studies

Subacute doses were identical to those used in the acute studies. Each subacute study animal was given the acute dosage once a day for each of five consecutive days (24 hours apart).

C. Mutagenicity Testing Protocols

1. Host-Mediated Assay

Flow Laboratories ICR random-bred male mice were used in this study. In the acute and subacute studies ten animals, 25-30 g each, were employed at each dose level. Solvent and positive controls were run at all times. The positive control (dimethyl nitrosamine) was run by the acute system only at a dose of 100 mg/kg for Salmonella. For yeast, ethyl methane sulfonate (EMS) intramuscularly injected at a dose of 350 mg/kg was used. The solvents used and the toxicity data are presented in the Results and Discussion Section of the report.

The indicator organisms used in this study were: (1) two histidine auxotrophs (his G-46, TA-1530) of Salmonella typhimurium, and (2) a diploid strain (D-3) of Saccharomyces cerevisiae. The induction of reverse mutation was determined with the Salmonella; mitotic recombination was determined with yeast. Chemicals were evaluated directly by in vitro bacterial and yeast studies prior to, or concurrent with, the studies in



mice. Only animals on the subacute studies were not fed the evening prior to compound administration. The Salmonella were carried in tryptone yeast extract gel, transferred weekly. They were transferred to tryptone yeast extract broth 48 hours before use: they were transferred a second time from broth to broth 24 hours prior to use, and again 8 hours before use. The mouse inoculum was prepared by transferring 4 ml of the 8-hour broth culture to 50 ml broth bottles which had been prewarmed at 37°C. Exponential log-phase organisms were inoculated intraperitoneally into the mice approximately 2-1/2 hours later when the appropriate density indicating 3.0×10^8 cells/ml was reached. The Saccharomyces was carried in yeast complete agar. The inoculum was prepared by harvesting the organisms from the surface of the plates with sterile saline. The cells were washed three times with sterile saline and suspended in a concentration of 5.0×10^8 cells/ml. Two ml of the suspension was inoculated into each mouse intraperitoneally. Total plate counts on Salmonella were on tryptone yeast extract and for Saccharomyces on yeast complete medium.

a. Acute study

Three dosage levels (usage, intermediate [determined as discussed previously], and LD_5) were administered orally by intubation to ten mice. Positive controls and negative vehicle controls were included in each study. All animals received 2 ml of the indicator organism intraperitoneally. Each ml contained 3.0×10^8 cells for Salmonella and 5.0×10^8 cells for Saccharomyces. Three hours later, each animal was killed and 2 ml of sterile saline was introduced intraperitoneally. As much fluid as possible was then aseptically removed from the peritoneal cavity. Dilution blanks for bacteria containing 4.5 ml of sterile saline were prepared in advance. Tenfold serial



BIONETICS

dilutions were made of each peritoneal exudate (0.5 ml exudate + 4.5 ml saline) yielding a concentration series from 10^0 (undiluted peritoneal exudate) through 10^{-7} . For enumeration of total bacterial counts, the 10^{-6} and 10^{-7} dilutions were plated on tryptone yeast extract agar, 3 plates/sample, 0.2 ml sample/plate. Each sample was spread over the surface of the plate using a bent glass rod immersed in 95% ethanol and flamed just prior to use. In plating for the total mutant counts on minimal agar, the 10^0 dilution was used, 0.2 ml being plated on each of 5 plates. The plating procedure was identical to that followed for the tryptone yeast extract agar plates. All plates were incubated at 37°C, tryptone yeast extract agar plates for 18 hours and minimal agar plates for 40 hours. For yeast mitotic recombination, dilution blanks containing 4.5 ml of sterile saline were prepared in advance. Tenfold serial dilutions were made of each sample yielding a series from 10^0 to 10^{-5} . Samples of 0.1 ml of the 10^{-5} , 10^{-4} , and 10^{-3} dilutions were removed and plated on complete medium (10 plates each). All plates were incubated at 30°C for 40 hours. The 10^{-5} dilutions were used to determine total populations and the 10^{-4} and 10^{-3} plates were examined after an additional 40 hours at 4°C for red sectors indicating a mutation. Bacterial scoring was calculated as follows:

Total mutants on 5 plates x appropriate exponent =
CFU/ml (CFU is Colony Forming Units) of sample plated
CFU/ml x one/dilution factor ($10^0 - 10^{-7}$) = CFU/ml in undiluted exudate. The mutation frequency (MF) calculated for each sample was:

$$MF = \frac{\text{total mutant cells}}{\text{total population}}$$

$$MF_t/MF_c = \frac{MF \text{ of experimental sample}}{MF \text{ of control sample}}$$

$$(MF_t/MF_c = 1.00 \text{ for control sample})$$



Yeast mitotic recombinants (presumptive ade 2, his 8 homozygotes) were seen as red colonies or as red sectors on a normally white yeast colony. The plates (from 10^{-4} and 10^{-3} dilutions) were scanned under the 10X lens of a dissecting scope to enumerate the red colonies and sectors. Population determinations were made from the 10^{-5} dilution plates. A recombinant frequency (RF) was calculated:

$$RF = \frac{\text{total recombinants counted}}{\text{total number colonies screened}}$$

b. Subacute study

Similar groups of animals at each dose level received five oral doses of the test compound 24 hours apart. Within 30 minutes after the last dosing, the animals were inoculated with the test organism and handled in the same fashion as those in the acute study.

c. In vitro study

Cultures of S. typhimurium histidine auxotrophs (G-46 and TA-1530) were plated on appropriate media. The test compound was then added to the plate, either in the form of a microdrop of solution (0.01 to 0.25 ml) applied to a small filter paper disc resting on the agar or a small crystal applied directly to the agar. Tenfold serial dilutions of the culture were employed and plated so as not to miss the optimum cell density for mutant growth. Mutant colonies were observed and scored. Strain D-3 Saccharomyces cells at proper dilutions were shaken with the test compound, diluted, and plated at 50% survival level or above (see HMA Supplementary Materials and Methods). Red sectors were then scored and the frequency calculated after suitable incubation. Negative and positive controls were run concurrently. The positive control was EMS for Salmonella and Saccharomyces. The in vitro Salmonella tests were reported



as (+) or (-) or questionable; the in vitro Saccharomyces tests were reported as sample concentrations, percent survival, and recombinants/ 10^5 survivors. For the Saccharomyces a 50% survival level, e.g., an arbitrary 5.0% w/v test level, was used when no LD₅₀ was determinable.

2. Cytogenetic Studies

a. In vivo study

Ten to twelve week old, male, albino rats obtained from a closed colony (random-bred) were used. A total of 59 animals in the acute study and 18 animals in the subacute study was used, as illustrated in the following protocol.

Number of Animals Used

Acute Study

<u>Treatment</u>	<u>Time Killed After Administration</u>		
	<u>6 Hours</u>	<u>24 Hours</u>	<u>48 Hours</u>
High Level	5	5	5
Intermediate Level	5	5	5
Low Level	5	5	5
Positive Control	0	0	5
Negative Control	3	3	3

Subacute Study

Five doses 24 hours apart; animals killed 6 hours after last dose.

<u>Treatment</u>	<u>Killed After Administration</u>
High Level	5
Intermediate Level	5
Low Level	5
Negative Control	3

All animals were dosed by gastric intubation.

Four hours after the last compound administration, and two hours prior to killing, each animal was given 4 mg/kg of colcemid intra-



BIONETICS

peritoneally in order to arrest the bone marrow cells in C-mitosis. Animals were killed by using CO₂, and the adhering muscle and epiphysis of one femur were removed. The marrow "plug" was removed with a tuberculin syringe and an 18 gauge needle, aspirated into 5 ml of Hanks' balanced salt solution (BSS) in a test tube and capped. The specimens were centrifuged at 1,500 RPM in a table-top centrifuge for 5 minutes, decanted, and 2 ml of hypotonic 0.5% KCl solution was added with gentle agitation to resuspended the cells. The specimens were then placed in a 37°C water bath for 20 minutes in order to swell the cells. Following centrifugation for 5 minutes at 1,500 RPM, the supernatant was decanted and 2 ml of fixative (3:1 absolute methanol:glacial acetic acid) was added. The cells were resuspended in the fixative with gentle agitation, capped, and placed at 4°C for 30 minutes. The specimens were again centrifuged, decanted, 2 ml of prepared fixative was added, and the cells were resuspended and placed at 4°C overnight.

The following day the specimens were again centrifuged, decanted and 0.3 - 0.6 ml of freshly prepared fixative was added to obtain a suitable density. The cells were resuspended and 2 - 3 drops of the suspension were allowed to drop onto a clean, dry slide held at 15° from the horizontal. As the suspension flowed to the edge of the slide, it was ignited by an alcohol burner and allowed to flame. Following ignition, the slides were allowed to dry at room temperature overnight. Duplicate slides were prepared. The slides were stained using a 5% Giemsa solution (Giemsa buffer pH 7.2) for 20 minutes, rinsed in acetone, 1:1 acetone:xylene, and placed in fresh xylene for 30 minutes. The slides were then mounted using Permount (Fisher Scientific) and 24 x 50 mm coverglasses. The coverglasses were selected to be 0.17 mm \pm 0.005 mm in thickness by use of a coverglass micrometer. The preparations

were examined using Leitz Ortholux I & II microscopes with brightfield optics and xenon light sources. These specimens were scanned with 10X and 24X objectives and suitable metaphase spreads that were countable were then examined critically using 40X, 63X or 100X oil immersion flatfield apochromatic objectives. Oculars were either 12X or 16X widefield periplanatics and the tube magnification either 1X or 1.25X. The filters used were either a didymium (BG20) or a Schott IL570 $m\mu$ interference filter.

The chromosomes of each cell were counted and only diploid cells were analyzed. They were scored for chromatid gaps and breaks, chromosome gaps and breaks, reunions, cells with greater than ten aberrations, polyploidy, pulverization, and any other chromosomal aberrations which were observed. They were recorded on the currently used forms and expressed as percentages on the summary sheets. Fifty metaphase spreads were scored per animal. Mitotic indices were obtained by counting at least 500 cells and the ratio of the number of cells in mitosis/the number of cells observed was expressed as the mitotic index.

Positive controls in the acute study consisted of animals which had been given the known mutagen Triethylene Melamine (TEM) administered intraperitoneally at a level of 0.30 mg/kg. Negative controls on the acute and subacute studies consisted of the vehicle in which the compound was administered. The dosage levels, solvents and toxicity data are included in the Results and Discussion Section of the report.

b. In vitro study

Human embryonic lung cultures (WI-38) which were negative for adventitious agents (viruses, mycoplasma) which may interfere

were used. These cells were employed at passage level 19. The cells had been transferred using 0.025% trypsin and planted in 32 oz. prescription bottles containing 40 ml of tissue culture medium. When growth was approximately 95% confluent the cells were removed from the glass using trypsin, centrifuged, and frozen in tissue culture medium containing dimethyl sulfoxide (DMSO). Cells were frozen in vials in the vapor phase of liquid nitrogen at a concentration of 2×10^6 cells/ml. When needed, the vials were removed from liquid nitrogen, quick-thawed in a 37°C water bath, washed free of DMSO, suspended in tissue culture medium (minimal essential medium [MEM] plus 1% glutamine, 200 units/ml of penicillin and 200 µg/ml of streptomycin and 15% fetal calf serum) and planted in milk dilution bottles at a concentration of 5×10^5 cells/ml. The test compound was added at three dose levels using three bottles for each level, 24 hours after planting. The dose levels required a preliminary determination of a tissue culture toxicity. This was accomplished by adding logarithmic doses of the compound in saline to a series of tubes containing 5×10^5 cells/ml which were almost confluent. The cells were examined at 24, 48, and 72 hours. Any cytopathic effect (CPE) or inhibition of mitoses was scored as toxicity. Five more closely spaced dose levels were employed within the two logarithmic dosages, the higher of which showed toxicity and the lower no effect. The solvents used and the range finding data are presented in the toxicity data report under Results and Discussion. The dose level below the lowest toxic level was employed as the high level. Logarithmic dose levels were employed for the medium and low levels.

Cells were incubated at 37°C and examined twice daily to determine when an adequate number of mitoses were present. Cells were harvested by shaking when sufficient mitoses were observed, usually 24 - 48

hours after planting, centrifuged, and fixed in absolute methanol:glacial acetic acid (3:1) for 30 minutes.

The specimens were centrifuged, decanted, and suspended in acetic acid-orcein stain (2.0%) and a drop of suspension placed on a clean dry slide. Selected coverglasses 0.17 mm in thickness were placed on the suspension and the excess stain gently expressed from the slide. The coverglasses were sealed with clear nail polish and examined immediately.

The microscopes, objectives, oculars, filters and light sources were enumerated under the metaphase description. Positive controls used were TEM (at a concentration of 0.1 mcg/ml dissolved in saline) and negative controls which consisted of the vehicle in which the test compound was dissolved, which was 0.85% saline. Data were reported on forms currently used and expressed as percentages on the anaphase summary sheets.

3. Dominant Lethal Assay

In this test, male and female random bred rats from a closed colony were employed. These animals were 10-12 weeks old at the time of use. Ten male rats were assigned to each of 5 groups; 3 dose levels selected as described above, a positive control (triethylene melamine) (TEM) and a negative control (solvent only). The positive control was administered intraperitoneally. Administration of the test compound was orally by intubation in both the acute study (1 dose) and in the subacute study (1 dose per day for 5 days). Following treatment, the males were sequentially mated to 2 females per week for 8 weeks (7 weeks in the subacute study). Two virgin female rats were housed with a male for 5 days (Monday through Friday). These two females were removed and housed in a cage until killed. The male was rested on Saturday and Sunday and two new females introduced to the cage on

Monday. It has been our experience that conception has taken place in more than 90% of the females by Friday and that the two day rest is beneficial to the male as regards subsequent weekly matings. Females were killed using CO₂ at 14 days after separating from the male, and at necropsy the uterus was examined for deciduomata (early deaths), late fetal deaths and total implantations.

Sufficient animals were provided in our experimental design to accommodate for any reduction in the number of conceptions. Each male was mated with two females per week, and this provided for an adequate number of implantations per group per week (200 minimum) for negative controls, even if there was a fourfold reduction in fertility of implantations. Results were analyzed according to the statistical procedures described in Supplementary Materials and Methods. Corpora lutea, early fetal deaths, late fetal deaths and total implantations per uterine horn were recorded on the raw data sheets, which are submitted separately.

D. Supplementary Materials and Methods

1. Host-Mediated Assay In Vitro and Formulae

a. Bacterial in vitro plate tests

This method has been published by Ames: The Detection of Chemical Mutagens with Enteric Bacteria, in Chemical Mutagens; Principles and Methods for Their Detection, Vol. 1, Chapter 9, pp. 267-282, A. Hollaender, Editor, Plenum Press, New York (1971).

b. In vitro for mitotic recombination

(1) Strain D-3 was grown to stationary phase on complete medium agar plates at 30°C (3-4 days). Cells were rinsed from the plates and washed twice in saline and cell concentration determined spectro-

photometrically. (A standard curve previously determined for colony forming units versus % transmittance at 545 mu was easily used.)

(2) Cells from the concentration suspension were diluted appropriately into 0.067 M Phosphate buffer pH 7.2 to provide 5×10^7 cells/ml in a total of 25 ml.

(3) The test chemical was first tested for 4 hours at 30°C, with shaking, at concentrations which permitted determination of the 50% survival level. Then, if not included in the first experiment, the compound was tested again only at the 50% survival level. If 50% survival level could not be determined, the arbitrary test level of 5% w/v was used.

(4) Following treatment, cells were diluted and plated on complete agar medium for determination of total population and red sectors. Total surviving population was conveniently measured on plates of 10^{-4} and 10^{-5} dilutions using 0.2 ml per plate (5 plates), and sectors determined on plates of 10^{-3} and 10^{-4} dilutions using 0.2 ml per plate (5 plates). Plates were incubated for 2 days at 30°C followed by a holding period of 2 days at 4°C to promote color development with limited enlargement of the colonies. Red sectors were scored by systematically scanning the plates with a dissecting microscope at 10X magnification.

(5) The frequency of red sectors can then be calculated and may be expressed conveniently as sectors per 10^5 survivors for comparison with untreated controls.

(6) Ethyl Methane Sulfonate (EMS) was employed as the positive control in both in vitro systems.

c. Minimal medium (bacteria):

Spizizen's Minimal Medium:



BIONETICS

4X Salt Solution:

$(\text{NH}_4) \text{SO}_4$	8.0 gm
K_2HPO_4	56.0 gm
KH_2PO_4	24.0 gm
Na Citrate	4.0 gm
Mg SO_4	0.8 gm
Biotin	0.004 gm
H_2O	qs to 1 liter
	Sterilize by autoclaving (121°C/15 min.)

Medium:

4X Salt Solution	:250 ml	
5.0% Glucose (sterile)	:100 ml	(If histidine is added at concentration of 30 mg/liter, this becomes a complete bacterial medium.)
1.5% Bacto-agar (sterile)	:650 ml	

d. Complete medium (bacteria):

Bacto-Tryptone	1.0 gm
Yeast-Extract	0.5 gm
Bacto-Agar	2.0 gm
Distilled H_2O	100.0 ml

Sterilize by autoclaving (121°C for 15 minutes).

e. Complete medium (yeast):

KH_2PO_4	1.5 gm
MgSO_4	0.5 gm
$(\text{NH}_4)_2\text{SO}_4$	4.5 gm



BIONETICS

Peptone	3.5 gm
Yeast-Extract	5.0 gm
Glucose	20.0 gm
Agar	20.0 gm
Distilled H ₂ O	1000.0 ml

Sterilize by autoclaving (121°C for 15 minutes).

2. Cytogenetics In Vitro Preparation of Anaphase Chromosomes (from Nichols, 1970)

"Anaphase preparations may be made by several methods. One convenient approach is to grow cells directly on coverslips in petri dishes. With human fibroblasts 400,000 cells added to a 22 x 44 mm coverslip in a 50 mm petri dish grown in a 5% CO₂ atmosphere in air has proved very satisfactory. When adequate numbers of mitoses are visualized directly utilizing an inverted microscope (usually 48 to 92 hours after planting) the coverslip is transferred to absolute ethanol for 15 minutes for fixation. They are then stained with any one of a number of suitable stains (Fuelgen, May-Grunwald-Giemse, orcein) and attached to a slide with mounting media for evaluation. Anaphase preparations may also be prepared on cells grown in suspension or cells from a monolayer that have been put into suspension. In this instance the cells are centrifuged and fixed with the squash fixative. They are then suspended in the stain and a drop of the suspension put on the slide and covered with a coverslip. However, in this case, only the excess stain is gently expressed from under the coverslip and no squashing is carried out. In anaphase preparations no pretreatment with colchicine or hypotonic expansion is used and no technique for spreading the cells is used, so that the spindle and normal relationships of the chromosomes are not disturbed."



BIONETICS

3. Statistical Analyses of Dominant Lethal Studies

The following statistical analyses were employed as a means of analyzing the results of the dominant lethal studies.

a. The fertility index

The number of pregnant females/number of mated females with the chi-square was used to compare each treatment to the control. Armitage's trend was used for linear proportions to test whether the fertility index was linearly related to arithmetic or log dose.

b. Total number of implantations

The t-test was used to determine significant differences between average number of implantations per pregnant female for each treatment compared to the control. Regression techniques were used to determine whether the average number of implantations per female was related to the arithmetic or log dose.

c. Total number of corpora lutea

The t-test was used to determine significant differences between average number of corpora lutea per pregnant female for each treatment compared to the control.

d. Preimplantation losses

Preimplantation losses were computed for each female by subtracting the number of implantations from the number of corpora lutea. Freeman-Tukey transformation was used on the preimplantation losses for each female and then the t-test was used to compare each treatment to control. Regression technique was used to determine whether the average number of preimplantation losses per female was related to the arithmetic or log dose.

e. Dead implants

Dead implants were treated the same as pre-implantation losses.

f. One or more dead implants

The proportion of females with one or more dead implants was computed, each treatment compared to control by chi-square test and Armitage's trend used for linear proportions to see if proportions were linearly related to either arithmetic or log dose. Also, probit regression analysis was used to determine whether the probit of the proportions was related to log dose.

g. Two or more dead implants

The proportion of females with two or more dead implants computed was treated same as above (f).

h. Dead implants per total implants

Dead implants per total implants were computed for each female and used Freeman-Tukey arc-sine transformation on data for each female; then used t-test to compare each treatment to control.

Historical control data was compiled on a continuous basis as studies were completed. In addition to comparing each treatment to control, as outlined above, each treatment was compared to a historical control.

In order to take variation between males into account, a nested model was used. An analysis of across weeks is also provided.

In addition to these tests, the distribution forms of the various parameters were tested in order to evaluate the appropriateness of some of the tests being used. Certain correlations between parameters may exist and were examined as one step to determine the appropriateness of models. If necessary, alternate test methods were implemented.



The results are presented in tabular form with the addition of historical control information. In addition to these tables, a written report of all findings is provided. As information became available from the on-going investigation of these data, it was reported and suggestions included for changes to the methods of analysis. The statistical reports give the level of significance using both a one-tailed and two-tailed test. Finally, a summary sheet for each study is provided.



MODEL

$$y_{ijk} = \mu + \alpha_i + c_{ij} + e_{ijk}$$

$i = 1, 2$ Group

$j = 1, 2, \dots, 10$ Males within each group

$k = 1, 2$ Females within Males within Groups

ASSUMPTIONS:

$$\alpha_1 + \alpha_2 = 0, \quad c_{ij} \sim \text{nid}(0, \sigma_c^2)$$

$$e_{ijk} \sim \text{nid}(0, \sigma^2)$$

Males are randomly drawn from infinite population

S.U.	d.f.	S.S.	MS	E(MS)	F
TOTAL	39	$\sum \sum \sum (y_{ijk} - \bar{y}_{...})^2$			
GROUPS	1	$20 \sum (\bar{y}_{i..} - \bar{y}_{...})^2$	S_1^2	$\sigma^2 + 20\sigma_c^2 + 20\sigma_e^2$	
MALES WITHIN GROUPS	18	$2 \sum \sum (\bar{y}_{ij.} - \bar{y}_{i..})^2$	S_2^2	$\sigma^2 + 2\sigma_c^2$	
REMAINDER	20	$\sum \sum \sum (y_{ijk} - \bar{y}_{ij.})^2$	S_3^2	σ^2	

E. References

1. Host-Mediated Assay

- a. Gabridge, M.G., Denunzio, A. and Legator, M.S.:
Nature, 221:68, 1969.
- b. Gabridge, M.G., Denunzio, A. and Legator, M.S.:
Science, 163:689, 1969.
- c. Gabridge, M.G. and Legator, M.S.: Proc. Soc.
Exptl. Biol. Med., 130:831, 1969.
- d. Gabridge, M.G., Oswald, E.J. and Legator, M.S.:
Mut. Res., 7:117, 1969.
- e. Legator, M.S. and Malling, H.V.: In, Environmental
Chemical Mutagens, A. Hollaender (Ed.), Plenum
Publishing Corp., New York, in press.

2. Cytogenetics

- a. Nichols, V.W.: Personal communication.
- b. Legator, M.S.: In, Laboratory Diagnosis of Diseases
Caused by Toxic Agents, F. W. Sunderman and F. W.
Sunderman (Ed.), Warren H. Green, Inc., St. Louis,
pp. 17-22, 1970.
- c. Hsu, T.C. and Patton, J.L.: Technical Addendum in,
Comparative Mammalian Cytogenetics, K. Benirschke
(Ed.), Springer-Verlag, New York, pp. 454-460, 1969.
- d. Legator, M.S. et al.: Cytogenetic studies in rats
of cyclohexylamine, a metabolite of cyclamate.
Science, 165:1139, 1969.



3. Dominant Lethal

- a. Bateman, A.J.: Genet. Res. Comb., 1:381, 1960.
- b. Bateman, A.J.: Nature, 210:205, 1966.
- c. Ehling, U.H., Cumming, R.B. and Malling, H.V.:
Mut. Res., 5:417, 1968.
- d. Epstein, S.S. and Shafner, H.: Nature, 219:
385, 1968.



F. Abbreviations

1. mu = micron
2. mcg = ug = microgram
3. g = gram
4. kg = kilogram
5. ml = milliliter
6. rpm = revolutions per minute
7. °C = degrees centigrade
8. pH = power of the hydrogen ion concentration to the base 10
9. M = molar solution
10. conc. = concentration
11. MTD = maximum tolerated dosage = High = LD₅ if determined
or else exceedingly high dose, such as 5 g/kg
12. INT = intermediate = medium level
13. USE = usage level if known = low level
14. BSS = balanced salt solution
15. C-metaphase = cells arrested in metaphase, using colchine
or colcemid
16. LD₅₀ = that dosage which produced 50% mortality in the
group of animals treated
17. LD₅ = that dosage which produced 5% mortality in the group
of animals treated
18. NC = negative control
19. PC = positive control
20. AU = acute usage level (low level)
21. AI = acute intermediate level (medium level)
22. AMTD = acute maximum tolerated dose level (LD₅ level,
high level)



- 23. SAU = subacute usage level (low level)
- 24. SAI = subacute intermediate level (medium level)
- 25. SA LD₅ = subacute LD₅ level (MTD level, high level)
- 26. CO₂ = carbon dioxide
- 27. DMN = Dimethyl nitrosamine
- 28. EMS = Ethyl methane sulfonate
- 29. TEM = Triethylene melamine
- 30. DMSO = Dimethyl sulfoxide
- 31. MEM = minimal essential medium (Eagle's)
- 32. CPE = cytopathic effect
- 33. his = histidine marker
- 34. D-3 = mitotic recombinant strain of Saccharomyces
- 35. mf = mean mutant frequency
- 36. MFt/MFc = mean mutant frequency of the test compound group compared to mean mutant frequency of the negative control group
- 37. CFU = colony forming units
- 38. WI-38 = code name for a strain of human embryonic lung tissue culture cells
- 39. Rec x 10⁵ = mitotic recombinants x 10⁵
- 40. Mean B/A = mean frequency
- 41. tot. scr. = total scored
- 42. tot. = total
- 43. χ^2 = a test of variation in the data from the computed regression line - tested in these studies at the 5% level
- 44. Aber. = aberrations
- 45. Frag. = fragment
- 46. HMA = host-mediated assay

