



BIONETICS

MUTAGENIC EVALUATION OF
COMPOUND FDA 73-3

SODIUM HEXAMETAPHOSPHATE

Mutagenic Evaluation of Compound FDA 73-3-Sodium Hexametaphosphate

7315 Wisconsin Avenue
Bethesda, Maryland
20014

LBI PROJECT # 2468

MUTAGENIC EVALUATION OF
COMPOUND FDA 73-3
SODIUM HEXAMETAPHOSPHATE

SUBMITTED TO
FOOD & DRUG ADMINISTRATION
DEPARTMENT OF HEALTH, EDUCATION AND WELFARE
ROCKVILLE, MARYLAND

SUBMITTED BY
LITTON BIONETICS, INC.
5516 NICHOLSON LANE
KENSINGTON, MARYLAND

JANUARY 10, 1975



BIONETICS

TABLE OF CONTENTS

	<u>Page No.</u>
Evaluation Summary	1
I. Objective	2
II. Materials	2
III. Methods	3
IV. Solubility Properties	7
V. Toxicity and Dosage Determinations	8
VI. Non-activation Plate Tests	9
VII. Activation Plate Tests	11
VIII. Non-activation Suspension Tests/ <u>Salmonella</u>	17
IX. Activation Suspension Tests/ <u>Salmonella</u>	19
X. Non-activation Suspension Tests/ <u>Saccharomyces</u>	27
XI. Activation Suspension Tests/ <u>Saccharomyces</u>	29
XII. Summary of Test Results	35
XIII. Interpretation and Conclusions	41
Appendix - Summary of Tests Evaluating DMSO for Genetic Activity in <u>Salmonella</u> and <u>Saccharomyces</u>	



BIONETICS

EVALUATION SUMMARY

Compound FDA 73-3, Sodium Hexametaphosphate, did not exhibit genetic activity in any of the in vitro tests employed in this evaluation.



DATE: January 10, 1975

SPONSOR: Food and Drug Administration

SUBJECT: Mutagenic Evaluation of Compound FDA 73-3

I. OBJECTIVE

The objective of this study was to assess the genetic activity of the test material in microbial assays with and without the addition of mammalian metabolic enzyme preparations.

II. MATERIALS

A. Test Material

Sodium Hexametaphosphate
Stauffer #0138

B. Tissue Homogenates and Supernatants

The tissue homogenates and 9,000 x g supernatants were prepared from liver, lung and testes of the following mammalian species: Mouse - ICR random bred adult males; rat - Sprague-Dawley adult males; and primate - Macaca mulatta adult males.

C. Indicator Organisms

The indicator organisms used for all tests are described below:

- Saccharomyces cerevisiae, strain D4: $\frac{\alpha}{a}$, $\frac{ade\ 2-2}{ade\ 2-1}$, $\frac{try\ 5-12}{try\ 5-27}$
- Salmonella typhimurium, strains:
TA-1535; hisG, uvrB, rfa (missense mutation)
TA-1537; hisC, uvrB, rfa (- frameshift mutation)
TA-1538; hisD, uvrB, rfa (+ frameshift mutation)

D. Reaction Mixture

The following reaction mixture was employed in the activation tests:

	<u>Component</u>	<u>Final Concentration/ml</u>
1.	TPN (sodium salt)	6 μ M
2.	Isocitric acid	49 μ M
3.	Tris buffer, pH 7.4	28 μ M
4.	MgCl ₂	1.7 μ M
5.	Isocitric dehydrogenase	1.0 Unit
6.	Tissue homogenate or cell fraction	72 mg

Components 1-4 were combined and frozen as a "core" reaction mixture to which the other components were added just prior to use.

E. Positive Control Compounds

Table 1 lists chemicals for positive controls in the direct and activation assays.

TABLE 1
POSITIVE CONTROLS USED IN DIRECT AND ACTIVATION ASSAYS

<u>ASSAY</u>	<u>CHEMICAL^a</u>	<u>SOLVENT</u>	<u>PROBABLE MUTAGENIC SPECIFICITY^b</u>
Non-activation	Ethylmethane sulfonate	Water or saline	BPS
	2-Nitrosofluorene	Dimethylsulfoxide ^c	FS
	Quinacrine or Quinacrine mustard	Water or saline	FS
Activation	Dimethylnitrosamine	Water or saline	BPS
	2-Acetylaminofluorene	Dimethylsulfoxide ^c	FS

^a Concentrations given in the Results Section.

^b BPS = base-pair substitution; FS = frameshift.

^c Previously shown to be non-mutagenic, see Appendix.

III. METHODS

A. Toxicity

The solubility, toxicity and doses for all chemicals were determined prior to screening.

Each chemical was tested for survival against strains TA-1537 and D4 over a range of doses to determine the 50% survival dose. Bacteria were tested in phosphate buffer, pH 7.4, for one hour at 37°C on a shaker. Yeasts were tested in phosphate buffer, pH 7.4, for four hours at 30°C on a shaker. The 50% survival dose was determined from the survival curve and the 1/4 and 1/2 50% doses calculated.



If no toxicity was obtained for a chemical with a given strain, then a maximum dose of 5% (w/v) was used against the strain.

Unless otherwise specified, the doses calculated for the tests in buffer were applied to the activation tests. The solubility of the test chemical under treatment conditions is stated in the Results Section.

B. Plate Tests

Only three bacteria strains were tested in qualitative plate tests. In the non-activation procedure, approximately 10^9 cells of a log phase culture of the bacterial indicator strains were spread over the surface of a minimal plate, and a measured amount of the test chemical was placed in the center of the test plate. In activation tests, the test chemical was added to the cells, and an aliquot of the mixture was spread on the surface of the test plate. The reaction mixture (0.1 ml) plus tissue extract was then spotted on the surface of the plate. Positive and solvent controls were included. All plates were incubated at 37°C for four days and then scored. Each compound (Test, Positive Control and Solvent Control) was done in duplicate. The results were scored as + or -. Concentrations of the positive control compounds are listed in the Results Section.

C. Suspension Tests

1. Non-activation

Log-phase bacteria and stationary-phase yeast cultures of the indicator organisms were grown in complete broth, washed and resuspended in 0.9% saline to densities of 1×10^9 cells/ml and 5×10^7 cells/ml, respectively. This constituted the working stock for tests of a group of test chemicals and their respective controls. Tests were conducted in 30 ml plastic tissue culture flasks. Cells plus appropriate volume(s) of the test chemical were added to the flasks to give a final volume of 2 ml. Solvent replaced the test chemical in the negative controls. Treatment was at 30°C for four hours for yeast tests and at 37°C for one hour for bacterial tests. All flasks were shaken during treatment. Following treatment, the flasks were set in ice. Aliquots of cells were removed, diluted in sterile saline (4°C) and plated on the appropriate complete media. Undiluted samples from flasks containing the bacteria were plated on minimal selective medium. Samples from a 10^{-1} dilution of treated cells were plated on the selected media for enumeration of gene conversion with strain D4. Bacterial plates were scored after incubation for 48 hours at 37°C. The yeast plates were incubated at 30°C for 3-5 days before scoring.

2. Activation

Bacteria and yeast cells were grown and prepared as described in the non-activation tests except that the cell densities were increased approximately five-fold for working stock suspensions. Measured amounts of the test and



control chemicals plus 0.25 ml of the stock cell suspension were added to a 30 ml plastic tissue homogenate. All flasks (bacteria and yeast) were incubated at 37°C with shaking. The treatment times as well as the dilutions, plating procedures and scoring of the plates were the same as described for non-activation tests.

D. Preparation of Tissue Homogenates and 9,000 x g Cell Fractions

1. Mice

Male mice (sufficient to provide the necessary quantities of organs for testes, lung and liver homogenates) were killed by cranial blow, decapitated and bled. The three organs were immediately dissected from the animal using aseptic techniques and placed in ice-cold 0.25 M sucrose buffered with Tris at pH of 7.4. Upon collection of the desired quantity of organs, they were washed twice with fresh buffered sucrose and completely homogenized with a motor-driven homogenizing unit at 4°C. The whole organ homogenate obtained from this step was divided into two samples. One sample was frozen at -80°C and the other was centrifuged for 20 minutes at 9,000 x g in a refrigerated centrifuge. The supernatant from the centrifuged sample was retained and frozen at -80°C. These two frozen samples were used for the activation studies.

2. Rats

The same procedures as described for mice were used for this mammal.

3. Primates

The liver, lungs and testes were aseptically removed from freshly killed adult male rhesus (*M. mulatta*) monkeys. Each organ was cut into a number of pieces each sufficient for one week's studies. The tissues were labeled and frozen at -80°C until needed. Tissue homogenates and 9,000 x g supernatants were prepared as described for mice.

E. Data Recording and Reporting

Following the specified incubation periods all population plates were scored by an automatic colony counter and the results from each plate of a set were recorded, in ink, in bound data books. Information necessary for identification of the specific experiment as well as the presence of any contaminant micro-organisms was recorded with each set of plate counts. All minimal or other types of selective media plates were hand scored and the results recorded along with the respective population data. For bacteria strains the number of colonies recorded from either the population or selective plates represents that number in 1 ml of test suspension plated. The numbers recorded for the yeast strain D4 represent the number in 0.5 ml of test suspension plated.



Frequencies were mechanically calculated and double checked. All data presented in the Results Section of this report consists of the actual sum of all raw data plate counts and only the frequencies are calculated figures.



IV. SOLUBILITY PROPERTIES OF THE TEST COMPOUND

1. NAME OR DESCRIPTION OF TEST COMPOUND:

Sodium Hexametaphosphate

2. TEST SOLVENT AND DESCRIPTION OF SOLUBILITY
OF THE TEST CHEMICAL UNDER TREATMENT
CONDITIONS:

This compound was soluble at the treatment concentration employed in this evaluation. All tests were conducted in an aqueous environment.

3. OTHER COMMENTS:

Project No. 2468



BIONETICS

V. TOXICITY AND DOSAGE DETERMINATIONSCOMPOUND FDA 73-3

	D4	TA-1537	
Range of concentrations of the test compound used to determine the 50% survival level	Dose No.	% Concentration	% Concentration
	1	1	0.1
	2	2	0.5
	3	3	1.0
	4	4	2.5
	5	5	5.0
Survival Results	Dose No.	% Survival	% Survival
	Control	100	100
	1	100	41
	2	100	24
	3	100	21
	4	100	16
	5	100	7
Test Date: <u>9-15-74</u>			
Concentrations of the test chemical required for mutagenicity tests	Dose	% Concentration	% Concentration
	Plate Test	-	0.035
	$\frac{1}{4}$ 50% Survival	-	0.018
	$\frac{1}{2}$ 50% Survival	-	0.035
	Other	L 2.5 H 5.0	

Project No. 2468

VI. NON-ACTIVATION PLATE TESTS

DATE: 11-20-74

Test	Compound	Concentration/plate	<u>TA-1535</u>		<u>TA-1537</u>		<u>TA-1538</u>	
			T-1	T-2	T-1	T-2	T-1	T-2
PC	EMS	0.05 ml undiluted chemical	>10 ³	>10 ³				
	QM	0.25 mg			>10 ²	>10 ²		
	NF	0.25 mg					>10 ²	>10 ²
SC	SALINE	-	2	1	2	4		
	DMSO	<10%					5	1

NOTE: PC = positive control
 SC = solvent control
 T-1 = trial 1
 T-2 = trial 2
 EMS = ethyl methanesulfonate
 QM = quinacrine mustard
 NF = nitrosofluorene
 DMSO = dimethyl sulfoxide
 (c) = contamination present

Project No. 2468

NON-ACTIVATION PLATE TESTS

DATE: 11-20-74

Test	Compound	Concentration	<u>TA-1535</u>		<u>TA-1537</u>		<u>TA-1538</u>	
			T-1	T-2	T-1	T-2	T-1	T-2
TC	FDA 73-3	0.035%	1	1	2	1	1	3

NOTE: TC = test compound
T-1 = trial 1
T-2 = trial 2
(c) = contamination present

Project No. 2468

VII. ACTIVATION PLATE TESTS

SPECIES: MOUSE

DATE: 11-20-74

Test	Organ	Compound	Concentration/plate	TA-1535		TA-1537		TA-1538	
				T-1	T-2	T-1	T-2	T-1	T-2
PC	Li	DMNA	25 μ moles	>10 ³	>10 ³				
		AAF	1.25 mg			44	43	>10 ²	>10 ²
	Lu	DMNA	25 μ moles	2	4				
		AAF	1.25 mg			9	3	13	8
	T	DMNA	25 μ moles	1	5				
		AAF	1.25 mg			6	10	3	3
SC	-	DMNA	25 μ moles	3	0				
	-	AAF	1.25 mg			10	5	1	0
	-	Saline	-	1	1				
	-	DMSO	<10%			12	10	6	7

NOTE: PC = positive control
 SC = solvent and chemical controls
 AAF = 2-acetylaminofluorene
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung

T = testes
 T-1 = trial 1
 T-2 = trial 2
 DMSO = dimethyl sulfoxide
 (c) = contamination present

Project No. 2468

ACTIVATION PLATE TESTS

SPECIES: MOUSE

DATE: 11-20-74

Test	Organ	Compound	Concentration	<u>TA-1535</u>		<u>TA-1537</u>		<u>TA-1538</u>	
				T-1	T-2	T-1	T-2	T-1	T-2
TC	Li	FDA 73-3	0.035%	1	0	8	9	(c)	22
	Lu	FDA 73-3	0.035%	2	0	8	4	0	5
	T	FDA 73-3	0.035%	1	4	10	7	5	5

NOTE: TC = test compound
 Li = liver
 Lu = lung
 T = testes
 T-1 = trial 1
 T-2 = trial 2
 (c) = contamination present

Project No. 2468

ACTIVATION PLATE TESTS

SPECIES: RAT

DATE: 11-20-74

Test	Organ	Compound	Concentration/plate	TA-1535		TA-1537		TA-1538	
				T-1	T-2	T-1	T-2	T-1	T-2
PC	Li	DMNA	25 μ moles	>10 ²	>10 ²				
		AAF	1.25 mg			41	30	>10 ²	>10 ²
	Lu	DMNA	25 μ moles	1	0				
		AAF	1.25 mg			7	10	5	0
	T	DMNA	25 μ moles	3	0				
		AAF	1.25 mg			14	17	10	3
SC	-	DMNA	25 μ moles	3	0				
	-	AAF	1.25 mg			10	5	1	0
	-	Saline	-	1	1				
	-	DMSO	<10%			12	10	6	7

NOTE:

PC = positive control
 SC = solvent and chemical controls
 AAF = 2-acetylaminofluorene
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung

T = testes
 T-1 = trial 1
 T-2 = trial 2
 DMSO = dimethyl sulfoxide
 (c) = contamination present

Project No. 2468

ACTIVATION PLATE TESTS

SPECIES: RAT

DATE: 11-20-74

Test	Organ	Compound	Concentration	<u>TA-1535</u>		<u>TA-1537</u>		<u>TA-1538</u>	
				T-1	T-2	T-1	T-2	T-1	T-2
TC	Li	FDA 73-3	0.035%	0	2	12	12	7	(c)
	Lu	FDA 73-3	0.035%	3	2	12	6	0	3
	T	FDA 73-3	0.035%	1	0	9	8	10	5

NOTE: TC = test compound
 Li = liver
 Lu = lung
 T = testes
 T-1 = trial 1
 T-2 = trial 2
 (c) = contamination present

Project No. 2468

ACTIVATION PLATE TESTS

SPECIES: MONKEY

DATE: 11-20-74

Test	Organ	Compound	Concentration/plate	TA-1535		TA-1537		TA-1538	
				T-1	T-2	T-1	T-2	T-1	T-2
PC	Li	DMNA	25 μ moles	>10 ²	>10 ²				
		AAF	1.25 mg			32	46	>10 ²	>10 ²
	Lu	DMNA	25 μ moles	0	4				
		AAF	1.25 mg			20	13	2	4
	T	DMNA	25 μ moles	1	1				
		AAF	1.25 mg			10	11	3	7
SC	-	DMNA	25 μ moles	3	0				
	-	AAF	1.25 mg			10	5	1	0
	-	Saline	-	1	1				
	-	DMSO	<10%			12	10	6	7

NOTE: PC = positive control
 SC = solvent and chemical controls
 AAF = 2-acetylaminofluorene
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung

T = testes
 T-1 = trial 1
 T-2 = trial 2
 DMSO = dimethyl sulfoxide
 (c) = contamination present

Project No. 2468

ACTIVATION PLATE TESTS

SPECIES: MONKEY

DATE: 11-20-74

Test	Organ	Compound	Concentration	<u>TA-1535</u>		<u>TA-1537</u>		<u>TA-1538</u>	
				T-1	T-2	T-1	T-2	T-1	T-2
TC	Li	FDA 73-3	0.035%	1	2			10	8
	Lu	FDA 73-3	0.035%	1	0			10	5
	T	FDA 73-3	0.035%	0	0			11	7

NOTE: TC = test compound
 Li = liver
 Lu = lung
 T = testes
 T-1 = trial 1
 T-2 = trial 2
 (c) = contamination present

Project No. 2468

VIII. NON-ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS:
POSITIVE AND SOLVENT CONTROL RESULTS

DATE: 10-9-74

Test	Indicator Strain	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	TA-1535	EMS	0.05 %	6.33	6990	1104.27
	TA-1537	QM	0.01 mg/ml	4.05	469	115.80
	TA-1538	NF	1.25 mg/ml	4.92	241	48.98
SC	TA-1535	SALINE	-	5.47	8	1.46
	TA-1537	SALINE	-	4.32	51	11.81
	TA-1538	DMSO	-	5.09	54	10.61

NOTE: PC = positive control
 SC = solvent control
 EMS = ethyl methanesulfonate
 QM = quinacrine mustard
 NF = nitrosofluorene
 DMSO = dimethyl sulfoxide
 (c) = contamination present

Project No. 2468



BIONETICS

NON-ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS

DATE: 10-9-74

Test	Indicator Strain	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
TC	TA-1535	FDA 73-3	H	5.67(104)	12	2.12
TC	TA-1535	FDA 73-3	L	5.56(102)	7	1.26
TC	TA-1537	FDA 73-3	H	4.18(97)	43	10.29
TC	TA-1537	FDA 73-3	L	4.94(114)	45	9.11
TC	TA-1538	FDA 73-3	H	5.14(101)	38	7.39
TC	TA-1538	FDA 73-3	L	7.12(140)	72	10.11

NOTE: TC = test compound
H = high dose
L = low dose
(c) = contamination present
() = percent survival

Project No. 2468



BIONETICS

**IX. ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS:
POSITIVE AND SOLVENT CONTROL RESULTS**

SPECIES: MOUSE

DATE: 10-7-74

Strain TA-1535

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	DMNA	100 μ moles/ml	3.00	2195	731.67
	Lu	DMNA	100 μ moles/ml	1.67	11	6.59
	T	DMNA	100 μ moles/ml	1.78	10	5.62
SC	-	DMNA	100 μ moles/ml	5.41	12	2.22
	-	SALINE	-	4.54	11(c)	2.42

DATE: 10-8-74

Strain TA-1537

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml	5.63	85	15.10
	Lu	AAF	1.25 mg/ml	5.86	29	4.95
	T	AAF	1.25 mg/ml	5.53	12	2.17
SC	-	AAF	1.25 mg/ml	4.24	36	8.49
	-	DMSO	-	5.74	38	6.62

DATE: 10-9-74

Strain TA-1538

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml	8.58	256	29.84
	Lu	AAF	1.25 mg/ml	7.74	55	7.11
	T	AAF	1.25 mg/ml	6.17	53	8.59
SC	-	AAF	1.25 mg/ml	6.79	48	7.07
	-	DMSO	-	7.90	46	5.82

NOTE: PC = positive control
 SC = solvent and chemical controls
 AAF = 2-acetylaminofluorene
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung
 T = testes
 DMSO = dimethyl sulfoxide

(c) = contamination present

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS

SPECIES: MOUSE

DATE: 10-7-74

Strain TA-1535

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
TC	Li	FDA 73-3	H	2.45(54)	4	1.63
		FDA 73-3	L	1.64(36)	3	1.83
	Lu	FDA 73-3	H	3.91(86)	6	1.54
		FDA 73-3	L	2.32(51)	4	1.72
	T	FDA 73-3	H	1.85(41)	8(c)	4.32
		FDA 73-3	L	1.63(36)	5	3.07

DATE: 10-8-74

Strain TA-1537

TC	Li	FDA 73-3	H	4.21(73)	33	7.84
		FDA 73-3	L	5.92(103)	50	8.45
	Lu	FDA 73-3	H	6.08(106)	33	5.43
		FDA 73-3	L	5.07(88)	29(c)	5.72
	T	FDA 73-3	H	5.63(98)	22	3.91
		FDA 73-3	L	5.83(102)	37	6.35

DATE: 10-9-74

Strain TA-1538

TC	Li	FDA 73-3	H	6.97(88)	46	6.60
		FDA 73-3	L	7.63(97)	35	4.59
	Lu	FDA 73-3	H	9.50(120)	42	4.42
		FDA 73-3	L	8.91(113)	32	3.59
	T	FDA 73-3	H	4.48(57)	17	3.80
		FDA 73-3	L	7.19(91)	20	2.78

NOTES: H = high dose
L = low dose
TC = test compound
Li = liver
Lu = lung
T = testes
(c) = contamination present
() = percent survival

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS:
POSITIVE AND SOLVENT CONTROL RESULTS

SPECIES: RAT

DATE: 10-11-74

Strain TA-1535

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	DMNA	100 μ moles/ml	8.06	1980	245.70
	Lu	DMNA	100 μ moles/ml	15.41	15	0.97
	T	DMNA	100 μ moles/ml	7.72	13	1.68
SC	-	DMNA	100 μ moles/ml	9.31	31	3.33
	-	SALINE	-	11.22	30	2.67

DATE: 10-31-74

Strain TA-1537

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml	3.13	92	29.39
	Lu	AAF	1.25 mg/ml	5.31	54	10.17
	T	AAF	1.25 mg/ml	2.83	25	8.83
SC	-	AAF	1.25 mg/ml	2.34	40	17.09
	-	DMSO	-	3.99	30	7.52

DATE: 10-23-74

Strain TA-1538

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml	2.29	98	42.80
	Lu	AAF	1.25 mg/ml	7.97	47	5.90
	T	AAF	1.25 mg/ml	6.41	39	6.08
SC	-	AAF	1.25 mg/ml	9.02	35	3.88
	-	DMSO	-	5.97	57	9.55

NOTE: PC = positive control
 SC = solvent and chemical controls (c) = contamination present
 AAF = 2-acetylaminofluorene
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung
 T = testes
 DMSO = dimethyl sulfoxide

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS

SPECIES: RAT

				Strain TA-1535		
DATE: 10-11-74						
Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
TC	Li	FDA 73-3	H	10.54(94)	14	1.33
		FDA 73-3	L	10.71(95)	22	2.05
	Lu	FDA 73-3	H	8.98(80)	11	1.23
		FDA 73-3	L	18.30(163)	13	0.71
	T	FDA 73-3	H	8.72(78)	13	1.49
		FDA 73-3	L	18.04(161)	10	0.55

				Strain TA-1537		
DATE: 10-31-74						
TC	Li	FDA 73-3	H	3.04(76)	26	8.55
		FDA 73-3	L	3.72(93)	32	8.60
	Lu	FDA 73-3	H	3.20(80)	37	11.56
		FDA 73-3	L	3.14(79)	37	11.78
	T	FDA 73-3	H	0.45(11)	11	24.44
		FDA 73-3	L	2.33(58)	8	3.43

				Strain TA-1538		
DATE: 10-23-74						
TC	Li	FDA 73-3	H	1.78(30)	37	20.79
		FDA 73-3	L	(c)2.15(36)	27	12.56
	Lu	FDA 73-3	H	2.09(35)	28	13.40
		FDA 73-3	L	2.95(49)	46	15.59
	T	FDA 73-3	H	1.03(17)	7	6.80
		FDA 73-3	L	2.97(50)	6	2.02

NOTES: H = high dose
L = low dose
TC = test compound
Li = liver
Lu = lung
T = testes
(c) = contamination present
() = percent survival

Project No. 2468



BIONETICS

**ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS:
POSITIVE AND SOLVENT CONTROL RESULTS**

SPECIES: RAT

DATE:

Strain TA-1535

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	DMNA	100 μ moles/ml			
	Lu	DMNA	100 μ moles/ml			
	T	DMNA	100 μ moles/ml			
SC	-	DMNA	100 μ moles/ml			
	-	SALINE	-			

DATE:

Strain TA-1537

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml			
	Lu	AAF	1.25 mg/ml			
	T	AAF	1.25 mg/ml			
SC	-	AAF	1.25 mg/ml			
	-	DMSO	-			

DATE: 11-25-74 (Repeated Doses)

Strain TA-1538

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml			
	Lu	AAF	1.25 mg/ml			
	T	AAF	1.25 mg/ml			
SC	-	AAF	1.25 mg/ml			
	-	DMSO	-	3.89	14	3.60

NOTE: PC = positive control
 SC = solvent and chemical controls
 AAF = 2-acetylaminofluorene
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung
 T = testes
 DMSO = dimethyl sulfoxide

(c) = contamination present

Project No. 2468



BIONETICS

Litton

ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS

SPECIES: RAT

DATE :				Strain TA-1535		
Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
TC	Li		H			
			L			
	Lu		H			
			L			
	T		H			
			L			

DATE :		Strain TA-1537	
TC	Li	H	
		L	
	Lu	H	
		L	
	T	H	
		L	

DATE: 11-25-74 (Repeated Doses)				Strain TA-1538		
TC	Li	FDA 73-3	H	3.93(101)	13	3.31
		FDA 73-3	L	3.62(93)	7	1.93
	Lu	FDA 73-3	H	2.32(60)	6	2.59
		FDA 73-3	L	3.62(93)	8	2.21
	T		H			
			L			

NOTES: H = high dose
 L = low dose
 TC = test compound
 Li = liver
 Lu = lung
 T = testes
 (c) = contamination present
 () = percent survival

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS:
POSITIVE AND SOLVENT CONTROL RESULTS

SPECIES: MONKEY

DATE: 10-25-74

Strain TA-1535

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	DMNA	100 μ moles/ml	5.24	1082	206.49
	Lu	DMNA	100 μ moles/ml	6.83	17	2.49
	T	DMNA	100 μ moles/ml	4.66	13	2.79
SC	-	DMNA	100 μ moles/ml	5.85	25	4.27
	-	SALINE	-	8.83	22	2.49

DATE: 10-24-74

Strain TA-1537

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml	3.78	101	26.72
	Lu	AAF	1.25 mg/ml	8.80	76	8.64
	T	AAF	1.25 mg/ml	8.68	60	6.91
SC	-	AAF	1.25 mg/ml	4.16	76	18.27
	-	DMSO	-	4.20	68	16.19

DATE: 10-30-74

Strain TA-1538

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
PC	Li	AAF	1.25 mg/ml	2.21	67	30.32
	Lu	AAF	1.25 mg/ml	2.03	28	13.79
	T	AAF	1.25 mg/ml	2.64	23	8.71
SC	-	AAF	1.25 mg/ml	2.94	27	9.18
	-	DMSO	-	3.25	45	13.85

NOTE: PC = positive control
SC = solvent and chemical controls
AAF = 2-acetylaminofluorene
DMNA = dimethylnitrosamine
Li = liver
Lu = lung
T = testes
DMSO = dimethyl sulfoxide

(c) = contamination present

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SALMONELLA INDICATOR STRAINS

SPECIES: MONKEY

DATE: 10-25-74

Strain TA-1535

Test	Organ	Compound	Concentration	Total Cells/ mlx10 ⁸	his+ Revertants/ ml	his+ Revertants/10 ⁸ Survivors
TC	Li	FDA 73-3	H	6.14(70)	10	1.63
		FDA 73-3	L	6.57(74)	27	4.11
	Lu	FDA 73-3	H	5.45(62)	10	1.84
		FDA 73-3	L	6.38(72)	26(c)	4.08
	T	FDA 73-3	H	4.23(48)	8	1.89
		FDA 73-3	L	3.48(39)	12	3.45

DATE: 10-24-74

Strain TA-1537

TC	Li	FDA 73-3	H	(c)3.06(73)	56	18.30
		FDA 73-3	L	4.53(108)	83	18.32
	Lu	FDA 73-3	H	5.69(135)	68	11.95
		FDA 73-3	L	(c)4.59(109)	63	13.73
	T	FDA 73-3	H	3.19(76)	29(c)	9.09
		FDA 73-3	L	3.90(93)	31(c)	7.95

DATE: 10-30-74

Strain TA-1538

TC	Li	FDA 73-3	H	2.82(87)	35	12.41
		FDA 73-3	L	2.83(87)	33	11.66
	Lu	FDA 73-3	H	3.18(98)	51	16.04
		FDA 73-3	L	3.33(102)	55	16.52
	T	FDA 73-3	H	2.71(83)	16	5.90
		FDA 73-3	L	4.08(126)	13(c)	3.19

NOTES: H = high dose
L = low dose
TC = test compound
Li = liver
Lu = lung
T = testes
(c) = contamination present
() = percent survival

Project No. 2468



BIONETICS

X. NON-ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4

DATE: 11-1-74

Strain D4							
Test	Compound	Concentration	Total Population Screened ^a	Number of Convertants ^b		Convertants Per 10 ⁵ Survivors	
				Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
PC	EMS	1.0 %	7.98	715	806	89.60	112.73
SC	Saline	-	9.62	63	42	6.55	4.37

NOTE: PC = positive control
 SC = solvent control
 EMS = ethyl methanesulfonate

a = number x 10⁵
 b = number at 10⁻¹ dilution
 (c) = contamination present

Project No. 2468



BIONETICS

NON-ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4

DATE: 11-1-74

Strain D4							
Test	Compound	Concentration	Total Population Screened ^a	Number Convertants ^b		Convertants Per 10 ⁵ Survivors	
				Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
TC	FDA 73-3	H	7.30(76)	47	36(c)	6.44	4.93
	FDA 73-3	L	(c)6.48(67)	38	25(c)	5.86	3.86

NOTE: TC = test compound
 H = high dose
 L = low dose
 a = number x 10⁵
 b = number at 10⁻¹ dilution
 (c) = contamination present
 () = percent survival

Project No. 2468



BIONETICS

**XI. ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4:
POSITIVE AND SOLVENT CONTROL RESULTS**

SPECIES: MOUSE

DATE: 11-26-74

Strain D4								
Test	Organ	Compound	Concentration	Total Population Screened ^a	Number of Convertants ^b		Convertants Per 10 ⁵ Survivors	
					Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
PC	Li	DMNA	150 μ moles/ml	7.93	65	70	8.20	8.83
	Lu	DMNA	150 μ moles/ml	7.22	44	33	6.09	4.57
	T	DMNA	150 μ moles/ml	7.39	27	40	3.65	5.41
SC	-	DMNA	150 μ moles/ml	8.76	54	31(c)	6.16	3.54
	-	SALINE	-	8.66	48	40	5.54	4.62

NOTE: PC = positive control
 SC = solvent and chemical controls
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung
 T = testes

 a = number x 10⁵
 b = number at 10⁻¹ dilution
 (c) = contamination present

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4

SPECIES: MOUSE

DATE: 11-26-74

				Strain D4				
Test	Organ	Compound	Concentration	Total Population Screened ^a	Number of Converstants ^b		Converstants Per 10 ⁵ Survivors	
					Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
TC	Li	FDA 73-3	H	6.53(75)	43	42	6.58	6.43
		FDA 73-3	L	9.89(114)	76(c)	46	7.68	4.65
	Lu	FDA 73-3	H	8.08(93)	48(c)	45(c)	5.94	5.57
		FDA 73-3	L	7.60(88)	46(c)	45	6.05	5.92
	T	FDA 73-3	H	8.37(97)	50	39	5.97	4.66
		FDA 73-3	L	8.02(93)	38	34	4.74	4.24

NOTE:

TC = test compound
H = high dose
L = low dose
Li = liver
Lu = lung
T = testes
a = number x 10⁵
b = number at 10⁻¹ dilution
(c) = contamination present
() = percent survival

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4:
POSITIVE AND SOLVENT CONTROL RESULTS

SPECIES: RAT

DATE: 10-25-74

Strain D4

Test	Organ	Compound	Concentration	Total Population Screened ^a	Number of Convertants ^b		Convertants Per 10 ⁵ Survivors	
					Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
PC	Li	DMNA	150 μ moles/ml	5.84(c)	55	60	9.42	10.27
	Lu	DMNA	150 μ moles/ml	8.18	32	29	3.91	3.55
	T	DMNA	150 μ moles/ml	4.45	1	9	0.23	2.02
SC	-	DMNA	150 μ moles/ml	9.61	26	24	2.71	2.50
	-	SALINE	-	12.35	42	36	3.40	2.92

NOTE: PC = positive control
 SC = solvent and chemical controls
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung
 T = testes

 a = number x 10⁵
 b = number at 10⁻¹ dilution
 (c) = contamination present

ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4

SPECIES: RAT

DATE: 10-25-74

Strain D4								
Test	Organ	Compound	Concentration	Total Population Screened ^a	Number of Converstants ^b		Converstants Per 10 ⁵ Survivors	
					Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
TC	Li	FDA 73-3	H	7.04(57)	0	25	0	3.55
		FDA 73-3	L	9.77(79)	20(c)	31(c)	2.05	3.17
	Lu	FDA 73-3	H	9.78(79)	34	36	3.48	3.68
		FDA 73-3	L	7.97(65)	30	28(c)	3.76	3.51
T		FDA 73-3	H	6.96(56)	29	16	4.17	2.30
		FDA 73-3	L	9.07(73)	33	54(c)	3.64	5.95

NOTE: TC = test compound
 H = high dose
 L = low dose
 Li = liver
 Lu = lung
 T = testes
 a = number x 10⁵
 b = number at 10⁻¹ dilution
 (c) = contamination present
 () = percent survival

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4:
POSITIVE AND SOLVENT CONTROL RESULTS

SPECIES: MONKEY

DATE: 11-7-74

Strain D4								
Test	Organ	Compound	Concentration	Total Population Screened ^a	Number of Convertants ^b		Convertants Per 10 ⁵ Survivors	
					Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
PC	Li	DMNA	150 μ moles/ml	7.44	73	21	9.81	2.82
	Lu	DMNA	150 μ moles/ml	9.39	37	34	3.94	3.62
	T	DMNA	150 μ moles/ml	7.04(c)	50	27	7.10	3.84
SC	-	DMNA	150 μ moles/ml	9.12(c)	37	13	4.06	1.43
	-	SALINE	-	9.79	28	13	2.86	1.33

NOTE: PC = positive control
 SC = solvent and chemical controls
 DMNA = dimethylnitrosamine
 Li = liver
 Lu = lung
 T = testes

a = number $\times 10^5$
 b = number at 10^{-1} dilution
 (c) = contamination present

Project No. 2468



BIONETICS

ACTIVATION SUSPENSION TESTS
WITH SACCHAROMYCES INDICATOR STRAIN D4

SPECIES: MONKEY

DATE: 11-7-74

Strain D4

Test	Organ	Compound	Concentration	Total Population Screened ^a	Number of Conyertants ^b		Conyertants Per 10 ⁵ Survivors	
					Ade ⁺	Try ⁺	Ade ⁺	Try ⁺
TC	Li	FDA 73-3	H	10.27(105)	45	15	4.38	1.46
		FDA 73-3	L	9.19(94)	25	43	2.72	4.68
	Lu	FDA 73-3	H	8.12(83)	30	20	3.69	2.46
		FDA 73-3	L	7.86(81)	48	25	6.11	3.18
	T	FDA 73-3	H	10.01(102)	15	35	1.50	3.50
		FDA 73-3	L	7.51(77)	64	20	8.64	2.66

NOTE: TC = test compound
H = high dose
L = low dose
Li = liver
Lu = lung
T = testes
a = number x 10⁵
b = number at 10⁻¹ dilution
(c) = contamination present
() = percent survival

Project No. 2468



BIONETICS

XII. SUMMARY OF TEST RESULTS

COMPOUND FDA 73-3

A. Suspension Tests

Test ^a	<u>Activation</u>		<u>Salmonella Reversion Frequencies (x 10⁻⁸)</u>			<u>Saccharomyces D4 Conversion Frequencies (x 10⁻⁵)</u>	
	Species ^b	Organ ^c	TA-1535	TA-1537	TA-1538	Ade ⁺	Try ⁺
NA-PC	-	-	1104.27	115.80	48.98	89.60	112.73
NA-NC	-	-	1.46	11.81	10.61	6.55	4.37
NA-H	-	-	2.12	10.29	7.39	6.44	4.93
NA-L	-	-	1.26	9.11	10.11	5.86	3.86
A-NC (-C)	-	-	2.42	6.62	5.82	5.54	4.62
A-NC (+C)	-	-	2.22	8.49	7.07	6.16	3.54
A-PC	M	Li	731.67	15.10	29.84	8.20	8.83
A-PC	M	Lu	6.59	4.95	7.11	6.20	4.57
A-PC	M	T	5.62	2.17	8.59	6.09	5.41
A-H	M	Li	1.63	7.84	6.60	6.58	6.43
A-L			1.83	8.45	4.59	7.68	4.65
A-H	M	Lu	1.54	5.43	4.42	5.94	5.57
A-L			1.72	5.72	3.59	6.05	5.92
A-H	M	T	4.32	3.91	3.80	5.97	4.66
A-L			3.07	6.35	2.78	4.74	4.24

^a NA = non activation
NC = negative control
PC = positive control
A = activation
H = high dose
L = low dose

^b M = mouse
Mo = monkey
R = rat

^c Li = liver
Lu = lung
T = testes

(-C) = solvent control
(+C) = chemical control

COMPOUND FDA 73-3

B. Plate Tests

Test ^a	<u>Activation</u>		<u>Salmonella Responses</u>		
	Species ^b	Organ ^c	TA-1535	TA-1537	TA-1538
NA-PC	-	-	+	+	+
NA-NC	-	-	-	-	-
NA-H	-	-	-	-	-
A-NC (-C)	-	-	-	-	-
A-NC (+C)	-	-	-	-	-
A-PC	M	Li	+	+	+
A-PC	M	Lu	-	-	-
A-PC	M	T	-	-	-
A-H	M	Li	-	-	-
A-H	M	Lu	-	-	-
A-H	M	T	-	-	-

^a NA = non activation
 NC = negative control
 PC = positive control
 A = activation
 H = high dose
 L = low dose

^b M = mouse
 Mo = monkey
 R = rat

^c Li = liver
 Lu = lung
 T = testes

(-C) = solvent control
 (+C) = chemical control

Project 2468

SUMMARY OF TEST RESULTSCOMPOUND FDA 73-3A. Suspension Tests

Test ^a	<u>Activation</u>		<u>Salmonella Reversion Frequencies (x 10⁻⁸)</u>			<u>Saccharomyces D4 Conversion Frequencies (x 10⁻⁵)</u>	
	Species ^b	Organ ^c	TA-1535	TA-1537	TA-1538	Ade ⁺	Try ⁺
NA-PC	-	-					
NA-NC	-	-					
NA-H	-	-					
NA-L	-	-					
A-NC (-C)	-	-	2.67	7.52	9.55	3.40	2.92
A-NC (+C)	-	-	3.33	17.09	3.88	2.71	2.50
A-PC	R	Li	245.70	29.39	42.80	9.42	10.27
A-PC	R	Lu	0.97	10.17	5.90	3.91	3.55
A-PC	R	T	1.68	8.83	6.08	0.23	2.02
A-H	R	Li	0.97	8.55	3.31*	--	3.55
A-L			1.69	8.60	1.93*	2.05	3.17
A-H	R	Lu	1.23	11.56	2.59*	3.48	3.68
A-L			2.38	11.78	2.21*	3.76	3.51
A-H	R	T	1.10	24.44	6.80	4.17	2.30
A-L			1.62	3.43	2.02	3.64	5.95

^a NA = non activation
 NC = negative control
 PC = positive control
 A = activation
 H = high dose
 L = low dose

^b M = mouse
 Mo = monkey
 R = rat

^c Li = liver
 Lu = lung
 T = testes

(-C) = solvent control
 (+C) = chemical control

* Date from repeat tests

COMPOUND FDA 73-3

B. Plate Tests

Test ^a	<u>Activation</u>		<u>Salmonella Responses</u>		
	Species ^b	Organ ^c	TA-1535	TA-1537	TA-1538
NA-PC	-	-			
NA-NC	-	-			
NA-H	-	-			
A-NC (-C)	-	-	-	-	-
A-NC (+C)	-	-	-	-	-
A-PC	R	Li	+	+	+
A-PC	R	Lu	-	-	-
A-PC	R	T	-	-	-
A-H	R	Li	-	-	-
A-H	R	Lu	-	-	-
A-H	R	T	-	-	-

^a NA = non activation
 NC = negative control
 PC = positive control
 A = activation
 H = high dose
 L = low dose

^b M = mouse
 Mo = monkey
 R = rat

^c Li = liver
 Lu = lung
 T = testes

(-C) = solvent control
 (+C) = chemical control

Project 2468



BIONETICS

SUMMARY OF TEST RESULTS

COMPOUND FDA 73-3

A. Suspension Tests

Test ^a	<u>Activation</u>		<u>Salmonella Reversion Frequencies (x 10⁻⁸)</u>			<u>Saccharomyces D4 Conversion Frequencies (x 10⁻⁵)</u>	
	Species ^b	Organ ^c	TA-1535	TA-1537	TA-1538	Ade ⁺	Try ⁺
NA-PC	-	-					
NA-NC	-	-					
NA-H	-	-					
NA-L	-	-					
A-NC (-C)	-	-	2.49	16.19	13.85	2.86	1.33
A-NC (+C)	-	-	4.27	18.27	9.18	4.06	1.43
A-PC	Mo	Li	206.49	26.72	30.32	9.81	2.82
A-PC	Mo	Lu	2.49	8.64	13.79	3.94	3.62
A-PC	Mo	T	2.79	6.91	8.71	7.10	3.84
A-H	Mo	Li	1.63	18.30	12.41	4.38	1.46
A-L			4.11	18.32	11.66	2.72	4.68
A-H	Mo	Lu	1.84	11.95	16.04	3.69	2.46
A-L			4.08	13.73	16.52	6.11	3.18
A-H	Mo	T	1.89	9.09	5.90	1.50	3.50
A-L			3.45	7.95	3.19	8.64	2.66

^a NA = non activation
NC = negative control
PC = positive control
A = activation
H = high dose
L = low dose

^b M = mouse
Mo = monkey
R = rat

^c Li = liver
Lu = lung
T = testes

(-C) = solvent control
(+C) = chemical control

Project 2468

COMPOUND FDA 73-3B. Plate Tests

Test ^a	<u>Activation</u>		<u>Salmonella Responses</u>		
	Species ^b	Organ ^c	TA-1535	TA-1537	TA-1538
NA-PC	-	-			
NA-NC	-	-			
NA-H	-	-			
A-NC (-C)	-	-	-	-	-
A-NC (+C)	-	-	-	-	-
A-PC	Mo	Li	+	+	+
A-PC	Mo	Lu	-	-	-
A-PC	Mo	T	-	-	-
A-H	Mo	Li	-	-	-
A-H	Mo	Lu	-	-	-
A-H	Mo	T	-	-	-

^a NA = non activation
 NC = negative control
 PC = positive control
 A = activation
 H = high dose
 L = low dose

^b M = mouse
 Mo = monkey
 R = rat

^c Li = liver
 Lu = lung
 T = testes

(-C) = solvent control
 (+C) = chemical control

Project 2468

XIII. INTERPRETATION AND CONCLUSIONS

Compound FDA 73-3, Sodium Hexametaphosphate, was evaluated for genetic activity in a series of in vitro microbial assays with and without metabolic activation. The following results were obtained:

A. Salmonella typhimurium

1. Plate Tests

At a concentration of 0.00025%, this compound was not mutagenic for TA-1535, TA-1537 or TA-1538 in direct or activation plate tests.

2. Non-activation Suspension Tests

These tests were negative.

3. Activation Suspension Tests

These tests were all negative. Several dose levels with rat tissues were repeated using strain TA-1538. The initial tests appeared unusually high, but the repeat data indicated that the original tests were aberrant and did not reflect true mutagenic activity.

The test with rat testes and TA-1537 at the high dose level of FDA 73-3 was unusually high, but no other data supported the response and it was considered to be a random test fluctuation or the result of contaminants.

B. Saccharomyces cerevisiae

1. Non-activation Suspension Tests

These tests were negative.

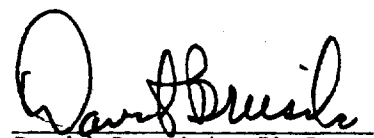
2. Activation Suspension Tests

These tests were negative.

C. Conclusions

Compound FDA 73-3, Sodium Hexametaphosphate, was not genetically active for bacterial and yeast indicator organisms under the conditions of this evaluation.

SUBMITTED BY:


David Brusick, Ph.D.
Director
Department of Genetics



BIONETICS

APPENDIX

SUMMARY OF TESTS EVALUATING DMSO FOR GENETIC
ACTIVITY IN SALMONELLA AND SACCHAROMYCES



BIONETICS

COMPOUND DIMETHYSULFOXIDEA. Suspension Tests

Test	Activation		Salmonella Reversion Frequencies ($\times 10^{-8}$)		Saccharomyces D4 Conversion Frequencies ($\times 10^{-5}$)	
	Species ^a	Organ ^b	TA-1535	TA-1538	Ade ⁺	Try ⁺
<u>Non-activation</u>						
Control (-C)	-	-	0.74	0.88	32.51	4.34
High Dose ^c	-	-	1.91	1.05	28.32	2.95
Low Dose ^d	-	-	0.53	1.37	40.73	0.49
<u>Activation</u>						
Control (+C)	-	-	1.80	0.36	38.27	2.47
Control (-C)	-	-	1.43	1.04	37.12	2.64
High Dose ^c	M	Li	0.34	1.07	47.77	2.75
	M	Lu	0.59	0.58	36.29	1.39
	M	T	0.62	0.30	34.35	1.94
Low Dose ^d	M	Li	-	0.87	34.02	1.18
	M	Lu	0.43	3.14	42.30	1.40
	M	T	0.11	0.39	45.95	2.32

Note: (-C) = solvent control and (+C) = test chemical control without homogenate

a M = mouse
Mo = monkey
R = rat

b Li = liver
Lu = lung
T = testes

c Bacteria = 3%
Yeast = 5%

d Bacteria = 1.5%
Yeast = 2.5%



BIONETICS

COMPOUND DIMETHYSULFOXIDE

B. Plate Tests

<u>Test</u>	<u>Activation</u>		<u>Salmonella Responses</u>		
	<u>Species^a</u>	<u>Organ^b</u>	TA-1535	TA-1537	TA-1538
<u>Non-activation</u>					
Control (-C)	-	-	-	-	-
Test compound (3%)	-	-	-	-	-
<u>Activation</u>					
Control (+C)	-	-	-	-	-
Control (-C)	-	-	-	-	-
Test compound (3%)	M	Li	-	-	-
	M	Lu	-	-	-
	M	T	-	-	-
	R	Li	-	-	-
	R	Lu	-	-	-
	R	T	-	-	-
	Mo	Li	-	-	-
	Mo	Lu	-	-	-
	Mo	T	-	-	-

Note: (-C) = solvent control and (+C) = chemical control without homogenate

a M = mouse
Mo = monkey
R = rat

b Li = liver
Lu = lung
T = testes