

BIONETICS

SUMMARY OF MUTAGENICITY
SCREENING STUDIES
HOST-MEDIATED ASSAY
CYTOGENETICS
DOMINANT LETHAL ASSAY
CONTRACT FDA 71-268
COMPOUND FDA 71-51

METHOCEL

E20

Summary of mutagenicity screening studies, host-mediated assay cytogenetics dominant
lethal assay-Contract FDA 71-268 & Compound FDA 71-51 (Methocel) 11/22/74

LBI PROJECT #2446

SUMMARY OF MUTAGENICITY
SCREENING STUDIES
HOST-MEDIATED ASSAY
CYTOGENETICS
DOMINANT LETHAL ASSAY
CONTRACT FDA 71-268
COMPOUND FDA 71-51
METHOCEL

SUBMITTED TO

FOOD & DRUG ADMINISTRATION
DEPARTMENT OF HEALTH, EDUCATION AND WELFARE
ROCKVILLE, MARYLAND

SUBMITTED BY

LITTON BIONETICS, INC.
5516 NICHOLSON LANE
KENSINGTON, MARYLAND

NOVEMBER 22, 1974



BIONETICS

TABLE OF CONTENTS

	Page No.
I. REPORT	1
A. Introduction	2
B. Objective	3
C. Compound	3
1. Test Material	3
2. Dosages	3
D. Methods	4
E. Summary	4
1. Host-Mediated Assay	4
2. Cytogenetics	4
a. <u>In vivo</u>	4
b. <u>In vitro</u>	5
3. Dominant Lethal	5
F. Results and Discussion	5
1. Toxicity Data - Test I	5
a. <u>In vivo</u>	5
b. <u>In vitro</u>	6
c. Toxicity data sheets	8
2. Host-Mediated Assay - Test I	12
a. Host-mediated assay summary sheets.	14
b. Host-mediated assay data sheets ...	16
3. Toxicity Data - Test II	45
a. Acute	45
b. Subacute	45
c. Toxicity data sheets	46
4. Host-Mediated Assay - Test II	48
a. Host-mediated assay summary sheets.	49
b. Host-mediated assay data sheets ...	52
5. Cytogenetics - Test I	71
a. <u>In vivo</u>	71
b. <u>In vitro</u>	71
c. Cytogenetic summary sheets	72
6. Cytogenetics - Test II	76
Cytogenetic summary sheets	77
7. Dominant Lethal Study - Test I	80
a. Acute study	80
b. Subacute study	80
c. Dominant lethal assay summary sheets	81
8. Dominant Lethal Study - Test II	98
Dominant lethal assay summary sheets	99
II. MATERIALS AND METHODS	116
A. Animal Husbandry	116
1. Animals (Rats and Mice)	116
2. Preparation of Diet	116
3. Husbandry	116



BIONETICS

TABLE OF CONTENTS (continued)

Page No.

II. MATERIALS AND METHODS (continued)

B.	Dosage Determination	116
1.	Acute LD ₅₀ and LD ₅ Determination	116
2.	Subacute Studies	118
C.	Mutagenicity Testing Protocols	118
1.	Host-Mediated Assay	118
a.	Acute study	119
b.	Subacute study	121
c.	<u>In vitro</u> study	121
2.	Cytogenetic Studies	122
a.	<u>In vivo</u> study	122
b.	<u>In vitro</u> study	124
3.	Dominant Lethal Assay	126
D.	Supplementary Materials and Methods	127
1.	Host-Mediated Assay <u>In Vitro</u> and Formulae.	127
a.	Bacterial <u>in vitro</u> plate tests	127
b.	<u>In vitro</u> for mitotic recombination.	127
c.	Minimal medium (bacteria)	128
d.	Complete medium (bacteria)	129
e.	Complete medium (yeast)	129
2.	Cytogenetics <u>In Vitro</u> Preparation of Anaphase Chromosomes	130
3.	Statistical Analyses of Dominant Lethal Studies	131
a.	The fertility index	131
b.	Total number of implantations	131
c.	Total number of <u>corpora lutea</u>	131
d.	Preimplantation losses	131
e.	Dead implants	132
f.	One or more dead implants	132
g.	Two or more dead implants	132
h.	Dead implants per total implants ..	132
E.	References	135
1.	Host-Mediated Assay	135
2.	Cytogenetics	135
3.	Dominant Lethal	136
F.	Abbreviations	137



BIONETICS

I. REPORT

A. Introduction

Litton Bionetics, Inc. (LBI) has investigated the possible mutagenicity of compounds selected and provided by the Food and Drug Administration under Contract 71-268. LBI's investigation utilized the three mammalian test systems herein described -- Host-Mediated Assay, Cytogenetic Studies and Dominant Lethal Assay. These tests provide information as to the types of genetic damage caused by environmental compounds -- pesticides, chemicals, food additives, drugs and cosmetics.

The Host-Mediated Assay is based upon the assumption that the action of a mutagen on the genetics of bacteria is similar to that in man. This is further strengthened by the use of an eukaryotic organism (Saccharomyces cerevisiae). Since the mutation frequencies are well established for the indicator organism, any deviation due to the action of the test compound is readily detectable. As some compounds are mutagenic in bacteria and not in the host animal, and vice versa, this test is able to differentiate an action which may have been due to hosts' ability to detoxify or potentiate a suspected mutagen. This action is dependent upon the ability of the compound to gain access to the peritoneal cavity. Coupled with the direct action of the compound on the indicator organism in vitro, the assay provides a clear insight into host-mediation of mutagenicity.

Cytogenetics provides a valuable tool for the direct observation of chromosomal damage in somatic cells. Alteration of the chromosome number and/or form in somatic cells may be an index of mutation. These studies utilized examination of bone marrow cells arrested in C-metaphase from rats exposed to the test compound as compared to positive and negative control animals. If mutational



changes occur, the types of damage expected due to the action of chemicals are structural rearrangements, breaks and other forms of damage to the chromosomal complement of the cells exposed.

For the in vitro cytogenetic studies, we have a more rapid and inexpensive means of determining chromosomal damage. This is accomplished by observing cells in anaphase. As the chromatids separate and move along the spindle, aberrations may occur. Chromatids which do not migrate to the daughter cells may lead to uneven distribution of parts or of entire chromatids (mitotic nondysjunction). These give rise to "side arm" bridges which have been interpreted as point stickiness or localized failures of chromosome duplication point errors. These aberrations (bridges, pseudochiasmata, multipolar cells, acentric fragments, etc.) are extremely sensitive indicators of genetic damage.

The Dominant Lethal Test is an accurate and sensitive measure of the amount and type of fetal wastage which may occur following administration of a potential mutagen. Dominant lethal mutations are indicators of lethal genetic lesions. The effects of mutagens on the chromosomal complement of the spermatozoa of treated males results in alterations of form and number of chromosomes. Structural rearrangements and aneuploidy may lead to the production of non-viable zygotes, early and late fetal deaths, abortions and congenital malformations. In addition, aberrations could lead to sterility or reduced reproductive capacity of the F_1 generation. The action of a mutagen on specific portions of spermatogenesis is also apparent in this test.

B. Objective

The purpose of these studies is to determine any mutagenic effect of the test compound by employing the Host-Mediated Assay, Cytogenetic Studies



and the Dominant Lethal Assay, both in vivo and in vitro tests are employed with the cytogenetic and microbial test systems. These tests and their descriptions are referenced in the Appendices A through F.

C. Compound

1. Test Material

Compound FDA 71-51, Methocel, MC Premium, Lot Number MM030728C, as supplied by the Food and Drug Administration.

2. Dosages

The animal employed, the determination of the dosage levels and the route of administration are contained in the technical discussion.

The dosage levels employed for compound FDA 71-51 are as follows for the Cytogenetic Studies in vivo in rats.

	<u>Test I⁺</u>	<u>Test II⁺</u>
Low Level	4.75 mg/kg	-----
Intermediate Level	47.5 mg/kg	-----
LD ₅	475.0 mg/kg	5000 mg/kg
Negative Control	Saline	Food Pellet
Positive Control (TEM*)	0.3 mg/kg	0.3 mg/kg

The dosage levels employed for compound FDA 71-51 are as follows for the Host-Mediated Assay in vivo in mice.

	<u>Test I⁺</u>	<u>Test II⁺</u>
Low Level	4.75 mg/kg	-----
Intermediate Level	47.5 mg/kg	-----
LD ₅	475.0 mg/kg	5000 mg/kg
Negative Control	Saline	Saline
Positive Control (EMS**)	350 mg/kg	350 mg/kg
(DMN***)	100 mg/kg	100 mg/kg

* Triethylene Melamine

** Ethyl Methane Sulfonate

*** Dimethyl Nitrosamine

+ These two tests were performed at different time intervals.



BIONETICS

The dosage levels employed for compound FDA 71-51 are as follows for the Dominant Lethal Assay in vivo in rats.

	<u>Test I</u> ⁺	<u>Test II</u> ⁺
Low Level	4.75 mg/kg	-----
Intermediate Level	47.5 mg/kg	-----
LD ₅	475.0 mg/kg	5000 mg/kg
Negative Control	Saline	Saline
Positive Control (TEM*)	0.3 mg/kg	0.3 mg/kg

The in vitro Cytogenetic Studies were performed employing three logarithmic dose levels.

Low Level	8.0 mcg/ml
Medium Level	80.0 mcg/ml
High Level	800.0 mcg/ml
Negative Control	Saline
Positive Control (TEM*)	0.1 mcg/ml

The discussion of this test is contained in the technical discussion.

D. Methods

The protocols employed are explained in Appendices C and D.

E. Summary

1. Host-Mediated Assay

This compound caused no significant increase in mutant or recombinant frequencies when tested in vitro and in vivo against Salmonella TA-1530 and G-46 and Saccharomyces D3, respectively.

2. Cytogenetics

a. In vivo

The compound produced no detectable significant aberration of the bone marrow metaphase chromosomes of rats when administered orally at the dosage levels employed in this study.

*Triethylene Melamine

+These two tests were performed at different time intervals.



BIONETICS

b. In vitro

The compound produced no significant aberration in the anaphase chromosomes of human tissue culture cells when tested at the dosage levels employed in this study.

3. Dominant Lethal

This compound was considered to be non-mutagenic in rats in the Dominant Lethal Assay when using the dosages employed in this study.

F. Results and Discussion

1. Toxicity Data - Test I

a. In vivo

Compound FDA 71-51 was suspended in 0.85% saline and administered to ten male rats by intubation. The average weight of the animals was 250 grams and each received a dose of 5000 mg/kg. All animals were found dead within 24 hours. Necropsy indicated a swollen stomach with reddened patches of the mucosa.

Dose levels of 100, 200, 500, 1000, 2000 and 3000 mg/kg were selected to determine an acute LD₅₀. The toxicity data is presented on the LD₅₀ reporting form using the Litchfield-Wilcoxon method.

The LD₅₀ was determined as 720 mg/kg. The LD₅ dose level was derived from the raw data LD₅₀ probit line (uncorrected). The LD₅₀ determined from both corrected probit line and uncorrected probit line was within confidence limits of each other. The acute doses used were LD₅ - 475.0 mg/kg, intermediate - 47.5 mg/kg and low - 4.75 mg/kg. The subacute dose levels used were the same as those for the acute. The data on the dose levels, numbers of animals and necropsy findings are presented in the toxicity data sheets.



b. In vitro

The compound was suspended in 0.85% saline at the concentrations listed above. It was introduced into tubes containing WI-38 cells in a logarithmic phase of growth. The cells were observed for cytopathic effect (CPE) and the presence of mitoses at 24 and 48 hours.

<u>Tube No.</u>	<u>No. of Cells</u>	<u>Conc. mcg/ml</u>	<u>CPE</u>	<u>Mitoses</u>
1	5 X 10 ⁵	1000	+	+
2	"	1000	+	+
3	"	500	-	+
4	"	500	-	+
5	"	200	-	+
6	"	200	-	+
7	"	100	-	+
8	"	100	-	+
9	"	10	-	+
10	"	10	-	+



Since an inhibition of mitoses was observed, a closer range of concentrations was employed, as follows.

1	5×10^5	1000	+	+
2	"	1000	+	+
3	"	800	-	+
4	"	800	-	+
5	"	600	-	+
6	"	600	-	+
7	"	500	-	+
8	"	500	-	+
9	"	400	-	+
10	"	400	-	+

The 800 mcg/ml concentration was used as the high level, 80 mcg/ml as the intermediate and 8 mcg/ml as the low.



c. TOXICITY DATA SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST I



TOXICITY DATA
COMPOUND FDA 71-51

Solvent: 0.85% saline

Dosage Form: Suspension

Animals: Male rats with an average body weight of 250 grams. All animals were observed for ten (10) days.

Range Finding:

	<u>Dose</u> <u>mg/kg</u>	<u># Dead</u> <u># Animals</u>	<u>Day of Death and Necropsy</u>
	5000	10/10	Day 1 (10): Swollen stomach; reddened mucosa with patches.
LD ₅₀ :	100	0/5	None
	200	0/5	None
	500	1/5	Day 3 (1): Swollen stomach; reddened mucosa with patches.
	1000	4/5	Day 1 (3) and Day 2 (1): Swollen stomach; reddened mucosa with patches.
	2000	5/5	Day 1 (5): Swollen stomach; reddened mucosa with patches.
	3000	5/5	Day 1 (5): Swollen stomach; reddened mucosa with patches.



LD₅₀ REPORTING FORM USING LITCHFIELD-WILCOXON METHOD

DOSE EFFECT CURVE FOR Compound FDA 71-51 Methocel

DOSE	PROPORTION	OBSERVED PERCENT	EXPECTED PERCENT	OBS-EXPT PERCENT	CONTRIB. TO (CHI) ²
100	0/5	0	0		
200	0/5	0	0		
500	1/5	20	18		
1000	4/5	80	82		
2000	5/5	100	100		

Total animals = 20

Total =

Number Doses, K = 4

(CHI)² = 0.32

Animals/Dose = 5

Degrees of Freedom, n=k-2= 2

(CHI)² for n of k-2 = 5.99

since 0.32 is less than 5.99,
therefore data not significantly
heterogeneous

LD₈₄ = 1050

LD₅₀ = 720

LD₁₆ = 500

$fLD_{50} = S \frac{2.77}{\sqrt{N!}} = \frac{1.45}{\sqrt{N!}} = \frac{2.77}{\sqrt{N!}} = \frac{1.45}{\sqrt{10}} = \frac{2.77}{\sqrt{10}} = (1.45)^{.876} = 1.39$


LD₅₀ × feD₅₀ = (720) (1.39) = 1000

$\frac{LD_{50}}{fLD_{50}} = \frac{(720)}{(1.39)} = 518$

LD₅₀ and 19/20 Confidence Limits = $P \{ 518 \leq LD_{50} \leq 1000 \} = .95$

Attached should be a plot of the dose-effect curve on log-probit paper.

#Animals Dead/ #Animals Used

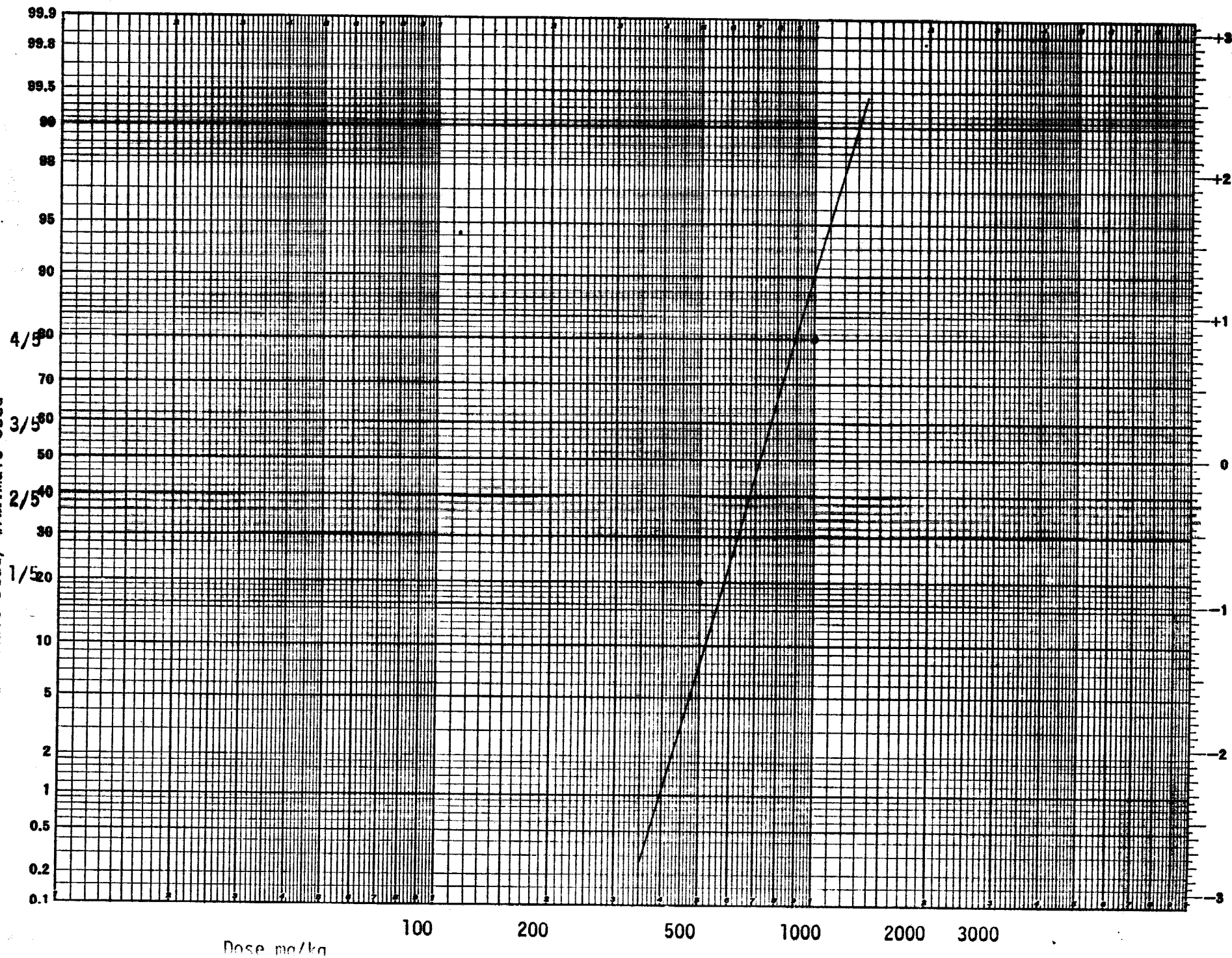

FDA Compound, 71-51 Methocel

2. Host-Mediated Assay - Test I

Compound FDA 71-51 caused no significant increase in mutant or recombinant frequencies when tested in the Host-Mediated Assay and in vitro tests against Salmonella TA-1530 and G-46 and Saccharomyces D3, respectively.



EVALUATION SHEET

Compound: FDA-71-51

Indicator Strain	In Vitro	In Vivo		
		Possible Low Recoveries	Controls	Other Comments
TA-1530	pos.	NC	NC OK	1. All doses negative 2. PC a little low
Acutes 11/6/72		PC		
		AL	PC OK	
S-acutes 11/9/72	neg.	AI		
		AH	SANC OK	
		SANC		
		SAL		
		SAI		
		SAH		
<hr/>				
G-46				
Acutes 11/27/72	pos.	NC	NC OK	1. All doses negative 2. PC a little low
		PC		
		AL	PC OK	
S-acutes 12/1/72	neg.	AI		
		AH	SANC OK	
		SANC		
		SAL		
		SAI		
		SAH		
<hr/>				
D3				
10/30/72 Acutes	pos.	NC	NC Low	1. Plating densities on Acute NC are too high reducing the recomb. frequencies too much.
		PC		
		AL	PC OK	
11/24/72 S-acutes	neg.	AI		2. Even though AI and AH doses appear positive against the acute NC they are not, against the S-acute negative control which is in fact an acute NC and more in line with expected NC frequencies.
		AH	SANC OK	
		SANC	SAPC OK	
		SAL		
		SAI		
		SAH		

Summary: The results for the bacterial tests are clearly negative and should be acceptable. The PC values for both TA-1530 and G-46 are somewhat low but the individual animal data all appears consistent. The results for D3 are more difficult to interpret because of the low value for the NC in the acute test. Because the S-acute negative control and other past negative control values have been around 5 recomb/10⁵ surv., I do not consider the acute doses to be positive. If this interpretation is not acceptable additional tests may be required to resolve the problem.

These should be acceptable results for D3

David Brunk

a. **HOST-MEDIATED ASSAY SUMMARY SHEETS**

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST I



HOST MEDIATED ASSAY

SUMMARY SHEET

COMPOUND: FDA 71-51

	SALMONELLA		SACCHAROMYCES D-3			
	TA1530		G-46			
	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MRF (X 10E-5)	MRT/MRC
ACUTE	11/4/72		11/27/72		10/30/72	
NC	.68		.60		2.86	
PC	7.54	11.09	16.23	27.05	60.17	21.04
AL	1.14	1.68	.90	1.50	4.61	1.61
AI	1.86	2.74	.91	1.52	6.94	2.43
LD5	1.94	2.85	1.02	1.70	10.32	3.61
SUBACUTE	11/4/72		12/1/72		11/24/72	
NC	.42		.59		4.94	
SL	.43	1.02	.93	1.58	5.86	1.19
SI	.88	2.10	.78	1.32	10.10	2.04
SLD5	.64	1.52	.77	1.31	8.93	1.81

IN VITRO	TA1530	G-46	% CONC	% SURVIVAL	R X 10E5
TCPD	-	-	10.0	100.0	1
NC	-	-	-	100.0	4
PC	+	+	0.5	50.2	387

STOP
SRU'S: .7

b. HOST-MEDIATED ASSAY DATA SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST I



BIONETICS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 6, 1972

ANIMAL NUMBER	RAW CFU X 10E7/0.0ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE. (C/B) X 10E-8
1	12.20	2.03	2.00	.93
2	16.60	2.77	1.00	.36
3	12.60	2.10	1.00	.48
4	34.00	5.67	2.00	.35
5	13.50	2.25	2.00	.69
6	13.70	2.62	2.00	.76
7	19.80	3.30	3.00	.91

NO. OF ANIMALS EQUALS 7
TOTAL CFU OUT OF RANGE EQUALS 0

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	2.96	1.86	.68
RANGE	3.63	2.00	.63
MAX	5.67	3.00	.98
MIN	2.03	1.00	.35

NO OUTLIERS

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: POSITIVE CONTROL - DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 6, 1972

ANIMAL NUMBER	RAW CFU X 10E7/0.6ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	27.10	4.52	30.00	6.64
2	26.70	4.45	18.00	4.04
3	16.10	2.68	20.00	7.45
4	22.10	3.68	28.00	7.60
5	20.20	3.37	48.00	14.26
6	21.80	3.53	29.00	7.98
7	27.40	4.57	22.00	4.82

NO. OF ANIMALS EQUALS 7

NO. OF CONTAMINATED EQUALS 3

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.84	27.86	7.54
RANGE	1.88	30.00	10.21
MAX	4.57	48.00	14.26
MIN	2.68	18.00	4.04

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.92	24.50	6.42
RANGE	1.88	12.00	3.94
MAX	4.57	30.00	7.98
MIN	2.68	18.00	4.04

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: LOW - 4.75 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 6, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	70.50	11.75	8.00	.68
2	38.40	6.40	11.00	1.72
3	24.50	4.08	9.00	2.20
4	35.50	5.92	1.00	.17
5	61.00	10.17	12.00	1.18
6	80.40	13.40	8.00	.60
7	16.70	2.78	4.00	1.44

NO. OF ANIMALS EQUALS 7
 NO. OF CONTAMINATED EQUALS 2
 TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.79	7.57	1.14
RANGE	10.62	11.00	2.04
MAX	13.40	12.00	2.20
MIN	2.78	1.00	.17

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: INTERMEDIATE - 47.5 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 6, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	37.50	6.25	10.00	1.60
2	23.40	3.90	4.00	1.03
3	48.10	8.02	5.00	.62
4	8.70	1.45	11.00	7.59
5	13.20	2.20	2.00	.91
6	60.70	10.12	4.00	.40
7	10.20	1.70	6.00	3.53
8	60.80	10.13	8.00	.79
9	19.00	3.17	1.00	.32

NO. OF ANIMALS EQUALS 9
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	5.21	5.67	1.86
RANGE	8.68	10.00	7.27
MAX	10.13	11.00	7.59
MIN	1.45	1.00	.32

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	5.69	5.00	1.15
RANGE	8.43	9.00	3.21
MAX	10.13	10.00	3.53
MIN	1.70	1.00	.32

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: LD5 - 475.0 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 6, 1972

ANIMAL NUMBER	A	B	C	D
	RAW CFU X 10E7/0.6ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	19.80	3.30	1.00	.30
2	18.30	3.05	14.00	4.59
3	21.00	3.50	2.00	.57
4	20.40	3.40	4.00	1.18
5	12.00	2.00	10.00	5.00
6	18.40	3.07	7.00	2.28
7	19.00	3.17	3.00	.95
8	19.50	3.25	2.00	.62

NO. OF ANIMALS EQUALS 8
 TOTAL CFU OUT OF RANGE EQUALS 1
 SAMPLES WITH ZERO MUTANTS EQUAL 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	3.09	5.38	1.94
RANGE	1.50	13.00	4.70
MAX	3.50	14.00	5.00
MIN	2.00	1.00	.30

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: NEGATIVE CONTROL - SALINE (SUBACUTE TRIALS)

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 9, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	31.00	5.17	2.00	.39
2	31.20	5.20	3.00	.58
3	64.70	10.78	4.00	.37
4	31.20	5.20	2.00	.38
5	42.30	7.05	4.00	.57
6	22.50	3.75	1.00	.27
7	49.50	8.25	6.00	.73
8	67.10	11.18	1.00	.09

NO. OF ANIMALS EQUALS 8
TOTAL CFU OUT OF RANGE EQUALS 1
SAMPLES WITH ZERO MUTANTS EQUAL 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.07	2.88	.42
RANGE	7.43	5.00	.64
MAX	11.18	6.00	.73
MIN	3.75	1.00	.09

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: LOW - 4.75 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: NOVEMBER 9, 1972

	A	B	C	D
ANIMAL NUMBER	RAW CFU X 10E7/0.6ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	40.20	6.70	1.00	.15
2	31.20	5.20	1.00	.19
3	32.20	5.37	2.00	.37
4	35.20	5.87	2.00	.34
5	31.90	5.32	4.00	.75
6	18.00	3.00	1.00	.33
7	33.50	5.58	2.00	.36
8	39.60	6.60	6.00	.91

NO. OF ANIMALS EQUALS 8

NO. OF CONTAMINATED EQUALS 1

TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	5.45	2.38	.43
RANGE	3.70	5.00	.76
MAX	6.70	6.00	.91
MIN	3.00	1.00	.15

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: INTERMEDIATE - 47.50 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE DATE STARTED: NOVEMBER 9, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	34.50	5.75	6.00	1.04
2	50.60	8.43	5.00	.59
3	45.20	7.53	6.00	.80
4	46.10	7.68	7.00	.91
5	48.90	8.15	8.00	.98
6	41.40	6.90	8.00	1.16
7	43.80	7.30	3.00	.41
8	31.90	5.32	6.00	1.13

NO. OF ANIMALS EQUALS 8
TOTAL CFU OUT OF RANGE EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.13	6.13	.88
RANGE	3.12	5.00	.75
MAX	8.43	8.00	1.16
MIN	5.32	3.00	.41

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: LD5 - 475.0 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: NOVEMBER 9, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	30.10	5.02	1.00	.20
2	33.10	5.52	1.00	.18
3	45.40	7.57	2.00	.26
4	31.90	5.32	3.00	.56
5	35.40	5.90	5.00	.85
6	48.10	8.02	7.00	.87
7	46.00	7.67	7.00	.91
8	31.50	5.25	7.00	1.33
9	32.00	5.33	3.00	.56

NO. OF ANIMALS EQUALS 9
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	6.18	4.00	.64
RANGE	3.00	6.00	1.15
MAX	8.02	7.00	1.33
MIN	5.02	1.00	.18

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	6.29	3.63	.55
RANGE	3.00	6.00	.73
MAX	8.02	7.00	.91
MIN	5.02	1.00	.18

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 27, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	42.10	7.02	3.00	.43
2	50.00	8.33	7.00	.84
3	35.20	5.87	4.00	.68
4	80.70	13.45	5.00	.37
5	81.20	13.53	7.00	.52
6	31.40	5.23	5.00	.96
7	91.00	15.17	4.00	.26
8	51.20	8.53	6.00	.70

NO. OF ANIMALS EQUALS 8
TOTAL CFU OUT OF RANGE EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	9.64	5.13	.60
RANGE	9.93	4.00	.69
MAX	15.17	7.00	.96
MIN	5.23	3.00	.26
NO OUTLIERS			

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: POSITIVE CONTROL - DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 27, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	54.20	9.03	107.00	11.84
2	46.60	7.77	72.00	9.27
3	37.70	6.28	210.00	33.42
4	61.90	10.32	130.00	12.60
5	41.20	6.87	97.00	14.13
6	57.20	9.53	121.00	12.69
7	52.20	8.70	142.00	16.32
8	62.20	10.37	117.00	11.29
9	53.40	8.90	218.00	24.49

NO. OF ANIMALS EQUALS 9

NO. OF CONTAMINATED EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.64	134.89	16.23
RANGE	4.08	146.00	24.15
MAX	10.37	218.00	33.42
MIN	6.28	72.00	9.27

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.94	125.50	14.08
RANGE	3.50	146.00	15.22
MAX	10.37	218.00	24.49
MIN	6.87	72.00	9.27

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LOW - 4.75 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 27, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	42.60	7.10	7.00	.99
2	37.00	6.17	6.00	.97
3	57.60	9.60	7.00	.73
4	36.70	6.12	9.00	1.47
5	65.10	10.85	8.00	.74
6	39.10	6.52	11.00	1.69
7	48.30	8.05	9.00	1.12
8	79.00	13.17	3.00	.23
9	41.20	6.87	1.00	.15

NO. OF ANIMALS EQUALS 9
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.27	6.78	.90
RANGE	7.05	10.00	1.54
MAX	13.17	11.00	1.69
MIN	6.12	1.00	.15

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: INTERMEDIATE - 47.50 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 27, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	32.10	5.35	6.00	1.12
2	46.00	7.67	13.00	1.70
3	42.30	7.05	5.00	.71
4	31.40	5.23	4.00	.76
5	54.50	9.08	8.00	.88
6	59.10	9.85	7.00	.71
7	43.00	7.17	5.00	.70
8	43.70	7.28	5.00	.69

NO. OF ANIMALS EQUALS 8
NO. OF CONTAMINATED EQUALS 1
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.34	6.63	.91
RANGE	4.62	9.00	1.01
MAX	9.85	13.00	1.70
MIN	5.23	4.00	.69

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.29	5.71	.80
RANGE	4.62	4.00	.43
MAX	9.85	8.00	1.12
MIN	5.23	4.00	.69

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LD5 - 475.0 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 27, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	30.60	5.10	9.00	1.76
2	30.20	5.03	5.00	.99
3	41.20	6.87	5.00	.73
4	58.10	9.68	5.00	.52
5	83.90	13.98	8.00	.57
6	37.30	6.22	5.00	.80
7	60.30	10.05	5.00	.50
8	56.20	9.37	9.00	.96
9	31.00	5.17	12.00	2.32

NO. OF ANIMALS EQUALS 9
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.94	7.00	1.02
RANGE	8.95	7.00	1.83
MAX	13.98	12.00	2.32
MIN	5.03	5.00	.50

NO. OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: NEGATIVE CONTROL - SALINE (SUBACUTES)

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: DECEMBER 1, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8	
1	31.80	5.30	7.00	1.32	*
2	57.00	9.50	4.00	.42	
3	68.70	11.45	7.00	.61	
4	91.40	15.23	8.00	.53	
5	81.40	13.57	3.00	.22	
6	83.30	13.88	9.00	.65	
7	137.00	22.83	6.00	.26	
8	60.60	10.10	7.00	.69	

NO. OF ANIMALS EQUALS 8
TOTAL CFU OUT OF RANGE EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	12.73	6.38	.59
RANGE	17.53	6.00	1.10
MAX	22.83	9.00	1.32
MIN	5.30	3.00	.22

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	13.80	6.29	.48
RANGE	13.33	6.00	.47
MAX	22.83	9.00	.69
MIN	9.50	3.00	.22

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LOW - 4.75 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: DECEMBER 1, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	55.00	9.17	5.00	.55
2	54.80	9.13	5.00	.55
3	42.40	7.07	9.00	1.27
4	63.20	10.53	3.00	.28
5	57.10	9.52	9.00	.95
6	38.00	6.33	13.00	2.05
7	71.70	11.95	8.00	.67
8	32.20	5.37	6.00	1.12

NO. OF ANIMALS EQUALS 8
NO. OF DEAD ANIMALS EQUALS 1
NO. OF CONTAMINATED EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.63	7.25	.93
RANGE	6.58	10.00	1.77
MAX	11.95	13.00	2.05
MIN	5.37	3.00	.28

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.96	6.43	.77
RANGE	6.58	6.00	.99
MAX	11.95	9.00	1.27
MIN	5.37	3.00	.28

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: INTERMEDIATE - 47.50 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: DECEMBER 1, 1972

	A	B	C	D
ANIMAL NUMBER	RAW CFU X 10E7/0.6ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	75.40	12.57	3.00	.24
2	50.10	8.35	5.00	.60
3	98.50	16.42	3.00	.18
4	54.10	9.02	7.00	.78
5	30.80	5.13	6.00	1.17
6	34.00	5.67	7.00	1.24
7	47.80	7.97	10.00	1.26

NO. OF ANIMALS EQUALS 7

NO. OF DEAD ANIMALS EQUALS 1

NO. OF CONTAMINATED EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	9.30	5.86	.78
RANGE	11.28	7.00	1.07
MAX	16.42	10.00	1.26
MIN	5.13	3.00	.18
NO. OUTLIERS			

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: LD5 - 475.00 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: DECEMBER 1, 1972

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	37.70	6.28	5.00	.80
2	59.50	9.92	4.00	.40
3	45.80	7.63	8.00	1.05
4	90.10	15.02	7.00	.47
5	69.00	11.50	11.00	.96
6	43.80	7.30	9.00	1.23
7	62.50	10.42	9.00	.86
8	61.30	10.22	4.00	.39

NO. OF ANIMALS EQUALS 8

NO. OF CONTAMINATED EQUALS 2

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	9.79	7.13	.77
RANGE	8.73	7.00	.64
MAX	15.02	11.00	1.23
MIN	6.28	4.00	.39

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: OCTOBER 30, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5	
1	1123.00	1.12	14.00	12.47	*
2	1300.00	1.30	2.00	1.54	
3	1700.00	1.70	0.	0.	
4	1400.00	1.40	2.00	1.43	
5	1100.00	1.10	3.00	2.73	
6	460.00	.46	0.	0.	
7	1800.00	1.80	6.00	3.33	
8	544.00	.54	0.	0.	
TOTAL		9.43	27.00		

NO. OF ANIMALS EQUALS 8
TOTAL SCREENED OUT OF RANGE EQUALS 2

MEAN C/MEAN B = 2.86

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.18	3.38	2.69
RANGE	1.34	14.00	12.47
MAX	1.80	14.00	12.47
MIN	.46	0.	0.

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 1.57

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.19	1.86	1.29
RANGE	1.34	6.00	3.33
MAX	1.80	6.00	3.33
MIN	.46	0.	0.

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: POSITIVE CONTROL - EMS - 350 MG/KG I.M.

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: OCTOBER 30, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	156.00	.16	24.00	153.85
2	784.00	.78	28.00	35.71
3	235.00	.24	17.00	72.34
4	204.00	.20	18.00	88.24
5	270.00	.27	14.00	51.85
6	395.00	.39	22.00	55.70
7	482.00	.48	29.00	60.17

TOTAL 2.53 152.00

NO. OF ANIMALS EQUALS 7
NO. OF CONTAMINATED EQUALS 1
TOTAL SCREENED OUT OF RANGE EQUALS 2

MEAN C/MEAN B = 60.17

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.36	21.71	73.98
RANGE	.63	15.00	118.13
MAX	.78	29.00	153.85
MIN	.16	14.00	35.71

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 54.01

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.39	21.33	60.67
RANGE	.58	15.00	52.52
MAX	.78	29.00	88.24
MIN	.20	14.00	35.71

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LOW - 4.75 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: OCTOBER 30, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	312.00	.31	8.00	25.64 *
2	1550.00	1.55	4.00	2.58
3	642.00	.64	5.00	7.79
4	1340.00	1.34	4.00	2.99
5	583.00	.58	4.00	6.86
6	1420.00	1.42	6.00	4.23
7	1560.00	1.56	5.00	3.21
8	1270.00	1.27	4.00	3.15

TOTAL 8.68 40.00

NO. OF ANIMALS EQUALS 8
NO. OF CONTAMINATED EQUALS 1
TOTAL SCREENED OUT OF RANGE EQUALS 1

MEAN C/MEAN B = 4.61

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.08	5.00	7.05
RANGE	1.25	4.00	23.06
MAX	1.56	8.00	25.64
MIN	.31	4.00	2.58

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 3.83

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.20	4.57	4.40
RANGE	.98	2.00	5.21
MAX	1.56	6.00	7.79
MIN	.58	4.00	2.58

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: INTERMEDIATE - 47.50 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: OCTOBER 30, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	460.00	.46	10.00	21.74
2	529.00	.53	8.00	15.12
3	1470.00	1.47	16.00	10.88
4	572.00	.57	5.00	8.74
5	1160.00	1.16	14.00	12.07
6	1860.00	1.86	3.00	1.61
7	1170.00	1.17	4.00	3.42
8	1300.00	1.30	4.00	3.08
9	1130.00	1.13	3.00	2.65

TOTAL

9.65

67.00

NO. OF ANIMALS EQUALS 9

TOTAL SCREENED OUT OF RANGE EQUALS 1

MEAN C/MEAN B = 6.94

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.07	7.44	8.81
RANGE	1.40	13.00	20.13
MAX	1.86	16.00	21.74
MIN	.46	3.00	1.61

NO OUTLIERS

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LD5 - 475.00 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: OCTOBER 30, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	1240.00	1.24	7.00	5.65
2	1780.00	1.78	8.00	4.49
3	747.00	.75	15.00	20.08
4	900.00	.90	3.00	3.33
5	840.00	.84	16.00	19.05
6	1070.00	1.07	22.00	20.56
7	1930.00	1.93	7.00	3.63
8	410.00	.41	14.00	34.15
TOTAL		8.92	92.00	

NO. OF ANIMALS EQUALS 8
TOTAL SCREENED OUT OF RANGE EQUALS 2

MEAN C/MEAN B = 10.32

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.11	11.50	13.87
RANGE	1.52	19.00	30.81
MAX	1.93	22.00	34.15
MIN	.41	3.00	3.33

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 9.17

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.22	11.14	10.97
RANGE	1.18	19.00	17.23
MAX	1.93	22.00	20.56
MIN	.75	3.00	3.33

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: NEGATIVE CONTROL - SALINE (SUBACUTES)

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 24, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	486.00	.49	3.00	6.17
2	557.00	.56	4.00	7.18
3	806.00	.81	0.	0.
4	333.00	.33	5.00	15.02
5	427.00	.43	2.00	4.68
6	307.00	.31	1.00	3.26
7	722.00	.72	3.00	4.16
8	814.00	.81	4.00	4.91
TOTAL		4.45	22.00	

NO. OF ANIMALS EQUALS 8
TOTAL SCREENED OUT OF RANGE EQUALS 2

MEAN C/MEAN B = 4.94

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.56	2.75	5.67
RANGE	.51	5.00	15.02
MAX	.81	5.00	15.02
MIN	.31	0.	0.

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 4.13

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.59	2.43	4.34
RANGE	.51	4.00	7.18
MAX	.81	4.00	7.18
MIN	.31	0.	0.

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: POSITIVE CONTROL - EMS - 350 MG/KG IM (SUBACUTES)

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: NOVEMBER 24, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	914.00	.91	32.00	35.01
2	595.00	.59	28.00	47.06
3	411.00	.41	40.00	97.32
4	910.00	.91	39.00	42.86
5	364.00	.36	27.00	74.18
6	542.00	.54	20.00	36.90
7	521.00	.52	35.00	67.18

TOTAL

4.26

221.00

NO. OF ANIMALS EQUALS 7

TOTAL SCREENED OUT OF RANGE EQUALS 3

MEAN C/MEAN B = 51.91

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.61	31.57	57.22
RANGE	.55	20.00	62.31
MAX	.91	40.00	97.32
MIN	.36	20.00	35.01

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 47.06

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.64	30.17	50.53
RANGE	.55	19.00	39.16
MAX	.91	39.00	74.18
MIN	.36	20.00	35.01

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LOW - 4.75 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: NOVEMBER 24, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	554.00	.55	4.00	7.22
2	1360.00	1.36	9.00	6.62
3	728.00	.73	4.00	5.49
4	870.00	.87	5.00	5.75
5	385.00	.38	0.	0.
6	1950.00	1.95	6.00	3.08
7	924.00	.92	7.00	7.58
8	911.00	.91	10.00	10.98
TOTAL		7.68	45.00	

NO. OF ANIMALS EQUALS 8
TOTAL SCREENED OUT OF RANGE EQUALS 2

MEAN C/MEAN B = 5.86

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.96	5.63	5.84
RANGE	1.57	10.00	10.98
MAX	1.95	10.00	10.98
MIN	.38	0.	0.

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 5.17

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.97	5.00	5.10
RANGE	1.57	9.00	7.58
MAX	1.95	9.00	7.58
MIN	.38	0.	0.

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: INTERMEDIATE - 47.5 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: NOVEMBER 24, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	488.00	.49	5.00	10.25
2	308.00	.31	1.00	3.25
3	507.00	.51	2.00	3.94
4	387.00	.39	0.	0.
5	355.00	.35	2.00	5.63
6	766.00	.77	4.00	5.22
7	622.00	.62	6.00	9.65
8	496.00	.50	12.00	24.19
9	625.00	.62	14.00	22.40

TOTAL 4.55 46.00

NO. OF ANIMALS EQUALS 9

NO. OF CONTAMINATED EQUALS 1

MEAN C/MEAN B = 10.10

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.51	5.11	9.39
RANGE	.46	14.00	24.19
MAX	.77	14.00	24.19
MIN	.31	0.	0.

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: LD5 - 475.00 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: NOVEMBER 24, 1972

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5	
1	462.00	.46	10.00	21.65	*
2	740.00	.74	7.00	9.46	
3	512.00	.51	0.	0.	
4	610.00	.61	6.00	9.84	
5	526.00	.53	6.00	11.41	
6	367.00	.37	2.00	5.45	
7	814.00	.81	5.00	6.14	
TOTAL		4.03	36.00		

NO. OF ANIMALS EQUALS 7
TOTAL SCREENED OUT OF RANGE EQUALS 3

MEAN C/MEAN B = 8.93

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.58	5.14	9.13
RANGE	.45	10.00	21.65
MAX	.81	10.00	21.65
MIN	.37	0.	0.

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 7.28

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.59	4.33	7.05
RANGE	.45	7.00	11.41
MAX	.81	7.00	11.41
MIN	.37	0.	0.

STOP

3. Toxicity Data - Test II

a. Acute

Compound FDA 71-51, Methocel, was prepared by adding 5000 mg/kg of body weight to a pellet of feed and allowing a group of six male rats to eat it as a single dose. The average body weight of the six rats was 208 grams.

No signs of toxicity or abnormal behavior was observed in the 7-day observation period. No deaths occurred. At termination all animals were killed and on necropsy no gross findings were observed.

The acute oral LD₅₀ in rats for compound FDA 71-51 is considered to be greater than 5000 mg/kg.

b. Subacute

Compound FDA 71-51, Methocel, was incorporated in a small quantity of food as in the acute study. The test substance was administered to a group of ten male rats (average body weight 176 grams), daily for five days at a dosage of 5000 mg/kg. No signs of toxicity or abnormal behavior were observed in the five-day period of compound administration or in the observation period which followed. The total period of observation was 14 days when the animals were terminated and gross necropsies performed. No abnormal gross findings were observed. The 14-day subacute oral LD₅₀ for compound FDA 71-51 is considered to be greater than 5000 mg/kg.



c. TOXICITY DATA SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST II



TOXICITY DATA
COMPOUND FDA 71-51

Solvent: None

Dosage Form: Food pellets

Animals: Acute - Male rats with an average body weight of 208 grams.

All animals were observed for seven days.

Subacute - Male rats with an average body weight of 176 grams.

All animals were observed for 14 days.


Oral Study: Acute - No deaths (0/6) following a single oral dose and 7-day observation period. No signs of toxicity and no gross necropsy findings. The acute oral LD₅₀ is greater than 5000 mg/kg.

Subacute - No deaths (0/10) following daily dose for five days and a total of 14 days observation period. No signs of toxicity and no gross necropsy findings. The five-day subacute oral LD₅₀ is greater than 5000 mg/kg.



4. Host-Mediated Assay - Test II

Compound FDA 71-51, Methocel, was not genetically active in any of the assays at a dose level of 5000 mg/kg.

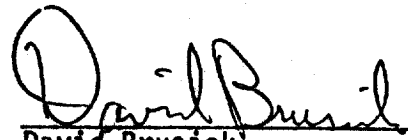
David Brusick



a. **HOST-MEDIATED ASSAY SUMMARY SHEETS**

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST II



BIONETICS

HOST MEDIATED ASSAY

SUMMARY SHEET

COMPOUND: FDA 71-51

	SALMONELLA		SACCHAROMYCES D-3			
	TA1530	G-46				
	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MRF (X 10E-5)	MRT/MRC
ACUTE						
NC	3.47		1.19		11.40	
PC	73.93	21.31	113.11	95.05	85.29	7.48
AL	0.	0.	0.	0.	0.	0.
AI	0.	0.	0.	0.	0.	0.
AH	3.69	1.06	.85	.71	18.27	1.60
SUBACUTE						
NC	1.00		1.00		1.00	
SL	0.	0.	0.	0.	0.	0.
SI	0.	0.	0.	0.	0.	0.
SH	0.	0.	0.	0.	0.	0.

IN VITRO

TA1530

G-46

% CONC D-3
% SURVIVAL R X 10E5

NC
PC

STOP
SRU'S: .4

HOST MEDIATED ASSAY

SUMMARY SHEET

COMPOUND: FDA 71-51

	SALMONELLA				SACCHAROMYCES D-3	
	TA1530		G-46			
	MMF (X 10E-8)	MFT/MFC	MMF (X 10E-8)	MFT/MFC	MRF (X 10E-5)	MRT/MRC
ACUTE						
NC	2.04		1.38		10.34	
PC	69.00	33.82	306.21	221.89	99.56	9.63
AL	0.	0.	0.	0.	0.	0.
AI	0.	0.	0.	0.	0.	0.
AH	0.	0.	0.	0.	0.	0.
SUBACUTE						
NC	2.04		1.38		10.34	
SL	0.	0.	0.	0.	0.	0.
SI	0.	0.	0.	0.	0.	0.
SH	2.76	1.35	.82	.59	11.63	1.12

IN VITRO TA1530 G-46 D-3
 % CONC % SURVIVAL R X 10E5

NC
 PC

STOP
 SRU'S: .4

b. HOST-MEDIATED ASSAY DATA SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST II



BIONETICS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 6, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	46.50	7.75	22.00	2.84
2	43.00	7.17	27.00	3.77
3	58.60	9.77	23.00	2.35
4	37.20	6.20	23.00	3.71
5	41.70	6.95	26.00	3.74
6	44.90	7.48	15.00	2.00
7	38.50	6.42	28.00	4.36
8	69.80	11.63	13.00	1.12
9	62.80	10.47	28.00	2.68
10	31.90	5.32	43.00	8.09

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.92	24.80	3.47
RANGE	6.32	30.00	6.97
MAX	11.63	43.00	8.09
MIN	5.32	13.00	1.12

SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.20	22.78	2.95
RANGE	5.43	15.00	3.25
MAX	11.63	28.00	4.36
MIN	6.20	13.00	1.12

STJP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1538

DOSE LEVEL: POSITIVE CONTROL - DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 6, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	50.30	8.38	1050.00	125.25
2	73.10	12.18	508.00	41.70
3	52.10	8.68	364.00	41.92
4	37.90	6.32	553.00	87.54
5	50.90	8.48	417.00	49.15
6	34.90	5.82	997.00	171.40
7	42.60	7.10	372.00	52.39
8	39.20	6.53	303.00	46.38
9	35.80	5.97	296.00	49.61

NO. OF ANIMALS EQUALS 9
NO. OF CONTAMINATED EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.72	540.00	73.93
RANGE	6.37	754.00	129.71
MAX	12.18	1050.00	171.40
MIN	5.82	296.00	41.70

* SUMMARY WITH OUT IERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.96	482.88	61.74
RANGE	6.22	754.00	83.55
MAX	12.18	1050.00	125.25
MIN	5.97	296.00	41.70

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: HIGH - 5000 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 6, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	47.60	7.93	45.00	5.67
2	67.20	11.20	29.00	2.59
3	69.00	11.50	25.00	2.17
4	37.70	6.28	27.00	4.30
5	66.50	11.08	34.00	3.07
6	40.30	6.72	27.00	4.02
7	67.30	11.22	27.00	2.41
8	57.70	9.62	45.00	4.68
9	46.00	7.67	33.00	4.30

NO. OF ANIMALS EQUALS 9

NO. OF DEAD ANIMALS EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	9.25	32.44	3.69
RANGE	5.22	20.00	3.50
MAX	11.50	45.00	5.67
MIN	6.28	25.00	2.17

NO OUTLIERS

TOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 29, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	54.80	9.13	13.00	1.42
2	74.50	12.42	18.00	1.45
3	31.30	5.22	16.00	3.07
4	53.80	8.97	8.00	.89
5	42.80	7.13	27.00	3.78
6	31.90	5.32	13.00	2.45
7	65.00	10.83	20.00	1.85
8	51.70	8.62	10.00	1.16
9	57.80	9.63	20.00	2.08
10	67.80	11.30	26.00	2.30

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.86	17.10	2.04
RANGE	7.20	19.00	2.89
MAX	12.42	27.00	3.78
MIN	5.22	8.00	.89

NO OUTLIERS

STOP

RU'S: .6

SWITCH INS: SL252

SAL

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: POSITIVE CONTROL - DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 29, 1974

	A	B	C	D
ANIMAL NUMBER	RAW CFU X 10E7/0.6ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	34.60	5.77	317.00	54.97
2	57.10	9.52	1062.00	111.59
3	48.30	8.05	494.00	61.37
4	32.80	5.47	518.00	94.75
5	54.50	9.08	658.00	72.44
6	50.10	8.35	207.00	24.79
7	49.50	8.25	349.00	42.30
8	37.60	6.27	562.00	89.68
9	49.20	8.20	390.00	47.56
10	45.80	7.63	691.00	90.52

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	7.66	524.80	69.00
RANGE	4.05	855.00	86.80
MAX	9.52	1062.00	111.59
MIN	5.47	207.00	24.79

NO OUTLIERS

STOP

RU'S: .6

ISWITCH IN\$: SL251

SAL

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA TA1530

DOSE LEVEL: HIGH - 5000 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 29, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	66.30	11.05	22.00	1.99
2	34.50	5.75	15.00	2.61
3	96.20	16.03	14.00	.87
4	62.30	10.38	24.00	2.31
5	37.50	6.25	23.00	3.68
6	32.50	5.42	21.00	3.88
7	39.40	6.57	26.00	3.96

NO. OF ANIMALS EQUALS 7

NO. OF CONTAMINATED EQUALS 3

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.78	20.71	2.76
RANGE	10.62	12.00	3.09
MAX	16.03	26.00	3.96
MIN	5.42	14.00	.87

NO OUTLIERS

STOP

SRU'S: .5

SWITCH IN: SL255

ISAL

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: JANUARY 16, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	53.10	8.85	13.00	1.47
2	46.60	7.77	12.00	1.55
3	61.00	10.17	13.00	1.28
4	36.50	6.08	7.00	1.15
5	50.30	8.38	11.00	1.31
6	61.90	10.32	10.00	.97
7	57.10	9.52	6.00	.63

NO. OF ANIMALS EQUALS 7
TOTAL CFU OUT OF RANGE EQUALS 3

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.73	10.29	1.19
RANGE	4.23	7.00	.91
MAX	10.32	13.00	1.55
MIN	6.08	6.00	.63

SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.59	11.00	1.29
RANGE	4.23	6.00	.58
MAX	10.32	13.00	1.55
MIN	6.08	7.00	.97

TOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: POSITIVE CONTROL - DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: JANUARY 16, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	80.20	13.37	869.00	65.01
2	67.70	11.28	1281.00	113.53
3	128.10	21.35	1209.00	56.63
4	69.80	11.63	1296.00	111.40
5	60.10	10.02	1632.00	162.93
6	56.20	9.37	909.00	97.04
7	68.50	11.42	1129.00	98.89
8	58.00	9.67	2107.00	217.96
9	69.60	11.60	1097.00	94.57

NO. OF ANIMALS EQUALS 9
TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	12.19	1281.00	113.11
RANGE	11.98	1238.00	161.33
MAX	21.35	2107.00	217.96
MIN	9.37	869.00	56.63

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	12.50	1177.75	100.00
RANGE	11.98	763.00	106.30
MAX	21.35	1632.00	162.93
MIN	9.37	869.00	56.63

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: HIGH - 5000 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: JANUARY 16, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.6ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	144.20	24.03	22.00	.92
2	119.00	19.83	15.00	.76
3	107.30	17.88	19.00	1.06
4	253.10	42.18	9.00	.21
5	109.20	18.20	19.00	1.04
6	160.40	26.73	9.00	.34
7	86.20	14.37	18.00	1.25
8	102.80	17.13	11.00	.64
9	107.60	17.93	33.00	1.84
10	142.20	23.70	11.00	.46

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	22.20	16.60	.85
RANGE	27.82	24.00	1.63
MAX	42.18	33.00	1.84
MIN	14.37	9.00	.21

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	22.67	14.78	.74
RANGE	27.82	13.00	1.04
MAX	42.18	22.00	1.25
MIN	14.37	9.00	.21

SEOP
FI 027 E ILLS I-FORMAT INPUT
M. 777 INS
10022 I LIST ITEM AT 100156. INPUT FIELD IS <R>
NA KBACK SEQUENCE
FI 027 ENTRY LINE ADDRESS CALLER LINE ADDRESS
SFIO SEIO 007817 55. 001824

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: JANUARY 18, 1974

	A	B	C	D
ANIMAL NUMBER	RAW CFU X 10E7/0.6ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	32.10	5.35	9.00	1.68
2	53.20	8.87	22.00	2.48
3	115.60	19.27	7.00	.36
4	40.20	6.70	7.00	1.04
5	45.20	7.53	11.00	1.46
6	39.40	6.57	13.00	1.98
7	39.30	6.55	8.00	1.22
8	50.10	8.35	7.00	.84

NO. OF ANIMALS EQUALS 8

NO. OF CONTAMINATED EQUALS 1

TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	8.65	10.50	1.38
RANGE	13.92	15.00	2.12
MAX	19.27	22.00	2.48
MIN	5.35	7.00	.36

NO OUTLIERS

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: POSITIVE CONTROL - DMN - 100 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: JANUARY 18, 1974

ANIMAL NUMBER	A RAW CFU X 10E7/0.5ML	B TOTAL CFU X 10E8/1.0ML	C TOTAL NO. MUTANTS X 10E0/1.0ML	D MUTATION FRE (C/B) X 10E-8
1	42.40	7.07	4376.00	619.23
2	30.20	5.03	1611.00	320.06
3	46.30	7.72	2907.00	376.71
4	36.20	6.03	2073.00	343.58
5	40.70	6.78	1404.00	206.97
6	46.30	7.72	1364.00	176.76
7	36.60	6.10	1376.00	225.57
8	53.20	8.87	1968.00	221.95
9	45.50	7.58	2270.00	299.33
10	39.40	6.57	1786.00	271.97

NO. OF ANIMALS EQUALS 10

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	6.95	2113.50	306.21
RANGE	3.83	3012.00	442.48
MAX	8.87	4376.00	619.23
MIN.	5.03	1364.00	176.76

* SUMMARY WITH OUT LIERS REMOVED

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	6.93	1862.11	271.43
RANGE	3.83	1543.00	199.95
MAX	8.87	2907.00	376.71
MIN	5.03	1364.00	176.76

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SALMONELLA G-46

DOSE LEVEL: HIGH - 5000 MG/KG

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: JANUARY 18, 1974

	A	B	C	D
ANIMAL NUMBER	RAW CFU X 10E7/0.5ML	TOTAL CFU X 10E8/1.0ML	TOTAL NO. MUTANTS X 10E0/1.0ML	MUTATION FRE (C/B) X 10E-8
1	40.80	6.80	9.00	1.32
2	25.80	4.30	6.00	1.40
3	74.80	12.47	3.00	.24
4	64.70	10.78	8.00	.74
5	121.20	20.20	14.00	.60
6	81.80	13.63	12.00	.88
7	109.00	18.17	9.00	.50

NO. OF ANIMALS EQUALS 7
 NO. OF DEAD ANIMALS EQUALS 1
 NO. OF CONTAMINATED EQUALS 1
 TOTAL CFU OUT OF RANGE EQUALS 1

	COL. B (X 10E8)	COL. C (X 10E0)	COL. D (X 10E-8)
MEAN	12.34	8.71	.82
RANGE	15.90	11.00	1.15
MAX	20.20	14.00	1.40
MIN	4.30	3.00	.24

NO OUTLIERS

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: NEGATIVE CONTROL - SALINE

TREATMENT: IN VIVO, ORAL, SUBACUTE

DATE STARTED: MARCH 13, 1974

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5	
1	519.00	.52	10.00	19.27	*
2	762.00	.76	7.00	9.19	
3	815.00	.81	11.00	13.50	
4	786.00	.79	11.00	13.99	
5	965.00	.96	9.00	9.33	
6	580.00	.58	9.00	15.52	
7	743.00	.74	6.00	8.08	
8	827.00	.83	8.00	9.67	
9	890.00	.89	7.00	7.87	
10	392.00	.39	5.00	12.76	
TOTAL		7.28	83.00		

NO. OF ANIMALS EQUALS 10

MEAN C/MEAN B = 11.40

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.73	8.30	11.92
RANGE	.57	6.00	11.40
MAX	.96	11.00	19.27
MIN	.39	5.00	7.87

* SUMMARY WITH OUTLIERS REMOVED

MEAN C/MEAN B = 10.80

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.75	8.11	11.10
RANGE	.57	6.00	7.65
MAX	.96	11.00	15.52
MIN	.39	5.00	7.87

STOP
RU'S:.7

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: POSITIVE CONTROL - EMS - 35 MG/KG I.M.

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARC 13, 1974

ANIMAL NUMBER	RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB. CFU SCREENED X 10E-5
1	794.00	.79	107.00	134.76
2	699.00	.70	66.00	94.42
3	635.00	.63	75.00	118.11
4	811.00	.81	47.00	57.95
5	804.00	.80	51.00	63.43
6	736.00	.74	74.00	100.54
7	756.00	.76	63.00	83.33
8	617.00	.62	23.00	37.28
9	1136.00	1.11	77.00	67.78
10	879.00	.88	88.00	100.11
TOTAL		7.87	671.00	

NO. OF ANIMALS EQUALS 10

MEAN C/MEAN B = 85.99

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	.79	67.1	87.77
RANGE	.52	84.00	97.48
MAX	1.14	107.00	134.78
MIN	.62	23.00	37.28

NO OUTLIERS

TOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND: FDA 71-51

ORGANISM: SACCHAROMYCES D-3

DOSE LEVEL: HIGH - 5000 MG/KG

TREATMENT: IN VIVO, ORAL, ACUTE

DATE STARTED: MARCH 13, 1974

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	1182.00	1.18	29.00	24.53
2	1939.00	1.94	27.00	13.92
3	741.00	.74	22.00	29.69
4	1099.00	1.10	23.00	20.93
5	995.00	.99	15.00	15.08
6	1643.00	1.64	17.00	10.35
7	1160.00	1.16	24.00	20.69
8	720.00	.72	13.00	18.06
9	1141.00	1.14	24.00	21.03

TOTAL

10.62

194.00

NO. OF ANIMALS EQUALS

TOTAL SCREENED OUT OF RANGE EQUALS

MEAN C/MEAN B =

18.27

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.18	21.56	19.36
RANGE	1.22	16.00	19.34
MAX	1.94	29.00	29.69
MIN	.72	13.00	10.35

NO OUTLIERS

STOP

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND : FDA 71-51

ORGANISM : SACCHAROMYCES D-3

DOSE LEVEL : NEGATIVE CONTROL - SALINE

TREATMENT : IN VIVO, ORAL, SUBACUTE

DATE STARTED : 7/15/74

ANIMAL NUMBER	A RAW CFU X 10E5/1.0ML	B TOTAL CFU SCREENED X 10E5/1.0ML	C TOTAL RECOMBINANTS /1.0ML	D RECOMB/CFU SCREENED X 10E-5
1	1979.00	1.98	26.00	13.14
2	1612.00	1.61	16.00	9.93
3	1768.00	1.77	14.00	7.92
4	1558.00	1.56	17.00	10.91
5	1574.00	1.57	11.00	6.99
6	1547.00	1.55	15.00	9.70
7	1567.00	1.57	21.00	13.40
8	1726.00	1.73	17.00	9.85
9	1248.00	1.25	9.00	7.21
10	1477.00	1.48	20.00	13.54

NO. OF ANIMALS EQUALS 10

NO. OF DEAD EQUALS 0

NO. OF CONTAMINATED EQUALS 0

NO. OF CFU OUT OF RANGE EQUALS 0

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.61	16.60	10.26
RANGE	0.73	17.00	6.55
MAX	1.98	26.00	13.54
MIN	1.25	9.00	6.99

MEAN C/MEAN B = 10.34

NO OUTLIERS

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND : FDA 71-51

ORGANISM : SACCHAROMYCES D-3

DOSE LEVEL : POSITIVE CONTROL - EMS-350 MG/KG

TREATMENT : IN VIVO, ORAL, ACUTE

DATE STARTED : 7/15/74

ANIMAL NUMBER	A	B	C	D	
	RAW CFU X 10E5/1.0ML	TOTAL CFU SCREENED X 10E5/1.0ML	TOTAL RECOMBINANTS /1.0ML	RECOMB/CFU SCREENED X 10E-5	
1	433.00	0.43	114.00	263.28	*
2	985.00	0.98	159.00	161.42	
3	1172.00	1.17	103.00	87.88	
4	1193.00	1.19	150.00	125.73	
5	1256.00	1.26	70.00	55.73	
6	1385.00	1.38	74.00	53.43	
7	563.00	0.56	86.00	152.75	
8	1621.00	1.62	101.00	62.31	

2 ANIMALS NOT ACCOUNTED FOR. VALUES COMPUTED FOR 8 ANIMALS

NO. OF DEAD EQUALS 0

NO. OF CONTAMINATED EQUALS 2

NO. OF CFU OUT OF RANGE EQUALS 0

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.08	107.13	120.32
RANGE	1.19	89.00	209.85
MAX	1.62	159.00	263.28
MIN	0.43	70.00	53.43

MEAN C/MEAN B = 99.56

* SUMMARY WITH OUTLIERS REMOVED

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.17	106.14	99.89
RANGE	1.06	89.00	107.99
MAX	1.62	159.00	161.42
MIN	0.56	70.00	53.43

MEAN C/MEAN B = 90.89

HOST MEDIATED ASSAY REPORT SHEET

COMPOUND : FDA 71-51 *

ORGANISM : SACCHAROMYCES D-3

DOSE LEVEL : HIGH - 5000. MG/KG

TREATMENT : IN VIVO, ORAL, SUBACUTE

DATE STARTED : 7/15/74

ANIMAL NUMBER	A	B	C	D
	RAW CFU X 10E5/1.0ML	TOTAL CFU SCREENED X 10E5/1.0ML	TOTAL RECOMBINANTS /1.0ML	RECOMB/CFU SCREENED X 10E-5
1	832.00	0.83	16.00	19.23
2	729.00	0.73	11.00	15.09
3	681.00	0.68	18.00	26.43
4	2773.00	2.77	23.00	8.29
5	512.00	0.51	6.00	11.72
6	1737.00	1.74	2.00	1.15
7	1940.00	1.94	21.00	10.82
8	920.00	0.92	23.00	25.00
9	1539.00	1.54	17.00	11.05
10	1234.00	1.23	13.00	10.53

NO. OF ANIMALS EQUALS 10

NO. OF DEAD EQUALS 0

NO. OF CONTAMINATED EQUALS 0

NO. OF CFU OUT OF RANGE EQUALS 0

	COL. B (X 10E5)	COL. C (X 10E0)	COL. D (X 10E-5)
MEAN	1.29	15.00	13.93
RANGE	2.26	21.00	25.28
MAX	2.77	23.00	26.43
MIN	0.51	2.00	1.15

MEAN C/MEAN B = 11.63

NO OUTLIERS

5. Cytogenetics - Test I

a. In vivo

(1) Acute study

There were no aberrations observed in the negative control group or in the test compound dosage level groups. A depression of the mitotic indices was observed in all three dosage level groups of the test compound at 48 hours. The expected severe chromosomal damage due to the positive control compound was observed.

(2) Subacute study

No aberrations were observed in the negative control group or in the three dosage level groups of the test compound. There was a marked depression of the mitotic indices of all three dosage level groups of the test compound. This was especially noticed in the LD₅ dosage level.

b. In vitro

The negative control group contained 2 cells with bridges, one of which contained a fragment. The low, medium and high levels of the test compound all contained 2% of the cells with bridges. The positive control group was within normal limits for the positive control compound effect.



c. CYTOGENETIC SUMMARY SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST I



METHOCEL
FDA 71-51
ACUTE STUDY
METAPHASE SUMMARY SHEET
TEST I

<u>Compound</u>	<u>Dosage (mg/kg)</u>	<u>Time*</u>	<u>No. of Animals</u>	<u>No. of Cells</u>	<u>Mitotic Index %***</u>	<u>% Cells with Breaks</u>	<u>% Cells with Reunion</u>	<u>% Cells Other Aber.**</u>	<u>% Cells with Aber.++</u>
Negative Control	Saline	6	3	150	6	0	0	0	0
		24	3	150	8	0	0	0	0
		48	3	150	6	0	0	0	0
Low Level	4.75	6	5	250	5	0	0	0	0
		24	5	250	6	0	0	0	0
		48	5	206	2	0	0	0	0
Intermediate Level	47.5	6	5	250	4	0	0	0	0
		24	5	250	4	0	0	0	0
		48	5	157	2	0	0	0	0
High Level	475	6	5	250	5	0	0	0	0
		24	5	250	4	0	0	0	0
		48	5	87	2	0	0	0	0
Positive Control TEM	0.3	48	5	232	4	3	21	14(a), 7.3(f), 1.3(pp)	35.3

* Time of kill after injection (hours).

** Cells that have polyploidy (P), pulverization (pp), or greater than 10 aberrations (a).

*** Percent of cells in mitosis: 500 cells observed/animal.

++ Duplicate aberrations in a single cell will cause this to be a % less than a summation of the % aberration seen.

METHOCEL
FDA 71-51
SUBACUTE STUDY
METAPHASE SUMMARY SHEET
TEST I

<u>Compound</u>	<u>Dosage (mg/kg)*</u>	<u>No. of Animals</u>	<u>No. of Cells</u>	<u>Mitotic Index %***</u>	<u>% Cells with Breaks</u>	<u>% Cells with Reunion</u>	<u>% Cells Other Aber.**</u>	<u>% Cells with Aber.</u>
Negative Control	Saline	3	150	4	0	0	0	0
Low Level	4.75	5	250	3	0	0	0	0
Intermediate Level	47.5	5	250	2	0	0	0	0
High Level	475	5	47	2	0	0	0	0

* Dosage 1X/day X 5 days.

** Cells that have polyploidy (P), pulverization (pp), fragments (f) or greater than 10 aberrations (a).

*** Percent of cells in mitosis: 500 cells observed/animal.

METHOCEL
FDA 71-51
ANAPHASE SUMMARY SHEET
TEST I

<u>Compound</u>	<u>Dosage (mcg/ml)</u>	<u>Mitotic Index**</u>	<u>No. of Cells</u>	<u>% Cells with Acentric Frag.</u>	<u>% Cells with Bridges</u>	<u>% Multipolar Cells</u>	<u>% Cells Other Aber.*</u>	<u>% Cells with Aber.++</u>
Low Level	8.0	1	100	0	2	0	0	2
Medium Level	80.0	1	100	0	2	0	0	2
High Level	800	1	100	0	2	0	0	2
Negative Control	Saline	1	100	1	2	0	0	2
Positive Control (TEM)	0.1	1	100	3	12	0	0	15

* Cells that have polyploidy (P), pulverization (pp), or greater than 10 aberrations (a).

** Percent of cells in mitosis: 200 cells observed/dose level.

++ Duplicate aberrations in a single cell will cause this to be a % less than a summation of the % aberration seen.

6. Cytogenetics - Test II

Compound FDA 71-51, Methocel, was administered to male rats with an average body weight of 300-350 grams. In the acute study (single dose) and in the subacute study (five doses) a dose of 5000 mg/kg was employed. Metaphase chromosome spreads were prepared from the bone marrow cells of these animals and scored for chromosomal aberrations. Neither the variety nor the number of these aberrations differed significantly from the negative controls; hence, compound FDA 71-51, Methocel, can be considered non-mutagenic as measured by the cytogenetic test.



CYTOGENETIC SUMMARY SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST II



BIONETICS

METHOCEL
FDA 71-51
ACUTE STUDY
METAPHASE SUMMARY SHEET
TEST II

<u>Compound</u>	<u>Dosage (mg/kg)</u>	<u>Time*</u>	<u>No. of Animals</u>	<u>No. of Cells</u>	<u>Mitotic Index %⁺⁺</u>	<u>No. of Cells w/ Breaks**</u>	<u>No. of Cells w/ Reunion**</u>	<u>No. of Cells With Other Aberrations**⁺</u>	<u>No. of Cells w/ Aber.**</u>
High Level	5000	6 hrs.	5	250	5.90	0	1(0.4)	0	1(0.4)
		24 hrs.	5	220	3.45	0	2(0.9)	0	2(0.9)
		48 hrs.	5	250	4.40	0	0	0	0
Negative Control	Food Pellet	6 hrs.	3	150	8.06	0	0	0	0
		24 hrs.	3	150	9.40	0	0	0	0
		48 hrs.	3	150	8.13	0	1(0.66)	0	1(0.66)
Positive Control TEM	0.3	24 hrs.	5	250	4.65	14(5.6)	34(13.6)	>12(4.8) 7f(2.8) 1pp(0.4) 1pu(0.4)	58(23.2)

* Time of kill after dosing.

** Numbers in () are percent aberrations per total cells counted.

+ Symbols: > = greater than 10 aberrations per cell; f = fragments; pp = polyploid; and pu = pulverization.

++ Based on a count of at least 500 cells per animal.

METHOCEL
FDA 71-51
SUBACUTE STUDY
METAPHASE SUMMARY SHEET
TEST II

<u>Compound</u>	<u>Dosage (mg/kg)</u>	<u>No. of Animals</u>	<u>No. of Cells</u>	<u>Mitotic Index %⁺⁺</u>	<u>No. of Cells w/ Breaks^{**}</u>	<u>No. of Cells w/ Reunion^{**}</u>	<u>No. of Cells w/ Other Aber.^{**}</u>	<u>No. of Cells w/ Aber.^{**}</u>
High Level	5000	5	250	4.06	0	3(1.2)	0	3(1.2)
Negative Control	Food Pellet	3	150	5.60	0	0	0	0

^{**} Numbers in () are percent aberrations per total cells counted.

⁺⁺ Based on a count of at least 500 cells per animal.

7. Dominant Lethal Study - Test I

a. Acute study

Significant differences between the negative control and experimental groups were shown in a few instances, but no strong indications of change were seen.

b. Subacute study

Significant differences between the negative control and experimental groups were shown in a few instances, but no strong indications of change were seen.



c. **DOMINANT LETHAL ASSAY SUMMARY SHEETS**

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST I

(Through error the computer had been programmed so that a double rounding off of numbers occurred at print out. In no way does this alter the statistics which are calculated on the full unrounded numbers.)



TABLE I
COMPOUND 51 STUDY ACUTE

FERTILITY INDEX

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
!		1	95/139=0.69	14/20=0.70	12/20=0.60	9/20=0.45*	10/20=0.50	14/20=0.70
!		2	103/139=0.75	16/20=0.80	16/20=0.80	11/19=0.58	12/20=0.60	15/20=0.75
		3	104/138=0.76	15/20=0.75	14/20=0.70	14/19=0.74	14/20=0.70	15/20=0.75
		4	118/140=0.85	18/20=0.90	17/20=0.85	17/20=0.85	15/20=0.75	14/20=0.70
		5	110/139=0.80	17/20=0.85	15/20=0.75	14/20=0.70	15/20=0.75	15/18=0.84
		6	109/139=0.79	19/20=0.95	14/20=0.70*	14/20=0.70*	18/20=0.90	17/20=0.85
		7	117/138=0.85	16/19=0.85	19/20=0.95	17/20=0.85	18/20=0.90	17/19=0.90
		8	116/140=0.83	17/20=0.85	17/20=0.85	18/20=0.90	17/20=0.85	15/20=0.75

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE II
COMPOUND 51 STUDY ACUTE

AVERAGE NUMBER OF IMPLANTATIONS PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
		1	1180/95=12.4	171/14=12.2	150/12=12.5	102/9=11.3	126/10=12.6	166/14=11.9
		2	1223/103=11.9	204/16=12.8	195/16=12.2	143/11=13.0 *@I	139/12=11.6@D	184/15=12.3
!		3	1276/104=12.3	159/15=10.6 *@D	154/14=11.0 @D	159/14=11.4	168/14=12.0@I	173/15=11.5
		4	1408/118=11.9	218/18=12.1	197/17=11.6	197/17=11.6	163/15=10.9	172/14=12.3
10		5	1290/110=11.7	176/17=10.4	175/15=11.7	167/14=11.9	179/15=11.9	172/15=11.5
		6	1292/109=11.9	220/19=11.6	163/14=11.6	173/14=12.4	221/18=12.3	204/17=12.0
		7	1436/117=12.3	190/16=11.9	226/19=11.9	204/17=12.0	218/18=12.1	209/17=12.3
		8	1353/116=11.7	198/17=11.7	212/17=12.5	199/18=11.1	199/17=11.7	186/15=12.4

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE III
COMPOUND 51 STUDY ACUTE

AVERAGE CORPORA LUTEA PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
		1	1322/95=13.9	182/14=13.0 @D	152/12=12.7	124/9=13.8 *@D	131/10=13.1	190/14=13.6
ε !	ε !	2	1359/103=13.2	229/16=14.3	221/16=13.8	154/11=14.0	148/12=12.3@D	202/15=13.5
		3	1364/104=13.1	201/15=13.4	176/14=12.6	175/14=12.5	182/14=13.0	222/15=14.8 *@I
		4	1532/118=13.0	252/18=14.0	223/17=13.1	217/17=12.8	187/15=12.5	189/14=13.5
ε !!	!	5	1428/110=13.0	220/17=12.9	215/15=14.3 *@I	195/14=13.9	215/15=14.3 @I	192/15=12.8
		6	1446/109=13.3	243/19=12.8	200/14=14.3@I	187/14=13.4	245/18=13.6	233/17=13.7
11		7	1543/117=13.2	224/16=14.0	256/19=13.5	226/17=13.3	235/18=13.1	233/17=13.7
		8	1599/116=13.8	224/17=13.2	228/17=13.4	229/18=12.7 @D	227/17=13.4	197/15=13.1

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ε AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, ε, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, ε, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

ε, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE IV
COMPOUND 51 STUDY ACUTE

AVERAGE PREIMPLANTATION LOSSES PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
		1	142/95 = 1.5	11/14 = 0.8	2/12 = 0.2 @D ** @D	22/9 = 2.4	5/10 = 0.5 * @D	24/14 = 1.7
		2	136/103 = 1.3	25/16 = 1.6	26/16 = 1.6	11/11 = 1.0	9/12 = 0.8	18/15 = 1.2
$\& \text{!!}$!	3	88/104 = 0.9	42/15 = 2.8 ** @DI	22/14 = 1.6 @I	16/14 = 1.1 @D	14/14 = 1.0* @D	49/15 = 3.3 ** @D
		4	124/118 = 1.1	34/18 = 1.9 * @I	26/17 = 1.5	20/17 = 1.2	24/15 = 1.6	17/14 = 1.2
$\& \text{!!}$		5	138/110 = 1.3	44/17 = 2.6	40/15 = 2.7 ** @DI	28/14 = 2.0 @I	36/15 = 2.4	20/15 = 1.3
1.2		6	154/109 = 1.4	23/19 = 1.2	37/14 = 2.6	14/14 = 1.0	24/18 = 1.3	29/17 = 1.7
$\& \text{!}$		7	107/117 = 0.9	34/16 = 2.1 ** @DI	30/19 = 1.6 @I	22/17 = 1.3	17/18 = 0.9* @D	24/17 = 1.4 @I
		8	246/116 = 2.1	26/17 = 1.5	16/17 = 0.9 * @D	30/18 = 1.7	28/17 = 1.7	11/15 = 0.7 * @D

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

$\&$ AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, $\&$, @ , * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, $\&$, @ , * = SIGNIFICANT AT P LESS THAN 0.01

* @ SIGNIFICANTLY DIFFERENT FROM CONTROL
 $\&$, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE V
COMPOUND 51 STUDY ACUTE

AVERAGE RESORPTIONS (DEAD IMPLANTS) PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
!		1	20/ 95=0.22	8/14=0.58	1/12=0.09@D	1/ 9=0.12	1/10=0.10@D	18/14=1.29 **@@I
		2	43/103=0.42	10/16=0.63	7/16=0.44	9/11=0.82	4/12=0.34	32/15=2.14*I **@@I
		3	53/104=0.51	8/15=0.54	9/14=0.65	6/14=0.43	3/14=0.22	22/15=1.47 @I
		4	53/118=0.45	9/18=0.50	8/17=0.48	8/17=0.48	10/15=0.67	11/14=0.79
!		5	60/110=0.55	14/17=0.83	11/15=0.74	4/14=0.29	5/15=0.34	18/15=1.20 @I
!		6	45/109=0.42	13/19=0.69	9/14=0.65	7/14=0.50	10/18=0.56	16/17=0.95
		7	53/117=0.46	12/16=0.75	9/19=0.48	10/17=0.59	15/18=0.84	8/17=0.48
		8	65/116=0.57	6/17=0.36	19/17=1.12	8/18=0.45	11/17=0.65	11/15=0.74

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VI
COMPOUND 51 STUDY ACUTE

PROPORTION OF FEMALES WITH ONE OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
		1	19/ 95=0.20	5/14=0.36	1/12=0.09	1/ 9=0.12	1/10=0.10	9/14=0.65 **
		2	32/103=0.32	6/16=0.38	4/16=0.25	6/11=0.55	4/12=0.34	11/15=0.74* **
		3	32/104=0.31	7/15=0.47	6/14=0.43	4/14=0.29	3/14=0.22	9/15=0.60 *
		4	39/118=0.34	7/18=0.39	4/17=0.24	6/17=0.36	4/15=0.27	4/14=0.29
		5	36/110=0.33	9/17=0.53	6/15=0.40	4/14=0.29	5/15=0.34	9/15=0.60 *
		6	36/109=0.34	8/19=0.43	5/14=0.36	4/14=0.29	6/18=0.34	8/17=0.48
		7	38/117=0.33	8/16=0.50	4/19=0.22	6/17=0.36	6/18=0.34	8/17=0.48
		8	44/116=0.38	6/17=0.36	8/17=0.48	3/18=0.17	10/17=0.59	7/15=0.47

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VII
COMPOUND 51 STUDY ACUTE

PORPORTION OF FEMALES WITH TWO OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
		1	1/ 95=0.02	2/14=0.15 **	0/12=0.0	0/ 9=0.0	0/10=0.0	5/14=0.36 **
		2	11/103=0.11	3/16=0.19	3/16=0.19	1/11=0.10	0/12=0.0	4/15=0.27
		3	16/104=0.16	1/15=0.07	1/14=0.08	2/14=0.15	0/14=0.0	3/15=0.20
		4	11/118=0.10	1/18=0.06	1/17=0.06	2/17=0.12	2/15=0.14	2/14=0.15
		5	16/110=0.15	2/17=0.12	4/15=0.27	0/14=0.0	0/15=0.0	3/15=0.20
15		6	9/109=0.09	4/19=0.22	3/14=0.22	3/14=0.22	2/18=0.12	3/17=0.18
		7	11/117=0.10	3/16=0.19	2/19=0.11	3/17=0.18	3/18=0.17	0/17=0.0
		8	18/116=0.16	0/17=0.0	3/17=0.18	2/18=0.12	1/17=0.06	3/15=0.20

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VIII
COMPOUND 51 STUDY ACUTE

DEAD IMPLANTS / TOTAL IMPLANTS

WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG	POSITIVE CONTROL
1	20/1180=0.02	8/171=0.05	1/150=0.01@d	1/102=0.01	1/126=0.01@d	18/166=0.11 **@
2	43/1223=0.04	10/204=0.05	7/195=0.04	9/143=0.07	4/139=0.03	32/184=0.18@I @I
3	53/1276=0.05	8/159=0.06	9/154=0.06	6/159=0.04	3/168=0.02* *@D	22/173=0.13
4	53/1408=0.04	9/218=0.05	8/197=0.05	8/197=0.05	10/163=0.07	11/172=0.07
5	60/1290=0.05	14/176=0.08	11/175=0.07	4/167=0.03@d @D	5/179=0.03* *@D	18/172=0.11
6	45/1292=0.04	13/220=0.06	9/163=0.06	7/173=0.05	10/221=0.05	16/204=0.08
7	53/1436=0.04	12/190=0.07	9/226=0.04	10/204=0.05	15/218=0.07	8/209=0.04
8	65/1353=0.05	6/198=0.04	19/212=0.09	8/199=0.05	11/199=0.06	11/186=0.06

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE HISTORICAL CONTROL GROUP

* = TWO-TAILED TEST

@ = ONE-TAILED TEST

ONE *,@ = SIGNIFICANT AT P LESS THAN 0.05

TWO *,@ = SIGNIFICANT AT P LESS THAN 0.01

*,@ SIGNIFICANTLY DIFFERENT FROM CONTROL

TABLE I
COMPOUND 51 STUDY SUBACUTE

FERTILITY INDEX

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
		1	92/139=0.67	12/20=0.60	10/20=0.50	15/20=0.75	12/20=0.60
		2	104/140=0.75	14/20=0.70	18/20=0.90	17/19=0.90	12/20=0.60
		3	101/139=0.73	18/20=0.90	16/20=0.80	15/20=0.75	13/20=0.65
		4	104/134=0.78	16/20=0.80	15/20=0.75	18/20=0.90	12/20=0.60
		5	108/139=0.78	14/18=0.78	18/19=0.95	17/20=0.85	14/20=0.70
		6	120/139=0.87	16/20=0.80	17/20=0.85	16/20=0.80	15/20=0.75
		7	117/135=0.87	18/20=0.90	18/20=0.90	18/20=0.90	16/20=0.80

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE II
COMPOUND 51 STUDY SUBACUTE

AVERAGE NUMBER OF IMPLANTATIONS PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
!		1	1084/ 92=11.8	147/12=12.3	126/10=12.6	192/15=12.8 @I	153/12=12.8
		2	1301/104=12.5	173/14=12.4	226/18=12.6	216/17=12.7	143/12=11.9
		3	1196/101=11.8	209/18=11.6	180/16=11.3	169/15=11.3	150/13=11.5
		4	1221/104=11.7	193/16=12.1	182/15=12.1	217/18=12.1	140/12=11.7
!	!	5	1299/108=12.0	163/14=11.6	219/18=12.2	200/17=11.8	155/14=11.1
		6	1437/120=12.0	189/16=11.8	178/17=10.5 *@@D	193/16=12.1	175/15=11.7
& !	!	7	1352/117=11.6	214/18=11.9	192/18=10.7	228/18=12.7 *@I	204/16=12.8 @I

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE III
COMPOUND 51 STUDY SUBACUTE

AVERAGE CORPORA LUTEA PER PREGNANT FEMALE

LOG ARITH DOSE DOSE WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
1	1218/92=13.2	167/12=13.9	137/10=13.7	223/15=14.9 *@I	161/12=13.4
2	1395/104=13.4	204/14=14.6	245/18=13.6	249/17=14.7	158/12=13.2
3	1290/101=12.8	245/18=13.6	201/16=12.6	194/15=12.9	177/13=13.6
4	1285/104=12.4	214/16=13.4	200/15=13.3 @I	231/18=12.8	153/12=12.8
5	1366/108=12.7	188/14=13.4	281/18=15.6@I **@@I	230/17=13.5	190/14=13.6
6	1580/120=13.2	229/16=14.3	226/17=13.3	208/16=13.0	187/15=12.5@d
7	1474/117=12.6	237/18=13.2	228/18=12.7	247/18=13.7 @I	211/16=13.2

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE IV
COMPOUND 51 STUDY SUBACUTE

AVERAGE PREIMPLANTATION LOSSES PER PREGNANT FEMALE

LOG ARITH DOSE DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
	1	134/ 92= 1.5	20/12= 1.7	11/10= 1.1	31/15= 2.1	8/12= 0.7 @D
	2	94/104= 0.9	31/14= 2.2	19/18= 1.1	33/17= 1.9	15/12= 1.3
E ! & !	3	94/101= 0.9	36/18= 2.0 @I	21/16= 1.3	25/15= 1.7	27/13= 2.1 *@I
E !	4	64/104= 0.6	21/16= 1.3	18/15= 1.2 *@@I	14/18= 0.8	13/12= 1.1
E&!! E&!!	5	67/108= 0.6	25/14= 1.8	62/18= 3.4 **@@I	30/17= 1.8 **@@I	35/14= 2.5 **@@I
!	6	143/120= 1.2	40/16= 2.5	48/17= 2.8 *@@I	15/16= 0.9	12/15= 0.8
! !	7	122/117= 1.0	23/18= 1.3	36/18= 2.0	19/18= 1.1	7/16= 0.4*@D @D

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE V
COMPOUND 51 STUDY SUBACUTE

AVERAGE RESORPTIONS (DEAD IMPLANTS) PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
		1	35/ 92=0.39	5/12=0.42	4/10=0.40	0/15=0.0 @D **@@D	4/12=0.34
		2	49/104=0.48	10/14=0.72	6/18=0.34	15/17=0.89 **@@I	5/12=0.42
		3	55/101=0.55	14/18=0.78	8/16=0.50	4/15=0.27	17/13=1.31
!&!!		4	61/104=0.59	5/16=0.32	3/15=0.20 *@D	1/18=0.06 **@@D	3/12=0.25
		5	71/108=0.66	7/14=0.50	9/18=0.50	10/17=0.59	11/14=0.79
!	!	6	47/120=0.40	15/16=0.94 *@@I	12/17=0.71	11/16=0.69 *@I	5/15=0.34*@D
21		7	59/117=0.51	11/18=0.62	15/18=0.84	13/18=0.73	10/16=0.63

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VI
COMPOUND 51 STUDY SUBACUTE

PROPORTION OF FEMALES WITH ONE OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000MG/KG
		1	28/ 92=0.31	3/12=0.25	4/10=0.40	0/15=0.0 *	3/12=0.25
		2	32/104=0.31	6/14=0.43	5/18=0.28	12/17=0.71 **	5/12=0.42
		3	34/101=0.34	8/18=0.45	5/16=0.32	3/15=0.20	6/13=0.47
!!		4	38/104=0.37	4/16=0.25	2/15=0.14	1/18=0.06 **	2/12=0.17
!!		5	49/108=0.46	5/14=0.36	4/18=0.23	8/17=0.48	6/14=0.43
22		6	33/120=0.28	10/16=0.63 **	6/17=0.36	10/16=0.63 **	3/15=0.20*
		7	34/117=0.30	8/18=0.45	9/18=0.50	6/18=0.34	7/16=0.44

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VII
COMPOUND 51 STUDY SUBACUTE

PROPORTION OF FEMALES WITH TWO OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
		1	6/ 92=0.07	2/12=0.17	0/10=0.0	0/15=0.0	1/12=0.09
		2	8/104=0.08	2/14=0.15	1/18=0.06	2/17=0.12	0/12=0.0
		3	14/101=0.14	3/18=0.17	3/16=0.19	1/15=0.07	3/13=0.24
		4	14/104=0.14	1/16=0.07	1/15=0.07	0/18=0.0	1/12=0.09
		5	18/108=0.17	1/14=0.08	1/18=0.06	2/17=0.12	2/14=0.15
		6	9/120=0.08	4/16=0.25*	3/17=0.18	1/16=0.07	1/15=0.07
		7	14/117=0.12	2/18=0.12	4/18=0.23	3/18=0.17	3/16=0.19

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VIII
COMPOUND 51 STUDY SUBACUTE

DEAD IMPLANTS / TOTAL IMPLANTS

WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 4.750 MG/KG	DOSE LEVEL 47.500 MG/KG	DOSE LEVEL 475.000 MG/KG
1	35/1084=0.04	5/147=0.04	4/126=0.04	0/192=0.0 @D **@@D	4/153=0.03
2	49/1301=0.04	10/173=0.06	6/226=0.03	15/216=0.07	5/143=0.04
3	55/1196=0.05	14/209=0.07	8/180=0.05	4/169=0.03	17/150=0.12
4	61/1221=0.05	5/193=0.03	3/182=0.02 **@@D	1/217=0.01 **@@D	3/140=0.03
5	71/1299=0.06	7/163=0.05	9/219=0.05	10/200=0.05	11/155=0.08
6	47/1437=0.04	15/189=0.08 @I	12/178=0.07	11/193=0.06	5/175=0.03
7	59/1352=0.05	11/214=0.06	15/192=0.08	13/228=0.06	10/204=0.05

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE HISTORICAL CONTROL GROUP

* = TWO-TAILED TEST
@ = ONE-TAILED TEST

ONE *,@ = SIGNIFICANT AT P LESS THAN 0.05
TWO *,@ = SIGNIFICANT AT P LESS THAN 0.01

*,@ SIGNIFICANTLY DIFFERENT FROM CONTROL

8. Dominant Lethal Study - Test II

Compound FDA 71-51, Methocel, was administered to ten male rats (400 grams) at a dose level of 5,000 mg/kg according to acute (single dose) and subacute (five doses). Each treated male rat was mated with two virgin female rats each week for seven (subacute) or eight (acute) weeks. Two weeks after mating, these female rats were sacrificed and the fertility index, preimplantation loss and lethal effects on the embryos were determined and compared with these same parameters calculated from negative (saline-dosed) and positive (0.3 mg/kg TEM-dosed) control animals.

The values calculated for these parameters from animals dosed with compound FDA 71-51, Methocel, did not significantly vary from those obtained from the negative controls. The significant increase in dead embryos per total implants registered for the week 8 acute sequence is primarily due to the fact that 16 dead embryos were present in a single female. If this female is treated as an outlier then the ratio obtained is more in line with the others.



DOMINANT LETHAL ASSAY SUMMARY SHEETS

CONTRACT FDA 71-268

COMPOUND FDA 71-51

METHOCEL

TEST II

(Through error the computer had been programmed so that a double rounding off of numbers occurred at print out. In no way does this alter the statistics which are calculated on the full unrounded numbers.)



BIONETICS

COMPOUND 51

TABLE I

STUDY ACUTE

FERTILITY INDEX

LOG ARITH DOSE DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
	1	138/199=0.70	16/ 20=0.80	14/ 20=0.70	10/ 20=0.50*
	2	154/199=0.78	15/ 20=0.75	17/ 20=0.85	12/ 20=0.60
	3	154/198=0.78	14/ 20=0.70	16/ 20=0.80	4/ 20=0.20**
	4	172/200=0.86	17/ 20=0.85	13/ 20=0.65*	6/ 20=0.30**
	5	160/199=0.81	15/ 20=0.75	16/ 20=0.80	15/ 20=0.75
	6	156/199=0.79	19/ 20=0.95	16/ 20=0.80	14/ 20=0.70*
	7	169/197=0.86	15/ 20=0.75	18/ 20=0.90	15/ 20=0.75
	8	166/200=0.83	16/ 20=0.80	18/ 20=0.90	14/ 20=0.70

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE II
COMPOUND 51 STUDY ACUTE

AVERAGE NUMBER OF IMPLANTATIONS PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
		1	1720/138=12.5	203/ 16=12.7	164/ 14=11.7	100/ 10=10.0**@D **@D
		2	1871/154=12.2	182/ 15=12.1	221/ 17=13.0 *@I	81/ 12= 6.8**@D **@D
		3	1867/154=12.1	173/ 14=12.4	203/ 16=12.7	25/ 4= 6.3**@D **@D
		4	2063/172=12.0	208/ 17=12.2	148/ 13=11.4	32/ 6= 5.3**@D **@D
		5	1906/160=11.9	193/ 15=12.9 *@I	195/ 16=12.2	149/ 15= 9.9**@D *@D
		6	1868/156=12.0	254/ 19=13.4 *@@I	195/ 16=12.2	185/ 14=13.2 *@I
		7	2082/169=12.3	183/ 15=12.2	217/ 18=12.1	184/ 15=12.3
		8	1978/166=11.9	199/ 16=12.4	218/ 18=12.1	182/ 14=13.0 *@@I

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE III
COMPOUND 51 STUDY ACUTE

AVERAGE CORPORA LUTEA PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
		1	1936/138=14.0	251/ 16=15.7 * ∂ I	202/ 14=14.4	126/ 10=12.6** ∂ D ∂ D
		2	2120/154=13.8	217/ 15=14.5	247/ 17=14.5	168/ 12=14.0
		3	2087/154=13.6	215/ 14=15.4 * ∂ ∂ I	248/ 16=15.5 ∂ I	54/ 4=13.5* ∂ D
		4	2299/172=13.4	232/ 17=13.7	188/ 13=14.5	77/ 6=12.8
		5	2132/160=13.3	220/ 15=14.7 * ∂ I	237/ 16=14.8 ∂ I	201/ 15=13.4
		6	2100/156=13.5	302/ 19=15.9 ** ∂ ∂ I	228/ 16=14.3	252/ 14=18.0 ** ∂ ∂ I
		7	2284/169=13.5	203/ 15=13.5	258/ 18=14.3	221/ 15=14.7 * ∂ I
		8	2330/166=14.0	217/ 16=13.6	260/ 18=14.4	212/ 14=15.1 ∂ I

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

$\&$ AND * = TWO-TAILED TEST
! AND ∂ = ONE-TAILED TEST

ONE !, $\&$, ∂ , * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, $\&$, ∂ , * = SIGNIFICANT AT P LESS THAN 0.01

* , ∂ SIGNIFICANTLY DIFFERENT FROM CONTROL
 $\&$, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE IV
COMPOUND 51 STUDY ACUTE

AVERAGE PREIMPLANTATION LOSSES PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
		1	216/138= 1.6	48/ 16= 3.0	38/ 14= 2.7	26/ 10= 2.6
		2	249/154= 1.6	35/ 15= 2.3 **@@I	26/ 17= 1.5	87/ 12= 7.3**@@I **@@I
		3	220/154= 1.4	42/ 14= 3.0 **@@I	45/ 16= 2.8 *@I	29/ 4= 7.3**@@I **@@I
		4	236/172= 1.4	24/ 17= 1.4	40/ 13= 3.1*@@I **@@I	45/ 6= 7.5**@@I **@@I
		5	226/160= 1.4	27/ 15= 1.8	42/ 16= 2.6	52/ 15= 3.5 **@@I
		6	232/156= 1.5	48/ 19= 2.5 *@@I	33/ 16= 2.1	67/ 14= 4.8 **@@I
		7	202/169= 1.2	20/ 15= 1.3	41/ 18= 2.3 *@@I	37/ 15= 2.5 @I
		8	352/166= 2.1	18/ 16= 1.1	42/ 18= 2.3@I	30/ 14= 2.1

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE V

STUDY ACUTE

AVERAGE RESORPTIONS (DEAD IMPLANTS) PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
		1	37/138=0.27	15/ 16=0.94 * ∂ I	6/ 14=0.43	88/ 10=8.80** $\partial\partial$ I ** $\partial\partial$ I
		2	72/154=0.47	13/ 15=0.87	9/ 17=0.53	69/ 12=5.75** $\partial\partial$ I ** $\partial\partial$ I
		3	89/154=0.58	3/ 14=0.22 ² * ∂ D	11/ 16=0.69 ∂ I	25/ 4=6.25** $\partial\partial$ I ** $\partial\partial$ I
		4	85/172=0.50	8/ 17=0.48	8/ 13=0.62	29/ 6=4.84** $\partial\partial$ I ** $\partial\partial$ I
		5	98/160=0.62	6/ 15=0.40	12/ 16=0.75	55/ 15=3.67** $\partial\partial$ I ** $\partial\partial$ I
		6	76/156=0.49	17/ 19=0.90	15/ 16=0.94	19/ 14=1.36 * ∂ I
		7	87/169=0.52	3/ 15=0.20 * ∂ D	6/ 18=0.34	19/ 15=1.27** $\partial\partial$ I ** $\partial\partial$ I
		8	87/166=0.53	9/ 16=0.57	33/ 18=1.84* ∂ I ** $\partial\partial$ I	15/ 14=1.08

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ϵ AND * = TWO-TAILED TEST

! AND ∂ = ONE-TAILED TEST

ONE !, ϵ , ∂ , * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, ϵ , ∂ , * = SIGNIFICANT AT P LESS THAN 0.01

*, ∂ SIGNIFICANTLY DIFFERENT FROM CONTROL

ϵ , ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE VI
COMPOUND 51 STUDY ACUTE

PROPORTION OF FEMALES WITH ONE OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
		1	31/138=0.23	8/ 16=0.50 * *	6/ 14=0.43	10/ 10=1.00** **
		2	52/154=0.34	6/ 15=0.40	8/ 17=0.48	12/ 12=1.00** **
		3	55/154=0.36	2/ 14=0.15	8/ 16=0.50*	4/ 4=1.00** **
		4	64/172=0.38	7/ 17=0.42	7/ 13=0.54	6/ 6=1.00* **
		5	60/160=0.38	6/ 15=0.40	10/ 16=0.63	15/ 15=1.00** **
		6	55/156=0.36	10/ 19=0.53	7/ 16=0.44	9/ 14=0.65 *
		7	62/169=0.37	3/ 15=0.20	5/ 18=0.28	11/ 15=0.74** **
		8	62/166=0.38	7/ 16=0.44	16/ 18=0.89** **	8/ 14=0.58

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE VII

STUDY ACUTE

PORPORTION OF FEMALES WITH TWO OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
		1	4/138=0.03	4/ 16=0.25 **	0/ 14=0.0 *	10/ 10=1.00** **
		2	18/154=0.12	4/ 15=0.27	1/ 17=0.06	12/ 12=1.00** **
		3	25/154=0.17	1/ 14=0.08	2/ 16=0.13	4/ 4=1.00** **
		4	15/172=0.09	1/ 17=0.06	1/ 13=0.08	5/ 6=0.84** **
		5	24/160=0.15	0/ 15=0.0	2/ 16=0.13	12/ 15=0.80** **
		6	17/156=0.11	4/ 19=0.22	4/ 16=0.25	4/ 14=0.29
		7	20/169=0.12	0/ 15=0.0	1/ 18=0.06	5/ 15=0.34* *
		8	20/166=0.13	2/ 16=0.13	3/ 18=0.17	3/ 14=0.22

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE VIII
STUDY ACUTE

DEAD IMPLANTS / TOTAL IMPLANTS

WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG	POSITIVE CONTROL
1	37/1720=0.02	15/203=0.07 *aI	6/164=0.04 aI	88/100=0.88**aI **aI
2	72/1871=0.04	13/182=0.07	9/221=0.04	69/ 81=0.85**aI **aI
3	89/1867=0.05	3/173=0.02 *aD	11/203=0.05aI	25/ 25=1.00**aI **aI
4	85/2063=0.04	8/208=0.04	8/148=0.05	29/ 32=0.91**aI **aI
5	98/1906=0.05	6/193=0.03	12/195=0.06	55/149=0.37**aI **aI
6	76/1868=0.04	17/254=0.07	15/195=0.08	19/185=0.10 aI
7	87/2082=0.04	3/183=0.02 *aD	6/217=0.03	19/184=0.10**aI *aI
8	87/1978=0.04	9/199=0.05	33/218=0.15*aI *aI	15/182=0.08

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE HISTORICAL CONTROL GROUP

* = TWO-TAILED TEST
a = ONE-TAILED TEST

ONE *,a = SIGNIFICANT AT P LESS THAN 0.05
TWO *,a = SIGNIFICANT AT P LESS THAN 0.01

*,a SIGNIFICANTLY DIFFERENT FROM CONTROL

COMPOUND 51

TABLE I

STUDY SUBACUTE

FERTILITY INDEX

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
		1	134/199=0.68	14/ 20=0.70	13/ 20=0.65
		2	152/200=0.76	14/ 20=0.70	17/ 19=0.90
		3	156/199=0.79	15/ 20=0.75	16/ 20=0.80
		4	153/194=0.79	17/ 20=0.85	15/ 20=0.75
		5	155/197=0.79	16/ 20=0.80	14/ 20=0.70
		6	168/199=0.85	18/ 20=0.90	17/ 20=0.85
		7	169/195=0.87	17/ 20=0.85	18/ 20=0.90

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE II
COMPOUND 51 STUDY SUBACUTE

AVERAGE NUMBER OF IMPLANTATIONS PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
		1	1630/134=12.2	162/ 14=11.6	160/ 13=12.3
		2	1904/152=12.5	182/ 14=13.0	236/ 17=13.9 **@@I
		3	1864/156=12.0	192/ 15=12.8	192/ 16=12.0
		4	1791/153=11.7	207/ 17=12.2	165/ 15=11.0
		5	1874/155=12.1	205/ 16=12.8	151/ 14=10.8*@D
		6	2021/168=12.0	216/ 18=12.0	203/ 17=11.9
		7	1954/169=11.6	219/ 17=12.9	217/ 18=12.1 **@@I

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ε AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, ε, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, ε, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
ε, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

TABLE III
COMPOUND 51 STUDY SUBACUTE

AVERAGE CORPORA LUTEA PER PREGNANT FEMALE

LOG ARITH DOSE DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
1	1855/134=13.8	195/ 14=13.9	187/ 13=14.4	
2	2110/152=13.9	232/ 14=16.6	283/ 17=16.7	
		**@@I	**@@I	
3	2085/156=13.4	244/ 15=16.3	238/ 16=14.9	
		**@@I	@I	
4	1953/153=12.8	244/ 17=14.4	190/ 15=12.7*@D	
		*@@I		
5	2079/155=13.4	254/ 16=15.9	203/ 14=14.5	
		**@@I		
6	2290/168=13.6	259/ 18=14.4	221/ 17=13.0@D	
7	2218/169=13.1	287/ 17=16.9	258/ 18=14.3**@@D	
		**@@I	*@@I	

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST
! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05
TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL
&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE IV

STUDY SUBACUTE

AVERAGE PREIMPLANTATION LOSSES PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	HISTORICAL CONTROL	WEEK	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
1		225/134= 1.7		33/ 14= 2.4	27/ 13= 2.1
2		206/152= 1.4		50/ 14= 3.6 **@@I	47/ 17= 2.8 **@@I
3		221/156= 1.4		52/ 15= 3.5 *@I	46/ 16= 2.9 *@@I
4		162/153= 1.1		37/ 17= 2.2 **@@I	25/ 15= 1.7
5		205/155= 1.3		49/ 16= 3.1 **@@I	52/ 14= 3.7 **@@I
6		269/168= 1.6		43/ 18= 2.4 @I	18/ 17= 1.1@D
7		264/169= 1.6		68/ 17= 4.0 **@@I	41/ 18= 2.3

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ε AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, ε, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, ε, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

ε, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE V

STUDY SUBACUTE

AVERAGE RESORPTIONS (DEAD IMPLANTS) PER PREGNANT FEMALE

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
		1	56/134=0.42	14/ 14=1.00	5/ 13=0.39
		2	77/152=0.51	18/ 14=1.29 @I	14/ 17=0.83
		3	90/156=0.58	16/ 15=1.07	12/ 16=0.75
		4	85/153=0.56	8/ 17=0.48	12/ 15=0.80
		5	96/155=0.62	14/ 16=0.88	21/ 14=1.50 *@@I
		6	84/168=0.50	17/ 18=0.95 *@@I	9/ 17=0.53@D
		7	107/169=0.64	20/ 17=1.18	17/ 18=0.95

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

& AND * = TWO-TAILED TEST

! AND @ = ONE-TAILED TEST

ONE !, &, @, * = SIGNIFICANT AT P LESS THAN 0.05

TWO !, &, @, * = SIGNIFICANT AT P LESS THAN 0.01

*, @ SIGNIFICANTLY DIFFERENT FROM CONTROL

&, ! SIGNIFICANT RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

PROPORTION OF FEMALES WITH ONE OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
		1	42/134=0.32	7/ 14=0.50	4/ 13=0.31
		2	50/152=0.33	8/ 14=0.58	6/ 17=0.36
		3	55/156=0.36	7/ 15=0.47	9/ 16=0.57
		4	58/153=0.38	6/ 17=0.36	8/ 15=0.54
		5	67/155=0.44	8/ 16=0.50	10/ 14=0.72 *
		6	58/168=0.35	12/ 18=0.67 **	5/ 17=0.30*
		7	61/169=0.37	7/ 17=0.42	10/ 18=0.56

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE VII

STUDY SUBACUTE

PORPORTION OF FEMALES WITH TWO OR MORE DEAD IMPLANTATIONS

LOG DOSE	ARITH DOSE	WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
		1	12/134=0.09	4/ 14=0.29 *	1/ 13=0.08
		2	14/152=0.10	4/ 14=0.29 *	5/ 17=0.30 *
		3	22/156=0.15	5/ 15=0.34	2/ 16=0.13
		4	18/153=0.12	2/ 17=0.12	2/ 15=0.14
		5	23/155=0.15	4/ 16=0.25	6/ 14=0.43 **
		6	20/168=0.12	4/ 18=0.23	3/ 17=0.18
		7	25/169=0.15	7/ 17=0.42 **	4/ 18=0.23

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT RELATIONSHIPS AND DIFFERENCES USING THE HISTORICAL CONTROL GROUP

ONE !,* = SIGNIFICANT AT P LESS THAN 0.05

TWO !,* = SIGNIFICANT AT P LESS THAN 0.01

* SIGNIFICANTLY DIFFERENT FROM CONTROL

! SIGNIFICANT LINEAR RELATIONSHIP WITH ARITH OR LOG DOSE (HEADING OF COLUMN)

COMPOUND 51

TABLE VIII
STUDY SUBACUTE

DEAD IMPLANTS / TOTAL IMPLANTS

WEEK	HISTORICAL CONTROL	NEGATIVE CONTROL	DOSE LEVEL 5000. MG/KG
1	56/1630=0.03	14/162=0.09	5/160=0.03
2	77/1904=0.04	18/182=0.10	14/236=0.06
3	90/1864=0.05	16/192=0.08	12/192=0.06
4	85/1791=0.05	8/207=0.04	12/165=0.07
5	96/1874=0.05	14/205=0.07	21/151=0.14 *aai
6	84/2021=0.04	17/216=0.08 *ai	9/203=0.04
7	107/1954=0.05	20/219=0.09	17/217=0.08

SYMBOLS ON FIRST LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE NEGATIVE CONTROL GROUP

SYMBOLS ON SECOND LINE DENOTE SIGNIFICANT DIFFERENCES USING
THE HISTORICAL CONTROL GROUP

* = TWO-TAILED TEST

a = ONE-TAILED TEST

ONE *,a = SIGNIFICANT AT P LESS THAN 0.05

TWO *,a = SIGNIFICANT AT P LESS THAN 0.01

*,a SIGNIFICANTLY DIFFERENT FROM CONTROL

APPENDICES

II. MATERIALS AND METHODS

A. Animal Husbandry

1. Animals (Rats and Mice)

Ten to twelve week old rats (280 to 350 g) and male mice (25 to 30 g) were fed a commercial 4% fat diet and water ad libitum until they were put on experiment. Flow Laboratories random-bred, closed colony, Sprague-Dawley CD strain rats were used in the cytogenetic studies. Flow Laboratories ICR male mice were employed in the Host-Mediated Assay.

2. Preparation of Diet

A commercial 4% fat diet was fed to all animals. Periodic tests to verify the absence of coliforms, Salmonella and Pseudomonas sp. were performed.

3. Husbandry

Animals were held in quarantine for 4-11 days. Mice were housed five to a cage and rats one to five to a cage. Animals were identified by ear punch. Sanitary cages and bedding were used, and changed two times per week, at which time water containers were cleaned, sanitized and filled. Once a week, cages were repositioned on racks; racks were repositioned within rooms monthly. Personnel handling animals or working within animal facilities wore head coverings and face masks, as well as suitable garments. Individuals with respiratory or other overt infections were excluded from the animal facilities.

B. Dosage Determination

1. Acute LD₅₀ and LD₅ Determination

Since the compounds proposed for testing are included in



the food additive regulations as "generally recognized as safe" (GRAS), it was expected that a large number of them would be sufficiently non-toxic so that determination of a LD_{50} or a LD_5 would be of no practical value. In fact, this has been our experience with previously tested compounds from this list. In the case of these relatively non-toxic compounds, attempts were made to assure that the amounts to be administered would not affect the animals by means (mechanical, physical, etc.) related to their bulk rather than to their toxicity. In the cases of certain compounds where a LD_{50} or a LD_5 could not be determined, an exceedingly high concentration, 5 g/kg, was employed and accepted as the LD_5 level. In cases where the toxicity was high enough to allow determination of a LD_5 , the following protocol was used.

Thirty rats of the strain chosen for studies described below and of approximately the age and weight specified were assigned at random to six groups. Each group was then given, using the chosen route of administration, one of a series of dosages of the test compound following a logarithmic dosage scheme. The series of dosages were derived from a consideration of whatever toxicity information was available for the particular test compound. The objective in selecting dosages was to choose values which would cause mortalities between 10% and 90%.

When information was inadequate to derive a suitable series of dosages, five rats were used to identify the proper range. Each of these was given one of a widely spaced (differing by 10X) series of doses. This was confidently expected to suffice for derivation of the series of dosages to be used in the LD_{50} determination.



The mortalities observed when the series of dosages were given to the 30 rats were then subjected to a probit analysis and calculation of LD_{50} , LD_5 , slope and confidence limits by the method of Litchfield and Wilcoxon. The highest dose level used was either a finite LD_5 or 5000 mg/kg. The intermediate level used was either 1/10 of the finite LD_5 or 2500 mg/kg. The low level used was either 1/100 of the finite LD_5 or 30 mg/kg.

2. Subacute Studies

Subacute doses were identical to those used in the acute studies. Each subacute study animal was given the acute dosage once a day for each of five consecutive days (24 hours apart).

C. Mutagenicity Testing Protocols

1. Host-Mediated Assay

Flow Laboratories ICR random-bred male mice were used in this study. In the acute and subacute studies ten animals, 25-30 g each, were employed at each dose level. Solvent and positive controls were run at all times. The positive control (dimethyl nitrosamine) was run by the acute system only at a dose of 100 mg/kg for Salmonella. For yeast, ethyl methane sulfonate (EMS) intramuscularly injected at a dose of 350 mg/kg was used. The solvents used and the toxicity data are presented in the Results and Discussion Section of the report.

The indicator organisms used in this study were: (1) two histidine auxotrophs (his G-46, TA-1530) of Salmonella typhimurium, and (2) a diploid strain (D-3) of Saccharomyces cerevisiae. The induction of reverse mutation was determined with the Salmonella; mitotic recombination was determined with yeast. Chemicals were evaluated directly by in vitro bacterial and yeast studies prior to, or concurrent with, the studies in

mice. Only animals on the subacute studies were not fed the evening prior to compound administration. The Salmonella were carried in tryptone yeast extract gel, transferred weekly. They were transferred to tryptone yeast extract broth 48 hours before use: they were transferred a second time from broth to broth 24 hours prior to use, and again 8 hours before use. The mouse inoculum was prepared by transferring 4 ml of the 8-hour broth culture to 50 ml broth bottles which had been prewarmed at 37°C. Exponential log-phase organisms were inoculated intraperitoneally into the mice approximately 2-1/2 hours later when the appropriate density indicating 3.0×10^8 cells/ml was reached. The Saccharomyces was carried in yeast complete agar. The inoculum was prepared by harvesting the organisms from the surface of the plates with sterile saline. The cells were washed three times with sterile saline and suspended in a concentration of 5.0×10^8 cells/ml. Two ml of the suspension was inoculated into each mouse intraperitoneally. Total plate counts on Salmonella were on tryptone yeast extract and for Saccharomyces on yeast complete medium.

a. Acute study

Three dosage levels (usage, intermediate [determined as discussed previously], and LD_5) were administered orally by intubation to ten mice. Positive controls and negative vehicle controls were included in each study. All animals received 2 ml of the indicator organism intraperitoneally. Each ml contained 3.0×10^8 cells for Salmonella and 5.0×10^8 cells for Saccharomyces. Three hours later, each animal was killed and 2 ml of sterile saline was introduced intraperitoneally. As much fluid as possible was then aseptically removed from the peritoneal cavity. Dilution blanks for bacteria containing 4.5 ml of sterile saline were prepared in advance. Tenfold serial



dilutions were made of each peritoneal exudate (0.5 ml exudate + 4.5 ml saline) yielding a concentration series from 10^0 (undiluted peritoneal exudate) through 10^{-7} . For enumeration of total bacterial counts, the 10^{-6} and 10^{-7} dilutions were plated on tryptone yeast extract agar, 3 plates/sample, 0.2 ml sample/plate. Each sample was spread over the surface of the plate using a bent glass rod immersed in 95% ethanol and flamed just prior to use. In plating for the total mutant counts on minimal agar, the 10^0 dilution was used, 0.2 ml being plated on each of 5 plates. The plating procedure was identical to that followed for the tryptone yeast extract agar plates. All plates were incubated at 37°C, tryptone yeast extract agar plates for 18 hours and minimal agar plates for 40 hours. For yeast mitotic recombination, dilution blanks containing 4.5 ml of sterile saline were prepared in advance. Tenfold serial dilutions were made of each sample yielding a series from 10^0 to 10^{-5} . Samples of 0.1 ml of the 10^{-5} , 10^{-4} , and 10^{-3} dilutions were removed and plated on complete medium (10 plates each). All plates were incubated at 30°C for 40 hours. The 10^{-5} dilutions were used to determine total populations and the 10^{-4} and 10^{-3} plates were examined after an additional 40 hours at 4°C for red sectors indicating a mutation. Bacterial scoring was calculated as follows:

Total mutants on 5 plates x appropriate exponent =
CFU/ml (CFU is Colony Forming Units) of sample plated
CFU/ml x one/dilution factor ($10^0 - 10^{-7}$) = CFU/ml in undiluted exudate. The mutation frequency (MF) calculated for each sample was:

$$MF = \frac{\text{total mutant cells}}{\text{total population}}$$

$$MF_t/MF_c = \frac{\text{MF of experimental sample}}{\text{MF of control sample}}$$

($MF_t/MF_c = 1.00$ for control sample)



Yeast mitotic recombinants (presumptive ade 2, his 8 homozygotes) were seen as red colonies or as red sectors on a normally white yeast colony. The plates (from 10^{-4} and 10^{-3} dilutions) were scanned under the 10X lens of a dissecting scope to enumerate the red colonies and sectors. Population determinations were made from the 10^{-5} dilution plates. A recombinant frequency (RF) was calculated:

$$RF = \frac{\text{total recombinants counted}}{\text{total number colonies screened}}$$

b. Subacute study

Similar groups of animals at each dose level received five oral doses of the test compound 24 hours apart. Within 30 minutes after the last dosing, the animals were inoculated with the test organism and handled in the same fashion as those in the acute study.

c. In vitro study

Cultures of S. typhimurium histidine auxotrophs (G-46 and TA-1530) were plated on appropriate media. The test compound was then added to the plate, either in the form of a microdrop of solution (0.01 to 0.25 ml) applied to a small filter paper disc resting on the agar or a small crystal applied directly to the agar. Tenfold serial dilutions of the culture were employed and plated so as not to miss the optimum cell density for mutant growth. Mutant colonies were observed and scored. Strain D-3 Saccharomyces cells at proper dilutions were shaken with the test compound, diluted, and plated at 50% survival level or above (see HMA Supplementary Materials and Methods). Red sectors were then scored and the frequency calculated after suitable incubation. Negative and positive controls were run concurrently. The positive control was EMS for Salmonella and Saccharomyces. The in vitro Salmonella tests were reported



as (+) or (-) or questionable; the in vitro Saccharomyces tests were reported as sample concentrations, percent survival, and recombinants/ 10^5 survivors. For the Saccharomyces a 50% survival level, e.g., an arbitrary 5.0% w/v test level, was used when no LD₅₀ was determinable.

2. Cytogenetic Studies

a. In vivo study

Ten to twelve week old, male, albino rats obtained from a closed colony (random-bred) were used. A total of 59 animals in the acute study and 18 animals in the subacute study was used, as illustrated in the following protocol.

Number of Animals Used

Acute Study

Treatment	Time Killed After Administration		
	6 Hours	24 Hours	48 Hours
High Level	5	5	5
Intermediate Level	5	5	5
Low Level	5	5	5
Positive Control	0	0	5
Negative Control	3	3	3

Subacute Study

Five doses 24 hours apart; animals killed 6 hours after last dose.

Treatment	Killed After Administration
High Level	5
Intermediate Level	5
Low Level	5
Negative Control	3

All animals were dosed by gastric intubation.

Four hours after the last compound administration, and two hours prior to killing, each animal was given 4 mg/kg of colcemid intra-



peritoneally in order to arrest the bone marrow cells in C-mitosis. Animals were killed by using CO₂, and the adhering muscle and epiphysis of one femur were removed. The marrow "plug" was removed with a tuberculin syringe and an 18 gauge needle, aspirated into 5 ml of Hanks' balanced salt solution (BSS) in a test tube and capped. The specimens were centrifuged at 1,500 RPM in a table-top centrifuge for 5 minutes, decanted, and 2 ml of hypotonic 0.5% KCl solution was added with gentle agitation to resuspended the cells. The specimens were then placed in a 37°C water bath for 20 minutes in order to swell the cells. Following centrifugation for 5 minutes at 1,500 RPM, the supernatant was decanted and 2 ml of fixative (3:1 absolute methanol:glacial acetic acid) was added. The cells were resuspended in the fixative with gentle agitation, capped, and placed at 4°C for 30 minutes. The specimens were again centrifuged, decanted, 2 ml of prepared fixative was added, and the cells were resuspended and placed at 4°C overnight.

The following day the specimens were again centrifuged, decanted and 0.3 - 0.6 ml of freshly prepared fixative was added to obtain a suitable density. The cells were resuspended and 2 - 3 drops of the suspension were allowed to drop onto a clean, dry slide held at 15° from the horizontal. As the suspension flowed to the edge of the slide, it was ignited by an alcohol burner and allowed to flame. Following ignition, the slides were allowed to dry at room temperature overnight. Duplicate slides were prepared. The slides were stained using a 5% Giemsa solution (Giemsa buffer pH 7.2) for 20 minutes, rinsed in acetone, 1:1 acetone:xylene, and placed in fresh xylene for 30 minutes. The slides were then mounted using Permount (Fisher Scientific) and 24 x 50 mm coverglasses. The coverglasses were selected to be 0.17 mm \pm 0.005 mm in thickness by use of a coverglass micrometer. The preparations

were examined using Leitz Ortholux I & II microscopes with brightfield optics and xenon light sources. These specimens were scanned with 10X and 24X objectives and suitable metaphase spreads that were countable were then examined critically using 40X, 63X or 100X oil immersion flatfield apochromatic objectives. Oculars were either 12X or 16X widefield periplanatics and the tube magnification either 1X or 1.25X. The filters used were either a didymium (BG20) or a Schott IL570 m μ interference filter.

The chromosomes of each cell were counted and only diploid cells were analyzed. They were scored for chromatid gaps and breaks, chromosome gaps and breaks, reunions, cells with greater than ten aberrations, polyploidy, pulverization, and any other chromosomal aberrations which were observed. They were recorded on the currently used forms and expressed as percentages on the summary sheets. Fifty metaphase spreads were scored per animal. Mitotic indices were obtained by counting at least 500 cells and the ratio of the number of cells in mitosis/the number of cells observed was expressed as the mitotic index.

Positive controls in the acute study consisted of animals which had been given the known mutagen Triethylene Melamine (TEM) administered intraperitoneally at a level of 0.30 mg/kg. Negative controls on the acute and subacute studies consisted of the vehicle in which the compound was administered. The dosage levels, solvents and toxicity data are included in the Results and Discussion Section of the report.

b. In vitro study

Human embryonic lung cultures (WI-38) which were negative for adventitious agents (viruses, mycoplasma) which may interfere

were used. These cells were employed at passage level 19. The cells had been transferred using 0.025% trypsin and planted in 32 oz. prescription bottles containing 40 ml of tissue culture medium. When growth was approximately 95% confluent the cells were removed from the glass using trypsin, centrifuged, and frozen in tissue culture medium containing dimethyl sulfoxide (DMSO). Cells were frozen in vials in the vapor phase of liquid nitrogen at a concentration of 2×10^6 cells/ml. When needed, the vials were removed from liquid nitrogen, quick-thawed in a 37°C water bath, washed free of DMSO, suspended in tissue culture medium (minimal essential medium [MEM] plus 1% glutamine, 200 units/ml of penicillin and 200 µg/ml of streptomycin and 15% fetal calf serum) and planted in milk dilution bottles at a concentration of 5×10^5 cells/ml. The test compound was added at three dose levels using three bottles for each level, 24 hours after planting. The dose levels required a preliminary determination of a tissue culture toxicity. This was accomplished by adding logarithmic doses of the compound in saline to a series of tubes containing 5×10^5 cells/ml which were almost confluent. The cells were examined at 24, 48, and 72 hours. Any cytopathic effect (CPE) or inhibition of mitoses was scored as toxicity. Five more closely spaced dose levels were employed within the two logarithmic dosages, the higher of which showed toxicity and the lower no effect. The solvents used and the range finding data are presented in the toxicity data report under Results and Discussion. The dose level below the lowest toxic level was employed as the high level. Logarithmic dose levels were employed for the medium and low levels.

Cells were incubated at 37°C and examined twice daily to determine when an adequate number of mitoses were present. Cells were harvested by shaking when sufficient mitoses were observed, usually 24 - 48

hours after planting, centrifuged, and fixed in absolute methanol:glacial acetic acid (3:1) for 30 minutes.

The specimens were centrifuged, decanted, and suspended in acetic acid-orcein stain (2.0%) and a drop of suspension placed on a clean dry slide. Selected coverglasses 0.17 mm in thickness were placed on the suspension and the excess stain gently expressed from the slide. The coverglasses were sealed with clear nail polish and examined immediately.

The microscopes, objectives, oculars, filters and light sources were enumerated under the metaphase description. Positive controls used were TEM (at a concentration of 0.1 mcg/ml dissolved in saline) and negative controls which consisted of the vehicle in which the test compound was dissolved, which was 0.85% saline. Data were reported on forms currently used and expressed as percentages on the anaphase summary sheets.

3. Dominant Lethal Assay

In this test, male and female random bred rats from a closed colony were employed. These animals were 10-12 weeks old at the time of use. Ten male rats were assigned to each of 5 groups; 3 dose levels selected as described above, a positive control (triethylene melamine) (TEM) and a negative control (solvent only). The positive control was administered intraperitoneally. Administration of the test compound was orally by intubation in both the acute study (1 dose) and in the subacute study (1 dose per day for 5 days). Following treatment, the males were sequentially mated to 2 females per week for 8 weeks (7 weeks in the subacute study). Two virgin female rats were housed with a male for 5 days (Monday through Friday). These two females were removed and housed in a cage until killed. The male was rested on Saturday and Sunday and two new females introduced to the cage on



Monday. It has been our experience that conception has taken place in more than 90% of the females by Friday and that the two day rest is beneficial to the male as regards subsequent weekly matings. Females were killed using CO₂ at 14 days after separating from the male, and at necropsy the uterus was examined for deciduomata (early deaths), late fetal deaths and total implantations.

Sufficient animals were provided in our experimental design to accommodate for any reduction in the number of conceptions. Each male was mated with two females per week, and this provided for an adequate number of implantations per group per week (200 minimum) for negative controls, even if there was a fourfold reduction in fertility of implantations. Results were analyzed according to the statistical procedures described in Supplementary Materials and Methods. Corpora lutea, early fetal deaths, late fetal deaths and total implantations per uterine horn were recorded on the raw data sheets, which are submitted separately.

D. Supplementary Materials and Methods

1. Host-Mediated Assay In Vitro and Formulae

a. Bacterial in vitro plate tests

This method has been published by Ames: The Detection of Chemical Mutagens with Enteric Bacteria, in Chemical Mutagens; Principles and Methods for Their Detection, Vol. 1, Chapter 9, pp. 267-282, A. Hollaender, Editor, Plenum Press, New York (1971).

b. In vitro for mitotic recombination

(1) Strain D-3 was grown to stationary phase on complete medium agar plates at 30°C (3-4 days). Cells were rinsed from the plates and washed twice in saline and cell concentration determined spectro-

photometrically. (A standard curve previously determined for colony forming units versus % transmittance at 545 mu was easily used.)

(2) Cells from the concentration suspension were diluted appropriately into 0.067 M Phosphate buffer pH 7.2 to provide 5×10^7 cells/ml in a total of 25 ml.

(3) The test chemical was first tested for 4 hours at 30°C, with shaking, at concentrations which permitted determination of the 50% survival level. Then, if not included in the first experiment, the compound was tested again only at the 50% survival level. If 50% survival level could not be determined, the arbitrary test level of 5% w/v was used.

(4) Following treatment, cells were diluted and plated on complete agar medium for determination of total population and red sectors. Total surviving population was conveniently measured on plates of 10^{-4} and 10^{-5} dilutions using 0.2 ml per plate (5 plates), and sectors determined on plates of 10^{-3} and 10^{-4} dilutions using 0.2 ml per plate (5 plates). Plates were incubated for 2 days at 30°C followed by a holding period of 2 days at 4°C to promote color development with limited enlargement of the colonies. Red sectors were scored by systematically scanning the plates with a dissecting microscope at 10X magnification.

(5) The frequency of red sectors can then be calculated and may be expressed conveniently as sectors per 10^5 survivors for comparison with untreated controls.

(6) Ethyl Methane Sulfonate (EMS) was employed as the positive control in both in vitro systems.

c. Minimal medium (bacteria):

Spizizen's Minimal Medium:



4X Salt Solution:

$(\text{NH}_4)_2\text{SO}_4$	8.0 gm
K_2HPO_4	56.0 gm
KH_2PO_4	24.0 gm
Na Citrate	4.0 gm
Mg SO_4	0.8 gm
Biotin	0.004 gm
H_2O	qs to 1 liter Sterilize by autoclaving (121°C/15 min.)

Medium:

4X Salt Solution	:250 ml
5.0% Glucose (sterile)	:100 ml (If histidine is added at concentration of 30 mg/liter, this becomes a complete bacterial medium.)
1.5% Bacto-agar (sterile)	:650 ml

d. Complete medium (bacteria):

Bacto-Tryptone	1.0 gm
Yeast-Extract	0.5 gm
Bacto-Agar	2.0 gm
Distilled H_2O	100.0 ml

Sterilize by autoclaving (121°C for 15 minutes).

e. Complete medium (yeast):

KH_2PO_4	1.5 gm
MgSO_4	0.5 gm
$(\text{NH}_4)_2\text{SO}_4$	4.5 gm



Peptone	3.5 gm
Yeast-Extract	5.0 gm
Glucose	20.0 gm
Agar	20.0 gm
Distilled H ₂ O	1000.0 ml

Sterilize by autoclaving (121°C for 15 minutes).

2. Cytogenetics In Vitro Preparation of Anaphase Chromosomes (from Nichols, 1970)

"Anaphase preparations may be made by several methods. One convenient approach is to grow cells directly on coverslips in petri dishes. With human fibroblasts 400,000 cells added to a 22 x 44 mm coverslip in a 50 mm petri dish grown in a 5% CO₂ atmosphere in air has proved very satisfactory. When adequate numbers of mitoses are visualized directly utilizing an inverted microscope (usually 48 to 92 hours after planting) the coverslip is transferred to absolute ethanol for 15 minutes for fixation. They are then stained with any one of a number of suitable stains (Fuelgen, May-Grunwald-Giemse, orcein) and attached to a slide with mounting media for evaluation. Anaphase preparations may also be prepared on cells grown in suspension or cells from a monolayer that have been put into suspension. In this instance the cells are centrifuged and fixed with the squash fixative. They are then suspended in the stain and a drop of the suspension put on the slide and covered with a coverslip. However, in this case, only the excess stain is gently expressed from under the coverslip and no squashing is carried out. In anaphase preparations no pretreatment with colchicine or hypotonic expansion is used and no technique for spreading the cells is used, so that the spindle and normal relationships of the chromosomes are not disturbed."



3. Statistical Analyses of Dominant Lethal Studies

The following statistical analyses were employed as a means of analyzing the results of the dominant lethal studies.

a. The fertility index

The number of pregnant females/number of mated females with the chi-square was used to compare each treatment to the control. Armitage's trend was used for linear proportions to test whether the fertility index was linearly related to arithmetic or log dose.

b. Total number of implantations

The t-test was used to determine significant differences between average number of implantations per pregnant female for each treatment compared to the control. Regression techniques were used to determine whether the average number of implantations per female was related to the arithmetic or log dose.

c. Total number of corpora lutea

The t-test was used to determine significant differences between average number of corpora lutea per pregnant female for each treatment compared to the control.

d. Preimplantation losses

Preimplantation losses were computed for each female by subtracting the number of implantations from the number of corpora lutea. Freeman-Tukey transformation was used on the preimplantation losses for each female and then the t-test was used to compare each treatment to control. Regression technique was used to determine whether the average number of preimplantation losses per female was related to the arithmetic or log dose.



e. Dead implants

Dead implants were treated the same as pre-implantation losses.

f. One or more dead implants

The proportion of females with one or more dead implants was computed, each treatment compared to control by chi-square test and Armitage's trend used for linear proportions to see if proportions were linearly related to either arithmetic or log dose. Also, probit regression analysis was used to determine whether the probit of the proportions was related to log dose.

g. Two or more dead implants

The proportion of females with two or more dead implants computed was treated same as above (f).

h. Dead implants per total implants

Dead implants per total implants were computed for each female and used Freeman-Tukey arc-sine transformation on data for each female; then used t-test to compare each treatment to control.

Historical control data was compiled on a continuous basis as studies were completed. In addition to comparing each treatment to control, as outlined above, each treatment was compared to a historical control.

In order to take variation between males into account, a nested model was used. An analysis of across weeks is also provided.

In addition to these tests, the distribution forms of the various parameters were tested in order to evaluate the appropriateness of some of the tests being used. Certain correlations between parameters may exist and were examined as one step to determine the appropriateness of models. If necessary, alternate test methods were implemented.



The results are presented in tabular form with the addition of historical control information. In addition to these tables, a written report of all findings is provided. As information became available from the on-going investigation of these data, it was reported and suggestions included for changes to the methods of analysis. The statistical reports give the level of significance using both a one-tailed and two-tailed test. Finally, a summary sheet for each study is provided.



MODEL

$$y_{ijk} = \mu + \alpha_i + c_{ij} + e_{ijk}$$

$i = 1, 2$ Group

$j = 1, 2, \dots, 10$ Males within each group

$k = 1, 2$ Females within Males within Groups

ASSUMPTIONS:

$$\alpha_1 + \alpha_2 = 0, \quad c_{ij} \sim \text{nid}(0, \sigma_c^2)$$

$$e_{ijk} \sim \text{nid}(0, \sigma^2)$$

Males are randomly drawn from infinite population

S.V.	d.f.	S.S.	MS	E(MS)	F
TOTAL	39	$\sum \sum \sum (y_{ijk} - \bar{y}_{...})^2$			
GROUPS	1	$20 \sum (\bar{y}_{i..} - \bar{y}_{...})^2$	S_1^2	$\sigma^2 + 20\sigma_c^2 + 20\sigma_e^2$	
MALES WITHIN GROUPS	18	$2 \sum \sum (\bar{y}_{ij.} - \bar{y}_{i..})^2$	S_2^2	$\sigma^2 + 2\sigma_c^2$	
REMAINDER	20	$\sum \sum \sum (y_{ijk} - \bar{y}_{ij.})^2$	S_3^2	σ^2	

E. References

1. Host-Mediated Assay

- a. Gabridge, M.G., Denunzio, A. and Legator, M.S.:
Nature, 221:68, 1969.
- b. Gabridge, M.G., Denunzio, A. and Legator, M.S.:
Science, 163:689, 1969.
- c. Gabridge, M.G. and Legator, M.S.: Proc. Soc.
Exptl. Biol. Med., 130:831, 1969.
- d. Gabridge, M.G., Oswald, E.J. and Legator, M.S.:
Mut. Res., 7:117, 1969.
- e. Legator, M.S. and Malling, H.V.: In, Environmental
Chemical Mutagens, A. Hollaender (Ed.), Plenum
Publishing Corp., New York, in press.

2. Cytogenetics

- a. Nichols, V.W.: Personal communication.
- b. Legator, M.S.: In, Laboratory Diagnosis of Diseases
Caused by Toxic Agents, F. W. Sunderman and F. W.
Sunderman (Ed.), Warren H. Green, Inc., St. Louis,
pp. 17-22, 1970.
- c. Hsu, T.C. and Patton, J.L.: Technical Addendum in,
Comparative Mammalian Cytogenetics, K. Benirschke
(Ed.), Springer-Verlag, New York, pp. 454-460, 1969.
- d. Legator, M.S. et al.: Cytogenetic studies in rats
of cyclohexylamine, a metabolite of cyclamate.
Science, 165:1139, 1969.



3. Dominant Lethal

- a. Bateman, A.J.: Genet. Res. Comb., 1:381, 1960.
- b. Bateman, A.J.: Nature, 210:205, 1966.
- c. Ehling, U.H., Cumming, R.B. and Malling, H.V.:
Mut. Res., 5:417, 1968.
- d. Epstein, S.S. and Shafner, H.: Nature, 219:
385, 1968.



F. Abbreviations

1. mu = micron
2. mcg = ug = microgram
3. g = gram
4. kg = kilogram
5. ml = milliliter
6. rpm = revolutions per minute
7. °C = degrees centigrade
8. pH = power of the hydrogen ion concentration to the base 10
9. M = molar solution
10. conc. = concentration
11. MTD = maximum tolerated dosage = High = LD₅ if determined
or else exceedingly high dose, such as 5 g/kg
12. INT = intermediate = medium level
13. USE = usage level if known = low level
14. BSS = balanced salt solution
15. C-metaphase = cells arrested in metaphase, using colchine
or colcemid
16. LD₅₀ = that dosage which produced 50% mortality in the
group of animals treated
17. LD₅ = that dosage which produced 5% mortality in the group
of animals treated
18. NC = negative control
19. PC = positive control
20. AU = acute usage level (low level)
21. AI = acute intermediate level (medium level)
22. AMTD = acute maximum tolerated dose level (LD₅ level,
high level)



23. SAU = subacute usage level (low level)
24. SAI = subacute intermediate level (medium level)
25. SA LD₅ = subacute LD₅ level (MTD level, high level)
26. CO₂ = carbon dioxide
27. DMN = Dimethyl nitrosamine
28. EMS = Ethyl methane sulfonate
29. TEM = Triethylene melamine
30. DMSO = Dimethyl sulfoxide
31. MEM = minimal essential medium (Eagle's)
32. CPE = cytopathic effect
33. his = histidine marker
34. D-3 = mitotic recombinant strain of Saccharomyces
35. mf = mean mutant frequency
36. MFt/MFc = mean mutant frequency of the test compound group
compared to mean mutant frequency of the negative control
group
37. CFU = colony forming units
38. WI-38 = code name for a strain of human embryonic lung
tissue culture cells
39. Rec x 10⁵ = mitotic recombinants x 10⁵
40. Mean B/A = mean frequency
41. tot. scr. = total scored
42. tot. = total
43. χ^2 = a test of variation in the data from the computed
regression line - tested in these studies at the 5% level
44. Aber. = aberrations
45. Frag. = fragment
46. HMA = host-mediated assay