



# BIONETICS

Mutagenic Evaluation of Compound FDA 75-92  
9/77

Potassium Sulfate

Final Report

MUTAGENICITY EVALUATION

OF

FDA 75-92  
POTASSIUM SULFATE

FINAL REPORT

G24

5516 Nicholson Lane  
Kensington, Maryland  
20795

MUTAGENICITY EVALUATION

OF

FDA 75-92  
POTASSIUM SULFATE

FINAL REPORT

SUBMITTED TO

GENETIC TOXICOLOGY BRANCH  
DIVISION OF TOXICOLOGY  
BUREAU OF FOODS  
U.S. FOOD AND DRUG ADMINISTRATION  
200 C STREET, S.W., ROOM 1066  
WASHINGTON, D.C.

SUBMITTED BY

LITTON BIONETICS, INC.  
5516 NICHOLSON LANE  
KENSINGTON, MARYLAND 20795

LBI PROJECT NO. 2672

SEPTEMBER 1977



BIONETICS

## TABLE OF CONTENTS

	Page No.
EVALUATION SUMMARY.....	1
I. <u>OBJECTIVE</u> .....	2
II. <u>MATERIALS</u> .....	2
A.   Test Compound.....	2
B.   Indicator Microorganisms.....	2
C.   Reaction Mixture.....	2
D.   Tissue Homogenates and Supernatants.....	3
E.   Positive Control Compounds.....	3
III. <u>METHODS</u> .....	3
A.   Toxicity.....	3
B.   Plate Tests.....	4
C.   Suspension Tests.....	4
D.   Preparation of Tissue Homogenates and 9,000 x g Cell Fractions.....	5
E.   Data Recording and Reporting.....	5
IV. <u>RESULTS SECTION</u>	
A.   Solubility Properties of the Test Compound.....	6
B.   Toxicity and Dosage Determinations for the Test Compound.....	6
C.   Plate Assay Results.....	7
D.   Suspension Assay Results.....	7
V. <u>INTERPRETATION OF RESULTS AND CONCLUSIONS</u> .....	15
VI. <u>EXPLANATION OF EVALUATION PROCEDURES FOR PLATE ASSAYS</u> ....	16
VII. <u>EXPLANATION OF EVALUATION PROCEDURES FOR SUSPENSION ASSAYS</u>	18
APPENDIX - Tabulation of Data.....	A-1



### EVALUATION SUMMARY

The test compound, FDA 75-92, Potassium Sulfate, did not exhibit mutagenic activity in any of the assays employed in these studies.



BIONETICS

DATE: July, 1977

SPONSOR: U.S. Food and Drug Administration

SUBJECT: Evaluation of Test Compound: FDA 75-92, Potassium Sulfate

I. OBJECTIVE

The objective of this study was to evaluate the test compound for genetic activity in microbial assays with and without the addition of mammalian metabolic activation preparations.

II. MATERIALS

A. Test Compound

1. Date Received: December 29, 1976

2. Description: White powder

B. Indicator Microorganisms

The following strains of indicator microorganisms were used in the evaluation:

Yeast Strain: Saccharomyces cerevisiae, strain D4

Bacteria Strains: Salmonella typhimurium, strains TA-1535  
TA-1537  
TA-1538  
TA-98  
TA-100

C. Reaction Mixture

The following reaction mixture was employed in the activation tests:

<u>Component</u>	<u>Final Concentration/ml</u>
1. TPN (sodium salt)	4 $\mu$ moles
2. Glucose-6-phosphate	5 $\mu$ moles
3. Sodium phosphate (dibasic)	100 $\mu$ moles
4. $MgCl_2$	8 $\mu$ moles
5. KCl	33 $\mu$ moles
6. Homogenate fraction equivalent to 25 mg of wet tissue.	

#### D. Tissue Homogenates and Supernatants

The tissue homogenates and 9,000 x g supernatants were prepared from tissues of the following mammalian species: Mouse - ICR random bred adult males; rat - Sprague-Dawley adult males; and monkey - Macaca mulatta adult males.

#### E. Positive Control Compounds

Table 1 lists chemicals for positive controls in the direct and activation assays.

TABLE 1  
POSITIVE CONTROLS USED IN DIRECT AND ACTIVATION ASSAYS

<u>Assay</u>	<u>Chemical</u> <sup>a</sup>	<u>Solvent</u>	<u>Probable Mutagenic Specificity</u>
Nonactivation	Methylnitrosoguanidine	Water or saline	BPS <sup>b</sup>
	Ethylmethanesulfonate	Water or saline	BPS <sup>b</sup>
	2-Nitrofluorene	Dimethylsulfoxide <sup>c</sup>	FS <sup>b</sup>
	Quinacrine mustard	Water or saline	FS <sup>b</sup>
Activation	Dimethylnitrosamine	Water or saline	BPS <sup>b</sup>
	2-Acetylaminofluorene	Dimethylsulfoxide <sup>c</sup>	FS <sup>b</sup>
	8-Aminoquinoline	Dimethylsulfoxide <sup>c</sup>	FS <sup>b</sup>
	2-Aminoanthracene	Dimethylsulfoxide <sup>c</sup>	BPS <sup>b</sup>

<sup>a</sup> Concentrations given in the Results Section

<sup>b</sup> BPS = base-pair substitution; FS = frameshift

<sup>c</sup> Previously shown to be non-mutagenic

### III. METHODS

#### A. Toxicity

The solubility, toxicity and doses for the test chemical were determined prior to screening.

The test chemical was tested for toxicity against specific indicator strains over a range of doses to determine the 50% survival dose. Bacteria were tested in phosphate buffer, pH 7.4, for one hour at 37°C on a shaker. Yeasts were tested in phosphate buffer, pH 7.4, for four hours at 30°C on a shaker. The 50% survival concentrations and the 1/4 and 1/2 50% doses calculated.

If no toxicity was obtained for the chemical with a given strain, then a maximum dose of 5% (w/v) was used.

Unless otherwise specified, the doses calculated for the tests in buffer were applied to the activation tests. The solubility of the test chemical under treatment conditions is stated in the Results Section.

## B. Plate Tests (Overlay Method)

Approximately  $10^8$  cells from an overnight culture of each indicator strain were added to test tubes containing 2.0 ml of molten agar supplemented with biotin and a trace of histidine. For nonactivation tests, the three dose levels of the test compound were added to the contents of the appropriate tubes and poured over the surfaces of selective agar plates. In activation tests 0.5 ml of a 9,000 x g tissue supernatant and required cofactors (core reaction mixture) were added to the overlay tubes. Three dose levels of the test chemical were added to the appropriate tubes, which were then mixed and the contents poured over the surface of a minimal agar (selective medium) plate and allowed to solidify. The plates were incubated for 48 to 72 hours at 37°C, and scored for the number of colonies growing on each plate. The concentrations of all chemicals are given in the Results Section. Positive and solvent controls using positive compounds that are active directly and those that require metabolic activation were run with each assay.

## C. Suspension Tests

### 1. Nonactivation

Bacteria and yeast cultures of the indicator organisms were grown in complete broth, washed and resuspended in 0.9% saline to densities of  $1 \times 10^{10}$  cells/ml and  $5 \times 10^9$  cells/ml, respectively. This constituted the working stock for tests of a group of test chemicals and their respective controls. Tests were conducted in plastic, 24-well tissue culture plates (Linbro). Cells plus appropriate volume(s) of the test chemical were added to the wells to give a final volume of 1.5 ml. The solvent replaced the test chemical in the negative controls. Treatment was at 30°C for four hours for yeast tests and at 37°C for one hour for bacterial tests. All flasks were shaken during treatment. Following treatment, the plates were set on ice. Aliquots of cells were removed, diluted in sterile saline (4°C) and plated on the appropriate complete media. Undiluted samples from flasks containing the bacteria were plated on minimal selective medium in reversion experiments. Samples from a  $10^{-1}$  dilution of treated cells were plated on the selected media for enumeration of gene conversion with strain D4. Bacterial plates were scored after incubation for 48 hours at 37°C. The yeast plates were incubated at 30°C for 3-5 days before scoring.

### 2. Activation

Bacteria and yeast cells were grown and prepared as described in the nonactivation tests. Measured amounts of the test and control chemicals plus 0.25 ml of the stock-cell suspension were added to wells of the Linbro plate containing the appropriate tissue fraction and reaction mixture. All flasks (bacteria and yeast) were incubated at 37°C with shaking. The treatment times as well as the dilutions, plating procedures and scoring of the plates were the same as described for nonactivation tests.

D. Preparation of Tissue Homogenates and 9,000 x g Cell Fractions

Male animals (except monkeys) sufficient to provide the necessary quantities of tissues were killed by cranial blow, decapitated and bled. Monkey tissues were obtained from freshly killed and bled male rhesus monkeys. Organs were immediately dissected from the animals using aseptic techniques and placed in ice-cold 0.15M KCl. Upon collection of the desired quantity of organs, they were washed twice with fresh KCl and completely homogenized with a motor-driven homogenizing unit at 4°C. The whole organ homogenate obtained from this step was divided into two samples. One sample was frozen at -80°C and the other was centrifuged for 20 minutes at 9,000 x g in a refrigerated centrifuge. The supernatant from the centrifuged sample was retained and frozen at -80°C. These two frozen samples were used for the activation studies. Protein and P-448 determinations were made for each lot of homogenate.

E. Data Recording and Reporting

1. Plate test assays

The numbers of colonies on each plate were counted and recorded on printed forms. These raw data were entered into a computer program designed to print out all data by test. The data are presented as revertants per plate for each indicator strain employed in the assay. The positive and solvent controls are provided as reference points.

2. Suspension assays

Following the specified incubation periods all population plates were scored by an automatic colony counter and the results from each plate of a set were recorded, in ink, on data processing forms. All minimal or other types of selective media plates were hand scored and the results recorded along with the respective population data. Other relevant experimental data were recorded on experimental definition forms. For bacteria strains the number of colonies recorded from either the population or selective plates represents that number in 1 ml of test suspension plated. The numbers recorded for the yeast strain D4 represent the number in 0.5 ml of test suspension plated. The data were then processed and printed from a computer program. All raw data sheets are dated and signed by the responsible technician.





IV. RESULTS SECTION

A. Solubility Properties of the Test Compound

1. Name or code designation of the test compound: FDA 75-92, Potassium Sulfate
2. Test solvent: \* Saline
3. Solubility of the test compound under treatment conditions: Soluble
4. Additional comments: White crystal

B. Toxicity and Dosage Determinations for the Test Compound

1. Test date for toxicity determination: April 4, 1977
2. The 50% survival level was determined for bacteria and yeast indicator organisms by conducting survival curves with the test compound at the following concentrations:

Percent Concentration (w/v or v/v)

5.0  
0.5  
0.05  
0.005  
0.0005

3. Concentrations of the test compound used in the mutagenicity tests:

<u>Test Doses</u>	<u>Percent Concentration</u>	
	<u>Bacteria</u>	<u>Yeast</u>
1/4 50% Survival	0.04	1.25
1/2 50% Survival	0.08	2.50
50% Survival	0.16	5.00

\*The concentration of solvent was equal to the highest volume of test material added.

C. Plate Test Results

The plate test results are summarized in the following table. The values presented in this table are the number of revertants per plate.

D. Suspension Assay Results

The suspension test results for the test compound are summarized in the tables following the plate test summary. The values presented in these tables are the calculated mutation frequencies for each control and experimental test point. The first table of the suspension set presents the results for the nonactivation assays, and the second table through the fourth table of the suspension set presents the results for the activation assays. A listing of computer codes and abbreviations is included for reference. Tabulation of all raw data is provided in the Appendix.

## SUMMARY OF TEST RESULTS

## PLATE TESTS

A. NAME OR CODE DESIGNATION OF THE TEST COMPOUND: 007778805

B. TEST DATE: MAY 18, 1977

TEST	SPECIES	ISSUE	B E V E R T A N I S P E R P L A T E									
			TA-1535		TA-1537		TA-1538		TA-98		TA-100	
			1	2	1	2	1	2	1	2	1	2
1. NON-ACTIVATION												
SOLVENT CONTROL*	---	---	28	21	22	35	17	16	34	27	148	143
POSITIVE CONTROL**	---	---	>1000	>1000	>1000	>1000	>1000	>1000	>1000	>1000	>1000	>1000
TEST 0.16000 %	---	---	12	11	14	13	17	18	29	22	114	112
0.08000 %	---	---	18	26	14	10	19	19	21	19	160	118
0.04000 %	---	---	16	16	11	17	12	17	30	37	133	149
2. ACTIVATION												
SOLVENT CONTROL*	MOUSE	LIVER	30	31	22	23	19	10	37	32	222	195
	RAT	LIVER	26	37	20	18	19	17	39	40	147	182
	MONKEY	LIVER	18	15	17	31	23	21	36	37	192	133
POSITIVE CONTROL***	MOUSE	LIVER	502	490	260	256	874	911	>1000	>1000	>1000	>1000
	RAT	LIVER	274	374	241	149	938	732	>1000	>1000	>1000	>1000
	MONKEY	LIVER	370	215	173	160	738	901	>1000	937	>1000	>1000
TEST 0.16000 %	MOUSE	LIVER	19	26	18	21	10	19	22	30	103	87
0.08000 %	MOUSE	LIVER	25	19	18	23	19	16	41	22	117	125
0.04000 %	MOUSE	LIVER	24	21	20	22	12	11	31	35	143	136
0.16000 %	RAT	LIVER	31	22	15	16	15	17	39	37	123	130
0.08000 %	RAT	LIVER	22	24	14	18	14	11	28	34	126	138
0.04000 %	RAT	LIVER	31	26	11	17	11	10	26	38	129	118
0.16000 %	MONKEY	LIVER	25	26	19	12	19	18	32	38	131	133
0.08000 %	MONKEY	LIVER	22	27	20	19	19	10	37	36	136	141
0.04000 %	MONKEY	LIVER	14	20	10	11	10	11	29	24	146	136

\* NON-ACTIVATION ASSAYS CONSIST OF THE CELLS PLUS THE TEST COMPOUND VEHICLE (SOLVENT). FOR ACTIVATION ASSAYS, THE OVERLAY CONTAINS THE ACTIVATION SYSTEM PLUS THE TEST COMPOUND VEHICLE.

** TA-1535	MNNG	2 UG/PLATE	*** TA-1535	ANTH	100 UG/PLATE
TA-1537	QM	20 UG/PLATE	TA-1537	AMQ	100 UG/PLATE
TA-1538	NF	100 UG/PLATE	TA-1538	AAF	100 UG/PLATE
TA-98	NF	100 UG/PLATE	TA-98	AAF	100 UG/PLATE
TA-100	MNNG	2 UG/PLATE	TA-100	ANTH	100 UG/PLATE

NOTE: CONCENTRATIONS ARE GIVEN IN MICROLITERS(UL) OR MICROGRAMS(UG) PER PLATE.

LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
REPORT EXR34

COMPOUND FREQUENCY SUMMARY REPORT 07/22/77

NONACTIVATION COMPOUND 007778805

TEST	ORG	TA100 HIS EX-8	TA1535 HIS EX-8	TA1537 HIS EX-8	TA1538 HIS EX-8	TA98 HIS EX-8	0000D4 ADE EX-5	0000D4 TRY EX-5	
NAN		85.71	3.58	11.59	5.46	13.83	19.91	6.89	CONTROLS
NAP		900.49	685.45	235.32	163.30	71.70	109.36	78.06	
<hr/>									
NA1		53.72	3.07	2.23	2.94	10.58	21.62	4.63	TEST DATA
NA2		64.36	3.13	5.71	3.98	9.85	14.88	5.14	
NA3		55.15	4.25	2.91	2.96	7.02	22.36	7.23	

LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
REPORT EXR34

COMPOUND FREQUENCY SUMMARY REPORT 07/22/77

SPECIES ICRFLO/MOUSE COMPOUND 007778805

TEST	ORG	TA100 HIS EX-8	TA1535 HIS EX-8	TA1537 HIS EX-8	TA1538 HIS EX-8	TA98 HIS EX-8	0000D4 ADE EX-5	0000D4 TRY EX-5	
ACT	A+C	84.43	7.09	5.62	9.43	9.46	6.93	15.63	NEGATIVE CONTROLS
ACT	A-C	55.84	4.18	5.48	8.48	8.79	6.90	22.67	
ACT	ALI	57.41	5.86	8.83	5.51	21.48	10.38	8.14	
ACT	ALU	58.69	8.09	8.92	7.42	14.70	5.15	15.00	
<hr/>									
ACT	PLI	182.07	78.46	97.15	149.24	87.28	53.36	91.82	POSITIVE CONTROLS
ACT	PLU	91.94	9.12	11.54	76.16	74.60	19.04	17.34	
<hr/>									
ACT	L11	31.09	4.73	3.54	2.78	25.04	9.54	10.37	TEST COMPOUND
ACT	L12	29.49	4.61	7.27	5.47	17.83	5.62	9.58	
ACT	L13	31.18	7.63	5.04	5.68	17.59	12.72	10.57	
ACT	LU1	84.48	5.80	8.63	4.60	18.85	11.25	9.20	
ACT	LU2	62.47	7.04	6.91	5.86	18.10	15.49	12.21	
ACT	LU3	70.41	5.34	6.63	5.00	10.85	12.27	11.83	

LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
REPORT EXR34

COMPOUND FREQUENCY SUMMARY REPORT 07/22/77

SPECIES SPRDAW/RAT

COMPOUND 007778805

TEST	ORG	TA100 HIS EX-8	TA1535 HIS EX-8	TA1537 HIS EX-8	TA1538 HIS EX-8	TA98 HIS EX-8	000004 ADE EX-5	000004 TRY EX-5
------	-----	----------------------	-----------------------	-----------------------	-----------------------	---------------------	-----------------------	-----------------------

ACT	A+C	20.79	6.55	31.96	10.56	13.11	15.60	10.18
ACT	A-C	18.64	40.87	3.14	5.06	15.56	12.85	7.44
ACT	ALI	83.48	12.32	12.17	10.82	36.48	15.01	10.34
ACT	ALI	64.68	4.05	6.98	10.33	17.83	15.67	10.27

NEGATIVE CONTROLS

---

ACT	PL1	247.96	196.42	55.58	88.95	297.10	110.32	71.61
ACT	PLU	52.63	7.32	17.21	160.05	124.87	17.35	7.69

---

POSITIVE CONTROLS

---

ACT	L11	91.00	1.96	3.06	6.48	26.44	16.60	16.28
ACT	L12	29.20	3.14	3.45	6.42	20.14	15.03	10.71
ACT	L13	39.62	5.26	3.13	4.38	25.22	19.53	15.35
ACT	LU1	52.46	4.55	9.16	7.94	11.02	16.68	8.84
ACT	LU2	44.52	6.06	9.20	7.76	13.61	12.53	10.97
ACT	LU3	51.08	6.01	6.62	14.63	8.94	18.42	9.16

TEST COMPOUND

LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
REPORT EXR34

COMPOUND FREQUENCY SUMMARY REPORT 07/22/77

SPECIES RHESUS/MONKEY COMPOUND 007778805

TEST	ORG	TA100 HIS EX-8	TA1535 HIS EX-8	TA1537 HIS EX-8	TA1538 HIS EX-8	TA98 HIS EX-8	000004 ADE EX-5	000004 TRY EX-5	
ACT	A+C	87.84	5.64	21.88	13.86	7.51	21.98	10.90	NEGATIVE CONTROLS
ACT	A-C	64.54	2.96	1.25	11.47	8.19	8.04	7.30	
ACT	ALI	80.09	5.54	2.09	12.37	14.55	16.25	6.21	
ACT	ALU	70.23	3.78	2.60	10.10	10.45	20.74	7.73	
<hr/>									
ACT	PLI	284.60	56.90	37.74	154.01	201.65	68.25	54.23	POSITIVE CONTROLS
ACT	PLU	69.20	4.38	13.64	7.74	8.22	20.46	7.34	
<hr/>									
ACT	LI1	77.80	4.37	4.73	16.72	13.18	9.98	4.21	TEST COMPOUND
ACT	LI2	77.87	6.80	5.66	20.16	12.33	7.73	3.83	
ACT	LI3	84.37	3.75	3.21	14.89	10.29	16.26	10.71	
ACT	LU1	77.49	5.17	3.60	8.52	9.53	7.21	5.17	
ACT	LU2	84.05	5.35	5.58	13.81	10.81	6.10	2.68	
ACT	LU3	84.23	6.16	7.19	9.06	10.22	14.26	4.21	

## DATA TABLE TERMS AND ABBREVIATIONS

ABBREVIATION OR TERM	DEFINITION OR EXPLANATION
COMPOUND	Client designated compound number appears in this column.
TEST CODES	<p> NAN = Nonactivation: Solvent Control  NAP = Nonactivation: Positive Control  NA1 = Nonactivation: Test Compound Dose 1  NA2, etc. = Reflects the other dose level(s) </p> <p> A+C = Negative Chemical Control for ACP  A-C = Activation: Solvent Control  ALI or A+T = Activation: Homogenate Control (Liver)  ALU = Activation: Homogenate Control (Lung)  ACP = Activation: Positive Control  ACT = Activation Test </p> <p> LI = Liver Tissue Activation Fraction  LU = Lung Tissue Activation Fraction  KI = Kidney Tissue Activation Fraction  TE = Testes Tissue Activation Fraction  1,2, etc. = Dose Levels </p>
CONCENTRATION	<p>All test compound dose levels are expressed as a whole number followed by an exponent (negative) identified by the appropriate units.</p> <p>Example: 0025-2PCT = 0.25 percent concentration</p>
POPU	Total number of viable cells in the plating sample raised to some exponent printed directly below the abbreviation (i.e., EP + 6 = $\times 10^6$ ).
MUT 1	Total number of mutants or convertants obtained from the sample plated raised to some exponent printed directly below the abbreviation (i.e., EP + 0 = $10^0$ ). For strain D4, MUT 1 represents the number of ADE+ convertants.
MUT 2	Only used for strain D4 and represents the number of TRY+ convertants in the plated sample.
FREQ 1	The calculated mutation or gene conversion frequency times the negative exponent written directly below. For strain D4, FREQ 1 represents the ADE+ value.
FREQ 2	Only used for strain D4 and represents the TRY+ conversion frequency.
CONTAM	Presence of contamination on any plates.



# DATA TABLE TERMS AND ABBREVIATIONS (continued)

ABBREVIATION OR TERM	DEFINITION OR EXPLANATION
AAF	2-Acetylaminofluorene
DMSO	Dimethylsulfoxide
DMN	Dimethylnitrosamine
EMS	Ethylmethanesulfonate
QM	Quinacrine Mustard
NF	Nitrofluorene
ANTH	2-Amino Anthracene
AMQ	8-Amino Quinoline
SPECIES	Animal Strains
SPRDAW	Sprague Dawley Rats
ICRFLO	Flow ICR Random Bred Mice
RHESUS	Rhesus Monkey ( <u>Macaca mulatta</u> )
MIXEDB	Dog, Mixed Breed
NEWZEA	New Zealand White Rabbit
UG	Microgram
UM	Micromole
ADE	Adenine
TRY	Tryptophan

V. INTERPRETATION OF RESULTS AND CONCLUSIONS

The test compound, FDA 75-92, Potassium Sulfate, was evaluated for genetic activity in a series of in vitro microbial assays with and without metabolic activation. The following results were obtained:

A. Salmonella typhimurium

1. Plate tests

The results of these tests were negative.

2. Nonactivation suspension tests

The results of these tests were negative.

3. Activation suspension tests

The results of these tests were negative.

B. Saccharomyces cerevisiae

1. Nonactivation suspension tests

The results of these tests were negative.

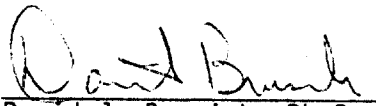
2. Activation suspension tests

The results of these tests were negative.

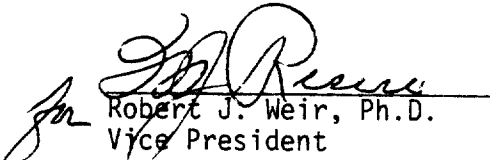
C. Conclusions

The test compound, FDA 75-92, Potassium Sulfate, did not exhibit mutagenic activity in any of the assays employed in these studies.

Submitted by:

  
David J. Brusick, Ph.D. 3/7/77  
Director Date  
Department of Molecular  
Toxicology

Reviewed by:

  
Robert J. Weir, Ph.D. 8/3/77  
Vice President Date

## VI. EXPLANATION OF EVALUATION PROCEDURES FOR PLATE ASSAYS

Plate test data consist of direct revertant colony counts obtained from a set of selective agar plates seeded with populations of mutant cells suspended in a semisolid overlay. Because the test chemical and cells are incubated in the overlay for 2-3 days, and a few cell divisions occur during the incubation period, the test is semiquantitative in nature. Although these features of the assay reduce the quantitation of results, they provide certain advantages not contained in a quantitative suspension test.

- The small number of cell divisions permits potential mutagens to act on replicating DNA which is often more sensitive than non-replicating DNA.
- The combined incubation of the compound and the cells in the overlay permit constant exposure of the indicator cells for 2-3 days.

### A. Surviving Populations

Plate test procedures do not permit exact quantitation of the number of cells surviving chemical treatment. At low concentrations of the test chemical, the surviving population on the treatment plates is essentially the same as the negative control plate. At high concentrations, the surviving population is usually reduced by some fraction. Our protocol normally employs dose levels that are selected such that the highest dose will show slight toxicity (as determined by subjective criteria) and several doses ranging down 1 to 2 logs lower.

### B. Dose Response Phenomena

The demonstration of dose-related increases in mutant counts is an important criterion in establishing mutagenicity. Factors which may modify dose response results for a mutagen would be the selection of doses that are too low (usually mutagenicity and toxicity are related). If the highest dose is far lower than a toxic concentration, no increases may be observed over the dose range selected. Conversely, if the lowest dose employed is highly cytotoxic, the test chemical may kill any mutants that are induced and the compound will not appear to be mutagenic.

### C. Control Tests

Positive and negative control assays are conducted with each experiment and consist of direct acting mutagens for nonactivation assays and mutagens that require metabolic biotransformation in activation assays. Negative controls consist of the test compound solvent in the overlay agar with the other essential components. The negative control plate for each strain gives a reference point to which the test data are compared. The positive control assay is conducted to demonstrate that the test systems are functional with known mutagens.

#### D. Evaluation Criteria for Ames Assay

Because the procedures used to evaluate the mutagenicity of the test chemical are semiquantitative, the criteria used to determine positive effects are inherently subjective and are based primarily on a historical data base. Most data sets are evaluated using the following criteria:

##### 1. Strains TA-1535, TA-1537, and TA-1538

If the solvent control value is within the normal range, a chemical that produces a positive dose response over three concentrations with the lowest increase equal to twice the solvent control value is considered to be mutagenic.

##### 2. Strains TA-98, TA-100, and D4

If the solvent control value is within the normal range, a chemical that produces a positive dose response over three concentrations with the highest increase equal to twice the solvent control value for TA-100 and two to three times the solvent control value for strains TA-98 and D4 is considered to be mutagenic. For these strains, the dose response increase should start at approximately the solvent control value.

##### 3. Pattern

Because TA-1535 and TA-100 were both derived from the same parental strain (G-46) and because TA-1538 and TA-98 were both derived from the same parental strain (D3052), there is a built-in redundancy in the microbial assay. In general the two strains of a set respond to the same mutagen and such a pattern is sought. It is also anticipated that if a given strain, e.g. TA-1537, responds to a mutagen in nonactivation tests it will generally do so in activation tests. (The converse of this relationship is not expected.) While similar response patterns are not required for all mutagens, they can be used to enhance the reliability of an evaluation decision.

##### 4. Reproducibility

If a chemical produces a response in a single test that cannot be reproduced in one or more additional runs, the initial positive test data loses significance.

The preceding criteria are not absolute and other extenuating factors may enter into a final evaluation decision. However, these criteria are applied to the majority of situations and are presented to aid those individuals not familiar with this procedure. As the data base is increased, the criteria for evaluation can be more firmly established.

## VII. EXPLANATION OF EVALUATION PROCEDURES FOR SUSPENSION ASSAYS

Data obtained from mutagenicity tests are evaluated on a test by test basis followed by an examination of the total response pattern using all the data. To facilitate this type of evaluation, we have prepared two separate formats in which data are processed. The first is the Compound Summary Backup Detail Sheet, which details the essential raw data from each experiment showing surviving population counts, total mutant or revertant counts, as well as, calculated mutation frequencies. This format permits close examination of each set of test data. The following considerations are part of any assessment.

### A. Surviving Population Counts

A certain level of chemically-induced toxicity is anticipated, but occasionally isolated tests or groups of tests show very low (<25%) survival compared to the tissue controls. Such isolated decreases may result from improper dilution procedures or defective growth media and decrease confidence in the calculated mutation frequencies especially if the total mutant counts appear unaffected. Data of this type are generally unacceptable and these experiments are routinely repeated at a lower dose level to reduce killing and increase confidence in the nature of the response.

### B. Total Mutant Counts

For nonmutagens, the mutant/surviving population ratio should be roughly equivalent for each test point in a given experiment. If the cell number drops in response to killing, the mutant number should decrease proportionately. A mutagenic chemical, however, will produce an altered mutant/surviving population ratio. Mutant numbers as well as calculated frequencies are compared to the negative control data. In certain instances, the mutant frequencies will increase with little or no change in the absolute number of mutants especially where the test chemical is toxic. Data of this type, although not necessarily aberrant, or even rare, must be viewed with special care to ensure that the increased frequencies were not the result of selective toxicity of the test chemical for the his cells. This phenomenon, referred to as selection, can lead to erroneous conclusions. Thus we attempt to keep the surviving population of cells high and look for positive responses that show increases in both numbers of mutants and mutation frequencies. Again, occasional isolated fluctuations in mutant counts are found that can be attributed to improper pipetting or media contamination. These fluctuations are usually easy to identify by inspection of the other data points in the experiment which will be negative.

#### C. Dose Response Phenomena

Dose-related increases in mutants and mutation frequencies are the most convincing data to have in assessing mutagenic activity of chemicals. In some cases, however, dose-related increases are not observed for mutagens. This depends considerably on the dose levels selected. The figure on the following page illustrates how one might obtain various types of dose-related responses by a mutagen based solely on dose selection. It also emphasizes the need to keep dose levels within a relatively low range of toxicity so that data are consistently on the uphill side of the hypothetical curve.

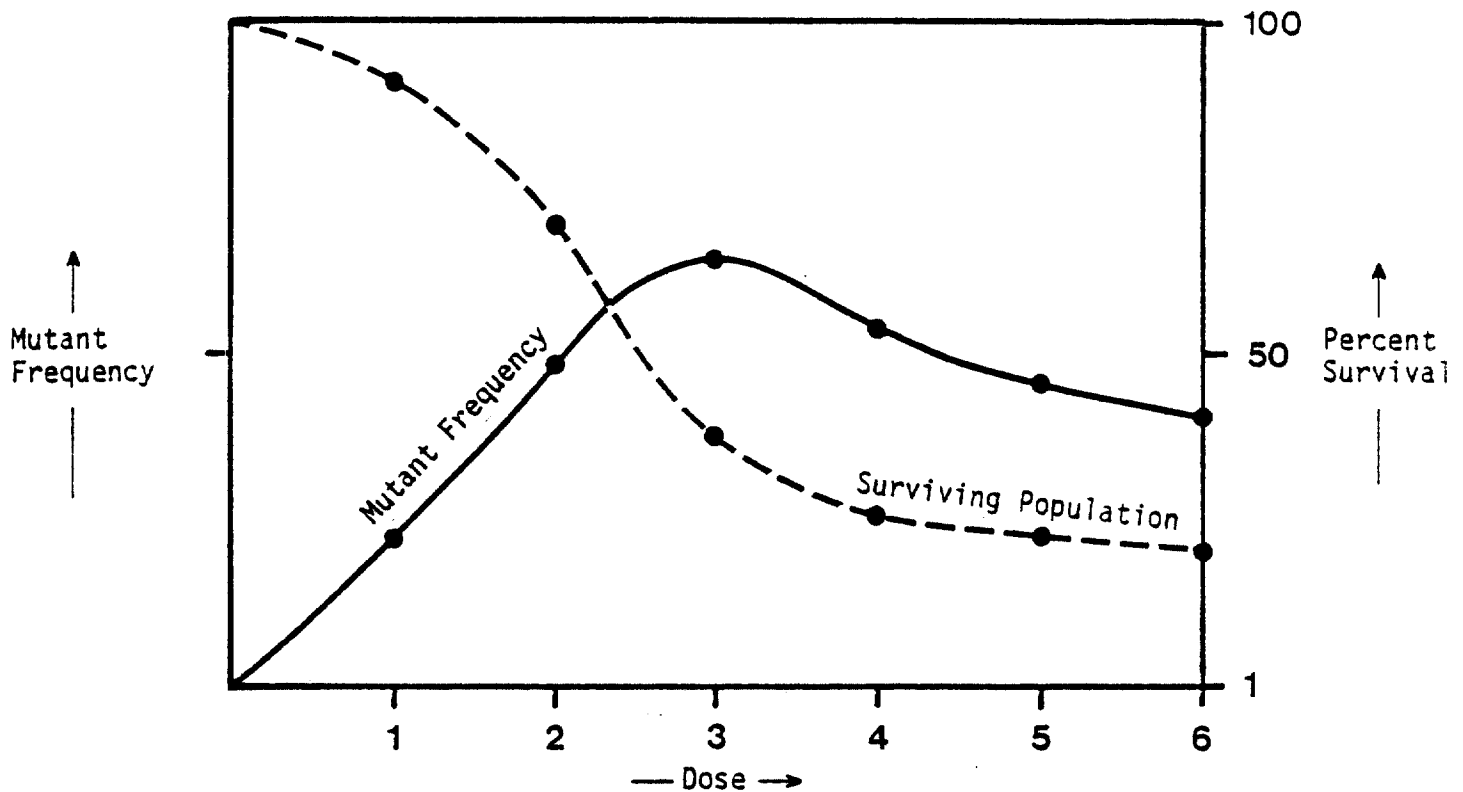
#### D. Control Tests

Positive and negative control tests are conducted with each experiment and consist of direct acting positive agents for nonactivation assays and chemicals that require metabolic transformation for activation assays. In nonactivation assays, the NAN control contain the test chemical solvent plus cells, but no chemical, and is used as a reference to assess the level of response obtained in the various tests. It is not possible at this time to put precise cut-off points where negative responses become positive responses. A statistical component for our computer program is under development and will be included when available. Positive controls are only used as relative reference points and to demonstrate that the system is functioning with known mutagens. In activation assays, three types of negative controls are run: (1) A solvent control minus the chemical and minus the activation system (A-C); (2) a control plus the positive control chemical minus the activation system (A+C); and (3) a control containing the activation system and the test chemical solvent (ALI or ALU). All three controls are used collectively to assess the level of response in the various activation tests. A chemical may appear positive when compared to an A-C control but not when compared to an A+T control. The value of each of the above controls with respect to their weight in evaluation is  $ALI \text{ or } ALU > A-C > A+C$ .

The other data format is the Compound Frequency Summary Report sheet in which all the calculated frequencies obtained for a given compound are displayed in a table. This format permits an overview of all data. The points form a matrix of information that should present a consistent pattern. Nonmutagens should produce a matrix with data frequencies clustered around the negative control values. Occasional random high or low fluctuations are not uncommon and seldom indicate true genetic activity. Mutagenic chemicals should, on the other hand, produce a set of consistent responses that demonstrate a logical pattern. The patterns depend on the mutagenic specificity of the chemical but can be easily recognized in the Compound Frequency Summary Report format.

These mutagenicity assays are designed to optimize the probability of recognizing mutagens from nonmutagens and, in most cases, they work well. Occasionally, the data points are such that a definitive conclusion cannot be made without additional data.

# HYPOTHETICAL MUTATION AND TOXICITY KINETICS



## HYPOTHETICAL EXPERIMENT

- (1) Dose levels 1, 2 & 3 were used
- (2) Dose levels 2, 3 & 4 were used
- (3) Dose levels 3, 4 & 5 were used

## OBSERVED DOSE RESPONSE

A typical positive dose response set of data would be obtained.

The intermediate dose level shows a higher mutation frequency than both the low dose and the high dose.

Here an inverted dose response would be observed with the highest dose level showing the lowest response.

APPENDIX  
Tabulation of Data



REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

EXPERIMENT	CONTRACT	223-76-2102			PROJECT	2672	
710253	DETECTOR	TA100	SPECIES	/			DATE - 07/22/77
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
		NAN	SOLVENT	0252	0216	85.71	0
		NAP	EMS 0.066%	0616	5547	900.49	0
007778805	NA1		0016-2 PCT.	0698	0375	53.72	0
007778805	NA2		0008-2 PCT.	0707	0455	64.36	0
007778805	NA3		0004-2 PCT.	0767	0423	55.15	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672		DATE - 07/22/77			
EXPERIMENT 710251		DETECTOR TA1535		SPECIES /			
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
		NAN	SOLVENT	0810	0029	3.58	0
		NAP	EMS 0.2%	0852	5840	685.45	0
007778805	NA1		0016-2 PCT.	0911	0028	3.07	0
007778805	NA2		0008-2 PCT.	0830	0026	3.13	0
007778805	NA3		0004-2 PCT.	0894	0038	4.25	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672		DATE - 07/22/77			
EXPERIMENT 710151	DETECTOR TA1537	SPECIES	/				
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	NAN		SOLVENT	0466	0054	11.59	0
	NAP		QM 13 UG/ML	0235	0553	235.32	0
007778805	NA1		0016-2 PCT.	1393	0031	2.23	0
007778805	NA2		0008-2 PCT.	0595	0034	5.71	0
007778805	NA3		0004-2 PCT.	1339	0039	2.91	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672		DATE - 07/22/77			
EXPERIMENT 705306	DETECTOR TA1538	SPECIES					
COMPOUND	TEST	ORG ID	CONCENTRATION	POPUS EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	NAN		SOLVENT	0403	0022	5.46	0
	NAP		NF 667 UG/ML	0376	0614	163.30	0
007778805	NA1		0016-2 PCT.	0544	0016	2.94	0
007778805	NA2		0008-2 PCT.	0503	0020	3.98	0
007778805	NA3		0004-2 PCT.	0507	0015	2.96	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
 COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672		DATE - 07/22/77			
EXPERIMENT 710252	DETECTOR TA98	SPECIES /					
COMPOUND	TEST	ORG ID	CONCENTRATION	POP1 EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	NAN		SOLVENT	0962	0133	13.83	0
	NAP		NF 667 UG/ML	0834	0598	71.70	0
007778805	NA1		0016-2 PCT.	1389	0147	10.58	0
007778805	NA2		0008-2 PCT.	1472	0145	9.85	0
007778805	NA3		0004-2 PCT.	1268	0089	7.02	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672							
EXPERIMENT 710952		DETECTOR 000004		SPECIES /		DATE - 07/22/77			
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+4	MUT1 EP+1	MUT2 EP+1	FREQ1 EP-5	FREQ2 EP-5	CONTAM
		NAN	SOLVENT	1175	0234	0081	19.91	6.89	1
		NAP	EMS 1.0 %	1463	1600	1142	109.36	78.06	0
007778805	NA1		0.5-0 PCT.	1619	0350	0075	21.62	4.63	0
007778805	NA2		0.25-1 PCT.	1693	0252	0087	14.88	5.14	0
007778805	NA3		0.125-2 PCT.	1176	0263	0085	22.36	7.23	1

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 710351	DETECTOR TA100	SPECIES ICRFLO/MOUSE	DATE - 07/22/77				
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		DMN 90 UM/ML	0989	0835	84.43	0
	A-C		SOLVENT	0899	0502	55.84	0
	ALI		TISSUE	1545	0887	57.41	0
	ALU		TISSUE	1869	1097	58.69	0
	ACP	LI	DMN 90 UM/ML	1160	2112	182.07	0
	ACP	LU	DMN 90 UM/ML	0496	0456	91.94	0
007778805	ACT	LI1	0016-2 PCT.	1518	0472	31.09	0
007778805	ACT	LI2	0008-2 PCT.	1577	0465	29.49	0
007778805	ACT	LI3	0004-2 PCT.	1411	0440	31.18	0
007778805	ACT	LU1	0016-2 PCT.	0612	0517	84.48	0
007778805	ACT	LU2	0008-2 PCT.	1860	1162	62.47	0
007778805	ACT	LU3	0004-2 PCT.	1355	0954	70.41	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 710353 DETECTOR TA1535 SPECIES ICRFLO/MOUSE DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		DMN 90 UM/ML	0649	0046	7.09	0
	A-C		SOLVENT	0431	0018	4.18	0
	ALI		TISSUE	0273	0016	5.86	0
	ALU		TISSUE	0346	0028	8.09	0
	ACP	LI	DMN 90 UM/ML	0687	0539	78.46	0
	ACP	LU	DMN 90 UM/ML	0570	0052	9.12	0
007778805	ACT	LI1	0016-2 PCT.	0465	0022	4.73	0
007778805	ACT	LI2	0008-2 PCT.	0456	0021	4.61	0
007778805	ACT	LI3	0004-2 PCT.	0367	0028	7.63	0
007778805	ACT	LU1	0016-2 PCT.	0379	0022	5.80	0
007778805	ACT	LU2	0008-2 PCT.	0355	0025	7.04	0
007778805	ACT	LU3	0004-2 PCT.	0637	0034	5.34	0



REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 710454 DETECTOR TA1537 SPECIES ICRFLO/MOUSE DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		AMQ 333 UG/ML	1334	0075	5.62	0
	A-C		SOLVENT	1680	0092	5.48	0
	ALI		TISSUE	0804	0071	8.83	0
	ALU		TISSUE	1020	0091	8.92	0
	ACP	LI	AMQ 333 UG/ML	1297	1260	97.15	0
	ACP	LU	AMQ 333 UG/ML	1144	0132	11.54	0
007778805	ACT	LI1	0016-2 PCT.	2345	0083	3.54	0
007778805	ACT	LI2	0008-2 PCT.	1514	0110	7.27	0
007778805	ACT	LI3	0004-2 PCT.	2025	0102	5.04	0
007778805	ACT	LU1	0016-2 PCT.	0997	0086	8.63	0
007778805	ACT	LU2	0008-2 PCT.	1173	0081	6.91	0
007778805	ACT	LU3	0004-2 PCT.	1131	0075	6.63	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 711251 DETECTOR TA1538 SPECIES ICRFLO/MOUSE DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-B	CONTAM
	A+C		ANTH 67 UG/ML	1219	0115	9.43	0
	A-C		SOLVENT	1226	0104	8.48	1
	ALI		TISSUE	1143	0063	5.51	0
	ALU		TISSUE	1145	0085	7.42	1
	ACP	LI	ANTH 67 UG/ML	1117	1667	149.24	1
	ACP	LU	ANTH 67 UG/ML	0302	0230	76.16	1
007778805	ACT	LI1	0016-2 PCT.	1187	0033	2.78	1
007778805	ACT	LI2	0008-2 PCT.	1096	0060	5.47	0
007778805	ACT	LI3	0004-2 PCT.	1127	0064	5.68	0
007778805	ACT	LU1	0016-2 PCT.	1087	0050	4.60	0
007778805	ACT	LU2	0008-2 PCT.	0990	0058	5.86	1
007778805	ACT	LU3	0004-2 PCT.	1620	0081	5.00	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 710354	DETECTOR TA98	SPECIES ICRFLO/MOUSE	DATE - 07/22/77				
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		ANTH 67 UG/ML	1797	0170	9.46	0
	A-C		SOLVENT	1513	0133	8.79	0
	ALI		TISSUE	0810	0174	21.48	0
	ALU		TISSUE	1095	0161	14.70	0
	ACP	LI	ANTH 67 UG/ML	0629	0549	87.28	0
	ACP	LU	ANTH 67 UG/ML	1134	0846	74.60	0
007778805	ACT	LI1	0016-2 PCT.	0703	0176	25.04	0
007778805	ACT	LI2	0008-2 PCT.	0819	0146	17.83	0
007778805	ACT	LI3	0004-2 PCT.	0921	0162	17.59	0
007778805	ACT	LU1	0016-2 PCT.	0854	0161	18.85	0
007778805	ACT	LU2	0008-2 PCT.	0917	0166	18.10	0
007778805	ACT	LU3	0004-2 PCT.	1318	0143	10.85	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 712952 DETECTOR 000004 SPECIES ICRFLO/MOUSE DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+4	MUT1 EP+1	MUT2 EP+1	FREQ1 EP-5	FREQ2 EP-5	CONTAM
	A+C		DMN 90 UM/ML	1184	0082	0185	6.93	15.63	0
	A-C		SOLVENT	1072	0074	0243	6.90	22.67	0
	ALI		TISSUE	1474	0153	0120	10.38	8.14	0
	ALU		TISSUE	1320	0068	0198	5.15	15.00	0
	ACP	LI	DMN 90 UM/ML	1222	0652	1122	53.36	91.82	0
	ACP	LU	DMN 90 UM/ML	1061	0202	0184	19.04	17.34	0
007778805	ACT	LI1	0005-0 PCT.	1582	0151	0164	9.54	10.37	0
007778805	ACT	LI2	0025-1 PCT.	1941	0109	0186	5.62	9.58	0
007778805	ACT	LI3	0125-2 PCT.	1580	0201	0167	12.72	10.57	0
007778805	ACT	LU1	0005-0 PCT.	1707	0192	0157	11.25	9.20	0
007778805	ACT	LU2	0025-1 PCT.	1827	0283	0223	15.49	12.21	0
007778805	ACT	LU3	0125-2 PCT.	1597	0196	0189	12.27	11.83	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 711053	DETECTOR TA100	SPECIES SPRDAW/RAT		DATE - 07/22/77			
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		DMN 90 UM/ML	0404	0084	20.79	0
	A-C		SOLVENT	0499	0093	18.64	0
	ALI		TISSUE	0230	0192	83.48	0
	ALU		TISSUE	0487	0315	64.68	0
	ACP	LI	DMN 90 UM/ML	0319	0791	247.96	0
	ACP	LU	DMN 90 UM/ML	0779	0410	52.63	0
007778805	ACT	LI1	0016-2 PCT.	0200	0182	91.00	0
007778805	ACT	LI2	0008-2 PCT.	0387	0113	29.20	0
007778805	ACT	LI3	0004-2 PCT.	0371	0147	39.62	0
007778805	ACT	LU1	0016-2 PCT.	0631	0331	52.46	0
007778805	ACT	LU2	0008-2 PCT.	0766	0341	44.52	0
007778805	ACT	LU3	0004-2 PCT.	0969	0495	51.08	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 711080 DETECTOR TA1535 SPECIES SPRDAW/RAT DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		DMN 90 UM/ML	0229	0015	6.55	0
	A-C		SOLVENT	0208	0085	40.87	0
	ALI		TISSUE	0763	0094	12.32	0
	ALU		TISSUE	0469	0019	4.05	0
	ACP	LI	DMN 90 UM/ML	0447	0878	196.42	0
	ACP	LU	DMN 90 UM/ML	0287	0021	7.32	0
007778805	ACT	LI1	0016-2 PCT.	0612	0012	1.96	0
007778805	ACT	LI2	0008-2 PCT.	0510	0016	3.14	0
007778805	ACT	LI3	0004-2 PCT.	0418	0022	5.26	0
007778805	ACT	LU1	0016-2 PCT.	0308	0014	4.55	0
007778805	ACT	LU2	0008-2 PCT.	0264	0016	6.06	0
007778805	ACT	LU3	0004-2 PCT.	0283	0017	6.01	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 710851 DETECTOR TA1537 SPECIES SPRDAW/RAT DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POP1 EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		AMQ 333 UG/ML	0341	0109	31.96	0
	A-C		SOLVENT	0446	0014	3.14	0
	ALI		TISSUE	0871	0106	12.17	0
	ALU		TISSUE	0573	0040	6.98	0
	ACP	LI	AMQ 333 UG/ML	0851	0473	55.58	0
	ACP	LU	AMQ 333 UG/ML	0552	0095	17.21	0
007778805	ACT	LI1	0016-2 PCT.	0360	0011	3.06	0
007778805	ACT	LI2	0008-2 PCT.	0493	0017	3.45	0
007778805	ACT	LI3	0004-2 PCT.	0831	0026	3.13	0
007778805	ACT	LU1	0016-2 PCT.	0502	0046	9.16	0
007778805	ACT	LU2	0008-2 PCT.	0522	0048	9.20	0
007778805	ACT	LU3	0004-2 PCT.	0604	0040	6.62	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 712551	DETECTOR TA1538	SPECIES SPRDAW/RAT	DATE - 07/22/77				
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		ANTH 67 UG/ML	0956	0101	10.56	0
	A-C		SOLVENT	1244	0063	5.06	2
	ALI		TISSUE	0536	0058	10.82	0
	ALU		TISSUE	0552	0057	10.33	0
	ACP	LI	ANTH 67 UG/ML	0977	0869	88.95	0
	ACP	LU	ANTH 67 UG/ML	0443	0709	160.05	0
007778805	ACT	LI1	0016-2 PCT.	0694	0045	6.48	2
007778805	ACT	LI2	0008-2 PCT.	0701	0045	6.42	2
007778805	ACT	LI3	0004-2 PCT.	0799	0035	4.38	2
007778805	ACT	LU1	0016-2 PCT.	0554	0044	7.94	2
007778805	ACT	LU2	0008-2 PCT.	0670	0052	7.76	2
007778805	ACT	LU3	0004-2 PCT.	0492	0072	14.63	2



REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 715151	DETECTOR TA98	SPECIES SPRDAW/RAT	DATE - 07/22/77				
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		ANTH 67 UG/ML	0915	0120	13.11	0
	A-C		SOLVENT	0405	0063	15.56	0
	ALI		TISSUE	0455	0166	36.48	0
	ALU		TISSUE	1127	0201	17.83	0
	ACP	LI	ANTH 67 UG/ML	0483	1435	297.10	0
	ACP	LU	ANTH 67 UG/ML	0961	1200	124.87	0
007778805	ACT	LI1	0016-2 PCT.	0537	0142	26.44	0
007778805	ACT	LI2	0008-2 PCT.	0879	0177	20.14	0
007778805	ACT	LI3	0004-2 PCT.	0559	0141	25.22	0
007778805	ACT	LU1	0016-2 PCT.	1143	0126	11.02	0
007778805	ACT	LU2	0008-2 PCT.	1234	0168	13.61	0
007778805	ACT	LU3	0004-2 PCT.	1778	0159	8.94	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 712901 DETECTOR 000004 SPECIES SPRDAW/RAT DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+4	MUT1 EP+1	MUT2 EP+1	FREQ1 EP-5	FREQ2 EP-5	CONTAM
	A+C		DMN 90 UM/ML	1051	0164	0107	15.60	10.18	0
	A-C		SOLVENT	1533	0197	0114	12.85	7.44	0
	ALI		TISSUE	1306	0196	0135	15.01	10.34	0
	ALU		TISSUE	1091	0171	0112	15.67	10.27	0
	ACP	LI	DMN 90 UM/ML	1173	1294	0840	110.32	71.61	0
	ACP	LU	DMN 90 UM/ML	1118	0194	0086	17.35	7.69	0
007778805	ACT	L11	0005-0 PCT.	1265	0210	0206	16.60	16.28	0
007778805	ACT	L12	0025-1 PCT.	1317	0198	0141	15.03	10.71	0
007778805	ACT	L13	0125-2 PCT.	1101	0215	0169	19.53	15.35	0
007778805	ACT	LU1	0005-0 PCT.	1211	0202	0107	16.68	8.84	0
007778805	ACT	LU2	0025-1 PCT.	1285	0161	0141	12.53	10.97	0
007778805	ACT	LU3	0125-2 PCT.	1026	0189	0094	18.42	9.16	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 710451	DETECTOR TA100	SPECIES RHESUS/MONKEY	DATE - 07/22/77				
COMPOUND	TEST	ORG ID	CONCENTRATION	PGPU EP+6	MUT1 EP+0	FREQ1 EP-B	CONTAM
	A+C		DMN 90 UM/ML	0839	0737	87.84	0
	A-C		SOLVENT	0863	0557	64.54	0
	ALI		TISSUE	0874	0700	80.09	0
	ALU		TISSUE	0823	0578	70.23	0
	ACP	LI	DMN 90 UM/ML	0617	1756	284.60	0
	ACP	LU	DMN 90 UM/ML	1013	0701	69.20	0
007778805	ACT	LI1	0016-2 PCT.	0991	0771	77.80	0
007778805	ACT	LI2	0008-2 PCT.	1044	0813	77.87	0
007778805	ACT	LI3	0004-2 PCT.	0883	0745	84.37	0
007778805	ACT	LU1	0016-2 PCT.	0924	0716	77.49	0
007778805	ACT	LU2	0008-2 PCT.	0953	0801	84.05	0
007778805	ACT	LU3	0004-2 PCT.	0742	0625	84.23	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 710452 DETECTOR TA1535 SPECIES RHESUS/MONKEY DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		DMN 90 UM/ML	0887	0050	5.64	0
	A-C		SOLVENT	0945	0028	2.96	0
	ALI		TISSUE	0903	0050	5.54	0
	ALU		TISSUE	0846	0032	3.78	0
	ACP	LI	DMN 90 UM/ML	1160	0660	56.90	0
	ACP	LU	DMN 90 UM/ML	1028	0045	4.38	0
007778805	ACT	LI1	0016-2 PCT.	1237	0054	4.37	0
007778805	ACT	LI2	0008-2 PCT.	1000	0068	6.80	0
007778805	ACT	LI3	0004-2 PCT.	1308	0049	3.75	0
007778805	ACT	LU1	0016-2 PCT.	0987	0051	5.17	0
007778805	ACT	LU2	0008-2 PCT.	1289	0069	5.35	0
007778805	ACT	LU3	0004-2 PCT.	0844	0052	6.16	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102 PROJECT 2672  
EXPERIMENT 710961 DETECTOR TA1537 SPECIES RHESUS/MONKEY DATE - 07/22/77

COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		AMQ 333 UG/ML	0352	0077	21.88	0
	A-C		SOLVENT	0641	0008	1.25	0
	ALI		TISSUE	0623	0013	2.09	0
	ALU		TISSUE	1229	0032	2.60	0
	ACP	LI	AMQ 333 UG/ML	1134	0428	37.74	0
	ACP	LU	AMQ 333 UG/ML	1158	0158	13.64	0
007778805	ACT	LI1	0016-2 PCT.	0528	0025	4.73	0
007778805	ACT	LI2	0008-2 PCT.	0760	0043	5.66	0
007778805	ACT	LI3	0004-2 PCT.	0841	0027	3.21	0
007778805	ACT	LU1	0016-2 PCT.	1084	0039	3.60	0
007778805	ACT	LU2	0008-2 PCT.	0735	0041	5.58	0
007778805	ACT	LU3	0004-2 PCT.	0612	0044	7.19	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 712553	DETECTOR TA1538	SPECIES RHESUS/MONKEY	DATE - 07/22/77				
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		ANTH 67 UG/ML	0700	0097	13.86	2
	A-C		SOLVENT	0689	0079	11.47	2
	ALI		TISSUE	0590	0073	12.37	2
	ALU		TISSUE	0802	0081	10.10	2
	ACP	LI	ANTH 67 UG/ML	0798	1229	154.01	2
	ACP	LU	ANTH 67 UG/ML	1046	0081	7.74	2
007778805	ACT	LI1	0016-2 PCT.	0598	0100	16.72	2
007778805	ACT	LI2	0008-2 PCT.	0511	0103	20.16	2
007778805	ACT	LI3	0004-2 PCT.	0591	0088	14.89	2
007778805	ACT	LU1	0016-2 PCT.	0610	0052	8.52	2
007778805	ACT	LU2	0008-2 PCT.	0420	0058	13.81	2
007778805	ACT	LU3	0004-2 PCT.	0872	0079	9.06	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672					
EXPERIMENT 710453	DETECTOR TA98	SPECIES RHESUS/MONKEY		DATE - 07/22/77			
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+6	MUT1 EP+0	FREQ1 EP-8	CONTAM
	A+C		ANTH 67 UG/ML	1852	0139	7.51	0
	A-C		SOLVENT	1660	0136	8.19	0
	ALI		TISSUE	0852	0124	14.55	0
	ALU		TISSUE	1167	0122	10.45	0
	ACP	LI	ANTH 67 UG/ML	1150	2319	201.65	0
	ACP	LU	ANTH 67 UG/ML	1789	0147	8.22	0
007778805	ACT	LI1	0016-2 PCT.	1085	0143	13.18	0
007778805	ACT	LI2	0008-2 PCT.	1192	0147	12.33	0
007778805	ACT	LI3	0004-2 PCT.	1545	0159	10.29	0
007778805	ACT	LU1	0016-2 PCT.	1417	0135	9.53	0
007778805	ACT	LU2	0008-2 PCT.	1416	0153	10.81	0
007778805	ACT	LU3	0004-2 PCT.	1604	0164	10.22	0

REPORT EXR33 LITTON BIONETICS MUTAGENIC ACTIVITY SYSTEM  
COMPOUND SUMMARY BACKUP DETAIL

CONTRACT 223-76-2102		PROJECT 2672							
EXPERIMENT 711302	DETECTOR 000004	SPECIES RHESUS/MONKEY	DATE - 07/22/77						
COMPOUND	TEST	ORG ID	CONCENTRATION	POPU EP+4	MUT1 EP+1	MUT2 EP+1	FREQ1 EP-5	FREQ2 EP-5	CONTAM
	A+C		DMN 90 UM/ML	1110	0244	0121	21.98	10.90	0
	A-C		SOLVENT	1343	0108	0098	8.04	7.30	0
	ALI		TISSUE	1385	0225	0086	16.25	6.21	0
	ALU		TISSUE	1268	0263	0098	20.74	7.73	0
	ACP	LI	DMN 90 UM/ML	1370	0935	0743	68.25	54.23	0
	ACP	LU	DMN 90 UM/ML	1212	0248	0089	20.46	7.34	1
007778805	ACT	LI1	0005-0 PCT.	1543	0154	0065	9.98	4.21	0
007778805	ACT	LI2	0025-1 PCT.	1565	0121	0060	7.73	3.83	0
007778805	ACT	LI3	0125-2 PCT.	1279	0208	0137	16.26	10.71	0
007778805	ACT	LU1	0005-0 PCT.	1179	0085	0061	7.21	5.17	0
007778805	ACT	LU2	0025-1 PCT.	1493	0091	0040	6.10	2.68	0
007778805	ACT	LU3	0125-2 PCT.	1283	0183	0054	14.26	4.21	0